

PRESSCARE G-SPORT³

PRESSOTHERAPY | 3 INFLATION PROGRAMS


DEVICE FOR PNEUMATIC COMPRESSION DEDICATED TO ATHLETES, REALIZED TO ACCELERATE MUSCULAR RECOVERY FROM POST-TRAINING/COMPETITION FATIGUE. LIGHT AND TRANSPORTABLE, EQUIPPED WITH A POWER BANK THAT ENABLES TREATMENTS ALSO RIGHT AFTER COMPETITION, DIRECTLY ON THE COURTSIDE.

TECHNICAL FEATURES

Outputs
2
Programs
3
Pressure
0-180 mmHg ±10%
Timer
up to 99 minutes
Power supply
- Mains, 100-240 V AC, 50-60 Hz - Rechargeable battery (Power Bank)
Size
180x190x135 mm
Weight
3 Kg

EQUIPMENT

· 1 Bag
· 1 Motor unit G-SPORT ³
· 2 Boots size "S" or size "M" or "L"
· 2 Connection cables
· 1 Power supply unit
· 1 Power Bank
· 1 Operating manual

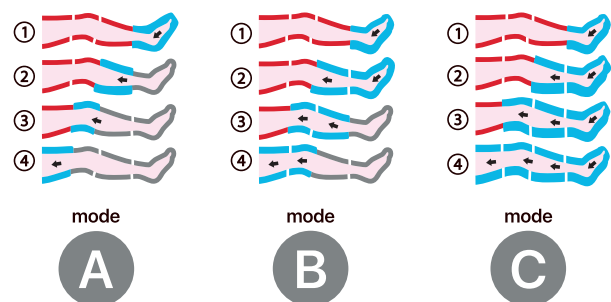

MODE A: the program sequence consists in the inflation and deflation of boot sectors one by one to execute an activating massage, ideal before the competition or training.

MODE B: the program sequence consists in the inflation and deflation of two chambers at a time to avoid fluid reflux downwards and make the treatment effective and very comfortable.

MODE C: the program sequence consists in the sequential inflation of all the chambers and their simultaneous deflation, ideal for a deeper massage.



PRESSCARE G-SPORT³



PRESSOTHERAPY | 3 INFLATION PROGRAMS

SIZE OF ACCESSORIES (4 SECTIONS)


Rev. En_01_2022

BOOTS SIZE "S"	(A) LENGTH GROIN - HEEL 72 cm (B) THIGH CIRCUMFERENCE 60 cm with EXTENSION 70 cm (C) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "M"	(A) LENGTH GROIN - HEEL 82 cm (B) THIGH CIRCUMFERENCE cm 67 con ESTENSIONE 77 cm (C) ANKLE CIRCUMFERENCE cm 37 con ESTENSIONE 42 cm
BOOTS SIZE "L"	(A) LENGTH GROIN - HEEL 92 cm (B) THIGH CIRCUMFERENCE cm 74 with EXTENSION 84 cm (C) ANKLE CIRCUMFERENCE cm 40 with EXTENSION 45 cm

SLEEVE	(A) 84 cm - (B) 63 cm (C) SHOULDER CIRCUMFERENCE 43 cm with EXTENSION 53 cm (D) WRIST CIRCUMFERENCE 31 cm with EXTENSION 36 cm
--------	--

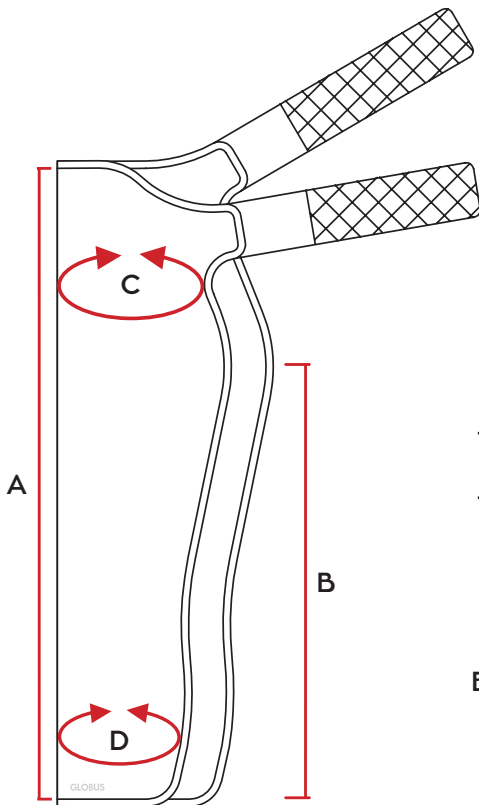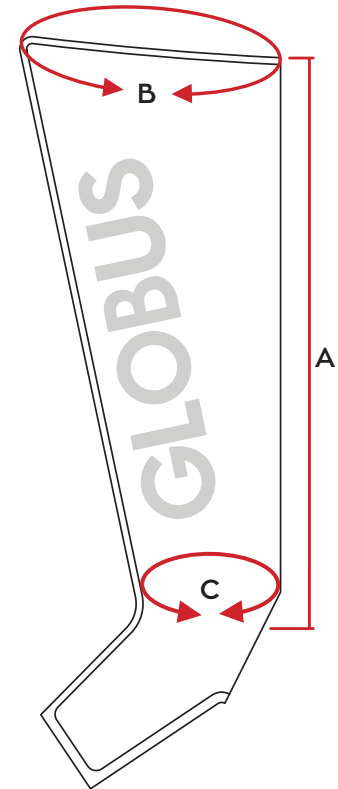

ABDOMINAL BAND	(A) 133 cm - (B) 38 cm ABDOMEN CIRCUMFERENCE from 82 cm to 114 cm
----------------	--

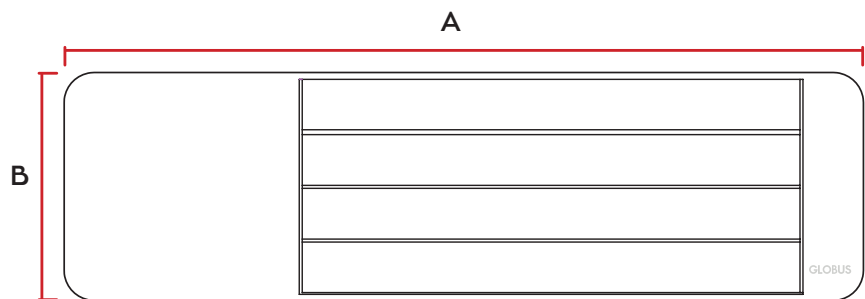

SPLITTER CABLE to use three accessories at the same time	SINGLE CABLE
--	--------------

N.B. All **EXTENSIONS** are to be considered an optional accessory available on request.

The **ACCESSORIES**, depending on the type of configuration, can be optional or included.