

# ACTIVA 700

## 4 CHANNELS ELECTROTHERAPY | 271 PROGRAMS

PRODUCT FOR BEAUTY AND FITNESS, DESIGNED WITH SPECIFIC PROGRAMS FOR FIRING, TONING AND REJUVENATION. THE SPECIAL HANDPIECE G-TRODE (SUPPLIED) MAKES ACTIVA 700 A TOP OF THE RANGE OF THE ELECTROSTIMULATION LINE.

60 BEAUTY | 58 FITNESS | 53 SPORT | 36 SERIAL SEQUENTIAL STIMULATION | 21 SPECIALS SPORT | 13 SKIN-FACE  
12 MICROCURRENTS | 12 TENS | 3 REHAB | 3 INCONTINENCE

### TECHNICAL FEATURES

#### Display

Visible area size 2,6"

#### Channels

4 independent (8 electrodes)

#### Frequency

0,3-150 Hz \*

#### Pulse amplitude

40-450 µs \*

#### Power

0-120 mA per channel

#### Power supply

Rechargeable batteries

#### Size

100x160x35 mm


#### Weight

454 gr

\* According to the type of current

### EQUIPMENT

- 1 Carrying case
- 1 Stimulator ACTIVA 700
- 1 Handpiece G-Trode
- 4 Cables for electrode connection
- 2 Cables for microcurrents and ionophoresis
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Gel flask
- 1 Charger
- 1 Operating manual

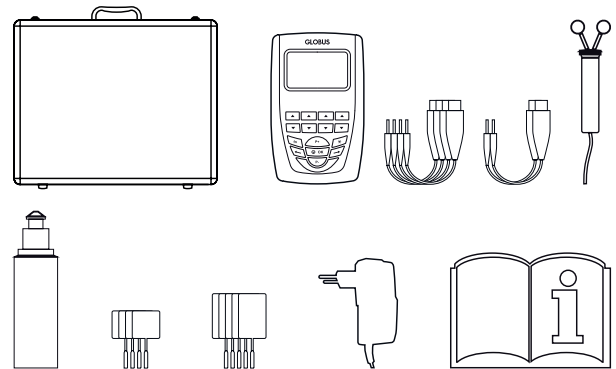

CE  
0476



TENS



MENS



3S



RUN TIME



LAST 10



FAVORITES



MULTI USER



PROGRAMMABLE



STIM LOCK



FUCTION 2+2



MY TRAINER



AUTO STIM



SYNCRO STIM



WORK TIME



RECHARGEABLE BATTERIES

# ACTIVA 700

## 4 CHANNELS ELECTROTHERAPY | 271 PROGRAMS

### PROGRAMS LIST

#### FITNESS

Firming  
Bio-Pulse firming  
Sculpting  
Bio-Pulse sculpting  
Toning  
Mass building  
Body sculpting  
Definition  
Jogging  
Anaerobic fitness  
Aerobic fitness  
Cramp prevention

#### BEAUTY

Drainage  
Bio-Pulse drainage  
Lipolysis  
Post-pregnancy lipolysis  
Bio-Pulse relaxing massage  
Toning massage  
Energizing massage  
Connective massage  
Swollen arms  
Face capillaries  
Definition  
Slim figure  
Skin tone improvement  
Post-pregnancy drainage  
Post-pregnancy firming  
Breast firming  
Breast sculpting

#### SPORT

Capillarization  
Warm-up  
Pre-competition warm-up  
Active recovery  
Maximum strength  
Endurance strength  
Explosive strength  
Reactivity  
Post-competition recovery  
Hypertrophy

#### SPECIAL SPORT

Running  
Skiing  
Golf

#### REHAB

Quadriceps atrophy (with knee prosthesis)  
Recovery after ACL surgery  
Shoulder subluxation prevention

#### TENS and PAIN

Conventional antalgic tens  
Modulated antalgic tens  
Endorphinic tens  
Menstrual pain  
Cronich pain  
Cervical pain  
Muscle pain  
Knee pain  
Scapulohumeral periarthritis  
Chronic lumbago  
Bursitis-tendinitis

#### G-PULSE

Breast microlifting  
Cleavage microlifting  
Face microlifting  
Skin elasticity  
Cellulite  
BioSkin collagen  
Wrinkles  
Eye wrinkles  
Expression wrinkles  
Stretch marks  
Nail strenghtening  
Bio peeling  
Hematoma

#### MICROCURRENTS

Epicondylitis  
Scapulohumeral periarthritis  
Contusion  
Sciatica  
Lumbago  
Acute pain  
Articular pain  
Stiff neck  
Whiplash  
Shoulder sprain  
Knee sprain  
Tendon inflammation

#### INCONTINENCE

Mixed incontinence  
Stress incontinence  
URGE incontinence

#### SERIAL SEQUENTIAL STIMULATION

The "3S" program list includes 36 parameter combinations.  
The "3S" programs are characterized by a delay activation of the channels 3 and 4 compared with the channels 1 and 2.  
The Serial Sequential Stimulation permits to stimulate the musculature in kinetic chain thanks to the differentiated activation times of the muscular groups involved.