

GLOBUS

ITALIAN EXCELLENCE

CATALOG **2022**

GLOBUS
ITALIAN EXCELLENCE

CONTENTS

ABOUT US	04
KINEO INTELLIGENT LOAD	06
REHAB LINE	24
- TECARTHERAPY Line	26
- MAGNETOTHERAPY Line	32
- ELECTROTHERAPY Line	44
- ULTRASOUND THERAPY Line	56
- LASER THERAPY Line for physiotherapy	60
- LASER THERAPY Line for podology	70
- PRESSOTHERAPY Line	76
- PREVENTION	80
BEAUTY & WELLNESS LINE	84
- TECAR BEAUTY Line	86
- ELECTROSTIMULATION Line	94
- RADIOFREQUENCY aesthetic line	106
- ULTRASOUND aesthetic line	112
- VIBRATION PLATFORM aesthetic line	120
SPORT LINE	122
- ELECTROSTIMULATION Line	124
- PRESSOTHERAPY	134
- BALL SHOOTING MACHINES Line	136
ACCESSORIES	144
- TECARTHERAPY AND TECAR BEAUTY	146
- MAGNETOTHERAPY	152
- ELECTROTHERAPY AND ELECTROSTIMULATION	155
- ULTRASOUND THERAPY	159
- LASER THERAPY	160
- RADIOFREQUENCY	161
- PRESSOTHERAPY	162
- BALL SHOOTING MACHINES	164
VETERINARY LINE	166
- TECARTHERAPY	168
- MAGNETOTHERAPY	170
- ELECTROTHERAPY	172
- ULTRASOUND THERAPY	176
- LASER THERAPY	178
- ACCESSORIES	182
KEY STRONG POINTS	186

ABOUT US

GLOBUS has been a world leading company in the production of portable electromedical devices for over 30 years. Our wide range includes complete lines of products for electrotherapy, diathermy, ultrasound therapy, magnetotherapy, laser therapy, vibration platforms and a functional evaluation system.

OUR MISSION

The mission of **GLOBUS** is to contribute to promoting health and well-being in the world by providing high quality products to rehabilitation and physiotherapy professionals as well as to private users who wish to improve their quality of life in total safety and benefit from the most modern technologies available on the market.

OUR HISTORY

GLOBUS was founded in 1984 in Codognè, in the Province of Treviso (Italy), by its current President and CEO, Pierpaolo Lucchetta, bronze in volleyball at the Los Angeles Olympics and award-winning professional until 1994.

The high-level experience of its President as an athlete has been fundamental for **GLOBUS** to design and develop different product lines for the treatment of several physical problems, muscle strengthening and the general improvement of sport performance. In order to develop safe, effective and highly technological products, **GLOBUS** collaborates with several universities and research institutes all over the world and scrupulously adheres to the most rigorous scientific guidelines.

Today, **GLOBUS** is a leading company in the rehabilitation and physiotherapy sector, exporting to more than 50 countries around the world. **GLOBUS** products are certified UNI CEI EN ISO 13485:2016 and Quality Assurance System in Europe, FDA in the United States and IEx in Brazil.





KINEO INTELLIGENT LOAD

CONTENTS

- KINEO MULTISTATION	08
- KINEO PULLEY & SQUAT	10
- KINEO LEG EXTENSION	12
- KINEO LEG PRO	14
- KINEO LEG PRESS	16
- KINEO CROSSOVER CABLE	18
- ACCESSORIES	20

KINEO MULTISTATION

TECHNICAL SPECIFICATIONS

· Code K0155M | Code K0155

Dimensions (WxDxH)

260x90x185 cm

Working area dimensions (WxDxH)

305x90x185 cm

Weight

342 Kg

Electric motor

Nominal power 1.6 kW

Power supply

Tension: 230 V - Frequency: 50/60 Hz

Maximal current: 8 A

Maximum cable extraction

4 m

Min/Max Leg Extension load

10-630 N/m

Min/Max Pulley load

1-110 Kg

Min/Max Squat load

1-220 Kg (up to 440 Kg with dedicated accessory)

Min/Max Isokinetic Speeds

0.1-4 m/s

Min/Max Leg Extension Isokinetic Speeds

12-778 °/s

FUNCTIONALITIES

- Wi-Fi connection and remote control
- 2 15.6" touch-screen displays
- 2 USB sockets
- Built-in memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. It provides a variable, controlled torque
- Pulley workstation
- Leg Extension workstation
- Squat workstation








ACCESSORIES

- Kit Leg Curl
- Kit Leg Press
- Lat Machine
- Wood squat platform

COMPLIANCE



SOFTWARE CONFIGURATIONS

TRAINING METHODS	VERSION 3.0	VERSION 5.0	VERSION 7.0
 ISOTONIC NIT (NO INERTIA)	✓	✓	✓
 ISOTONIC WIT (WITH INERTIA)	✓	✓	✓
 ELASTIC		✓	✓
 SMART VISCOUS		✓	✓
 VARIABLE LOAD			✓
 ISOMETRIC			✓
 ISOKINETIC			✓

TEST	VERSION 3.0	VERSION 5.0	VERSION 7.0
MAX CONCENTRIC STRENGTH TEST	✓	✓	✓
MAX CONCENTRIC STRENGTH BALANCE TEST	✓	✓	✓
POWER TEST	✓	✓	✓
MAX ECCENTRIC STRENGTH TEST		✓	✓
MAX ECCENTRIC STRENGTH BALANCE TEST		✓	✓
JUMP POWER TEST		✓	✓
MAX ISOMETRIC STRENGTH TEST		✓	✓
MAX ISOMETRIC STRENGTH BALANCE TEST		✓	✓
MAX ISOKINETIC STRENGTH TEST			✓
MAX ISOKINETIC STRENGTH BALANCE TEST			✓
ISOKINETIC CONCENTRIC/ECCENTRIC RATIO TEST			✓

KINEO PULLEY & SQUAT

TECHNICAL SPECIFICATIONS

· Code K0156M | · Code K0156

Dimensions (WxDxH)

185x90x185 cm

Working area dimensions (WxDxH)

185x90x185 cm

Weight

274 Kg

Electric motor

Nominal power 1.6 kW

Power supply

Tension: 230 V - Frequency: 50/60 Hz

Maximal current: 8 A

Maximum cable extraction

4 m

Min/Max Pulley load

1-110 Kg

Min/Max Squat load

1-220 Kg (up to 440 Kg with dedicated accessory)

Min/Max Isokinetic Speeds

0.1-4 m/s

FUNCTIONALITIES

- Wi-Fi connection and remote control
- 2 15.6" touch-screen displays
- 2 USB sockets
- Built-in memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. It provides a variable, controlled torque
- Pulley workstation
- Squat workstation








ACCESSORIES

- Kit Leg Curl
- Kit Leg Press
- Lat Machine
- Wood squat platform

COMPLIANCE



SOFTWARE CONFIGURATIONS

TRAINING METHODS	VERSION 3.0	VERSION 5.0	VERSION 7.0
 ISOTONIC NIT (NO INERTIA)	✓	✓	✓
 ISOTONIC WIT (WITH INERTIA)	✓	✓	✓
 ELASTIC		✓	✓
 SMART VISCOUS		✓	✓
 VARIABLE LOAD			✓
 ISOMETRIC			✓
 ISOKINETIC			✓

TEST	VERSION 3.0	VERSION 5.0	VERSION 7.0
MAX CONCENTRIC STRENGTH TEST	✓	✓	✓
MAX CONCENTRIC STRENGTH BALANCE TEST	✓	✓	✓
POWER TEST	✓	✓	✓
MAX ECCENTRIC STRENGTH TEST		✓	✓
MAX ECCENTRIC STRENGTH BALANCE TEST		✓	✓
JUMP POWER TEST		✓	✓
MAX ISOMETRIC STRENGTH TEST		✓	✓
MAX ISOMETRIC STRENGTH BALANCE TEST		✓	✓
MAX ISOKINETIC STRENGTH TEST			✓
MAX ISOKINETIC STRENGTH BALANCE TEST			✓
ISOKINETIC CONCENTRIC/ECCENTRIC RATIO TEST			✓

KINEO LEG EXTENSION

TECHNICAL SPECIFICATIONS

· Code K0249M | · Code K0249

Dimensions (WxDxH)

145x90x185 cm

Working area dimensions (WxDxH)

200x90x185 cm

Weigth

290 Kg

Electric motor

Nominal power 1.6 kW

Power supply

Tension: 230 V - Frequency: 50/60 Hz

Maximal current: 8 A

Min/Max Leg Extension load

10-630 nm

Min/Max Leg Extension Isokinetic Speeds

12-778 °/s

FUNCTIONALITIES

- Wi-Fi connection and remote control
- 2 15.6" touch-screen displays
- 2 USB sockets
- Built-in memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. It provides a variable, controlled torque
- Leg Extension workstation

ACCESSORIES

- Kit Leg Press
- Kit Squat
- Wood squat platform

COMPLIANCE



CE
0476

CE



SOFTWARE CONFIGURATIONS

TRAINING METHODS	VERSION 3.0	VERSION 5.0	VERSION 7.0
 ISOTONIC NIT (NO INERTIA)	✓	✓	✓
 ISOTONIC WIT (WITH INERTIA)	✓	✓	✓
 ELASTIC		✓	✓
 SMART VISCOUS		✓	✓
 VARIABLE LOAD			✓
 ISOMETRIC			✓
 ISOKINETIC			✓

TEST	VERSION 3.0	VERSION 5.0	VERSION 7.0
MAX CONCENTRIC STRENGTH TEST	✓	✓	✓
MAX CONCENTRIC STRENGTH BALANCE TEST	✓	✓	✓
POWER TEST	✓	✓	✓
MAX ECCENTRIC STRENGTH TEST		✓	✓
MAX ECCENTRIC STRENGTH BALANCE TEST		✓	✓
JUMP POWER TEST		✓	✓
MAX ISOMETRIC STRENGTH TEST		✓	✓
MAX ISOMETRIC STRENGTH BALANCE TEST		✓	✓
MAX ISOKINETIC STRENGTH TEST			✓
MAX ISOKINETIC STRENGTH BALANCE TEST			✓
ISOKINETIC CONCENTRIC/ECCENTRIC RATIO TEST			✓

KINEO LEG PRO

TECHNICAL SPECIFICATIONS

· Code K0417M | · Code K0417

Dimensions (WxDxH)

260x90x185 cm

Working area dimensions (WxDxH)

310x90x185 cm

Weight

365 Kg

Electric motor

Nominal power 1.6 kW

Power supply

Tension: 230 V - Frequency: 50/60 Hz

Maximal current: 8 A

Min/Max Leg Extension load

10-630 nm

Min/Max Leg Curl load

1-110 Kg

Min/Max Isokinetic Speeds

0.1-4 m/s

Min/Max Leg Extension Isokinetic Speeds

12-778 °/s

FUNCTIONALITIES

- Wi-Fi connection and remote control
- 2 15.6" touch-screen displays
- 2 USB sockets
- Built-in memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. It provides a variable, controlled torque
- Leg Extension workstation
- Leg Curl workstation

ACCESSORIES

- Kit Leg Press
- Kit Squat
- Wood squat platform

COMPLIANCE










CE
0476

CE



SOFTWARE CONFIGURATIONS

TRAINING METHODS	VERSION 3.0	VERSION 5.0	VERSION 7.0
 ISOTONIC NIT (NO INERTIA)	✓	✓	✓
 ISOTONIC WIT (WITH INERTIA)	✓	✓	✓
 ELASTIC		✓	✓
 SMART VISCOUS		✓	✓
 VARIABLE LOAD			✓
 ISOMETRIC			✓
 ISOKINETIC			✓

TEST	VERSION 3.0	VERSION 5.0	VERSION 7.0
MAX CONCENTRIC STRENGTH TEST	✓	✓	✓
MAX CONCENTRIC STRENGTH BALANCE TEST	✓	✓	✓
POWER TEST	✓	✓	✓
MAX ECCENTRIC STRENGTH TEST		✓	✓
MAX ECCENTRIC STRENGTH BALANCE TEST		✓	✓
JUMP POWER TEST		✓	✓
MAX ISOMETRIC STRENGTH TEST		✓	✓
MAX ISOMETRIC STRENGTH BALANCE TEST		✓	✓
MAX ISOKINETIC STRENGTH TEST			✓
MAX ISOKINETIC STRENGTH BALANCE TEST			✓
ISOKINETIC CONCENTRIC/ECCENTRIC RATIO TEST			✓

KINEO LEG PRESS

TECHNICAL SPECIFICATIONS

· Code K1350M | · Code K1350

Dimensions (WxDxH)

320x75x170 cm

Working area dimensions (WxDxH)

320x75x170 cm

Weigth

500 Kg

Electric motor

Nominal power 1.6 kW

Power supply

Tension: 230 V - Frequency: 50/60 Hz

Maximal current: 8 A

Min/Max load

10-440 Kg

Min/Max isokinetic speeds

0.1-4 m/s

FUNCTIONALITIES

- Wi-Fi connection and remote control
- 2 15.6" touch-screen displays
- 2 USB sockets
- Built-in memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. It provides a variable, controlled torque
- Leg Press workstation








ACCESSORIES

- Proprioceptive cushion with easy magnetic hooking

COMPLIANCE



SOFTWARE CONFIGURATIONS

TRAINING METHODS	VERSION 3.0	VERSION 5.0	VERSION 7.0
 ISOTONIC NIT (NO INERTIA)	✓	✓	✓
 ISOTONIC WIT (WITH INERTIA)	✓	✓	✓
 ELASTIC		✓	✓
 SMART VISCOUS		✓	✓
 VARIABLE LOAD			✓
 ISOMETRIC			✓
 ISOKINETIC			✓

TEST	VERSION 3.0	VERSION 5.0	VERSION 7.0
MAX CONCENTRIC STRENGTH TEST	✓	✓	✓
MAX CONCENTRIC STRENGTH BALANCE TEST	✓	✓	✓
POWER TEST	✓	✓	✓
MAX ECCENTRIC STRENGTH TEST		✓	✓
MAX ECCENTRIC STRENGTH BALANCE TEST		✓	✓
JUMP POWER TEST		✓	✓
MAX ISOMETRIC STRENGTH TEST		✓	✓
MAX ISOMETRIC STRENGTH BALANCE TEST		✓	✓
MAX ISOKINETIC STRENGTH TEST			✓
MAX ISOKINETIC STRENGTH BALANCE TEST			✓
ISOKINETIC CONCENTRIC/ECCENTRIC RATIO TEST			✓

KINEO CROSSOVER CABLE



KINEO CROSSOVER CABLE IS THE UNION OF TWO KINEO DEVICES EQUIPPED WITH THE LAT MACHINE ACCESSORY, WHICH ARE CONNECTED BY A CENTRAL BRIDGE. THIS GIVES NEW POSSIBILITIES AND NEW GEOMETRIES TO PERFORM EXERCISES WITH A FREE CABLE, ALLOWING TO USE A CABLE FOR EACH LIMB. THE TWO KINEO DEVICES WORK INDEPENDENTLY OF EACH OTHER, WITH THE TRADITIONAL SOFTWARE FUNCTIONS INCLUDED.



ACCESSORIES



PH 01

PH 01. LAT MACHINE ARM

KO495	LAT MACHINE ARM
-------	-----------------

K1640	KIT LAT MACHINE, COMPOSED OF LAT MACHINE ARM AND LAT MACHINE BENCH
-------	--



PH 02

PH 02. EXTENSION FOR LAT MACHINE ARM

KO314	EXTENSION FOR LAT MACHINE ARM
-------	-------------------------------



PH 03

PH 03. LAT MACHINE BENCH

KO378	LAT MACHINE BENCH
-------	-------------------



PH 04

PH 04. KIT LEG CURL

KO133	KIT LEG CURL
-------	--------------

ACCESSORIES



PH 05. JUMP PLATFORM

KO393 | JUMP PLATFORM

PH 05



PH 06. SQUAT PLATFORM

KO449 | SQUAT PLATFORM

PH 06



PH 07. STEP

KO034 | STEP

PH 07



PH 08. KINEO ACCESSORIES CART

KO762 | KINEO ACCESSORIES CART

PH 08

ACCESSORIES

KIT LEG PRESS

TECHNICAL SPECIFICATIONS

· Code K1161

Dimensions (WxDxH)

220x75x170 cm

Working area dimensions (WxDxH)

220x75x170 cm

Weight

230 Kg

Maximum seat weight

135 Kg

Min/Max load*

10-440 Kg

Min/Max isokinetic speed*

0.1-4 m/s

* The technical data refer to the KINEO device when connected.

KINEO LEG PRESS Accessory cannot work autonomously.

FUNCTIONALITIES

· Adjustable depth of the seat and 7 seatback positions

· High-quality Italian leather

HARDWARE

· Leg Press and Rowing workstation with Rowing bar

· Wheels with lock system

ACCESSORIES

· Proprioceptive cushion with easy magnetic hooking

COMPLIANCE





REHAB Line

CONTENTS

- TECARTHERAPY Line	26
- MAGNETOTHERAPY Line	32
- ELECTROTHERAPY Line	44
- ULTRASOUND THERAPY Line	56
- LASER THERAPY Line for physiotherapy	60
- LASER THERAPY Line for podology	70
- PRESSOTHERAPY Line	76
- PREVENTION	80

TECARTHERAPY

Line



DIACARE 7000

TECARTHERAPY

DIACARE 7000 IS THE MOST POWERFUL, VERSATILE AND COMPLETE DEVICE IN THE GLOBUS TECARTHERAPY RANGE. IT USES BOTH THE RESISTIVE AND THE CAPACITIVE MODE FOR THE TREATMENT OF ACUTE AND CHRONIC PATHOLOGIES OF THE MUSCULAR AND SKELETAL SYSTEM. THE DEVICE CAN BE USED BY ACTING ON THE PATIENT BOTH IN STATIC MANUAL MODE AND COMBINING PASSIVE MOVEMENTS OR, THANKS TO THE HANDS-FREE KIT, IT CAN ALSO BE USED DURING ACTIVE MOVEMENTS OF THE SUBJECT. DIACARE 7000 IS READY TO USE WITH APPROPRIATE ERGONOMIC TOOLS DESIGNED FOR THE MOBILIZATION AND THE ADVANCED TREATMENT OF ALL SOFT TISSUE DISEASES (FASCIA TOOLS).

TECHNICAL SPECIFICATIONS

Display

Size 7" color touch screen

Nominal Power (NP)

350 W $\pm 10\%$

Emission frequencies

400 kHz, 448 kHz, 470 kHz, 700 kHz, 1000 kHz, 1200 kHz, $\pm 10\%$

Processing methods

Capacitive and resistive

Capacitive electrodes

3 size \varnothing (30, 50, 70 mm)

Resistive electrodes

3 size \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1,8 A

Size

515x340x170 mm

Weight

11.5 Kg

STANDARD EQUIPMENT

· Cod. G5210

· 1 DIACARE 7000 tecarthrapy unit

· 1 Counter electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Bipolar handpiece \varnothing (50 mm)

· 1 Hands-free treatment kit

· 1 1000 ml conductive cream bottle

· 1 Power supply cable

· 1 Operating manual



STRONG POINTS



DIACARE 5000 | DIACARE 5000 RE

TECARTHERAPY | 14 PROGRAMS

THE DIACARE 5000 MODEL PROVIDES ALL THE QUALITY OF THE GLOBUS TECARTHERAPY IN A MACHINE OF SMALL SIZE AND WEIGHT, TO GUARANTEE MAXIMUM PORTABILITY AND KEEP ITS EFFECTIVENESS AND PROFESSIONALISM UNCHANGED.

TECHNICAL SPECIFICATIONS

Display

Size 4.3" color touch screen

Nominal Power (NP)

180 W $\pm 10\%$

Emission frequencies

400 kHz, 448 kHz, 470 kHz, 700 kHz, 1000 kHz, $\pm 10\%$

Processing methods

Capacitive and resistive

Capacitive electrodes

3 size \varnothing (30, 50, 70 mm)

Resistive electrodes

3 size \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1,8 A

Size


170x220x60 mm

Weight

1.04 Kg

STANDARD EQUIPMENT

- 1 Suitcase
- 1 DIACARE 5000 tecartherapy unit
- 1 Counter electrode (flexible plate)
- 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 250 ml conductive cream bottle
- 1 Power supply
- 1 Operating manual

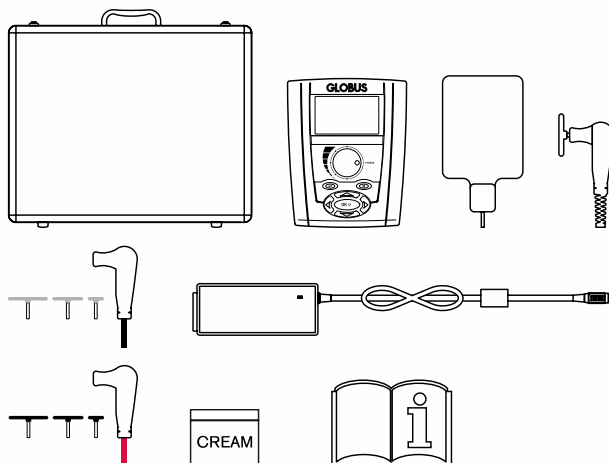

PURCHASE OPTIONS


- Code **G3488** DIACARE 5000
- Code **G3996** DIACARE 5000 RE *

* RECHARGEABLE TECARTHERAPY:

- 5 HOURS OF TREATMENT INCLUDED IN THE COST OF THE DEVICE
- CHEAP CHARGING
- GUARANTEED SUPPORT FOR 5 YEARS AND FORKLIFT
- CHARGING VIA INTERNET



STRONG POINTS

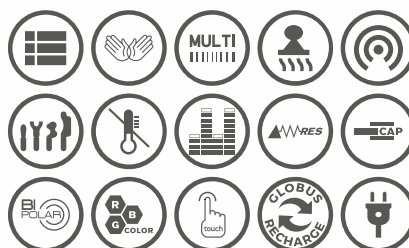

FASCIA TOOLS

FOR THE DIACARE 7000 AND 5000 LINE

GLOBUS has successfully combined the DIACARE technology with instrument assisted soft tissues mobilization (IASTM), bringing together the advantages of these two extraordinary therapeutic possibilities. The GLOBUS FASCIA TOOLS are special ergonomic electrodes designed for the mobilization and the advanced treatment of all soft tissue diseases.

Thanks to the joint application of tecartherapy and the myofascial electrodes, it is possible to:

- **ELIMINATE FASCIAL RESTRICTIONS AND SCAR ADHESIONS**
- **INCREASE THE EXTENSIBILITY OF THE CONNECTIVE TISSUE**
- **IMPROVE TISSUE VASCULARIZATION**
- **IMPROVE MUSCLE RECRUITMENT AND MOTOR CONTROL**



Code **G5651** KIT FASCIA TOOLS FOR MANUAL TREATMENTS



DiaBlade Large: tool for the treatment of the Lower Limb and the Trunk.



DiaBlade Small: tool for the treatment of the Upper Limb and the Cervical Rachis.



Dia Trigger Point: tool for the treatment of Trigger Points.



Dia Spine: tool for Back treatment and Vertebral Mobilization.



MAGNETOTHERAPY

Line



MAGNUM 3500 PRO

PROFESSIONAL MAGNETOTHERAPY | 70 PROGRAMS

THE ONLY GLOBUS MAGNETOTHERAPY DEVICE WITH FOUR CHANNELS THAT ALLOWS TO TREAT DIFFERENT BODY PARTS OR TWO SUBJECTS AT THE SAME TIME. MAGNUM 3500 PRO IS EQUIPPED WITH A WIDE RANGE OF ACCESSORIES FOR ALL TREATMENT NEEDS.

26 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS
6 INFLAMMATIONS | 6 ORTHOPEDICS | 3 PAIN | 3 TRAUMAS/EDEMAS/WOUNDS | 3 NIGHT PROGRAMS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4

Intensity can be emitted with flexible diffuser

- 200 Gauss average per channel
- 800 Gauss on average on the whole (for 4 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Free programming

Possibility of adding up to 10 custom programs

2+2 Function

Possibility of running 2 different programs at the same time, using 2 channels each

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions | Weight

170x220x60 mm | 820 gr

STANDARD EQUIPMENT

· Cod. G5434

· 1 Suitcase

· 1 MAGNUM 3500 PRO magnetotherapy unit

· 2 Flexible solenoids 300x100 mm

· 2 Soft solenoids 120x120 mm


· 2 "Pocket Pro" solenoids 155x95 mm

· 1 "Soft 4" mat 250x250 mm

· 1 Magnum Test

· 1 Power supply

· 1 Operating manual

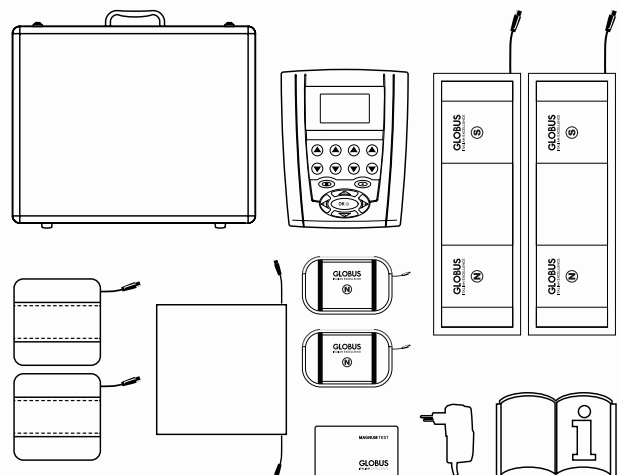

CE
0476



STRONG POINTS



MAGNUM 3000 PRO

PROFESSIONAL MAGNETOTHERAPY | 70 PROGRAMS

POWERFUL PROFESSIONAL MAGNETOTHERAPY EQUIPPED WITH MANY SPECIAL PROGRAMS AND FUNCTIONS THAT MAKE IT EFFECTIVE FOR THE TREATMENT OF SEVERAL OSTEOARTICULAR PATHOLOGIES AND FOR THE RECOVERY AFTER ORTHOPEDIC SURGERY. POSSIBILITY OF LOCAL AND TOTAL-BODY TREATMENTS.

26 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS
6 INFLAMMATIONS | 6 ORTHOPEDICS | 3 PAIN | 3 TRAUMAS/EDEMAS/WOUNDS | 3 NIGHT PROGRAMS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2

Intensity can be emitted with flexible diffuser

- 200 Gauss average per channel
- 400 Gauss on average on the whole (for 2 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Free programming

Possibility of adding up to 10 custom programs

1+1 Function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply


- Mains
- Rechargeable battery (not included)

Dimensions | Weight

170x220x60 mm | 820 gr

STANDARD EQUIPMENT

- 1 Bag
- 1 MAGNUM 3000 PRO magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 1 Magnum Test
- 1 Power supply
- 1 Operating manual

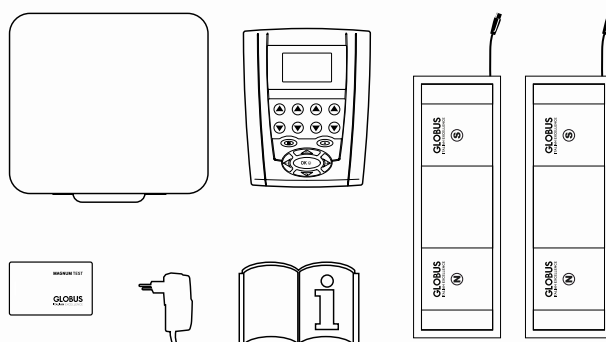

PURCHASE OPTIONS

- Code **G6034** MAGNUM 3000 PRO, 2 "Pocket Pro" solenoids
- Code **G5335** MAGNUM 3000 PRO, 2 flexible solenoids
- Code **G5643** MAGNUM 3000 PRO, 2 soft solenoids



CE
0476

STRONG POINTS



MAGNUM XL PRO

MAGNETOTHERAPY | 46 PROGRAMS

MAGNETOTHERAPY DEVICE WITH A COMPLETE RANGE OF PRESET PROGRAMS TO TREAT THE MOST COMMON OSTEOARTICULAR DISORDERS. THE DEVICE IS AMONG THE MOST POWERFUL ON THE MARKET AND IS EQUIPPED WITH TWO SOLENOIDS.

32 HEALTH AND REHAB | 14 AESTHETICS AND FITNESS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2

Intensity can be emitted with flexible diffuser

- 200 Gauss average per channel
- 400 Gauss on average on the whole (for 2 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


100x160x35 mm

Weight

280 gr

STANDARD EQUIPMENT

- 1 Bag
- 1 MAGNUM XL PRO magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 1 Magnum Test
- 1 Power supply
- 1 Operating manual

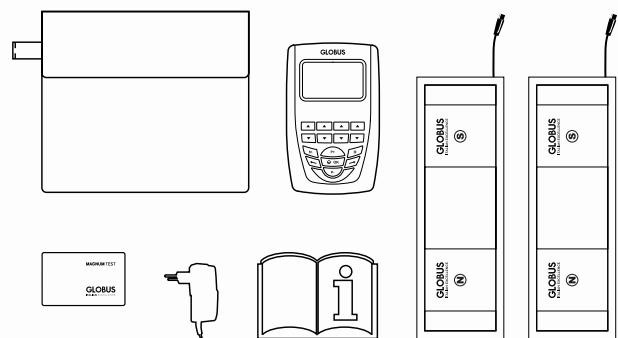

PURCHASE OPTIONS

- Code **G6033** MAGNUM XL PRO, 2 "Pocket Pro" solenoids
- Code **G3956** MAGNUM XL PRO, 2 flexible solenoids
- Code **G4278** MAGNUM XL PRO, 2 soft solenoids



STRONG POINTS



MAGNUM 2500

MAGNETOTHERAPY | 52 PROGRAMS

PORTABLE MAGNETOTHERAPY, IDEAL FOR HOME TREATMENTS. MAGNUM 2500 IS INDICATED FOR A BETTER RECOVERY FROM TRAUMAS AND FRACTURES, BUT ALSO TO PROVIDE RELIEF AND HELP ARTICULAR REGENERATION IN PATHOLOGIES LIKE ARTHROSIS AND ARTHRITIS.

16 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 14 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS
5 INFLAMMATIONS | 4 ORTHOPEDICS | 3 PAIN | 2 TRAUMAS/EDEMAS/WOUNDS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2

Intensity can be emitted with flexible diffuser

- 160 Gauss average per channel
- 320 Gauss on average on the whole (for 2 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Free programming

Possibility of adding up to 5 custom programs

1+1 Function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply

- Mains
- Rechargeable battery (not included)

Size


100x160x35 mm

Weight

280 gr

STANDARD EQUIPMENT

- 1 Bag
- 1 MAGNUM 2500 magnetotherapy unit
- 1 Flexible solenoid 300x100 mm
- 1 Magnum Test
- 1 Power supply
- 1 Operating manual

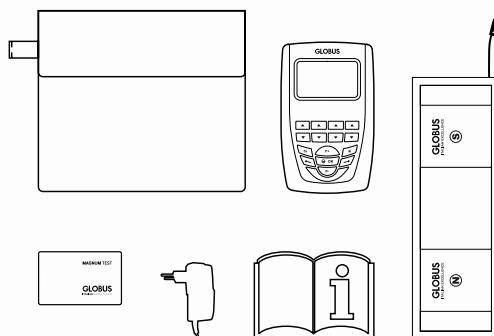

PURCHASE OPTIONS

- Code **G6191** MAGNUM 2500, 2 "Pocket Pro" solenoids
- Code **G5438** MAGNUM 2500, 1 flexible solenoid
- Code **G5645** MAGNUM 2500, 2 soft solenoids



STRONG POINTS



MAGNUM XL

MAGNETOTHERAPY | 26 PROGRAMS

DEVICE FOR HOME MAGNETOTHERAPY SUITABLE FOR TREATING THE MAIN OSTEOARTICULAR DISORDERS. EVEN EASIER TO USE WITH THE NEW EASY PROGRAM FUNCTION.

16 HEALTH AND REHAB | 10 AESTHETICS AND FITNESS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2

Intensity can be emitted with flexible diffuser

- 140 Gauss average per channel
- 280 Gauss on average on the whole (for 2 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


100x160x35 mm

Weight

260 gr

STANDARD EQUIPMENT

- 1 Bag
- 1 MAGNUM XL magnetotherapy unit
- 1 Flexible solenoid 300x100 mm
- 1 Power supply
- 1 Operating manual

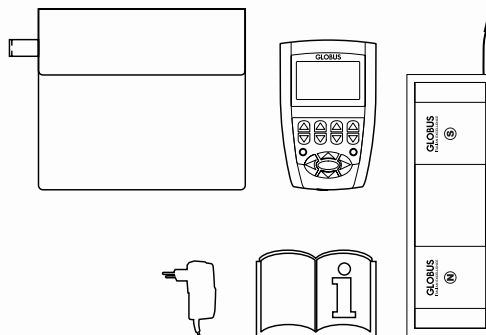

PURCHASE OPTIONS

- Code G3216 MAGNUM XL, 1 flexible solenoid
- Code G4279 MAGNUM XL, 2 soft solenoids



STRONG POINTS



MAGNUM L

MAGNETOTHERAPY | 8 PROGRAMS

VERY SIMPLE MAGNETOTHERAPY DEVICE, EQUIPPED WITH 8 PRE-SET PROGRAMS TO TREAT MAINLY FRACTURES AND ARTHROSIS.

8 HEALTH AND REHAB

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

1

Intensity can be emitted with flexible diffuser

120 Gauss on average

Frequency

5-200 Hz

Duty Cycle

5-50%

Power supply

Mains

Dimensions

100x160x35 mm

Weight

260 gr

STANDARD EQUIPMENT

· Code **G3947**


· 1 Bag

· 1 MAGNUM L magnetotherapy unit

· 1 Flexible solenoid 300x100 mm

· 1 Power supply

· 1 Operating manual

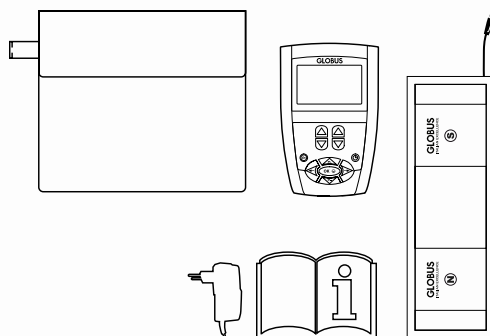

CE
0476



STRONG POINTS






PROGRAMS TABLE

PROGRAM	3500 PRO	3000 PRO	XL PRO	2500	XL	L
Algodystrophy / Sudeck disease	✓	✓	/	✓	/	/
Anti-inflammatory	✓	✓	/	✓	/	/
Arthritis	✓	✓	✓	✓	/	✓
Arthrosis - Osteoarthritis	✓	✓	✓	✓	✓	✓
Ankle sprain	✓	✓	✓	✓	/	/
Joint pain	✓	✓	✓	✓	/	/
Bone edemas	✓	✓	/	✓	/	/
Edemas-contusions	✓	✓	✓	✓	/	/
Epicondylitis / Epitrochleitis	✓	✓	✓	✓	/	✓
Fracture with plaster	✓	✓	*	✓	*	*
Fractures	✓	✓	✓	✓	✓	✓
Spinal fusion	✓	✓	/	/	/	/
Cartilage lesion	✓	✓	✓	✓	/	/
Back pain	✓	✓	✓	✓	✓	✓
Osteonecrosis	✓	✓	/	✓	/	/
Osteoporosis	✓	✓	✓	✓	✓	✓
Osteotomy	✓	✓	/	/	/	/
After ACL surgery	✓	✓	✓	✓	/	/
Hip - knee prosthesis	✓	✓	✓	✓	/	/
Pseudoarthritis	✓	✓	/	/	/	/
Delayed Fracture Consolidation	✓	✓	✓	✓	✓	/
Shoulder tendinitis / Rotator cuff	✓	✓	✓	✓	✓	/
Tendinitis	✓	✓	✓	/	/	/
Carpal tunnel	✓	✓	✓	✓	✓	/
Ulcers	✓	✓	/	/	/	/
Multifrequency program	✓	✓	/	1	/	/
Pre-set frequency programs	25	25	/	15	/	/
Night programs	3	3	/	/	/	/
Aesthetic programs	/	/	14	/	10	/

* Use the “Fractures” program.

Some programs are divided into body areas. For further information, see the technical sheets of the products.

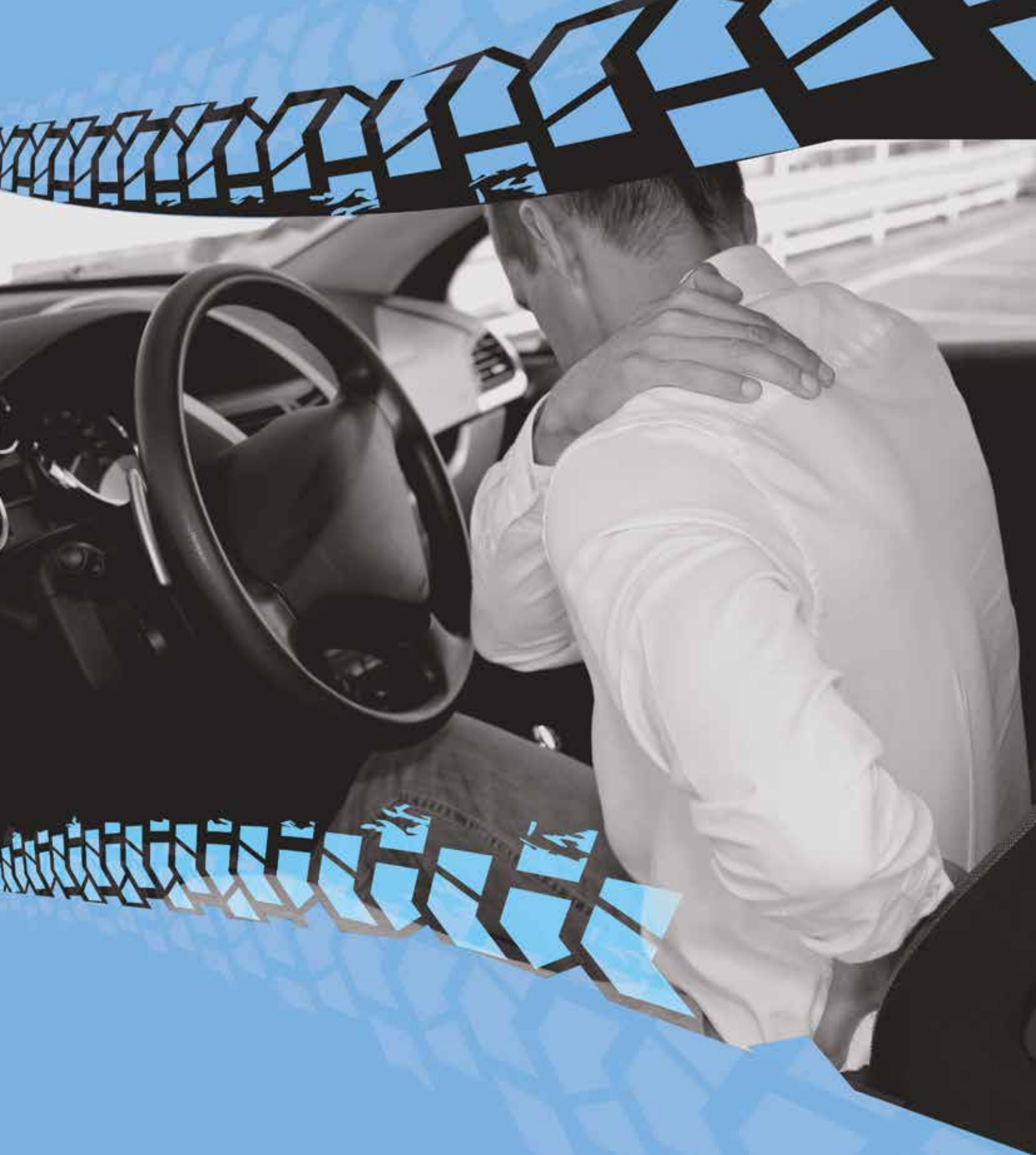

BACK PAIN?

with the **new MAGNUM 2 Pro Drive**
driving will be a **pleasure!**

MAGNUM 2 PRO DRIVE

MAGNETOTHERAPY | 41 PROGRAMS

MAGNUM 2 PRO DRIVE IS DESIGNED TO OFFER A SOLUTION TO THOSE WHO SUFFER FROM PAIN IN THE BACK OR OTHER PARTS OF THE BODY DUE TO TOO MUCH TIME SPENT DRIVING. SEDENTARY LIFESTYLE AND MANY HOURS BEHIND THE WHEEL, SUBJECTED TO THE VIBRATIONS OF THE VEHICLE, ARE SURELY AMONG THE MOST FREQUENT CAUSES OF MUSCULOSKELETAL PAIN. MAGNUM 2 PRO DRIVE IS AN EXCELLENT INSTRUMENTAL THERAPY THAT CAN HELP TO REDUCE THE INFLAMMATION AND PAIN OF THIS CATEGORY OF PEOPLE.

15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 6 INFLAMMATIONS | 6 ORTHOPEDICS
3 PAIN | 3 TRAUMAS/EDEMAS/WOUNDS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2

Intensity can be emitted with flexible diffuser

- 200 Gauss average per channel
- 400 Gauss on average on the whole (for 2 channels)

Frequency

5-200 Hz

Duty Cycle

5-50%

Free programming

Possibility of adding up to 5 custom programs

1+1 Function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

- Code G5429
- 1 Bag
- 1 MAGNUM 2 PRO DRIVE magnetotherapy unit
- 1 "Soft 4" mat 250x250 mm
- 1 Flexible solenoid 300x100 mm
- 1 Magnum Test
- 1 Power supply
- 1 Operating manual

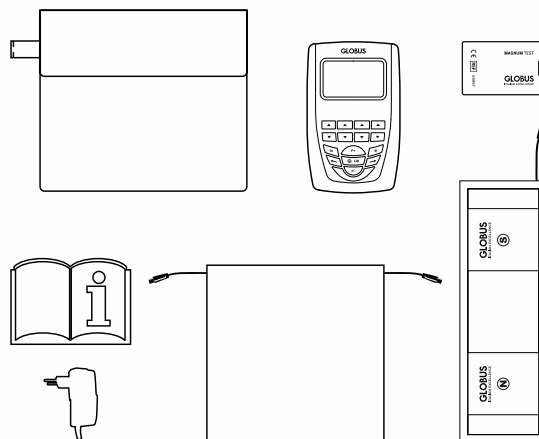

STRONG POINTS



ELECTROTHERAPY

Line



GENESY 3000

4-CHANNEL ELECTROTHERAPY | 423 PROGRAMS

TOP OF THE RANGE PROFESSIONAL PRODUCT WITH FOUR CHANNELS USED IN THE REHABILITATION SECTOR FOR THE TREATMENT OF THE MOST COMMON PATHOLOGIES OF THE MUSCULOSKELETAL SYSTEM.

236 ELECTROTHERAPY MENU (EMS/RUSSIAN) | 53 SPORT | 33 PAIN | 29 FITNESS | 23 MICROCURRENTS
16 BEAUTY | 11 REHAB | 9 DENERVATED | 9 INTERFERENTIALS | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

- Mains
- Rechargeable batteries

Size

170x220x60 mm

Weight

1.12 Kg

* Based on the type of current

STANDARD EQUIPMENT

· Code **G1034**

· 1 Bag

· 1 GENESY 3000 stimulator

· 4 Cables for electrode connection

· 2 Cables for microcurrents and ionophoresis


· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes


· 2 Ionophoresis electrodes + 2 bands

· 1 Power supply

· 1 Operating manual

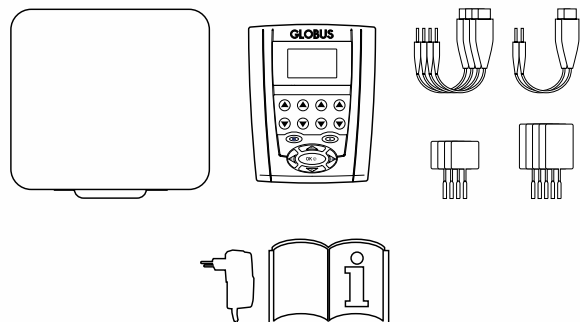

STRONG POINTS

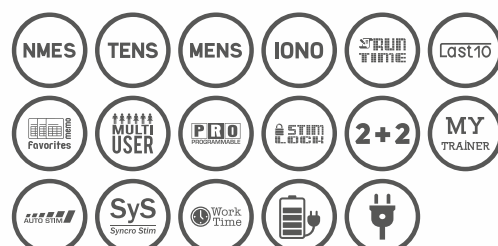

GENESY 1500

4-CHANNEL ELECTROTHERAPY | 414 PROGRAMS

FOUR-CHANNEL PROFESSIONAL PRODUCT USED IN THE MEDICAL AND REHABILITATION SECTOR. IT HAS ALL THE CURRENTS USED FOR REHABILITATION, PHYSIOTHERAPY AND MOTOR REHABILITATION THERAPY.

92 ELECTROTHERAPY MENU (EMS/RUSSIAN | 84 ACTION NOW | 54 SERIAL SEQUENTIAL STIMULATION | 53 SPORT
33 PAIN | 29 FITNESS | 23 MICROCURRENTS | 16 BEAUTY | 11 REHAB | 9 DENERVATED | 6 INTERFERENTIALS
3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

- Mains
- Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code **G3554**

· 1 Bag

· 1 GENESY 1500 stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Power supply


· 1 Operating manual

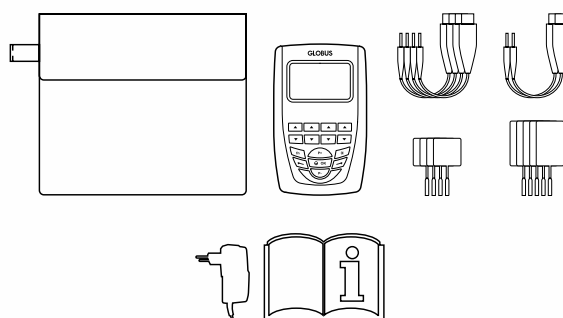

CE
0476

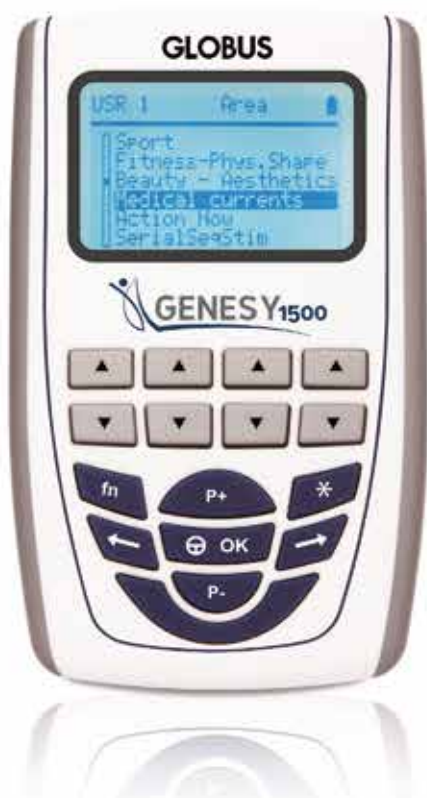

STRONG POINTS

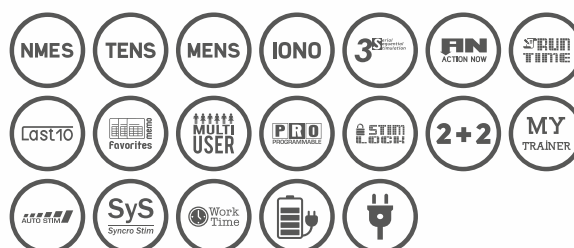

GENESY 600

4-CHANNEL ELECTROTHERAPY | 149 PROGRAMS

FOUR-CHANNEL PROFESSIONAL PRODUCT FOR DEMANDING USERS AND REHABILITATION PROFESSIONALS.

53 SPORT | 29 FITNESS | 20 PAIN | 12 MICROCURRENTS | 12 SERIAL SEQUENTIAL STIMULATION | 11 REHAB
7 ACTION NOW | 3 DENERVATED | 1 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm


Weight

450 gr


* Based on the type of current

STANDARD EQUIPMENT

- Code G3553
- 1 Bag
- 1 GENESY 600 stimulator
- 4 Cables for electrode connection
- 2 Cables for microcurrents and ionophoresis
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Battery charger
- 1 Operating manual

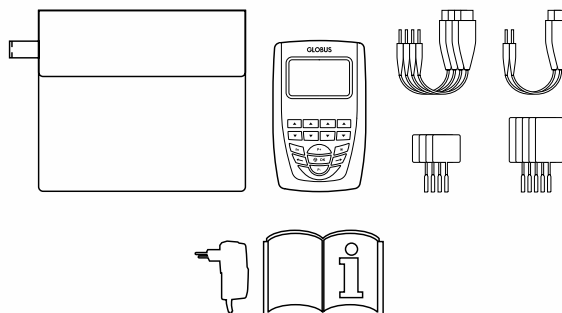

STRONG POINTS

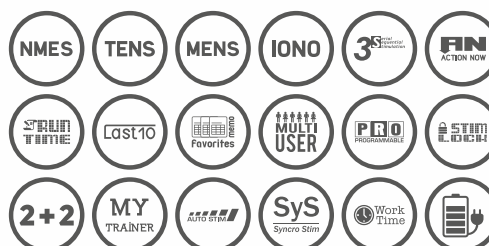

GENESY 300 PRO

4-CHANNEL ELECTROTHERAPY | 91 PROGRAMS

FOUR-CHANNEL PRODUCT WITH A WIDE RANGE OF PROGRAMS FOR PHYSIOTHERAPY.

29 FITNESS, BEAUTY, WELLNESS | 20 PAIN | 22 SPORT | 12 MICROCURRENTS | 6 REHAB
1 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm


Peso

440 gr

* Based on the type of current

STANDARD EQUIPMENT


- Code G3222
- 1 Bag
- 1 GENESY 300 PRO stimulator
- 4 Cables for electrode connection
- 2 Cables for microcurrents and ionophoresis
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Battery charger
- 1 Operating manual

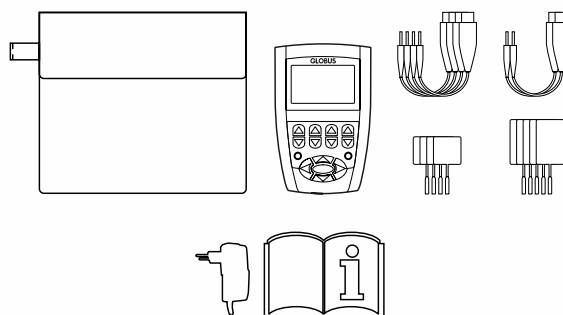

CE
0476



STRONG POINTS

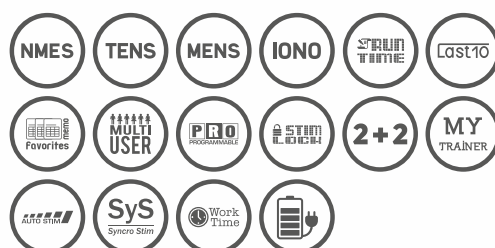

GENESY S2

2-CHANNEL ELECTROTHERAPY | 60 PROGRAMS

THE TWO-CHANNEL PRODUCT FOR PHYSIOTHERAPY.

16 TENS AND PAIN | 12 SPORT | 11 REHAB | 9 PREVENTION | 8 FITNESS, BEAUTY, WELLNESS
3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2 independent (4 electrodes)

Frequency

1-150 Hz

Pulse width

50-400 μ s

Intensity

0-100 mA per channel

Power supply

Rechargeable batteries

Size


80x130x25 mm

Weight

220 gr

STANDARD EQUIPMENT

- Code G3725
- 1 Bag
- 1 GENESY SII stimulator
- 2 Cables for electrode connection
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Battery charger
- 1 Operating manual

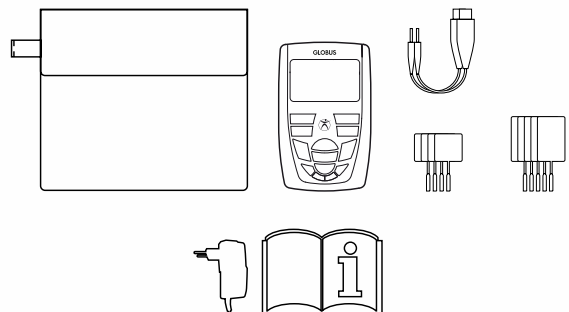

STRONG POINTS



DUO TENS

2-CHANNEL ELECTROTHERAPY | 20 PROGRAMS

THE IDEAL PRODUCT TO TREAT THE MOST COMMON PAINFUL SYMPTOMS.

16 TENS AND PAIN | 3 SPORT | 1 FITNESS, BEAUTY, WELLNESS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2 independent (4 electrodes)

Frequency

1-150 Hz

Pulse width

50-400 μ s

Intensity

0-100 mA per channel

Power supply

Rechargeable batteries

Size


80x130x25 mm

Weight

220 gr

STANDARD EQUIPMENT

- Code G3729
- 1 Bag
- 1 DUO TENS stimulators
- 2 Cables for electrode connection
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Battery charger
- 1 Operating manual

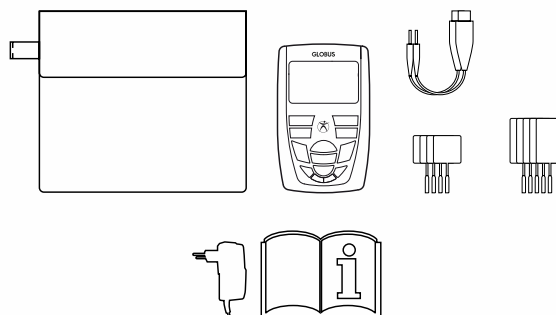

CE
0476



STRONG POINTS



INCONTINENCE of **women** and **men**

GET RID OF IT NOW
*with the electrostimulators of the **GENESY** line
for specific training of the pelvic muscles*



REHABILITATION OF THE PELVIC FLOOR

Pelvic floor rehabilitation aims at improving or solving problems such as male and female urinary incontinence in a non-invasive and painless way, using appropriate physical therapies in combination with manual therapy and specific exercises.

In particular, using electrostimulation and the special vaginal or anal probe, it is possible to recruit the hypotonic perineal muscles or reduce the detrusor reactivity.

This type of treatment is also used for pain relief to reduce chronic pelvic floor pain. Our **GENESY Line** devices include programs for the treatment of urge, stress and mixed incontinence. The treatments can be carried out by using the vaginal probe, for women, or the anal probe, for men.

OUR PROBES



VAGINAL PROBE APPLICATION



Code **G1188** SINGLE-PATIENT VAGINAL PROBE

ANAL PROBE APPLICATION




Code **G0757** SINGLE-PATIENT ANAL PROBE

Rehabilitation and **DYSPHAGIA**

Dysphagia is a difficulty in carrying out the act of swallowing.

For those affected, it becomes very complicated, or even impossible, to swallow both saliva and liquid and/or solid food.

In most cases, dysphagia is due to neurological pathologies, but it frequently affects elderly people who suffer from a physiological weakening of the muscles involved in swallowing.

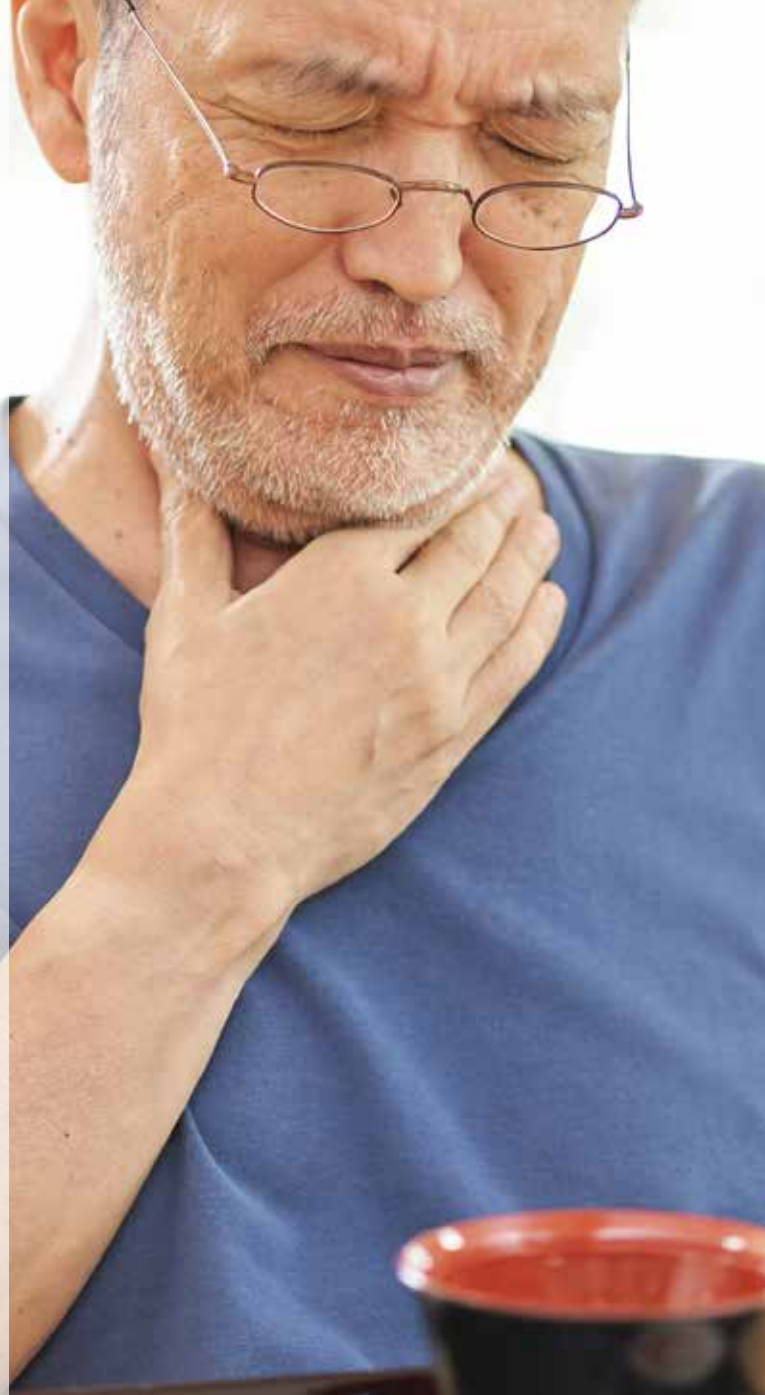

OUR DEVICES

AMI PRO and **AMI LIFE** are support devices for the rehabilitation of these patients and they act by promoting muscle trophism and improving proprioceptive neuromuscular sensitivity with the aim of re-educating the muscles involved in swallowing.

The devices intervene by inducing a non-invasive neuromuscular stimulation (NMES) that stimulates the muscles through electrodes positioned on the neck.

The protocols comply with the guidelines of the main scientific evidence currently available in the literature.

AMI PRO PRO is the top-of-the-range device, fully equipped with different stimulation programs with motor, sensorimotor and sensory goals, with different parameters to enable the most appropriate choice of the work cycle. It also features various facial rehabilitation programs as well as programs that can be edited to suit the professional's therapeutic needs.

AMI PRO has the exclusive special **"ACTION NOW"** function, to perform the stimulation only at the request of the operator or the patient.

By pressing the **"ACTION NOW"** button on the device, it is possible to coordinate muscle activation with the patient's voluntary swallowing, creating an actual self-feedback system.

AMI LIFE has 3 stimulation programs with motor, sensorimotor and sensory goals and a facial rehabilitation program.

It is not possible to edit the pre-set programs and for this reason it is suitable for patients who need to continue the therapy at home (always under the supervision of the doctor).

AMI PRO | AMI LIFE

4-CHANNEL ELECTROTHERAPY FOR DYSPHAGIA

TECHNICAL SPECIFICATION

Display

Size of the visible area 2.6

Channels

4 independent

Frequency

1-80 Hz

Pulse width

80-260 μ s

Intensity

max 20 mA per channel

Power supply

Rechargeable batteries

Dimensions | Weight


100x160x35 mm | 450 gr



CE
0476

STANDARD EQUIPMENT

- 1 Bag
- 1 Stimulator
- 4 Cables for electrode connection
- 4 Butterfly electrodes
- 1 Battery charger
- 1 Operating manual

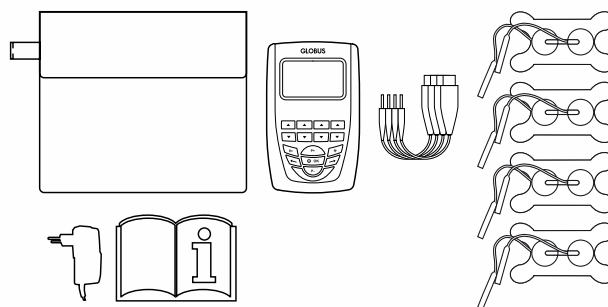

PURCHASE OPTIONS

- Code G5573 AMI PRO
- Code G5852 AMI LIFE

STRONG POINTS



STRONG POINTS



The use of these devices is only intended for doctors, authorized health professionals and their patients only under medical prescription and under the continuous supervision of the professional.

ULTRASOUND THERAPY

Line



MEDISOUND 3000

ULTRASOUND THERAPY | 49 PROGRAMS

TOP-OF-THE-RANGE PRODUCT SUITABLE FOR PHYSIOTHERAPY CENTERS AND HOSPITALS. BEING COMPACT, IT CAN BE EASILY TRANSPORTED IN ALL OPERATIONAL SITUATIONS.

31 AESTHETICS | 18 HEALTH

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

3 W/cm² ± 20%

Contact management

Yes

Emitted frequency

3 MHz and 1 MHz

Emission

Continuous and pulsed

Duty Cycle

10%-100%

Power suppl

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

820 gr

STANDARD EQUIPMENT

- Code G1033
- 1 Bag
- 1 MEDISOUND 3000 ultrasound therapy unit
- 1 Handpiece
- 1 Gel bottle
- 1 Power supply
- 1 Operating manual

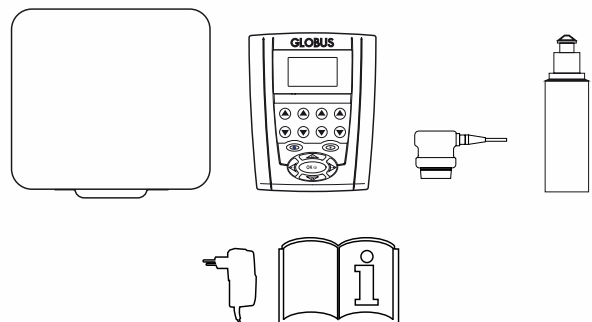

STRONG POINTS



HT 922

ULTRASOUND THERAPY

HOME CARE DEVICE USED FOR BOTH MEDICAL AND BEAUTY TREATMENTS WITH A COMBINED THERAPY THAT MAKES A SYNERGIC USE OF THE PAIN-RELIEVING AND ANTI-INFLAMMATORY EFFECTS OF ULTRASOUND THERAPY AND MICROCURRENTS.

9 PROGRAMS WITH MICROCURRENTS AND NEGATIVE IONS WITH DIFFERENT MODULATIONS
1 CUSTOMIZABLE FOR ULTRASOUND

TECHNICAL SPECIFICATIONS

Maximum power

5 W, 9 intensity levels

Emission frequency

1 MHz $\pm 5\%$

Head diameter

30 mm

Transducer

Integrated

Power supply

Mains

Time settings

1' to 30' with steps of 1

Size

50x190x50 mm

Weight

160 gr

STANDARD EQUIPMENT

· Code G5936

· 1 Bag

· 1 HT 922 devices


· 1 Cable for electrode connection

· 2 Self-adhesive square electrode

· 1 Gel bottle

· 1 Power supply

· 1 User manual

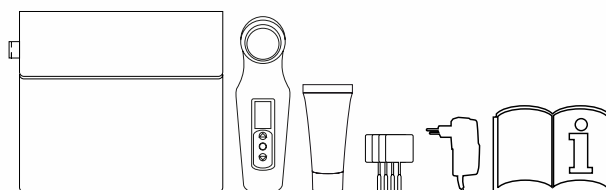

STRONG POINTS



LASER THERAPY

Line for physiotherapy



PHYSIOLASER 12.0 PRO

LASER THERAPY

THE "PRO" VERSION OF PHYSIOLASER 12.0 IS DESIGNED FOR THE PROFESSIONAL WHO USES THE DEVICE FREQUENTLY AND PREFERS TO KEEP IT ALWAYS NEARBY ON A CART. THE 12 W OF MAXIMUM POWER ALLOW RAPID AND EFFECTIVE TREATMENTS EVEN ON THE LARGEST BODY AREAS, USING THE 22 MM OPTIC INCLUDED. THE WIDE LIBRARY OF PROGRAMS ALSO GUARANTEES THE PROFESSIONAL TO TREAT ANY MUSCULOSKELETAL DISORDER AND SOME PROBLEMS RELATING TO WOUNDS AND ULCERS OF THE SKIN.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12 W $\pm 20\%$

Compatible with a wavelength of

808 nm

980 nm and 1064 nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


220x170x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

- Code **G5803**
- 1 Suitcase
- 1 PHYSIOLASER 12.0 PRO laser therapy unit
- 1 Handpiece HP 808 nm with Ø 11 mm and Ø 22 mm optics
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

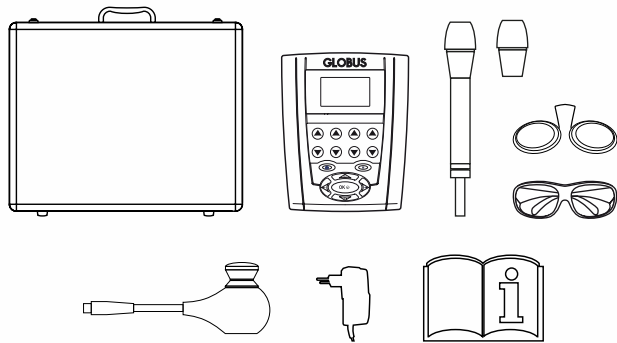

PURCHASE OPTIONS

- Code **G5995** HP 980 nm handpiece
- Code **G5996** 1064 nm handpiece



STRONG POINTS



PHYSIOLASER 12.0

LASER THERAPY

PHYSIOLASER 12.0 IS THE MOST USER-FRIENDLY AND POWERFUL LASER THERAPY DEVICE IN THE GLOBUS LINE. THE WIDE LIBRARY OF PROGRAMS ALSO ALLOWS THE PROFESSIONAL TO TREAT ANY MUSCULOSKELETAL DISORDER AND SOME PROBLEMS RELATING TO WOUNDS AND ULCERS OF THE SKIN. DEFINITELY ONE OF THE DEVICES WITH THE BEST PRICE-QUALITY RATIO IN THE ENTIRE RANGE.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12 W $\pm 20\%$

Compatible with a wavelength of

808 nm

980 nm and 1064 nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

- Code **G5741**
- 1 Suitcase
- 1 PHYSIOLASER 12. laser therapy unit
- 1 Handpiece HP 808 nm with Ø 11 mm e Ø 22 mm optics
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

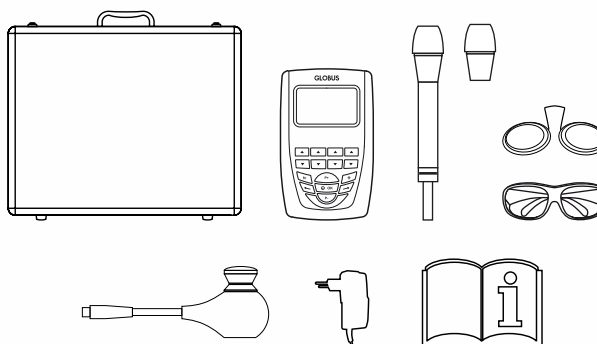

PURCHASE OPTIONS

- Code **G5995** HP 980 nm handpiece
- Code **G5996** 1064 nm handpiece



STRONG POINTS



PHYSIOLASER 6.0 PRO

LASER THERAPY

PHYSIOLASER 6.0 PRO IS A DEVICE OF SMALL WEIGHT AND SIZE, SUITABLE BOTH TO BE PLACED ON A STUDIO CART AND TO BE TRANSPORTED, IF NECESSARY, TO OTHER LOCATIONS. WITH ITS MAXIMUM POWER OF 6W, IT IS IDEAL FOR TREATING OSTEOARTICULAR AND MUSCULAR PATHOLOGIES, QUICKLY REDUCING PAIN AND EDEMA AND ACCELERATING THE RECOVERY OF INJURED AND INFLAMED TISSUES.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6 W $\pm 20\%$

Compatible with a wavelength of

808 nm

980 nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000 Hz

Power supply

- Mains

- Rechargeable battery

Size

220x170x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

· Code **G5764**

· 1 Suitcase

· 1 PHYSIOLASER 6.0 PRO laser therapy unit


· 1 Handpiece MP 808 nm with Ø 11 mm optic

· 2 Protection glasses (patient and operator)

· 1 Emergency button

· 1 Power supply

· 1 Operating manual

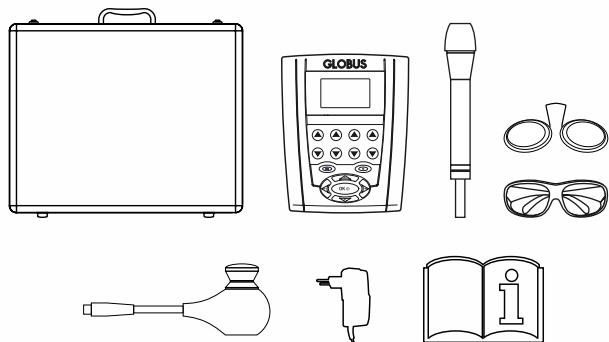

PURCHASE OPTIONS

· Code **G5994** MP 980 nm handpiece

· Code **G5996** 1064 nm handpiece

· Code **G6092** Ø 22 mm optic



STRONG POINTS



PHYSIOLASER 6.0

LASER THERAPY

PHYSIOLASER 6.0 IS A LIGHTWEIGHT PORTABLE DEVICE, SUITABLE TO ALWAYS BE CARRIED ALONG TO PERFORM TREATMENTS IN MULTIPLE CLINICS. WITH ITS MAXIMUM POWER OF 6 W, IT IS IDEAL FOR TREATING OSTEOARTICULAR AND MUSCULAR PATHOLOGIES, QUICKLY REDUCING PAIN AND EDEMA AND ACCELERATING THE RECOVERY OF INJURED AND INFLAMED TISSUES.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6 W $\pm 20\%$

Compatible with a wavelength of

808 nm
980 nm and 1064 nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic
Ø 22 mm optic (optional)

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

- Code **G5740**
- 1 Suitcase
- 1 PHYSIOLASER 6.0 laser therapy unit
- 1 Handpiece MP 808 nm with Ø 11 mm optic
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

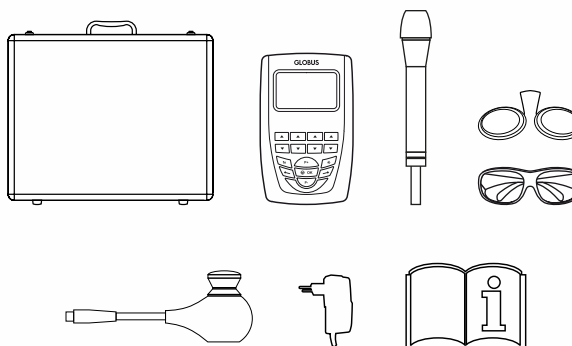

PURCHASE OPTIONS

- Cod. **G5994** MP 980 nm handpiece
- Cod. **G5996** 1064 nm handpiece
- Cod. **G6092** Ø 22 mm optic



STRONG POINTS



PHYSIOLASER 2.0

LASER THERAPY

LIGHT AND PORTABLE DEVICE, EQUIPPED WITH ALL THE MAIN PROGRAMS TO TREAT CHRONIC AND ACUTE PATHOLOGIES OF THE MUSCULOSKELETAL SYSTEM.

TECHNICAL SPECIFICATION

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

2 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code **G5739**

· 1 Suitcase

· 1 PHYSIOLASER 2.0 laser therapy unit


· 1 Handpiece Black 808 nm with Ø 11 mm optic

· 2 Protection glasses (patient and operator)

· 1 Emergency button

· 1 Power supply

· 1 Operating manual

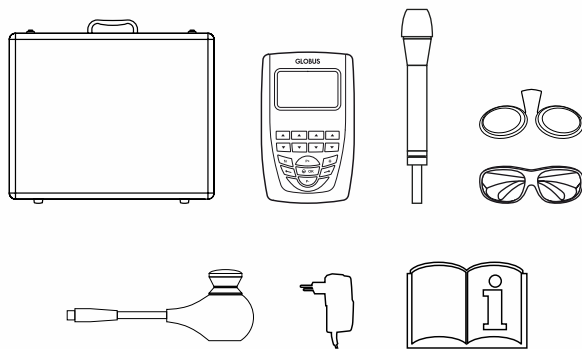

STRONG POINTS



PHYSIOLASER 1000

LASER THERAPY

PHYSIOLASER 1000 IS A PRODUCT OF THE PROFESSIONAL LINE. POWERFUL AND PORTABLE, THE PRODUCT IS INDICATED TO ACCELERATE THE RECOVERY OF INFLAMED TISSUES, REGENERATE CELLS AND IMPROVE MICROCIRCULATION.

TECHNICAL SPECIFICATION

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

1 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code G3579

· 1 Suitcase

· 1 PHYSIOLASER 1000 laser therapy unit


· 1 Handpiece Green 808 nm with Ø 11 mm optic

· 2 Protection glasses (patient and operator)

· 1 Emergency button

· 1 Power supply

· 1 Operating manual

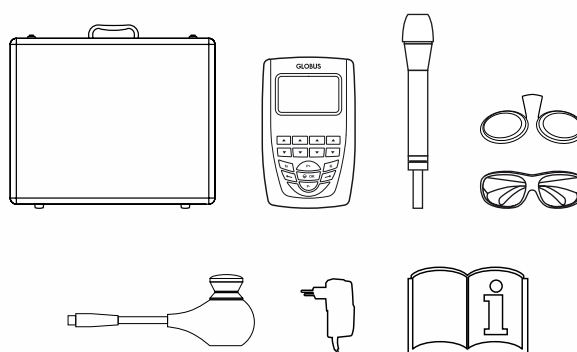

STRONG POINTS



PHYSIOLASER 500

LASER THERAPY

THE PRODUCT IS INDICATED TO ACCELERATE THE RECOVERY OF INFLAMED TISSUES, REGENERATE CELLS AND IMPROVE MICROCIRCULATION. PHYSIOLASER 500 PROMOTES THE NATURAL HEALING PROCESSES.

TECHNICAL SPECIFICATION

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

0.5 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code **G3786**

· 1 Suitcase

· 1 PHYSIOLASER 500 laser therapy unit


· 1 Handpiece White 808 nm with Ø 11 mm optic

· 2 Protection glasses (patient and operator)

· 1 Emergency button


· 1 Power supply

· 1 Operating manual

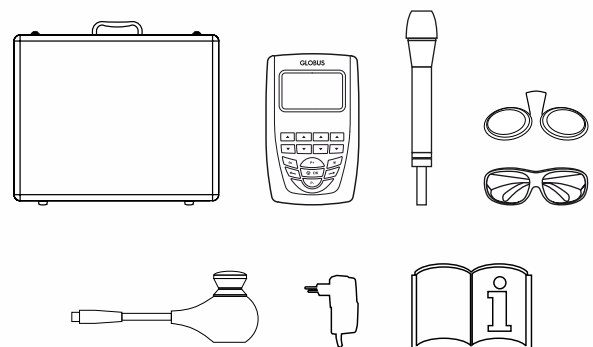

STRONG POINTS



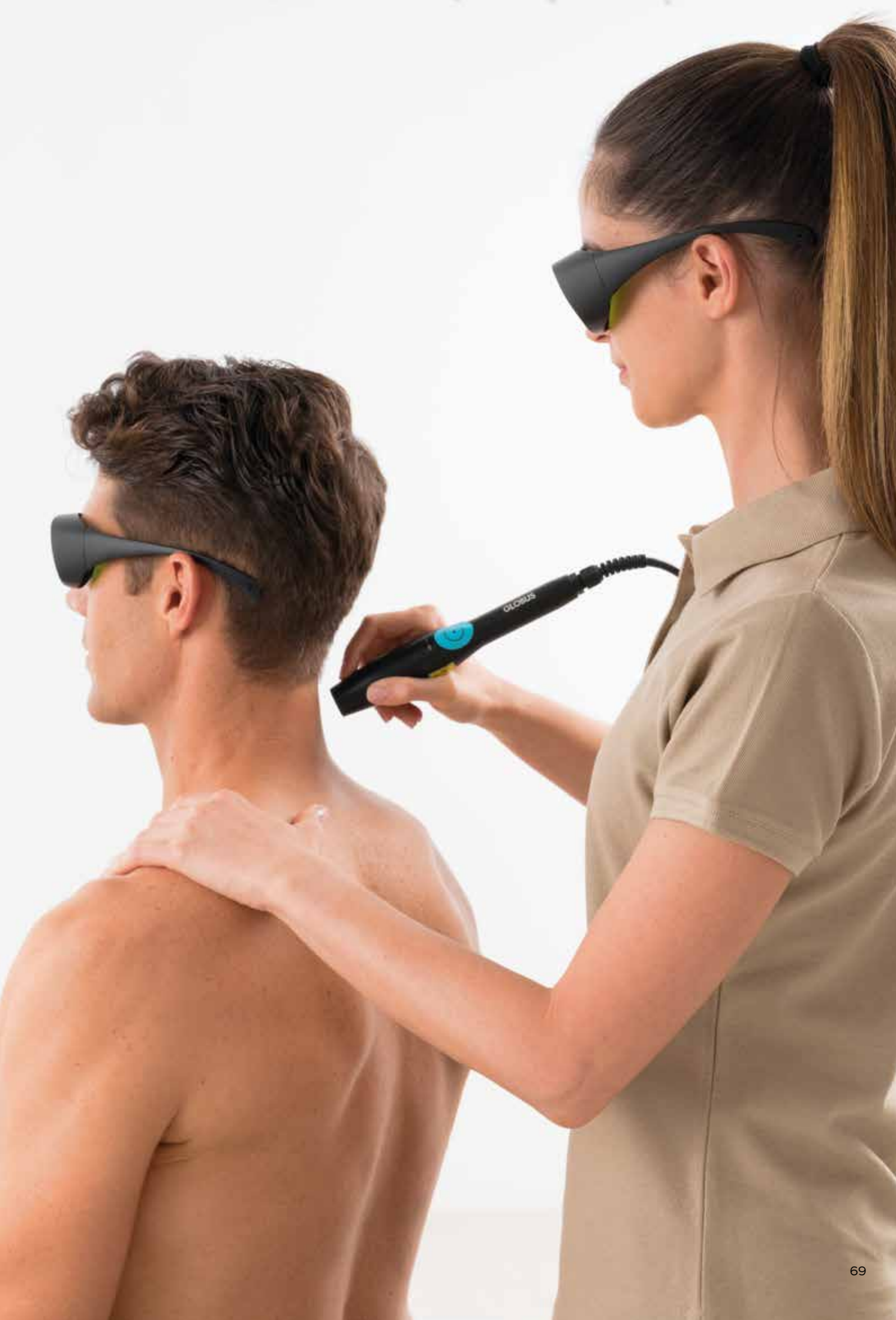

LASER THERAPY

Line for podology



PODCARE 12.0 PRO

LASER THERAPY

TOP OF THE RANGE AMONG THE LASER DEVICES FOR PODOLOGY. WITH ITS 12 W OF MAXIMUM POWER, IT ALLOWS TO PERFORM DIFFERENT PODOLOGICAL AND MUSCULOSKELETAL TREATMENTS OPTIMIZING TIME AND GUARANTEEING A DEEP ACTION.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

220x170x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

· Code **G5800**

· 1 Suitcase

· 1 PODCARE 12.0 PRO laser therapy unit


· 1 Handpiece HP 808 nm with Ø 11 mm and Ø 22 mm optics

· 2 Protection glasses (patient and operator)

· 1 Emergency button

· 1 Power supply

· 1 Operating manual

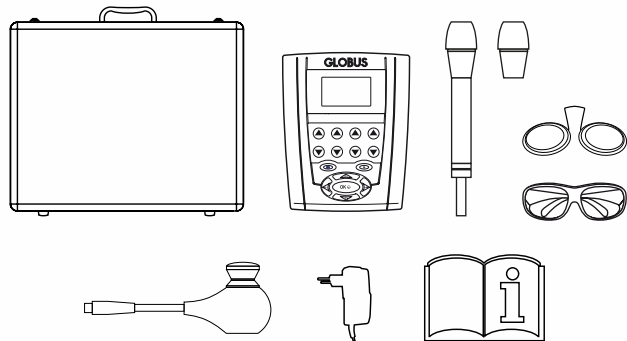

PURCHASE OPTIONS

· Code **G5995** HP 980 nm handpiece

· Code **G5996** 1064 nm handpiece



STRONG POINTS



PODCARE 6.0 PRO

LASER THERAPY

PODCARE 6.0 PRO IS A DEVICE OF SMALL WEIGHT AND SIZE, SUITABLE BOTH TO BE PLACED ON A STUDIO CART AND TO BE TRANSPORTED, IF NECESSARY, TO OTHER LOCATIONS. THE DEVICE IS IDEAL FOR PODOLOGICAL APPLICATIONS ON THE FOOT AND THE ANKLE, BUT IT ALSO HAS A SERIES OF SPECIFIC PROGRAMS FOR ULCERS, WOUNDS AND MUSCULOSKELETAL DISORDERS.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ottica, Ø 11 mm optic

Ottica, Ø 22 mm optic (optional)

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


220x170x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

- Code G5762
- 1 Suitcase
- 1 PODCARE 6.0 PRO laser therapy unit
- 1 Handpiece MP 808 nm with Ø 11 mm optic
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

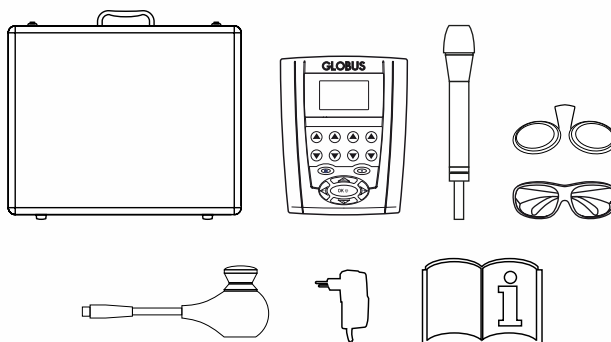

OPZIONI DI ACQUISTO

- Code G5994 MP 980 nm handpiece
- Code G5996 1064 nm handpiece
- Code G6092 Ø 22 mm optic



STRONG POINTS



PODCARE 6.0

LASER THERAPY

PODCARE 6.0 IS A LIGHTWEIGHT AND PORTABLE DEVICE, SUITABLE TO ALWAYS BE CARRIED ALONG TO PERFORM TREATMENTS IN MULTIPLE CLINICS. IT IS EQUIPPED WITH A 808NM DIODE WITH A MAXIMUM POWER OF 6W. THE DEVICE IS IDEAL FOR PODOLOGICAL APPLICATIONS ON THE FOOT AND THE ANKLE, BUT IT ALSO HAS A SERIES OF SPECIFIC PROGRAMS FOR ULCERS, WOUNDS AND MUSCULOSKELETAL DISORDERS.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic ((optional))

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

- Code **G5737**
- 1 Suitcase
- 1 PODCARE 6.0 laser therapy unit
- 1 Handpiece MP 808 nm with Ø 11 mm optic
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

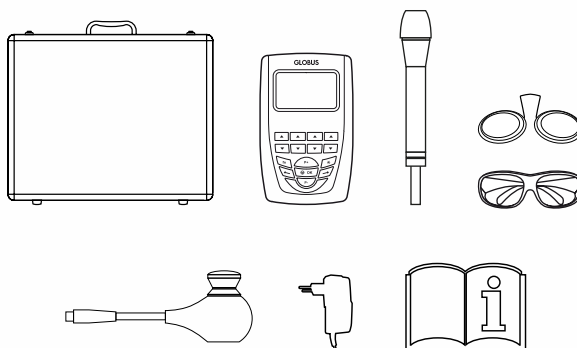

PURCHASE OPTIONS

- Code **G5994** MP 980 nm handpiece
- Code **G5996** 1064 nm handpiece
- Code **G6092** Ø 22 mm optic



CE
0476



STRONG POINTS



PODCARE 2.0

LASER THERAPY

PODCARE 2.0 IS A PORTABLE DEVICE WITH PROTOCOLS TO TREAT DIFFERENT MUSCULOSKELETAL PATHOLOGIES AND WITH SPECIAL PODOLOGICAL APPLICATIONS SUCH AS ONICOMYCOSIS, ONICOCRIPTOSIS AND OTHER TYPICAL INFLAMMATORY PATHOLOGIES OF THE FOOT AND THE LEG.

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

2 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic ((optional))

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

- Code **G5736**
- 1 Suitcase
- 1 PODCARE 2.0 laser therapy unit
- 1 Handpiece Black 808 nm with Ø 11 mm
- 2 Protection glasses (patient and operator)
- 1 Emergency button
- 1 Power supply
- 1 Operating manual

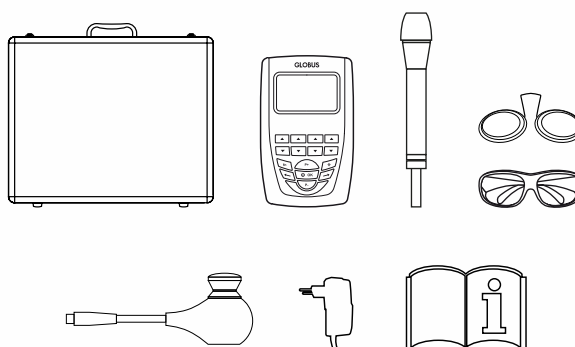

PURCHASE OPTIONS

- Code **G5547** Ø 3 mm Small optic



STRONG POINTS



PRESSOTHERAPY

Line



PRESSCARE G300M

PRESSOTHERAPY | 2 OUTPUTS - 3 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, AESTHETICS AND RELAXATION NEEDS. SUITABLE FOR LYMPHEDEMA, VENOUS CIRCULATION PROBLEMS, HEAVY LEGS, EDEMA, SWOLLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, MASSAGES, HEMATOMAS, VENOUS STASIS. IT CAN ALSO BE USED WITH THREE ACCESSORIES SIMULTANEOUSLY. IT IS RECOMMENDED TO CONSULT YOUR DOCTOR BEFORE USE.

TECHNICAL SPECIFICATIONS

Outputs

2

Programs

3

Pressure

0-180 mmHg $\pm 10\%$

Timer

Up to 99 minutes

Mains power supply

230 V AC, 50 Hz

Size

185x180x230 mm

Weight

3 Kg




CE
0598

STANDARD EQUIPMENT

· Depending on the version

· 1 Power supply

· 1 Operating manual

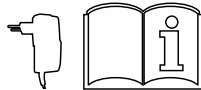

PURCHASE OPTIONS

· Code **G5411** PRESSCARE G300M-1, with 1 boot size "S"

· Code **G5412** PRESSCARE G300M-1, with 1 boot size "M"

· Code **G5413** PRESSCARE G300M-1, with 1 boot size "L"

· Code **G5414** PRESSCARE G300M-2, with 2 boots size "S"

· Code **G5415** PRESSCARE G300M-2, with 2 boots size "M"


· Code **G5416** PRESSCARE G300M-2, with 2 boots size "L"

· Code **G5417** PRESSCARE G300M-3, with 2 boots size "S",
1 abs/glutes band

· Code **G5418** PRESSCARE G300M-3, with 2 boots size "M",
1 abs/glutes band

· Code **G5419** PRESSCARE G300M-3, with 2 boots size "L",
1 abs/glutes band

· Code **G5420** PRESSCARE G300M-1B, , with 1 sleeve

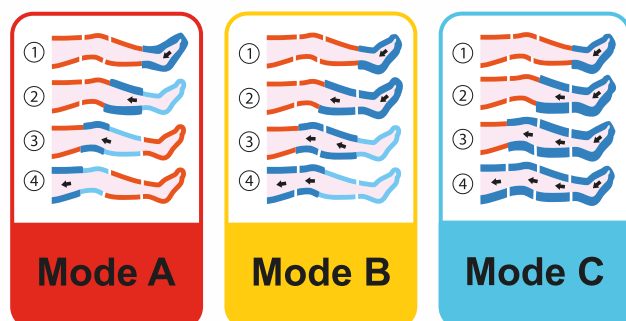

MODE A: the program sequence inflates and deflates the boot one sector at a time, from the most distal one (foot or hand) to the most proximal one.

MODE B: the program sequence inflates two chambers at a time, simulating the function of the boots with overlapping chambers.

MODE C: the program sequence carries out the pressure starting from the foot and keeping the previous chambers inflated until the entire part (leg, arm or abdomen) is inflated.

PRESSCARE G200M

PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, AESTHETICS AND RELAXATION NEEDS. SUITABLE FOR LYMPHEDEMA, VENOUS CIRCULATION PROBLEMS, HEAVY LEGS, EDEMA, SWOLLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, MASSAGES, HEMATOMAS, VENOUS STASIS. IT CAN ALSO BE USED WITH THREE ACCESSORIES SIMULTANEOUSLY. IT IS RECOMMENDED TO CONSULT YOUR DOCTOR BEFORE USE.

TECHNICAL SPECIFICATIONS

Outputs
2
Programs
2
Pressure
0-180 mmHg $\pm 10\%$
Timer
Up to 99 minutes
Mains power supply
230 V AC, 50 Hz
Size
185x180x230 mm
Weight
3 Kg



CE
0598


STANDARD EQUIPMENT

- Depending on the version
- 1 Power supply
- 1 Operating manual



PURCHASE OPTIONS

- Code **G5401** PRESSCARE G200M-1, with 1 boot size "S"
- Code **G5402** PRESSCARE G200M-1, with 1 boot size "M"
- Code **G5403** PRESSCARE G200M-1, with 1 boot size "L"
- Code **G5404** PRESSCARE G200M-2, with 2 boots size "S"
- Code **G5405** PRESSCARE G200M-2, with 2 boots size "M"
- Code **G5406** PRESSCARE G200M-2, with 2 boots size "L"
- Code **G5407** PRESSCARE G200M-3, with 2 boots size "S", 1 abs/glutes band
- Code **G5408** PRESSCARE G200M-3, with 2 boots size "M", 1 abs/glutes band
- Code **G5409** PRESSCARE G200M-3, with 2 boots size "L", 1 abs/glutes band
- Code **G5410** PRESSCARE G200M-1B, with 1 sleeve

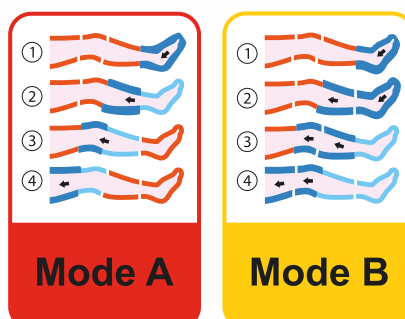

MODE A: the program sequence inflates and deflates the boot one sector at a time, from the most distal one (foot or hand) to the most proximal one.

MODE B: the program sequence inflates two chambers at a time, simulating the function of the boots with overlapping chambers.

PREVENTION



SANI-CARE OZONATORS

THE GLOBUS SANI-CARE LINE INCLUDES THREE MODELS OF OZONE GENERATORS FOR THE SANITIZATION AND PURIFICATION OF THE AIR IN CLOSED ENVIRONMENTS, TO INACTIVATE BACTERIA AND MAKE THE ENVIRONMENT HEALTHIER AND FRESHER. THE DEVICES ARE PARTICULARLY SUITABLE FOR ELIMINATING ODORS CAUSED BY PETS, WHICH CAN INEVITABLY PERVADE THE PLACES WHERE THE PETS LIVE AND ARE TRANSPORTED (FOR EXAMPLE CAR SEATS). THE RANGE INCLUDES THE SANI-CARE 3 "HOME-CAR" MODEL, RECOMMENDED FOR MOTOR VEHICLES AND SMALL AND MEDIUM SIZED CLOSED ENVIRONMENTS, AND THE TWO MODELS SANI-CARE 10 AND SANI-CARE 20, SUITABLE FOR MEDIUM AND LARGE CLOSED ENVIRONMENTS.

SANI-CARE 20

This ozone generator is suitable for sanitizing large environments (200-300 m²).

The device is suitable for sanitizing environments from the microorganisms in the air and on surfaces as well as to eliminate dust and bad smell.

The control panel allows to set the working time in order to program the sanitization based on the size of the room and a convenient function allows to manage the delayed start.

TECHNICAL SPECIFICATIONS

Code G6030

Ozone production

20 g/H

Surface covered

200-300 m²

Power

330 W

Power supply

Mains

Size

430(L) x 230(D) x 230(H) mm

Weight

7 Kg

SANI-CARE 10

Ozone generator suitable for rooms with a surface between 100-150 m².

The device is suitable for sanitizing environments from the microorganisms in the air and on surfaces as well as to eliminate dust and bad smell.

The control panel allows to set the working time in order to program the sanitization based on the size of the room and a convenient function allows to manage the delayed start.

TECHNICAL SPECIFICATIONS

Code G6031

Ozone production

10 g/H

Surface covered

100-150 m²

Power

100 W

Power supply

Mains

Size

260(L) x 170(D) x 300(H) mm

Weight

3.3 Kg

SANI-CARE 3

This device is an ozone and negative ion generator ideal for small environments such as cars, vans, buses and RVs.

The double function ozone + negative ions is suitable for the sanitization from the microorganisms in the air and on surfaces as well as to eliminate bad smell like cigarette smoke.

TECHNICAL SPECIFICATIONS

Code G6032

Ozone production

3 g/H

Surface covered

50 m²

Power

78 W

Power supply

DC 12 V and mains

Size

260(L) x 150(D) x 210(H) mm

Weight

2.6 Kg



PULSE OXIMETER - YM201

WITH ALARM FOR SATURATION AND PULSES

THE "YM201" PULSE OXIMETER ENABLES AN EASY EVALUATION OF THE OXYGEN SATURATION IN THE BLOOD AND OF THE HEART RATE OF A SUBJECT. IT IS A USER-FRIENDLY, PAINLESS AND ABSOLUTELY NON-INVASIVE TOOL, SINCE THE TEST IS MADE ON A FINGER BY MEANS OF AN INFRARED RADIATION.

USING THE PULSE OXIMETER CAN BE USEFUL TO EVALUATE AND MONITOR BREATHING FUNCTIONS IN HOSPITALS, CLINICS AND AT HOME.

THE TOOL IS PARTICULARLY USEFUL FOR PATIENTS SUFFERING FROM HEART FAILURE OR RESPIRATORY DISEASES LIKE: PNEUMONIA - BRONCHIAL ASTHMA - CHRONIC BRONCHITIS - COPD - OTHER PULMONARY AND PLEURAL DISEASES.

TECHNICAL SPECIFICATION

Display

Color OLED 0.96"

Oxygen saturation (SpO2)

Measurement range: 0~99%

Accuracy: $\pm 2\%$ (70%~100%)

Pulse rate

Measurement range: 25~250 bpm

Accuracy: ± 3 bpm

Perfusion index (PI)

0.1~20%

Power supply

2 Alkaline batteries AAA (LR03) included

Size

57x30x31 mm

Weight

54 gr

STANDARD EQUIPMENT

· Code **G6102**

· 1 YM201 PULSE OXIMETER unit

· 1 Holding cord

· 2 Alkaline batteries AAA (LR03)

· 1 Operating manual



CERTIFIED MEDICAL DEVICE

DEVICE FUNCTIONS

Detection and indication of the oxygen saturation in the blood (SpO2)

Pulse rate detection (heart rate)

Indication of the perfusion index (PI)

Accuracy and high data stability

Beat histogram

Battery charge indicator

Suitable for adults and children



BEAUTY & WELLNESS

Line

CONTENTS

- TECAR BEAUTY Line	86
- ELECTROSTIMULATION Line	94
- RADIOFREQUENCY Aesthetic line	106
- ULTRASOUND Aesthetic line	112
- VIBRATION PLATFORM Aesthetic line	120

TECAR BEAUTY

Line



TECAR BEAUTY 7000 MED

AESTHETIC TECAR THERAPY

THE BEST DEVICE WITH MEDICAL CERTIFICATION IN THE GLOBUS TECAR BEAUTY LINE. THE DEVICE IS EQUIPPED WITH A WIDE RANGE OF PROGRAMS INCLUDING BOTH AESTHETIC TREATMENTS AND PROTOCOLS FOR DIFFERENT MUSCULOSKELETAL PATHOLOGIES. THE TOUCHSCREEN NAVIGATION MENU ON A COLOR SCREEN INCLUDES AN ATLAS WITH SEVERAL PICTURES AND A PRACTICAL DESCRIPTION OF HOW TO PERFORM THE TREATMENTS. MOREOVER, TECAR BEAUTY 7000 MED ALLOWS EXPERT USERS TO CREATE AND STORE NEW PROGRAMS IN JUST A FEW STEPS.

TECHNICAL SPECIFICATIONS

Display

Size 7" color touch screen

Nominal Power (NP)

350 W $\pm 10\%$

Emission frequencies

400 kHz, 448 kHz, 470 kHz, 700 kHz, 1000 kHz, 1200 kHz, $\pm 10\%$

Processing methods

Capacitive and Resistive

Capacitive Electrodes

3 \varnothing (30, 50, 70 mm)

Resistive Electrodes

3 \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1.8 A

Size

515x340x170 mm

Weight

11.5 Kg

STANDARD EQUIPMENT

· Code **G6244**

· 1 TECAR BEAUTY 7000 MED unit

· 1 Counter electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Bipolar handpiece \varnothing (30 mm)

· 1 Hands-free treatment kit

· 1 1000 ml conductive cream bottle

· 1 Power supply cable

· 1 Operating manual



STRONG POINTS



TECAR BEAUTY 7000

AESTHETIC TECAR THERAPY

POWERFUL AND VERSATILE DEVICE FROM THE TECAR BEAUTY RANGE. NAVIGATION SYSTEM WITH COLOR TOUCHSCREEN TO MAKE TREATMENT SEARCH AND PARAMETER MODIFICATION EASIER, AN ATLAS THAT DESCRIBES THE TREATMENTS DIVIDED INTO DIFFERENT PARTS OF THE BODY WITH COLOR PHOTOS AND A SECTION WHERE THE OPERATOR CAN STRUCTURE AND STORE PERSONALIZED TREATMENTS.

THE BEST EXPRESSION OF TECAR BEAUTY FOR BEAUTY TREATMENTS!

TECHNICAL SPECIFICATIONS

Display

Size 7" color touch screen

Power

350 W $\pm 10\%$

Emission frequencies

400 kHz, 448 kHz, 470 kHz, 700 kHz, 1000 kHz, 1200 kHz, $\pm 10\%$

Processing methods

Capacitive and Resistive

Capacitive Electrodes

3 \varnothing (30, 50, 70 mm)

Resistive Electrodes

3 \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1.8 A

Size

515x340x170 mm

Weight

11.5 Kg

STANDARD EQUIPMENT

· Cod. G6248

· 1 TECAR BEAUTY 7000 unit

· 1 Counter electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Bipolar handpiece \varnothing (30 mm)

· 1 Hands-free treatment kit

· 1 500 ml hyaluronic acid conductive face cream


· 1 1000 ml conductive cream for lymphatic drainage with eascin, fucus and caffeine

· 1 Power supply cable

· 1 Operating manual



STRONG POINTS

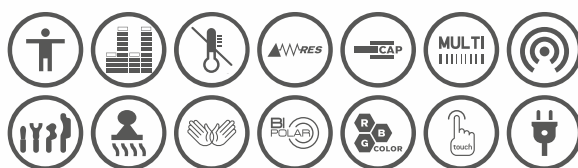

TECAR BEAUTY ELITE

AESTHETIC TECAR THERAPY

TECAR BEAUTY IS A NON-INVASIVE METHOD THAT FIGHTS AGING AND SKIN LAXITY. THE ELECTROMAGNETIC WAVES PRODUCED BY THE DEVICE PENETRATE INTO THE DERMIS AND RELEASE HEAT IN THE TISSUES. THIS STIMULATES THE SYNTHESIS OF NEW COLLAGEN AND ELASTIN, VISIBLY IMPROVING THE COMPACTNESS OF THE SKIN. IN ADDITION, THE HEAT FACILITATES MICROCIRCULATION, PROMOTING A BETTER OXYGENATION OF THE TISSUES. PROFESSIONAL DEVICE SUITABLE FOR AESTHETIC CENTERS.

TECHNICAL SPECIFICATIONS

Display

Size 7" color touch screen

Power

50 W $\pm 10\%$ Capacitive

25 W $\pm 10\%$ Resistive

Emission frequencies

400 kHz, 448 kHz, 470 kHz, 700 kHz,
1000 kHz, 1200 kHz, $\pm 10\%$

Processing methods

Capacitive and Resistive

Capacitive Electrodes

3 \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1,8 A

Size

515x340x170 mm

Weight

11.5 Kg

STANDARD EQUIPMENT

· Code **G6243**

· 1 RF BEAUTY ELITE unit

· 1 Counter electrode (flexible plate)

· 1 Monopolar capacitive handpiece with 3 electrodes
 \varnothing (30, 50, 70 mm)

· 1 Bipolar handpiece \varnothing (30 mm)

· 1 500 ml hyaluronic acid conductive face cream

· 1 1000 ml conductive cream for lymphatic drainage
with eascin, fucus and caffeine

· 1 Power supply cable

· 1 Operating manual



In compliance with the Interministerial Decree n. 206/2015, the device is suitable for use in beauty centers.

STRONG POINTS



TECAR BEAUTY 6000 MED

TECAR BEAUTY 6000 MED RE

AESTHETIC TECAR THERAPY

TECAR BEAUTY 6000 MED AND TECAR BEAUTY 6000 MED RE ARE PORTABLE DEVICES FOR CAPACITIVE AND RESISTIVE TRANSFER TECAR THERAPY. EQUIPPED WITH DIFFERENT PRE-SET PROGRAMS TO TREAT THE MAIN SKIN AND SUBCUTANEOUS TISSUE IMPERFECTIONS, THEY ALSO INCLUDE PROGRAMS TO TREAT THE MAIN MUSCULOSKELETAL PATHOLOGIES.

27 BEAUTY | 5 MEDICAL

TECHNICAL SPECIFICATIONS

Display

Size 4.3" color touch screen

Nominal Power (NP)

180 W $\pm 10\%$

Emission frequencies

400 kHz, 448 kHz, 470 kHz,
700 kHz, 1000 kHz, $\pm 10\%$

Processing methods

Capacitive and Resistive

Capacitive Electrodes

3 \varnothing (30, 50, 70 mm)

Resistive Electrodes

3 \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1.8 A

Size


170x220x60 mm

Weight

1.04 Kg

STANDARD EQUIPMENT

- 1 Suitcase
- 1 TECAR BEAUTY 6000 MED unit
- 1 Counter electrode (flexible plate)
- 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 Bipolar handpiece \varnothing (30 mm)
- 1 250 ml conductive cream bottle
- 1 Power supply
- 1 Operating manual

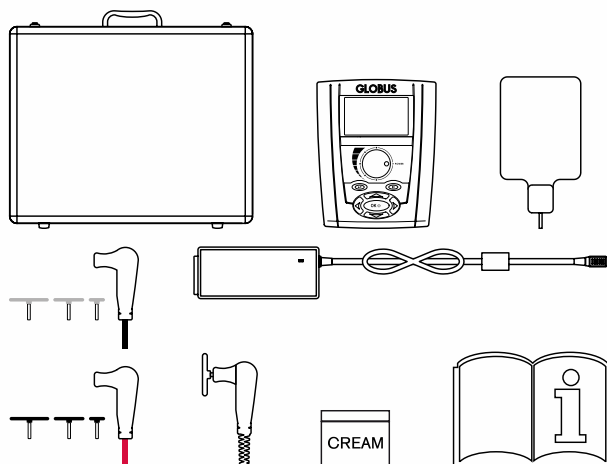

PURCHASE OPTIONS

- Code **G6246** TECAR BEAUTY 6000 MED
- Code **G6247** TECAR BEAUTY 6000 MED RE*

* RECHARGEABLE TECAR BEAUTY:

- 5 HOURS OF TREATMENT INCLUDED IN THE COST OF THE DEVICE
- CHEAP CHARGING
- GUARANTEED SUPPORT FOR 5 YEARS AND FORKLIFT
- CHARGING VIA INTERNET



CE
0476

STRONG POINTS



FASCIA TOOLS

FOR TECAR BEAUTY 7000 MED, 7000, ELITE AND 6000 MED

The anatomical shape of the special resistive handpieces makes the treatment particularly effective on thighs, gluteus, hips and arms. The combination of the manual treatment with the myofascial handpieces, together with the energy of TECAR BEAUTY, offers an extraordinary and new treatment possibility.

TECAR BEAUTY with the special resistive handpieces allows to:

- **MOBILIZE THE DENSIFICATION OF THE CONNECTIVE TISSUE**
- **VASCULARIZE AND TREAT FIBROUS AND SCLEROTIC NODULES**
- **WORK WITH MAXIMUM ERGONOMICS AND TREATMENT PRECISION**



Code **G6041** KIT FASCIA TOOLS BEAUTY



DiaBlade Large: handpiece for the treatment of the Lower limb, in particular the Thighs and the Gluteus.



DiaBlade Small: handpiece for the treatment of the Upper limb and the Neck area.



ELECTROSTIMULATION

Line





SHAPE & BEAUTY

4-CHANNEL ELECTROSTIMULATION

SHAPE & BEAUTY IS A PROFESSIONAL DEVICE FOR TONING AND FIRMING THE BODY TO SHAPE THE CURVES AND IMPROVE THE SHAPE OF THE BODY IN JUST A FEW SESSIONS.

THE PRESET PROTOCOLS INSIDE THE DEVICE WILL GUIDE THE CHOICE OF THE MOST SUITABLE APPLICATION ACCORDING TO THE GOAL, FOLLOWING A WORK PROGRESSION THAT WILL ALLOW TO OBTAIN VISIBLE RESULTS AFTER A FEW WEEKS.

THE DEVICE ALSO HAS RELAXING AND DECONTRACTING MASSAGE PROGRAMS AND VASCULARIZATION AND CAPILLARIZATION PROGRAMS, TO IMPROVE THE DRAINAGE OF EXCESS LIQUIDS.

THE G-TRODE HANDPIECE INCLUDED HAS A VERY POWERFUL SKIN REGENERATING AND BIOSTIMULATING EFFECT, ACTING ON EXPRESSION WRINKLES, STRETCH MARKS AND SMALL ACNE SCARS. SHAPE & BEAUTY IS THE PERFECT ALLY FOR YOUR BEAUTY CLINIC!

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

170x220x60 mm

Weight

1.12 Kg

* Based on the type of current

STANDARD EQUIPMENT

· Code G6230

· 1 Suitcase

· 1 SHAPE & BEAUTY stimulator

· 1 G-Trode handpiece - Skin regeneration

· 4 Cables for electrode connection


· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Gel bottle

· 1 Battery charger

· 1 Operating manual

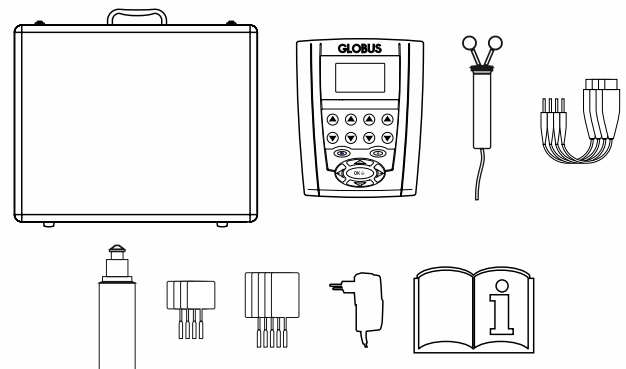

CE



STRONG POINTS



ACTIVA 700

4-CHANNEL ELECTROSTIMULATION | 271 PROGRAMS

PRODUCT FOR BEAUTY AND PHYSICAL SHAPE, DESIGNED WITH SPECIFIC PROGRAMS FOR FIRMING, TONING AND REJUVENATION. THE SPECIAL G-TRODE HANDPIECE (INCLUDED) MAKES ACTIVA 700 A TOP-OF-THE-RANGE IN THE ELECTROSTIMULATION LINE.

60 BEAUTY | 58 FITNESS | 53 SPORT | 36 SERIAL SEQUENTIAL STIMULATION | 21 SPECIAL SPORTS | 13 WRINKLES-SKIN
12 MICROCURRENTS | 12 TENS | 3 REHAB | 3 INCONTINENCE

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G3550

· 1 Suitcase

· 1 ACTIVA 700 stimulator

· 1 G-Trode handpiece

· 4 Cables for electrode connection

· 2 Cables for microcurrents and ionophoresis


· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Gel bottle

· 1 Battery charger


· 1 Operating manual

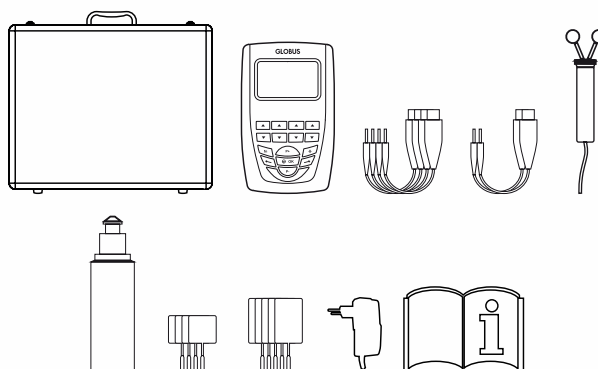

CE
0476



STRONG POINTS

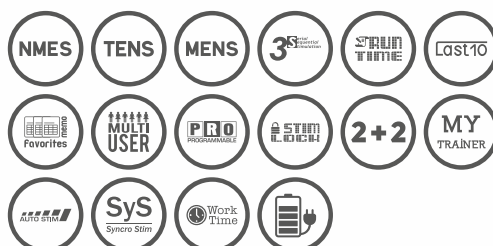

PREMIUM 400

4-CHANNEL ELECTROSTIMULATION | 258 PROGRAMS

MORE THAN 250 PROGRAMS ON THE PALM OF ONE HAND. PREMIUM 400 IS THE IDEAL ELECTROSTIMULATOR FOR SPORT AND FOR THE CARE OF YOUR BEAUTY, FOR HIM AND FOR HER.

58 FITNESS | 53 SPORT | 44 SPECIAL SPORTS | 36 BEAUTY | 23 MICROCURRENTS | 18 SERIAL SEQUENTIAL STIMULATION
12 TENS AND PAIN | 7 ACTION NOW | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G3551

· 1 Bag

· 1 PREMIUM 400 stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

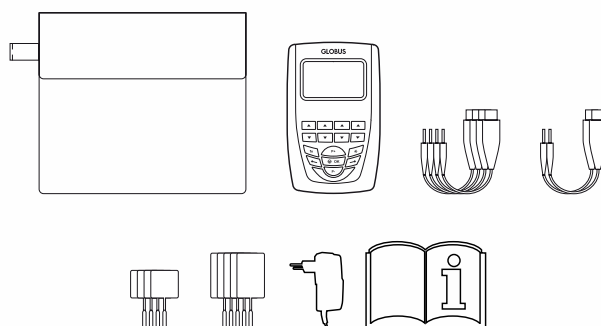

STRONG POINTS

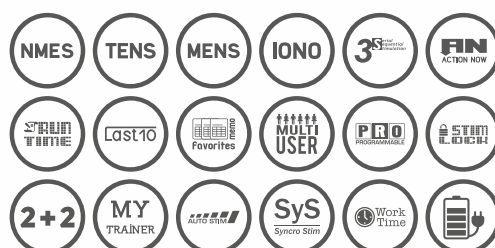

ELITE 150

4-CHANNEL ELECTROSTIMULATION | 150 PROGRAMS

PRODUCT DESIGNED FOR SPORT AND FITNESS, WITH SPECIAL ATTENTION TO BEAUTY TREATMENTS THANKS TO THE NEW SEQUENTIAL STIMULATION CURRENTS, EXCELLENT FOR THEIR DRAINING EFFECTS.

36 BEAUTY | 36 SERIAL SEQUENTIAL STIMULATION | 32 FITNESS | 31 SPORT | 12 TENS AND PAIN | 3 REHAB

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G1416

· 1 Bag

· 1 ELITE 150 stimulator



· 4 Cables for electrode connection

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

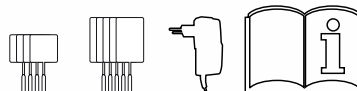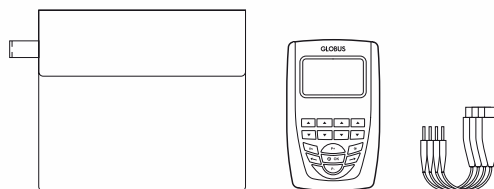

STRONG POINTS



ELITE

4-CHANNEL ELECTROSTIMULATION | 98 PROGRAMS

FOUR-CHANNEL PRODUCT FOR EVERY TREATMENT NEED.

32 SPORT | 29 BEAUTY | 26 FITNESS | 9 TENS | 2 REHAB

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

440 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code **G4300**

· 1 Bag

· 1 ELITE stimulator


· 4 Cables for electrode connection

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

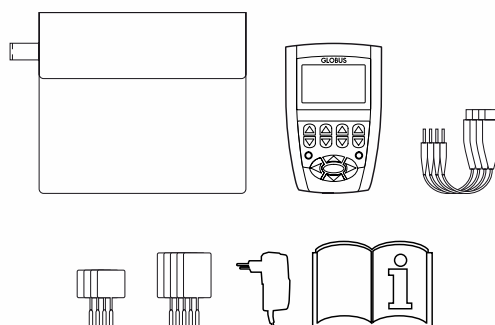

CE
0476



STRONG POINTS



ELITE S2

2-CHANNEL ELECTROSTIMULATION | 100 PROGRAMS

THE TWO-CHANNEL PRODUCT FOR EVERY TREATMENT NEED.

44 FITNESS, BEAUTY, WELLNESS | 30 SPORT | 16 TENS | 6 PREVENTION | 3 REHAB | 1 INCONTINENCE

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

2 independent (4 electrodes)

Frequency

1-150 Hz

Pulse width

50-400 μ s

Intensity

0-100 mA per channel

Power supply

Rechargeable batteries

Size

80x130x25 mm

Weight

220 gr

STANDARD EQUIPMENT

· Code G3552

· 1 Bag

· 1 ELITE SII stimulator



· 2 Cables for electrode connection

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

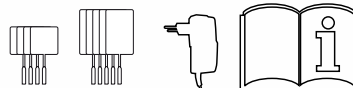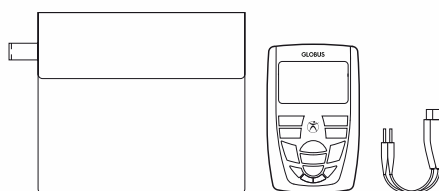

STRONG POINTS





FAST BAND

ABDOMINAL BAND

FAST BAND IS A BAND THAT IS EASILY CONNECTED TO ALL 2- AND 4-CHANNEL GLOBUS ELECTROSTIMULATORS. WITH THE FAST BAND, WITHOUT USING THE CLASSIC ELECTRODES, YOU CAN STIMULATE THE ABDOMINAL AND OBLIQUE MUSCLES, THE LUMBAR AREA AND THE GLUTES SIMPLY AND EFFECTIVELY.

TECHNICAL SPECIFICATIONS

Size

980x150 mm

Weight


214 gr

STANDARD EQUIPMENT

· Code **G0487**


· 1 FAST BAND

· 4 Clip adapters

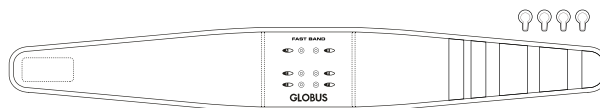

PURCHASE OPTIONS

· Code **G0488** FAST BODY KIT: "FAST BAND" + "FAST PAD"

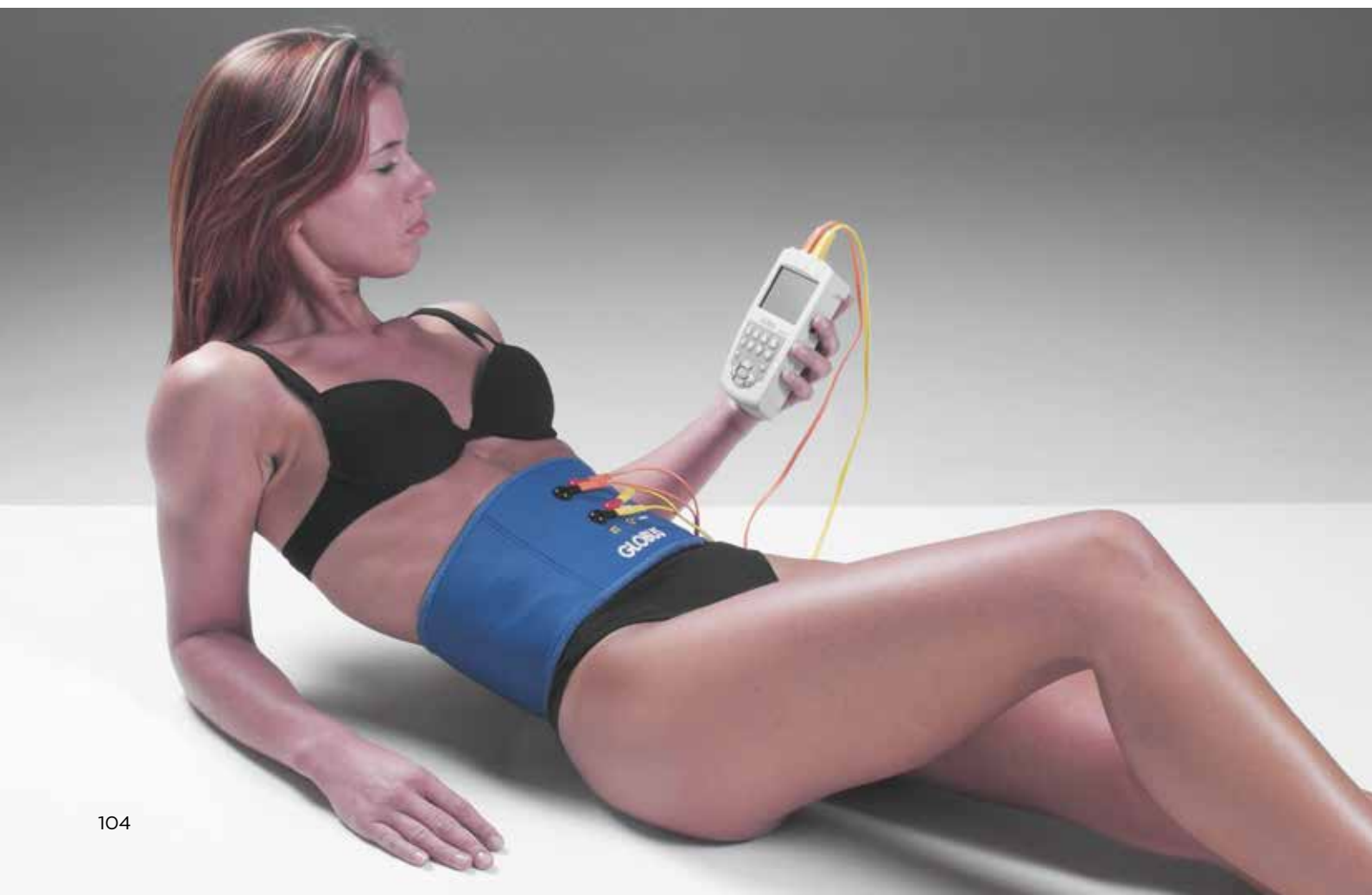

FAST PAD

LEG BAND

FAST PADS ARE BANDS THAT ARE EASILY CONNECTED TO ALL 2- AND 4-CHANNEL GLOBUS ELECTROSTIMULATORS. THEY ARE IDEAL TO PERFORM SPECIFIC STIMULATION TREATMENTS FOR THE LEGS IN A COMFORTABLE AND EFFECTIVE WAY. THEY ARE ALSO CHARACTERIZED BY AN OVAL SHAPE WITH LARGE CONTACT SURFACE AND ARE CONNECTED TO THE ELECTROSTIMULATOR THROUGH CLIP ADAPTERS.

TECHNICAL SPECIFICATION

Size

250x150 mm

Weight

60 gr (single band)


STANDARD EQUIPMENT

· Code **G0489**

· 2 FAST PADS

· 4 Clip adapters

· Gel

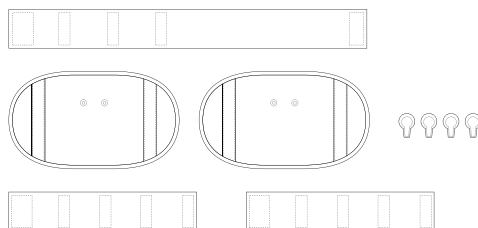

PURCHASE OPTIONS

· Cod. **G0488** FAST BODY KIT: "FAST BAND" + "FAST PAD"



Aesthetic line **RADIOFREQUENCY**



RF CLINIC PRO

RADIOFREQUENCY | 30 PROGRAMS

RF CLINIC PRO IS AN INNOVATIVE PROFESSIONAL DEVICE THAT GENERATES RADIOFREQUENCY WAVES FOR AESTHETIC APPLICATIONS, AGAINST THE AGING OF THE SKIN OF BOTH THE FACE AND THE WHOLE BODY. THANKS TO THE 50 MM BIPOLAR HANDPIECE IT IS POSSIBLE TO PERFORM FIRMING AND CELLULITE PROGRAMS.

13 SKIN | 10 CELLULITE | 7 FACE

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

- 1 RF Clinic 100, Ø 50 mm handpiece
- 1 RF Clinic 100 2-sphere handpiece

Power

200 W (50 Ω)

Emitted frequency

470 KHz to 2 MHz

Processing methods

Resistive

Power supply

- Mains
- Rechargeable battery (not included)

Size

170x220x60 mm

Weight

920 gr

STANDARD EQUIPMENT

· Code G1506

· 1 Suitcase

· 1 RF CLINIC PRO radiofrequency unit

· 1 RF Clinic 500 handpiece


· 1 RF Clinic 100 handpiece

· 1 50 ml hyaluronic acid conductive face cream

· 1 50 ml conductive cream for lymphatic drainage with eascin, fucus and caffeine

· 1 Battery charger

· 1 Operating manual

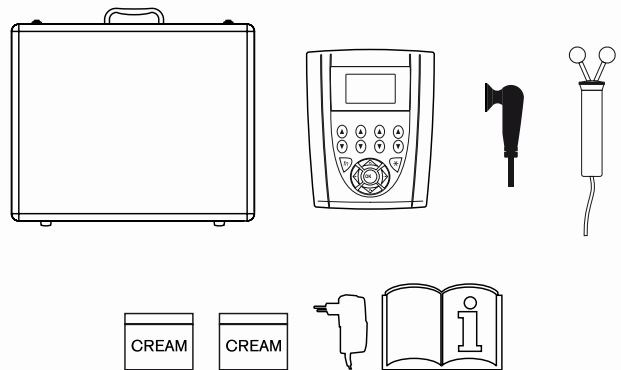

CE



STRONG POINTS



RF CLINIC BODY

RADIOFREQUENCY | 10 PROGRAMS

RF CLINIC BODY IS AN INNOVATIVE DEVICE THAT GENERATES RADIOFREQUENCY WAVES FOR AESTHETIC APPLICATIONS, AGAINST THE AGING OF THE SKIN OF BOTH THE FACE AND THE WHOLE BODY. IT ALSO FACILITATES THE VASCULARIZATION OF THE SKIN TISSUE, ELIMINATING THE IMPERFECTIONS CAUSED BY THE EXPOSURE TO SUN RAYS AND THE SKIN AGING.

6 SKIN | 4 FACE

TECHNICAL SPECIFICATIONSE

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

1 RF Clinic 100 2-sphere handpiece

Power

200 W (50 Ω)

Emitted frequency

470 KHz to 2 MHz

Processing methods

Resistive

Power supply

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

920 gr

STANDARD EQUIPMENT

- Code G1507
- 1 Suitcase
- 1 RF CLINIC BODY radiofrequency unit
- 1 RF Clinic 100 handpiece
- 1 50 ml hyaluronic acid conductive face cream
- 1 50 ml conductive cream for lymphatic drainage with eascin, fucus and caffeine
- 1 Power supply
- 1 Operating manual

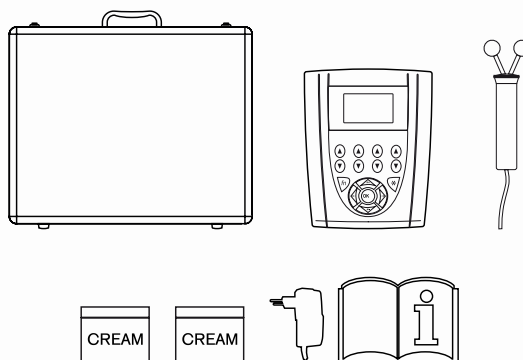

CE



STRONG POINTS



RF BEAUTY MINI

RADIOFREQUENCY | 3 PROGRAMS

ELEGANT AND ERGONOMIC, RF BEAUTY MINI IS THE PORTABLE RADIOFREQUENCY DEVICE FOR AESTHETIC APPLICATIONS AND AGAINST THE AGING OF THE SKIN OF THE FACE. EASY AND READY TO USE, IT IS INDICATED FOR THE REDUCTION OF WRINKLES AND FOR THE REJUVENATION OF THE SKIN AND OF THE EYE CONTOUR.

“SRT” EYE CONTOUR | “WR” WRINKLE REDUCTION | “SRT” SKIN REJUVENATION

TECHNICAL SPECIFICATIONS

Power

24 W

Emission frequency

1.15+0.1 MHz

Power supply

Mains

Size


55x155x255 mm

Weight

128 gr

STANDARD EQUIPMENT

- Code **G4259**
- 1 Suitcase
- 1 RF BEAUTY MINI radiofrequency unit
- 1 50 ml hyaluronic acid conductive face cream
- 1 50 ml conductive cream for lymphatic drainage with eascin, fucus and caffeine
- 1 Power supply
- 1 Operating manual

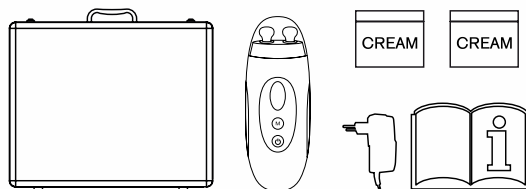

CE





Aesthetic line **ULTRASOUND**



LIPOZERO G39

AESTHETIC ULTRASOUND | 49 PROGRAMS

LIPOZERO G39 IS THE PROFESSIONAL LIPOLYTIC ULTRASOUND THAT FIGHTS CELLULITE AND ELIMINATES EXCESS FAT THANKS TO PROGRAMS CONCEIVED TO INTERVENE AT DIFFERENT FREQUENCIES AND ENSURE MAXIMUM EFFECTIVENESS IN EVERY BEAUTY TREATMENT.

10 MF-PLUS | 10 BEAUTY | 10 SKIN | 6 CAVIMAX 150 kHz CAVITATION
5 ADIPOSITY MOBILIZATION NorAd-526 | 4 MALE ADIPOSITY | 4 ULTRACAV 39 kHz

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

- 1 Ø 50 mm - IP68 handpiece
- 1 Ø 42 mm - IP68 handpiece

Power

3 W/cm² ± 20%

Emitted frequency

3 MHz - 1 MHz - 526 KHz - 150 KHz - 39 KHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

920 gr

STANDARD EQUIPMENT

- Code **G1404**
- 1 Suitcase
- 1 Ultrasound LIPOZERO G39
- 1 Multifrequency handpiece, Ø 42 mm - IP68
- 1 Handpiece to use with 39 KHz, Ø 50 mm - IP68
- 1 Gel bottle
- 1 Battery charger
- 1 Operating manual

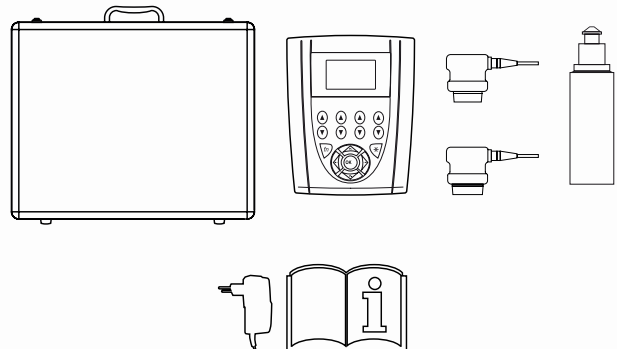

CE



STRONG POINTS



LIPOZERO G800

AESTHETIC ULTRASOUND | 30 PROGRAMS

DEVICE FOR PROFESSIONAL AESTHETICS, SUITABLE FOR FIGHTING THE IMPERFECTIONS OF CELLULITE AND LOCALIZED ADIPOSITY. AND TO IMPROVE SKIN APPEARANCE.
THE POSSIBILITY OF PROVIDING DIFFERENT FREQUENCIES WILL GUARANTEE A SPECIFIC AND EFFECTIVE TREATMENT FOR EVERY GOAL.

10 MF-PLUS | 10 BEAUTY | 10 SKIN

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

- 3 W/cm² ± 20% to 1.3 MHz and 3 MHz
- 1.5 W/cm² ± 20% to 1 MHz and 800 KHz

Emitted frequency

3 MHz - 1.3 MHz - 1 MHz - 800 KHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

920 gr

STANDARD EQUIPMENT

- Code G6232
- 1 Suitcase
- 1 Ultrasound LIPOZERO G800
- 1 Multifrequency handpiece, Ø 42 mm - IP68
- 1 Gel bottle
- 1 Battery charger
- 1 Operating manual

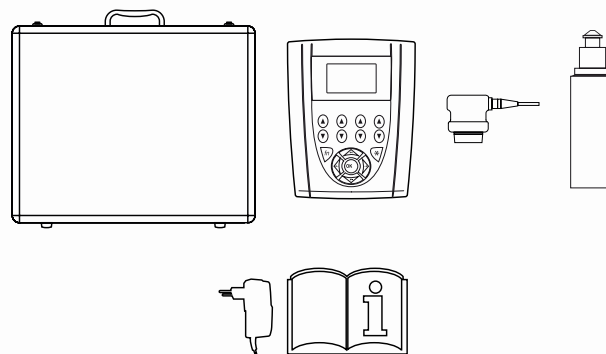

CE

In compliance with the Interministerial Decree n. 206/2015, the device is suitable for use in beauty centers.



STRONG POINTS



LIPOZERO G150

AESTHETIC ULTRASOUND | 42 PROGRAMS

ULTRASONIC CAVITATION IS THE MOST EFFECTIVE WORK METHOD TO FIGHT CELLULITE AND ELIMINATE EXCESS FAT. WITH LIPOZERO G150 YOU WILL BE ABLE TO BENEFIT FROM THE PROFESSIONAL TREATMENTS OF A BEAUTY CENTER COMFORTABLY AT HOME.

10 MF-PLUS | 9 CAVIMAX CAVITATION 150 kHz | 8 BEAUTY | 8 SKIN
5 ADIPOSITY MOBILIZATION NorAd-526 | 2 MALE ADIPOSITY

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

3 W/cm² ± 20%

Emitted frequency

3 MHz - 1 MHz - 526 KHz - 150 KHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

920 gr

STANDARD EQUIPMENT

- Code G1398
- 1 Suitcase
- 1 Ultrasound LIPOZERO G150
- 1 Multifrequency handpiece, Ø 42 mm - IP68
- 1 Gel bottle
- 1 Battery charger
- 1 Operating manual

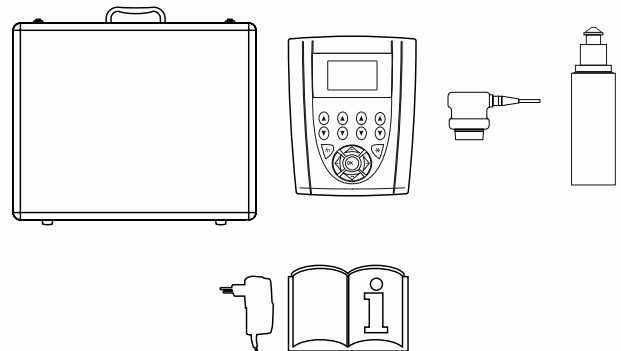

CE



STRONG POINTS





LIPOZERO EXCEL

AESTHETIC ULTRASOUND | 7 PROGRAMS

LIPOZERO EXCEL IS PART OF THE GLOBUS HOME CARE LINE AND HAS BEEN DESIGNED TO CARRY OUT THE PROFESSIONAL TREATMENTS OF AESTHETIC CENTERS AT HOME. THE PORTABLE DESIGN AND THE GRAPHIC MENU HAVE BEEN SPECIFICALLY DESIGNED FOR THE HOME USE AND MAKE THE DEVICE EASY TO USE EVEN BY LITTLE EXPERIENCED USERS. LIPOZERO EXCEL ALSO ALLOWS THE TREATMENT OF THE MOST COMMON IMPERFECTIONS RELATED TO CELLULITE OR FAT ACCUMULATION.

5 BEAUTY | 2 CELLULITE

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

1.2 W/cm² ± 20%

Emitted frequency

3 MHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Size

100x160x35 mm

Weight

270 gr

STANDARD EQUIPMENT

· Code G1297

· 1 Suitcase


· 1 Ultrasound LIPOZERO EXCEL

· 1 3 MHz, Ø 42 mm - IP68 handpiece

· 1 Gel bottle

· 1 Battery charger

· 1 Operating manual

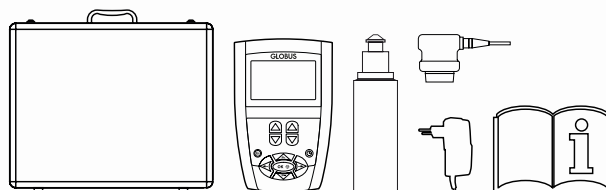

CE



STRONG POINTS



HT 906

AESTHETIC ULTRASOUND

HOME ULTRASOUND, IDEAL TO TREAT INFLAMMATIONS, JOINT AND MUSCLE PAIN. THE DEVICE ALSO ALLOWS TO TREAT IMPERFECTIONS LIKE LOCALIZED FAT AND CELLULITE.

1 CUSTOMIZABLE PROGRAM FOR ULTRASOUND

SPECIFICHE TECNICHE

Total peak power

2.5 W totali $\pm 20\%$

Average power per cm²

0.5 W/cm² $\pm 20\%$

Head diameter

30 mm

Emission frequency

1 MHz $\pm 5\%$

Emission

Continuous and pulsed

Mode

3 levels ("L" low - "M" medium - "H" high)

Timer

3 levels (5' - 10' - 15')

Power supply

Mains

Size

45x175x45 mm

Weight

110 gr

DOTAZIONE DI SERIE

· Code G5935


· 1 Bag

· 1 HT 906 device

· 1 Gel bottle

· 1 Power supply

· 1 Operating manual

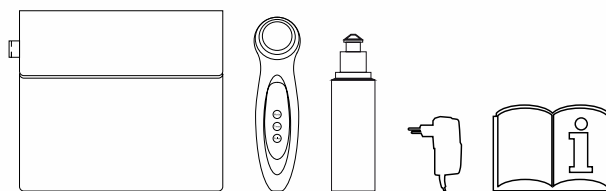

STRONG POINTS





Aesthetic line

VIBRATION PLATFORM

PHYSIOPLATE MY GOLD

VIBRATION PLATFORM

PHYSIOPLATE MY GOLD IS THE IDEAL VIBRATION PLATFORM FOR SPORT AND FITNESS LOVERS. IT CONTAINS MORE THAN 100 PROTOCOLS DIVIDED BY LEVELS AND WITH SPECIFIC EXERCISES FOR THE VARIOUS PARTS OF THE BODY.

TECHNICAL SPECIFICATIONS

Display

Size 7" color touch screen

Power

18 G

Vibration type

Vertical

Frequency range

15-60 Hz, with 1 Hz variation

Vibration width

Low/High

Voltage

230 V

Mains power supply

200-230 V AC, 50 Hz

Total amount of space

930x1550x760 mm

Plate size

930x600 mm

Maximum bearable weight

120 kg

Weight

170 kg

STANDARD EQUIPMENT

· Code **G1061**

· 1 PHYSIOPLATE MY GOLD vibration platform

· 1 Exercise mat

· 2 Exercise straps

· 1 Power supply

· 1 Operating manual

CE

ACCESSORIES



G1130

PHYSIOPLATE RED MAT

G1129

PHYSIOPLATE BLUE MAT



G1090

EXERCISE STRAP



SPORT
Line

CONTENTS

- ELECTROSTIMULATION SPORT Line	124
- PRESSOTHERAPY SPORT Line	134
- BALL SHOOTING MACHINE Line	136

The background is a dark, abstract composition featuring a complex network of thin, light-colored lines and geometric shapes. These elements create a sense of depth and movement, resembling a stylized architectural structure or a network of fibers. The overall tone is high-contrast and modern.

ELECTROSTIMULATION

SPORT Line



THE CHAMPION

4-CHANNEL ELECTROSTIMULATION | 441 PROGRAMS

IDEAL DEVICE FOR SPORTS LOVERS. THE DEVICE INCLUDES 12 SPORTS DISCIPLINES, EACH OF WHICH PROVIDES SPECIFIC PROGRAMS TO BUILD MUSCLES AND IMPROVE PERFORMANCE.

THE PROGRAMS ARE SUITABLE BOTH FOR THOSE WHO PRACTICE SPORTS AT AMATEUR LEVEL AND FOR PROFESSIONALS WHO CARRY OUT STRENGTHENING ACTIVITIES THAT THEY WANT TO INTEGRATE WITH TARGETED ELECTRO-STIMULATION.

90 SPECIAL SPORTS | 84 ACTION NOW | 70 BEAUTY | 58 FITNESS | 54 SERIAL SEQUENTIAL STIMULATION | 53 SPORT
23 MICROCURRENTS | 12 TENS AND PAIN | 3 REHAB | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G5015

· 1 Bag

· 1 THE CHAMPION stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

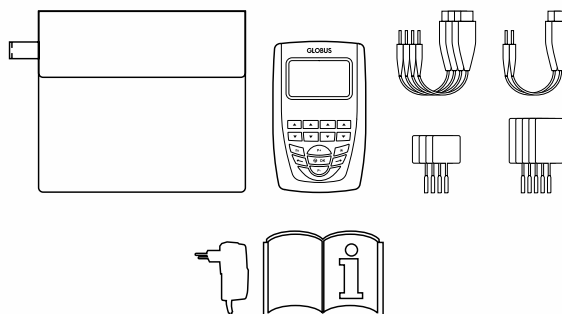

STRONG POINTS

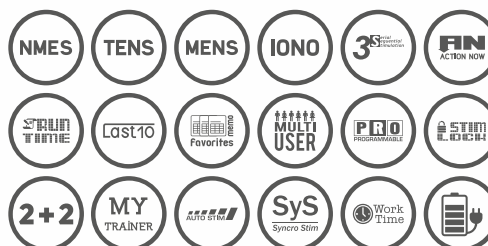

TRIATHLON PRO

4-CHANNEL ELECTROSTIMULATION | 424 PROGRAMS

DEVICE THAT, IN ADDITION TO A WIDE RANGE OF SPORT, FITNESS, BEAUTY AND REHAB PROGRAMS, INCLUDES A SERIES OF SPECIFIC PROGRAMS AIMED AT IMPROVING THE PHYSICAL PERFORMANCE OF ATHLETES. THE SPECIFIC PROGRAMS FOR TRIATHLON PRO ARE DIVIDED INTO AREAS CORRESPONDING TO THE THREE SPORTS DISCIPLINES (SWIMMING, CYCLING, RUNNING). A PRACTICAL USER GUIDE WILL HELP YOU TO CHOOSE THE MOST SUITABLE PROGRAMS ACCORDING TO THE DISTANCES YOU PRACTICE.

84 ACTION NOW | 60 SPECIAL TRIATHLON | 60 BEAUTY | 58 FITNESS | 54 SERIAL SEQUENTIAL STIMULATION
53 SPORT | 23 MICROCURRENTS | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 3 REHAB | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm


Weight

450 gr


* Based on the type of current

STANDARD EQUIPMENT

- Code **G3666**
- 1 Bag
- 1 TRIATHLON PRO stimulator
- 4 Cables for electrode connection
- 2 Cables for microcurrents and ionophoresis
- 4 Self-adhesive square electrodes
- 4 Self-adhesive rectangular electrodes
- 1 Battery charger
- 1 Operating manual

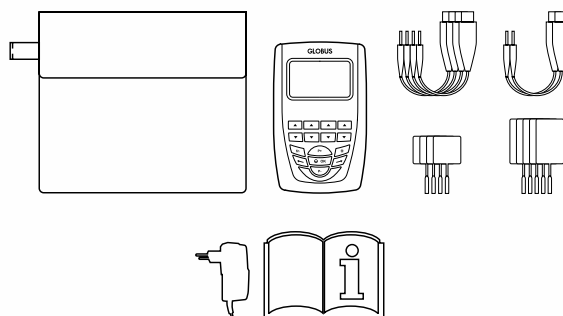

STRONG POINTS

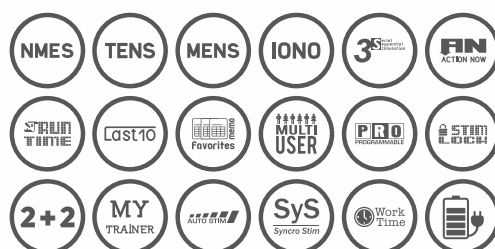

SWING PRO

4-CHANNEL ELECTROSTIMULATION | 400 PROGRAMS

IDEAL FOR PROFESSIONAL AND AMATEUR GOLFERS.

84 ACTION NOW | 60 BEAUTY | 58 FITNESS | 53 SPORT | 54 SERIAL SEQUENTIAL STIMULATION | 36 SPECIAL SPORTS
23 MICROCURRENTS | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 3 REHAB | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G3747

· 1 Bag

· 1 SWING PRO stimulator

· 1 G-Trode handpiece

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger


· 1 Operating manual

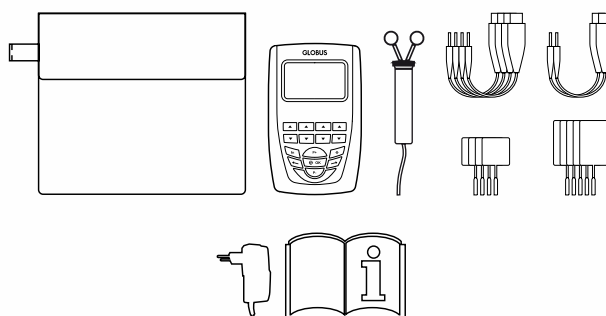

CE
0476



STRONG POINTS

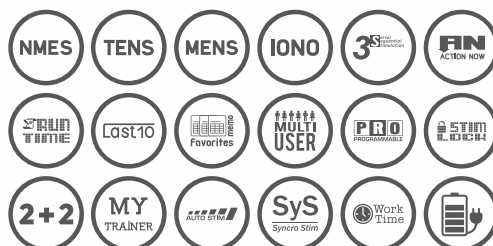

CYCLING PRO

4-CHANNEL ELECTROSTIMULATION | 270 PROGRAMS

DEVICE FOR THE CYCLISTS WHO WANT TO BE GOOD AND FEEL FIT DURING THEIR BICYCLE ACTIVITIES. THE SPECIFIC PROGRAMS FOR CYCLING ARE DIVIDED INTO MACRO-AREAS: PHYSICAL TRAINING, PREVENTION, PAIN-INJURY AND RECOVERY.

58 FITNESS | 53 SPORT | 46 SPECIAL CYCLING | 36 BEAUTY | 23 MICROCURRENTS
18 SERIAL SEQUENTIAL STIMULATION | 13 G-PULSE | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G4230

· 1 Bag

· 1 CYCLING PRO stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

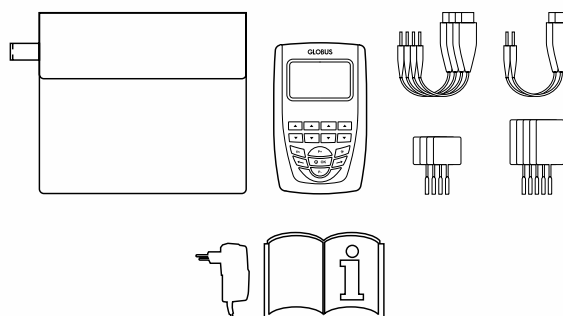

STRONG POINTS

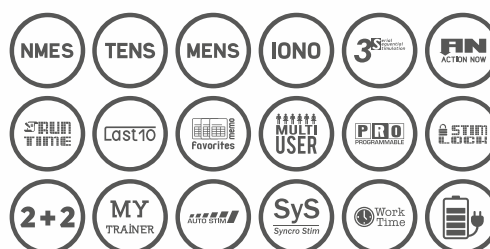

SAIL PRO

4-CHANNEL ELECTROSTIMULATION | 260 PROGRAMS

DEVICE WITH SPECIFIC PROGRAMS AIMED AT IMPROVING THE TRAINING STATE OF SAILBOAT LOVERS AND OF THOSE WHO PRACTICE WATER SPORTS, SUCH AS WINDSURF OR KITESURF.

THE AIM OF SAIL PRO IS TO HELP MAKE YOUR DAYS ON BOARD MORE PLEASANT AND SAFE AND, AT THE SAME TIME, IMPROVE PERFORMANCE IF YOU LIKE TO RACE.

58 FITNESS | 53 SPORT | 36 SPECIAL SAIL | 36 BEAUTY | 23 MICROCURRENTS
18 SERIAL SEQUENTIAL STIMULATION | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code **G4232**

· 1 Bag

· 1 SAIL PRO stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

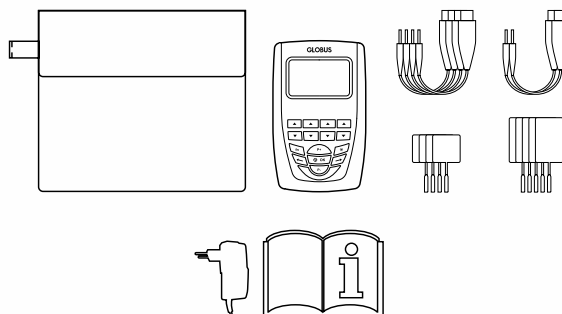

STRONG POINTS

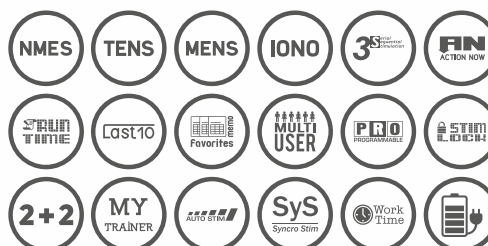

SOCCER PRO

4-CHANNEL ELECTROSTIMULATION | 258 PROGRAMS

SOCCER PRO IS THE ONLY ELECTROSTIMULATOR ON THE MARKET WITH SPECIFIC PROGRAMS TO IMPROVE THE PHYSICAL PERFORMANCE OF SOCCER PLAYERS. THE PROGRAMS ARE DIVIDED INTO FOUR LARGE AREAS CORRESPONDING TO THE MAIN FIELDS OF ATHLETIC TRAINING IN THE GAME OF SOCCER.

58 FITNESS | 53 SPORT | 36 BEAUTY | 34 SPECIAL SOCCER | 23 MICROCURRENTS

18 SERIAL SEQUENTIAL STIMULATION | 13 G-PULSE | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G4228

· 1 Bag

· 1 SOCCER PRO stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger


· 1 Operating manual

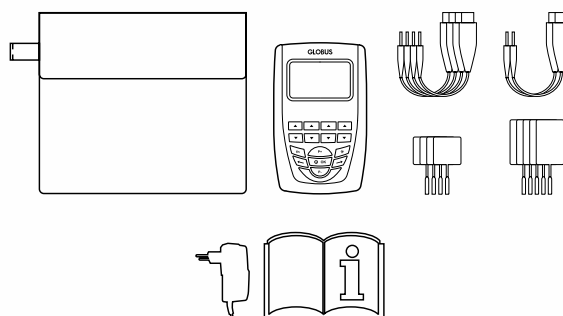

CE
0476



STRONG POINTS

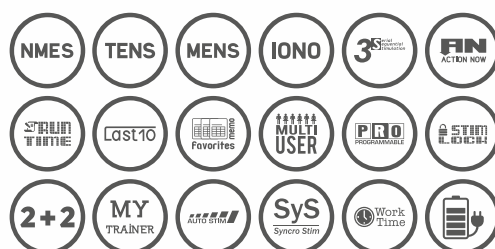

MOTO PRO

4-CHANNEL ELECTROSTIMULATION | 256 PROGRAMS

MOTO PRO ELECTROSTIMULATOR IS THE NEW TRAVEL COMPANION FOR ALL THE MOTORCYCLE ENTHUSIASTS. MOTO PRO WILL TAKE CARE OF ALL THE ASPECTS OF THE MOTORCYCLIST'S PHYSICAL SHAPE, DIFFERING PROPOSALS OF TRAINING WHICH WILL MAKE MORE PLEASANT HIS/HER RIDES. MOTO PRO IS REALLY A PERSONAL TRAINER! THE SPECIFIC PROGRAMS FOR MOTORCYCLE RIDERS ARE DIVIDED INTO THREE MACRO-AREAS.

58 FITNESS | 53 SPORT | 36 BEAUTY | 32 SPECIAL MOTORCYCLING | 23 MICROCURRENTS
18 SERIAL SEQUENTIAL STIMULATION | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code G4231

· 1 Bag

· 1 MOTO PRO stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes


· 4 Self-adhesive rectangular electrodes

· 1 Battery charger

· 1 Operating manual

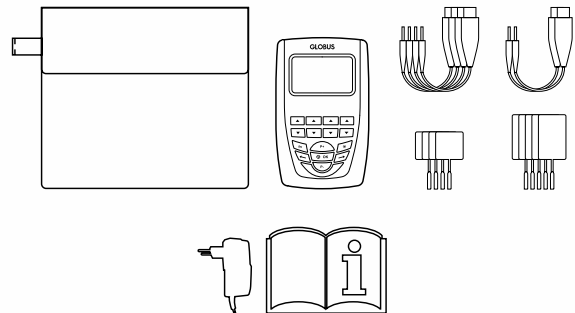

STRONG POINTS

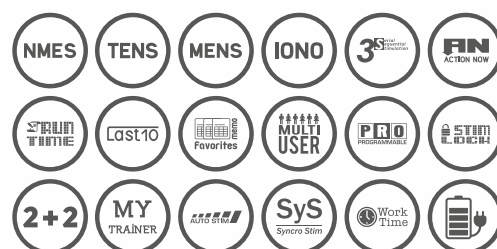

RUNNER PRO

4-CHANNEL ELECTROSTIMULATION | 254 PROGRAMS

RUNNER PRO IS A FUNDAMENTAL TOOL FOR THE MODERN RUNNER WHO WANTS TO OBTAIN IMPROVEMENTS IN PERFORMANCE. RUNNER PRO IS EQUIPPED WITH PROGRAMS SUITABLE TO STIMULATE THE MUSCLES OF THOSE WHO RUN AND EVERY RUNNER WILL BE ABLE TO PROGRAM AND CUSTOMIZE THE TRAINING SESSIONS ACCORDING TO PERSONAL NEEDS. THE PROGRAMS, DIVIDED INTO FOUR MACRO-GOALS, ARE SUITABLE BOTH FOR THE AMATEUR RUNNER AND FOR THE PROFESSIONAL.

58 FITNESS | 53 SPORT | 36 BEAUTY | 30 SPECIAL RUNNING | 23 MICROCURRENTS

18 SERIAL SEQUENTIAL STIMULATION | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Size of the visible area 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150 Hz *

Pulse width

40-450 μ s *

Intensity

0-120 mA per channel

Power supply

Rechargeable batteries

Size

100x160x35 mm

Weight

450 gr

* Based on the type of current

STANDARD EQUIPMENT

· Code **G4231**

· 1 Bag

· 1 RUNNER PRO stimulator

· 4 Cables for electrode connection


· 2 Cables for microcurrents and ionophoresis

· 4 Self-adhesive square electrodes

· 4 Self-adhesive rectangular electrodes

· 1 Battery charger


· 1 Operating manual

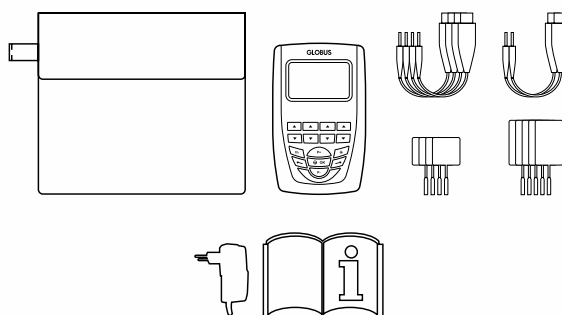

CE
0476



STRONG POINTS

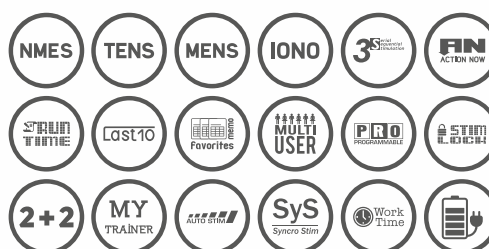



PRESSCARE G-SPORT³

PRESSOTHERAPY | 2 OUTPUTS - 3 INFLATION PROGRAMS

DEVICE FOR PNEUMATIC COMPRESSION DEDICATED TO ATHLETES, TO ACCELERATE MUSCULAR RECOVERY FROM POST-TRAINING/COMPETITION FATIGUE. LIGHT AND TRANSPORTABLE, EQUIPPED WITH A POWER BANK THAT ENABLES TREATMENTS ALSO RIGHT AFTER COMPETITION, DIRECTLY ON THE COURTSIDE.

TECHNICAL SPECIFICATIONS


Outputs
2
Programs
3
Pressure
0-180 mmHg $\pm 10\%$
Timer
Up to 99 minutes
Power supply
- Mains, 100-240 V AC, 50-60 Hz
- Rechargeable battery (Power Bank)
Size
180x190x135 mm
Weight
1.10 Kg

STANDARD EQUIPMENT

· Code G5380
· 1 Bag
· 1 G-SPORT ³ motor unit
· 2 Boots size "S"
· 2 Connection cables
· 1 Power supply
· 1 Power Bank
· 1 Operating manual

PURCHASE OPTIONS

· Code G5381 2 Boots size "M" included
· Code G5382 2 Boots size "L" included

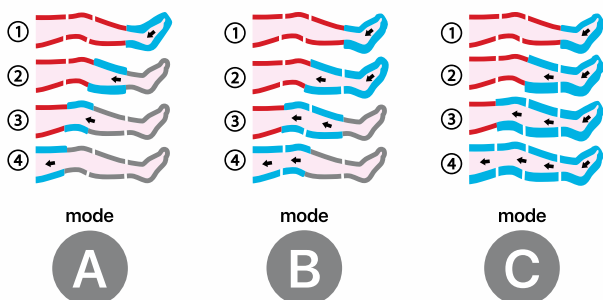

MODE A: the program sequence consists in the inflation and deflation of boot sectors one by one to execute an activating massage, ideal before the competition or training.


MODE B: the program sequence consists in the inflation and deflation of two chambers at a time to avoid downward fluid reflux and make the treatment effective and very comfortable.

MODE C: the program sequence consists in the sequential inflation of all the chambers and their simultaneous deflation, ideal for a deeper massage.



BALL SHOOTING MACHINES

Line

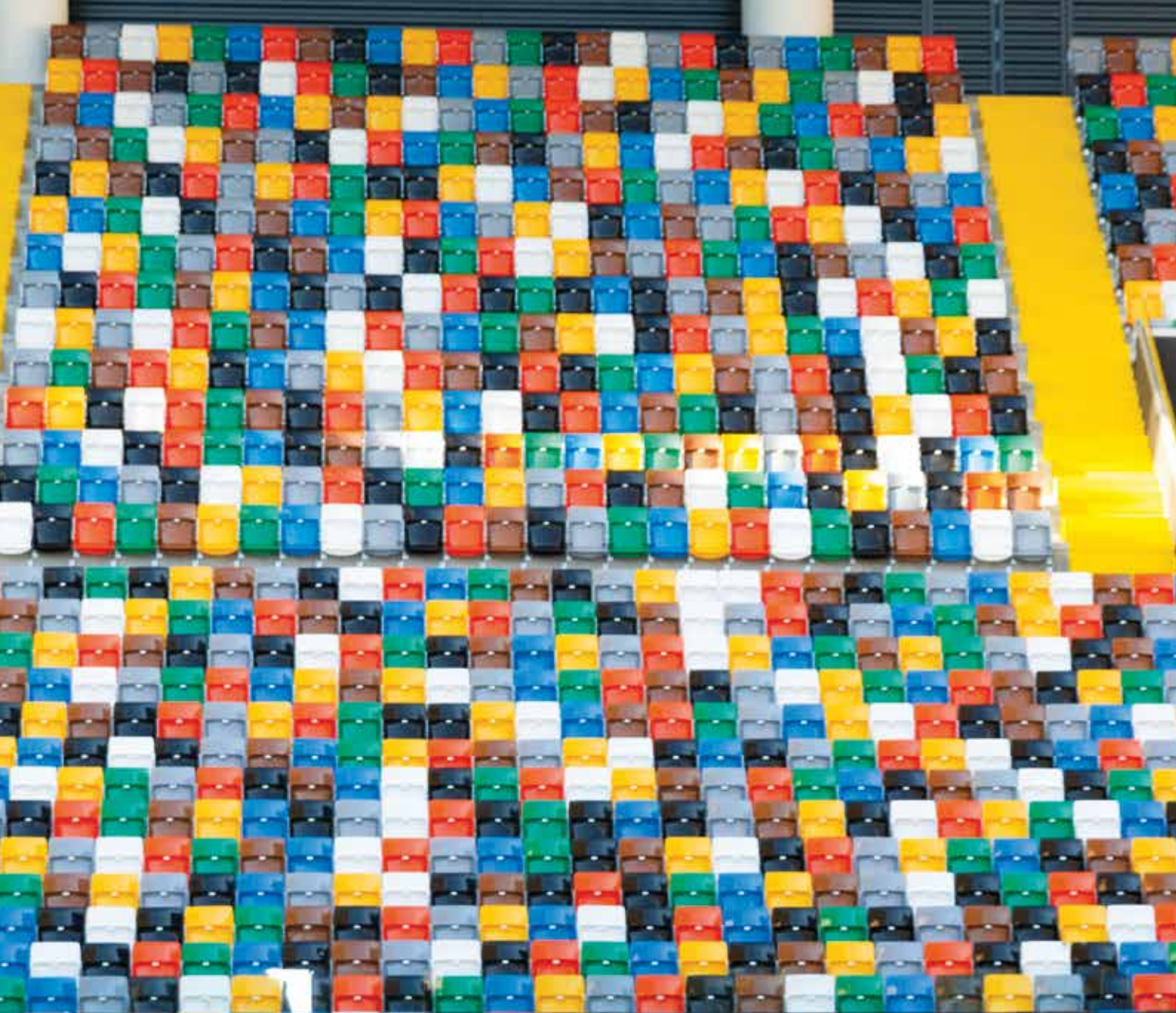

GLOBUS GLOBUS GLOBUS GLOBUS



EUROGOAL 1500

SOCCER BALL SHOOTING MACHINES

GLOBUS EUROGOAL 1500 IS A PROFESSIONAL BALL SHOOTING MACHINE FOR TRAINING, WHICH ALLOWS TO CARRY OUT WITH EXTREME REALISM CROSSES, SHOTS AND KICKS TYPICAL OF THE SOCCER GAME. IT IS CHARACTERIZED BY A VERY HIGH POWER THAT ALLOWS TO LAUNCH THE BALL UP TO 140 KM/H AND ESPECIALLY BY THE POSSIBILITY OF CREATING INFINITE TRAJECTORIES WITH ALL THE DESIRED EFFECTS.

THANKS TO THE PATENTED GLOBUS TRIAXIAL SYSTEM, THE COACH CAN ADJUST THE SHOT ON THE VERTICAL PLANE (UPWARD OR DOWNWARD), ON THE HORIZONTAL PLANE (RIGHT OR LEFT) AND, ABOVE ALL, TILT THE DEVICE ON THE OBLIQUE PLANE TO ENHANCE THE CURVE OF THE TRAJECTORIES. GLOBUS EUROGOAL 1500 IS, WITH NO DOUBT, THE IDEAL PARTNER OF THE COACH FOR HUNDREDS OF TRAINING EXERCISES, WITH HIGH RHYTHM AND MAXIMUM SHOOTING PRECISION.

TECHNICAL SPECIFICATIONS

· Code **G3305**

Dimensions (WxDxH)

115x70x60 cm

Weight

100 Kg

Shooting wheels

2, super-grip, wear-resistant rubber

Maximum ball shooting height

35 m

Lateral movement of the head

+30°/-30°

Vertical inclination

+40°/-5°

Shooting distance

100 m

Shooting speed

10 to 140 Km/h

Power supply

Double high-capacity battery
(5 actual hours or 1500 shots) with battery
charger included

Batteries

2, 12 V - 55 Ah

Shooting effects

Knuckleball - Magnus - Top Spin - Back Spin

Electric actuator for vertical elevation

Oblique inclination of the head

Lateral shot range lock

Electronic control unit with 5 shooting memories

Simultaneous and digital electronic
adjustment of the two wheels

STRONG POINTS

· Triaxial trajectory adjustment - GLOBUS patent

· No ball damage on the shooting wheels - GLOBUS patent

· Powerful kicks and spin shots, both in-swinging and out-swinging, from the right and from the left

· Digital display

· 5 buttons to store 5 shots with different speed and effects (spin)

· Ball speed control system that indicates when both wheels have reached the set speed (green LED)

· TOP SPIN: super curve of the trajectory thanks to the oblique movement of the triaxial system

· Compatible with different ball sizes

· Easy to move on the wheels

· Ball insertion guide to ensure the maximum repeatability of the launch

· Recommended for professional and semi-professional teams



EUROGOAL 600

SOCCER BALL SHOOTING MACHINES

GLOBUS EUROGOAL 600 IS A PROFESSIONAL BALL SHOOTING MACHINE FOR TRAINING, WHICH ALLOWS TO CARRY OUT EXTREMELY PRECISE SHOTS AND KICKS TYPICAL OF THE SOCCER GAME, WITH A POWER UP TO 110 KM/H.

THANKS TO A SPECIAL DEVICE, THE COACH CAN ADJUST THE SHOT ON THE VERTICAL PLANE (UPWARDS OR DOWNWARDS), TO MAKE GROUND SHOTS, HALF-HEIGHT SHOTS OR TO SHOOT AT THE CROSSING OF THE POSTS.

GLOBUS EUROGOAL 600 IS A TECHNICAL TRAINING TOOL THAT ALLOWS TO PERFORM MANY EXERCISES BOTH FOR GOALKEEPERS AND OTHER PLAYERS. IT CAN BE USED IN GLOBAL GAME SITUATIONS, SUCH AS CONSTRUCTIONS ON CROSSES AND CORNER KICKS.

TECHNICAL SPECIFICATIONS

· Code **G1503**

Dimensions (WxDxH)

95x70x60 cm

Weight

90 Kg

Shooting wheels

2, super-grip, wear-resistant rubber

Maximum ball shooting height

25 m

Lateral movement of the head

+30°/-30°

Vertical inclination

+40°/-5°

Shooting distance

60 m

Shooting speed

10 to 120 Km/h

Power supply

Double high-capacity battery
(5 actual hours or 1500 shots) with battery
charger included

Batteries

2, 12 V - 55 Ah

Shooting effects

Knuckleball - Magnus

Electric actuator for vertical elevation

Electronic control unit

Manual adjustment by means of two independent potentiometers

STRONG POINTS

· No ball damage on the shooting wheels - GLOBUS patent

· Powerful kicks and spin shots, both in-swinging and out-swinging, from the right and from the left

· Independent and micrometric adjustment of the wheel speed (999 levels)

· Compatible with different ball sizes

· Easy to move on the wheels

· Ball insertion guide to ensure the maximum repeatability of the launch

· Also recommended for youth sections and soccer schools



WINSHOT 1500

VOLLEYBALL SHOOTING MACHINES

WINSHOT 1500 IS THE TOP OF THE RANGE OF THE GLOBUS BALL SHOOTER LINE. IT IS PARTICULARLY SUITABLE FOR BOTH MALE AND FEMALE PROFESSIONAL TEAMS AS IT ALLOWS TO MAKE HUNDREDS OF SERVES, SPIKES AND SETS WITH MAXIMUM PRECISION AND WITH A PAUSE OF ONLY 1.5 SECONDS BETWEEN THE SHOTS. THIS MODEL CAN STORE 5 TYPES OF DIFFERENT SHOTS, BOTH IN TERMS OF POWER AND EFFECT, THANKS TO A SPECIAL ELECTRONIC BOARD. WINSHOT 1500 ALSO HAS THE "SPEED CONTROL SYSTEM" WHICH, THROUGH TWO RED AND GREEN LEDS IN THE MIDDLE OF THE CONTROL PANEL, ALLOWS TO CONTROL THE ACTUAL LAUNCH SPEED.

TECHNICAL SPECIFICATIONS

· Code **G3504**

Size

Minimum height with head off 190 cm

Minimum height with head on 260 cm

Maximum total height 370 cm

Operator platform height 150 cm

Maximum width without baskets 80 cm

Maximum width with one basket 120 cm

Maximum width with two baskets 160 cm

Machine depth

150 cm

Weight

220 Kg

Shooting height

230 cm to 340 cm

Shooting speed

10 to 140 Km/h

Power supply

2 12 V - 55 Ah rechargeable batteries with charger included

STRONG POINTS

- No ball damage on the shooting wheels - GLOBUS patent
- Digital display with simultaneous adjustment of the two shooting wheels
- Speed Control System which indicates with a green LED when the speed of the ball is the set one
- Possibility to store five types of shots with different power and effects
- Articulated, light and adjustable shooting head to use variable trajectories
- Easy to move around the court in all areas for team and individual work
- Silent and easy to use, with a targeted ball insertion that allows maximum shooting accuracy



WINSHOT 800

VOLLEYBALL SHOOTING MACHINES

WINSHOT 800 IS THE GLOBUS BALL SHOOTING MACHINE ESPECIALLY SUITABLE FOR PROFESSIONAL FEMALE TEAMS AND MID-LEVEL MALE TEAMS, AS IT ALLOWS TO SHOOT HUNDREDS OF SERVES, SPIKES AND SETS WITH MAXIMUM PRECISION. THIS MACHINE CAN LAUNCH BALLS, TO SERVE AND ATTACK, AT THE IMPRESSIVE SPEED OF 130 KM PER HOUR AT A HEIGHT OF 305 CM. THIS FEATURE MAKES WINSHOT 800 VERY USEFUL AND EFFECTIVE IN TRAINING, AS SERVES AND SPIKES IN COMPETITION HAVE REACHED A BALL SPEED THAT FORCES THE RECEPTION AND THE DEFENSE TO TRAIN WITH GREAT FREQUENCY AND ACCURACY, WITH SERVES HARDLY PERFORMED BY THE HUMAN ARM, SAFEGUARDING THE SHOULDERS OF ATHLETES AND COACHES.

TECHNICAL SPECIFICATIONS

· Code **G1334**

Size

Minimum height with head off 175 cm

Minimum height with head on 255 cm

Maximum total height 350 cm

Operator platform height 110 cm

Maximum width without baskets 80 cm

Maximum width with one basket 120 cm

Maximum width with two baskets 160 cm

Machine depth

130 cm

Weight

- No Battery kit 150 Kg

- With Battery kit 190 Kg

Shooting length

200 cm to 800 cm

Shooting speed

10 to 130 Km/h

Power supply

- Mains, with power supply included

- 2 12 V - 55 Ah rechargeable batteries with charger included (not included)

STRONG POINTS

· No ball damage on the shooting wheels - GLOBUS patent

· Shooting wheels adjusted separately with 999 speed levels

· Articulated, light and adjustable shooting head to use variable trajectories

· Mains power supply included or rechargeable batteries with 10-hour capacity (optional)

· Easy to move around the court in all areas for team and individual work

· Silent and easy to use, with a targeted ball insertion that allows maximum shooting accuracy



WINSHOT 500

VOLLEYBALL SHOOTING MACHINES

WINSHOT 500 IS THE GLOBUS BALL SHOOTING MACHINE ESPECIALLY SUITABLE FOR YOUNG MALE AND FEMALE TEAMS AND FOR MID-LEVEL MALE TEAMS, AS IT ALLOWS TO SHOOT HUNDREDS OF SERVES, SPIKES AND SETS WITH MAXIMUM PRECISION. THIS MACHINE CAN LAUNCH BALLS, TO SERVE AND ATTACK, AT THE IMPRESSIVE SPEED OF 120 KM PER HOUR AT A MAXIMUM HEIGHT OF 270 CM. WINSHOT 500 IS EASY TO DISASSEMBLE AND TO MOVE TO DIFFERENT PLACES. IT IS VERY FLEXIBLE AND LIGHT IN THE GYM, TO CARRY OUT SPECIFIC TECHNICAL SESSIONS FOR THE RECEPTION AND THE DEFENSE, BOTH INDIVIDUALLY AND IN TEAM. THESE FEATURES MAKE WINSHOT 500 VERY USEFUL AND EFFECTIVE DURING TRAINING TO IMPROVE THE TECHNIQUE OF YOUNG VOLLEYBALLERS, THROUGH A SIGNIFICANT NUMBER OF HIGH-QUALITY REPETITIONS.

TECHNICAL SPECIFICATIONS

· Code **G1333**

Size

Minimum height with head off 175 cm

Minimum height with head on 240 cm

Maximum total height 300 cm

Maximum width 80 cm

Machine depth

100 cm

Weight

- No Battery kit 70 Kg

- With Battery kit 110 Kg

Shooting length

210 cm to 270 cm

Shooting speed

10 to 120 Km/h

Power supply

- Mains, with power supply included

- 2 12 V - 55 Ah rechargeable batteries with charger included (not included)

STRONG POINTS

· No ball damage on the shooting wheels - GLOBUS patent

· Shooting wheels adjusted separately with 999 speed levels

· Articulated, light and adjustable shooting head to use variable trajectories

· Mains power supply included or rechargeable batteries with 10-hour capacity (optional)

· Easy to move around the court in all areas for team and individual work

· Silent and easy to use, with a targeted ball insertion that allows maximum shooting accuracy

· Quick and easy to disassemble for transport by car to various training locations



RUGBYKICK

RUGBY SHOOTING MACHINES

RUGBYKICK IS THE TOP OF THE RANGE AMONG THE GLOBUS TRAINING MACHINES. IT IS PARTICULARLY SUITABLE FOR PROFESSIONAL TEAMS, AS IT IS ABLE TO REPEAT HUNDREDS OF SERVICES, ATTACKING AND SETTING ACTIONS WITH MAXIMUM PRECISION, AT INTERVALS OF ONLY 1.5 SECONDS. THIS MODEL CAN STORE 5 SHOTS WITH DIFFERENT POWER AND EFFECT PARAMETERS THANKS TO A SPECIAL ELECTRONIC BOARD. THE RUGBYKICK SPEED CONTROL SYSTEM ALLOWS TO CONTROL THE SPEED OF THE BALL THANKS TO TWO LEDS (RED AND GREEN) IN THE MIDDLE OF THE CONTROL PANEL.

TECHNICAL SPECIFICATIONS

· Code **G5371**

Size

Maximum height 145 cm

Minimum height 130

Width 70 cm

Length 100 cm

Weight

95 Kg

Shooting height

230 cm to 80 cm

Shooting speed

10 to 120 Km/h

Power supply

2 rechargeable batteries 12 V - 55 Ah

STRONG POINTS

- No ball damage on the shooting wheels - GLOBUS patent
- Digital display with simultaneous adjustment of the two shooting wheels
- Speed Control System which indicates with a green LED when the speed of the ball is the set one
- Possibility to store five types of shots with different power and effects
- Articulated, light and adjustable shooting head to use variable trajectories
- Easy to move around the court in all areas for team and individual work
- Silent and easy to use, with a targeted ball insertion that allows maximum shooting accuracy



ACCESSORIES

CONTENTS

- TECARTHERAPY AND TECAR BEAUTY	146
- MAGNETOTHERAPY	152
- ELECTROTHERAPY	155
- ELECTROSTIMULATION	155
- ULTRASOUND THERAPY	159
- LASER THERAPY	160
- RADIOFREQUENCY	161
- PRESSOTHERAPY	162
- BALL SHOOTING MACHINES	164

ACCESSORIES FOR DIACARE 7000 - TECAR BEUTY 7000 MED TECAR BEAUTY 7000 - TECAR BEAUTY ELITE



PH 01

PH 01. BIPOLAR HANDPIECE FOR RESISTIVE TREATMENTS

G5128	BIPOLAR HANDPIECE Ø 30 mm
G5127	BIPOLAR HANDPIECE Ø 50 mm



PH 02

PH 02. MONOPOLAR RESISTIVE HANDPIECE AND ELECTRODES

G5129	MONOPOLAR RESISTIVE HANDPIECE
G5337	RESISTIVE ELECTRODE, Ø 30 mm HEAD
G5338	RESISTIVE ELECTRODE, Ø 50 mm HEAD
G5339	RESISTIVE ELECTRODE, Ø 70 mm HEAD



PH 03

PH 03. MONOPOLAR CAPACITIVE HANDPIECE AND ELECTRODES

G5130	MONOPOLAR CAPACITIVE HANDPIECE
G3306	CAPACITIVE ELECTRODE, Ø 30 mm HEAD
G3308	CAPACITIVE ELECTRODE, DA Ø 50 mm HEAD
G3310	CAPACITIVE ELECTRODE, Ø 70 mm HEAD



PH 04

PH 04. ACTIVE NEUTRAL HANDPIECE

G5133	RESISTIVE HANDPIECE TO USE INSTEAD OF THE NEUTRAL PLATE
-------	---



PH 05 (above) / PH 06 (below)

PH 05. SET OF NEUTRAL CABLE

G5123	NEUTRAL "A" CABLE, 4 mm, BLACK
G5022	ADAPTER FOR NEUTRAL CABLES, BLACK
G6020	"C" NEUTRAL BLACK CABLE WITH CLIP

PH 06. SET OF RESISTIVE CABLES

G5126	RESISTIVE "A" CABLE, 4 mm, BLUE
G5018	ADAPTER FOR RESISTIVE CABLES, BLUE
G6021	"C" RESISTIVE BLUE CABLE WITH CLIP

ACCESSORIES FOR DIACARE 5000 AND 5000 RE TECAR BEUTY 6000 MED AND 6000 MED RE



PH 01

PH 01. BIPOLAR HANDPIECE FOR RESISTIVE TREATMENTS

G3486	BIPOLAR HANDPIECE Ø 30 mm
G3487	BIPOLAR HANDPIECE Ø 50 mm



PH 02

PH 02. MONOPOLAR RESISTIVE HANDPIECE AND ELECTRODES

G3485	MONOPOLAR RESISTIVE HANDPIECE
G5337	RESISTIVE ELECTRODE, Ø 30 mm HEAD
G5338	RESISTIVE ELECTRODE, Ø 50 mm HEAD
G5339	RESISTIVE ELECTRODE, Ø 70 mm HEAD



PH 03

PH 03. MONOPOLAR CAPACITIVE HANDPIECE AND ELECTRODES


G3483	MONOPOLAR CAPACITIVE HANDPIECE
G3306	CAPACITIVE ELECTRODE, Ø 30 mm HEAD
G3308	CAPACITIVE ELECTRODE, Ø 50 mm HEAD
G3310	CAPACITIVE ELECTRODE, Ø 70 mm HEAD



PH 04

PH 04. ACTIVE NEUTRAL HANDPIECE

G5035	RESISTIVE HANDPIECE TO USE INSTEAD OF THE NEUTRAL PLATE
-------	---

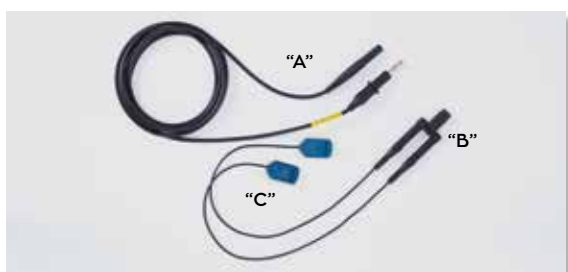

PH 05

PH 05. SET OF NEUTRAL CABLE

G5144	NEUTRAL "A" CABLE, 4 mm, BLACK
G5022	"B" ADAPTER FOR NEUTRAL CABLES, BLACK
G6020	"C" NEUTRAL BLACK CABLE WITH CLIP




PH 06

PH 06. SET OF RESISTIVE CABLES

G5017	RESISTIVE "A" CABLE, 4 mm, BLUE
G5018	"B" ADAPTER FOR RESISTIVE CABLES, BLUE
G6021	"C" RESISTIVE BLUE CABLE WITH CLIP

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE

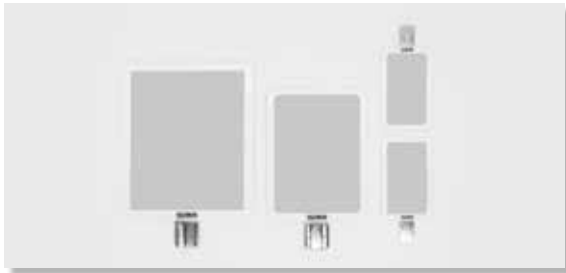

PH 01

PH 01. FLEXIBLE NEUTRAL PLATES FOR ALL THE DEVICES IN THE "DIACARE" AND "TECAR BEAUTY" LINE

G5393	FLEXIBLE STEEL LARGE PLATE, 200X250 mm
G5728	FLEXIBLE STEEL MEDIUM PLATE, 150X210 mm
G5985	FLEXIBLE STEEL SMALL PLATE, 75X130 mm
G5983	KIT FLEXIBLE STEEL SMALL PLATES, 75X130 mm, 2 PC



PH 02

PH 02. HANDS-FREE TREATMENT ELECTRODES FOR ALL THE DEVICES IN THE "DIACARE" AND "TECAR BEAUTY" LINE

G5986	ADHESIVE TECAR ELECTRODE, SMALL 74X88 mm, 1 PC
G5977	ADHESIVE TECAR ELECTRODE, MEDIUM 90X148 mm, 1 PC
G5987	ADHESIVE TECAR ELECTRODE, LARGE 101X202 mm, 1 PC



PH 03

PH 03. HANDS-FREE KIT PLUS FOR THE FOLLOWING DEVICES: "DIACARE 7000" - "TECAR BEAUTY 7000 MED" "TECAR BEAUTY 7000" AND "TECAR BEAUTY ELITE"

G5976	1 BAG 1 RESISTIVE CABLE 2 Y CABLES (1 BLUE, 1 BLACK) 4 CABLES WITH CROCODILES (2 BLUE, 2 BLACK) 4 MEDIUM ADHESIVE ELECTRODES 2 ELASTIC BANDS
-------	---



PH 04

PH 04. HANDS-FREE KIT PLUS FOR THE FOLLOWING DEVICES: "DIACARE 5000" AND "DIACARE 5000 RE" "TECAR BEAUTY 6000 MED" AND "TECAR BEAUTY 6000 MED RE"

G5981	1 BAG 1 RESISTIVE CABLE 2 Y CABLES (1 BLUE, 1 BLACK) 4 CABLES WITH CROCODILES (2 BLUE, 2 BLACK) 4 MEDIUM ADHESIVE ELECTRODES 2 ELASTIC BANDS
-------	---

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



PH 05

PH 05. NEUTRAL PLATE CABLE FOR THE FOLLOWING DEVICES:
"DIACARE 7000" - "TECAR BEAUTY 7000 MED"
"TECAR BEAUTY 7000" AND "TECAR BEAUTY ELITE"

G5123 | NEUTRAL PLATE CABLE, 4 mm, BLACK



PH 06

PH 06. PLATE CABLE FOR THE FOLLOWING DEVICES:
"DIACARE 5000" AND "DIACARE 5000 RE"
"TECAR BEAUTY 6000 MED" AND "TECAR BEAUTY 6000 MED RE"

G3503 | PLATE CABLE, 4 mm, BLACK



PH 07

PH 07. POWER-SUPPLY CABLE FOR THE FOLLOWING DEVICES:
"DIACARE 7000" - "TECAR BEAUTY 7000 MED"
"TECAR BEAUTY 7000" AND "TECAR BEAUTY ELITE"

G4349 | POWER-SUPPLY CABLE, 230 V- 50/60 Hz, SHUKO



PH 08

PH 08. POWER-SUPPLY FOR THE FOLLOWING DEVICES:
"DIACARE 5000" AND "DIACARE 5000 RE"
"TECAR BEAUTY 6000 MED" AND "TECAR BEAUTY 6000 MED RE"

G3517 | POWER-SUPPLY 48 V



PH 09

PH 09. TECAR CONDUCTIVE CREAMS FOR THE FOLLOWING DEVICES:
"DIACARE 7000" - "DIACARE 5000" AND "DIACARE 5000 RE"

G3368 | CONDUCTIVE TECAR CREAM, 250 ml

G3369 | CONDUCTIVE TECAR CREAM, 1000 ml

G5396 | CONDUCTIVE TECAR CREAM WITH ARNICA, 1000 ml



PH 10

PH 10. TECAR BEAUTY CONDUCTIVE CREAMS FOR THE FOLLOWING DEVICES:
"TECAR BEAUTY 7000 MED" - "TECAR BEAUTY 7000" - "TECAR BEAUTY ELITE"
"TECAR BEAUTY 6000 MED" AND "TECAR BEAUTY 6000 MED RE"

G6333 | CONDUCTIVE TECAR BEAUTY CREAM, 250 ml

G6332 | CONDUCTIVE TECAR BEAUTY CREAM, 1000 ml

G6226 | TECAR BEAUTY CREAM WITH HYALURONIC ACID, 500 ml

G6227 | TECAR BEAUTY CREAM FOR LYMPHATIC DRAINAGE, 1000 ml

OTHER ACCESSORIES FOR DIACARE AND TECAR BEAUTY LINE



G3371	GRAY ALUMINUM SUITCASE FOR "DIACARE 5000" AND "5000 RE" - "TECAR BEAUTY 6000 MED" AND "6000 MED RE"
G5591	CONVENIENT CART WITH WHEELS



G5520	TROLLEY "7000" AND "ELITE" LINE
G5521	TROLLEY WITH MULTIPLE SPACES TO TRANSPORT THE DEVICES

OTHER ACCESSORIES FOR DIACARE AND TECAR BEAUTY LINE



G6342 | HOLDER FOR DEVICES AND ACCESSORIES ONLY FOR "DIACARE 5000" AND "5000 RE" - "TECAR BEAUTY 6000 MED" AND 6000 MED RE"



G6341 | HOLDER FOR ACCESSORIES FOR "DIACARE 5000" AND "5000 RE" - "TECAR BEAUTY 6000 MED" AND "6000 MED RE"

MAGNETOTHERAPY ACCESSORIES



PH 01

PH 01. SINGLE "SOFT" SOLENOID

G4275 | SINGLE "SOFT" SOLENOID, 120X120 mm, 1 COIL, WASHABLE AND SANITIZABLE FABRIC



PH 02

PH 02. "SOFT 4" SOLENOID

G4276 | "SOFT 4" MAT, 250X250 mm, 4 COILS, 2 TOTAL CHANNELS, WASHABLE AND SANITIZABLE FABRIC



PH 03

PH 03. FLEXIBLE SOLENOID

G1408 | FLEXIBLE SOLENOID, 400X150 mm, 4 COILS



PH 04

PH 04. FLEXIBLE SOLENOID

G1407 | FLEXIBLE SOLENOID, 300X100 mm, 2 COILS



PH 05

PH 05. "POCKET PRO" SOLENOID

G5988 | "POCKET PRO" SOLENOID, 155x95 mm, 2 COILS. THE CODE REFERS TO 1 SOLENOID BUT THE PURCHASE IS MADE IN PAIRS

G6050 | POCKET PRO CASE



PH 06

PH 06. MAGNUM TEST

G5867 | INDICATOR OF CORRECT EMISSION OF THE MAGNETIC FIELD, 85X55 mm



PH 07

PH 07. ELASTIC BAND

G4280 | VELCRO ELASTIC BAND, 80X800 mm, BLACK

MAGNETOTHERAPY ACCESSORIES



PH 08

PH 08. SOLENOID COVER

G5878	FLEXIBLE SOLENOID COVER KIT (10 PIECES)
G5879	"SOFT" AND RIGID SOLENOID COVER KIT (10 PIECES)



PH 09



PH 10

PH 09. POWER SUPPLY


G3566	MAGNETO LINE POWER SUPPLY
-------	---------------------------

PH 10. POWER SUPPLY

G4048	MAGNETO LINE POWER SUPPLY
-------	---------------------------



PH 11

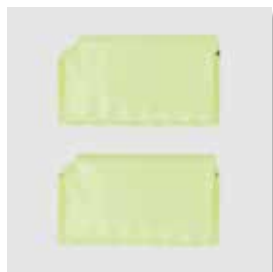

PH 12

PH 11. POWER SUPPLY

G4147	"MAGNUM 3500 PRO" POWER SUPPLY
-------	--------------------------------

PH 12. BATTERY PACK

G1120	BATTERY PACK - 12 V AA NI-MH FOR "MAGNUM 3500 PRO"
-------	--



PH 13



PH 14

PH 13. PACCO BATTERIA

G0699	BATTERY PACK - 7,2 V NI-MH FOR "MAGNUM 3000 PRO" - "MAGNUM 2500" - "MAGNUM XL PRO" "MAGNUM XL" AND "MAGNUM L"
-------	---

PH 14. Y CABLE

G4337	Y CABLE, TO USE 2 SINGLE SOFT SOLENOIDS ON THE SAME CHANNEL
-------	---



PH 15

PH 15. EXTENSION

G5063	EXTENSION FOR ACCESSORIES LENGTH 140 cm
-------	---


MAGNETOTHERAPY ACCESSORIES

MEMORY FOAM SINGLE MAT



DIMENSIONS: LENGTH 440 mm, WIDTH 700 mm

G5944	"MAT 100 MEMORY FOAM" MAT, 4 SOLENOIDS,
G5945	"E-MAT 100" MAT, ACTIVE EXTENSION, MEMORY FOAM
G5946	"N-MAT 100" MAT, NEUTRAL EXTENSION, MEMORY FOAM
G6048	EXTERNAL COVER FOR "MAT 100 MEMORY FOAM"
G6049	INTERNAL COVER FOR "MAT 100 MEMORY FOAM"

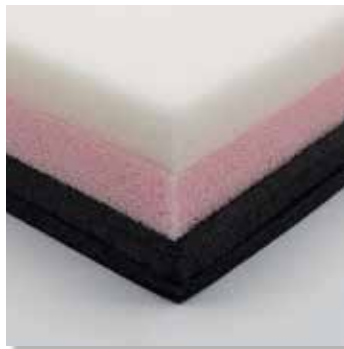

INTERNAL SECTION THAT HIGHLIGHTS THE DIFFERENT LAYERS OF THE MAT AND MATTRESS

"TOTAL BODY" MEMORY FOAM MAT



DIMENSIONS: LENGTH 1760 mm, WIDTH 700 mm, TO CONNECT TO THE "GLOBUS" MAGNETOTHERAPY DEVICES

G5947	"TOTAL BODY 200 MEMORY FOAM", 2 ACTIVE MATS OUT OF 4, 8 TOTAL SOLENOIDS
G5948	"TOTAL BODY 400 MEMORY FOAM", 4 ACTIVE MATS OUT OF 4, 16 TOTAL SOLENOIDS
G5949	"TOTAL BODY 400 XP MEMORY FOAM", 4 ACTIVE MATS OUT OF 4, 16 TOTAL SOLENOIDS ONLY COMPATIBLE WITH MAGNUM 3500 PRO

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 01

PH 01. ELECTRODES - 4 PCs., 50X50 mm

GO136	MYOTRODE "PLATINUM"
GO464	MYOTRODE "PLUS"
GO466	MYOTRODE "PREMIUM"



PH 02

PH 02. ELECTRODES - 4 PCs., 50X90 mm

GO137	MYOTRODE "PLATINUM"
GO465	MYOTRODE "PLUS"
GO467	MYOTRODE "PREMIUM"



PH 03

PH 03. DOUBLE INLET ELECTRODES - 2 PCs., 50X100 mm

GO458	MYOTRODE "PLATINUM"
-------	---------------------



PH 04

PH 04. FACE ELECTRODES - Ø 32 mm WITH CABLE

GO451	4 FACE ELECTRODES
-------	-------------------



PH 05

PH 05. 30 FACE ELECTRODES - Ø 26 mm WITH CLIP

GO871	30 FACE CLIP ELECTRODES
-------	-------------------------



PH 06



PH 07

PH 06. CLIP SQUARE ELECTRODES - 4 PCs., 46X47 mm

GO650	4 CLIP ELECTRODES
-------	-------------------

PH 07. CLIP RECTANGULAR ELECTRODES - 4 PCs., 45X98 mm

GO651	4 CLIP ELECTRODES
-------	-------------------

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES




PH 08

PH 08. BUCKSKIN ELECTRODES

GO452	MEDIUM BUCKSKIN ELECTRODES - 4 PCs, 80X60 mm
GO453	LARGE BUCKSKIN ELECTRODES - 4 PCs, 80X120 mm



PH 09

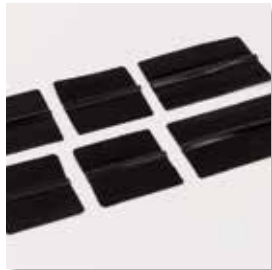

PH 10

PH 09. / PH 10. CONDUCTIVE SILICONE ELECTRODES

GO461	CONDUCTIVE SILICONE ELECTRODE PER METER
GO463	KIT 6 CONDUCTIVE SILICONE ELECTRODES



PH 11

PH 11. / PH 12. CONDUCTIVE ELASTIC BANDS

GO479	KIT 4 THIGH BANDS
GO480	KIT 8 "FITNESS TOP" LEG AND THIGH BANDS



PH 12

MORE ELASTIC BANDS

G1225	KIT 2 ELASTIC BANDS "S", 80X400 mm, BLUE
G1226	KIT 2 ELASTIC BANDS "M", 80X800 mm, BLUE



PH 13

PH 13. "EXTENSION BELT" BAND

G1601	"EXTENSION BELT" BAND, 450 mm (extension for "WINTEC SMART" and "FAST BAND")
-------	---



PH 14

PH 14. IONOPHORESIS ELECTRODES

GO890	MEDIUM CARBON ELECTRODE + POUCH, 50x50 mm		
GO885	LARGE CARBON ELECTRODE + POUCH, 60x85 mm		
G1210	POUCH, 50x50 mm	GO155	POUCH, 60x85 mm

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 15

PH 15. G-TRODE HANDPIECE

G1309 | G-TRODE BIPOLAR HANDPIECE FOR "SHAPE & BEAUTY" AND "ACTIVA 700"



PH 16

PH 16. COLORED CABLES FOR 4 CHANNEL DEVICES

G1839 | KIT 4 COLORED CABLES



PH 17

PH 17. Y CABLES FOR 4 CHANNEL DEVICES

GO441 | KIT 2 Y CABLES FOR 4 CHANNELS



PH 18



PH 19

PH 18. CABLES FOR 4 CHANNEL DEVICES

GO732 | GLOBUS SPARE CABLE FOR 4 CHANNELS

PH 19. MICROCURRENT CABLES

GO440 | KIT OF 2 CABLES FOR MICROCURRENTS



PH 20



PH 21

PH 20. CABLES FOR 2 CHANNEL DEVICES, ROUND PLUG

G3330 | "WINTEC" GRAY SPARE CABLE

PH 21. CABLES FOR 2 CHANNEL DEVICES, RECTANGULAR PLUG

G1171 | "WINTEC" BLU SPARE CABLE



PH 22



PH 23

PH 22. ADAPTERS

GO727 | 4 ADAPTERS FOR CLIP ELECTRODES, Ø 2 mm

PH 23. CABLES FOR 2 CHANNEL DEVICES

GO439 | KIT 2 Y CABLES (universal)

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 24



PH 25

PH 24. BATTERY PACK

G3870 | BATTERY PACK FOR 2 CHANNEL DEVICES

PH 25. BATTERY PACK

G1126 | BATTERY PACK FOR "ELITE"



PH 26



PH 27

PH 26. BATTERY PACK

G0699 | BATTERY PACK FOR 4 CHANNEL DEVICE

PH 27. BATTERY PACK

G1127 | BATTERY PACK FOR "GENESY 3000"



PH 28



PH 29

PH 28. POWER SUPPLY

G4048 | "GENESY 1500" POWER SUPPLY

PH 29. POWER SUPPLY

G4146 | "GENESY 3000" POWER SUPPLY



PH 30



PH 31

PH 30. BATTERY CHARGER

G5279 | BATTERY CHARGER FOR 2 CHANNEL DEVICES

PH 31. BATTERY CHARGER

G4154 | BATTERY CHARGER FOR 4 CHANNEL DEVICES



PH 32



PH 33

PH 32. BATTERY CHARGER

G4067 | BATTERY CHARGER FOR "GENESY 300 PRO" AND "ELITE"

PH 33. POWER SUPPLY

G3915 | CHARGER FOR AA/AAA BATTERIES ONLY FOR "DUO PRO"



PH 34



PH 35

PH 34. MOTOR POINT PEN

G1156 | MOTOR POINT PEN

PH 35. GEL

G0484 | 260 ml GEL BOTTLE, 10 PCs

G0869 | 260 ml GEL BOTTLE, 1 PC

G4034 | 1000 ml GEL CAN

G3609 | 5000 ml GEL CAN

ULTRASOUND THERAPY ACCESSORIES



PH 01

PH 01. Ø 16 mm HANDPIECES

G5516	Ø 16 mm HEAD HANDPIECE FOR "MEDISOUND 3000" - SD5 BOARD
G4060	Ø 16 mm HEAD HANDPIECE FOR "MEDISOUND 3000" - SD4 BOARD



PH 02



PH 03

PH 02. Ø 42 mm HANDPIECES

G5374	Ø 42 mm HEAD HANDPIECE FOR "MEDISOUND 3000" AND "LIPOZERO G150"
G3152	Ø 42 mm HANDPIECE WHITE HEAD FOR "MEDISOUND PRO 2" AND "LIPOZERO EXCEL"
G1402	Ø 42 mm HANDPIECE PER DISPOSITIVO "LIPOZERO G39" - SCHEDA SD4

PH 03. Ø 50 mm HANDPIECE

G1401	Ø 50 mm HANDPIECE, 39 kHz
-------	---------------------------



PH 04



PH 05

PH 04. POWER SUPPLY

G1362	"LIPOZERO G39" POWER SUPPLY
-------	-----------------------------

PH 05. POWER SUPPLY

G5216	POWER SUPPLY FOR "MEDISOUND 3000", LIPOZERO G800" AND "LIPOZERO G150"
-------	--



PH 06




PH 07

PH 06. POWER SUPPLY

G3566	POWER SUPPLY FOR "LIPOZERO EXCEL"
-------	-----------------------------------

PH 07. BATTERY PACK

G1120	BATTERY PACK FOR "MEDISOUND 3000", LIPOZERO G800" AND "LIPOZERO G150"
-------	--

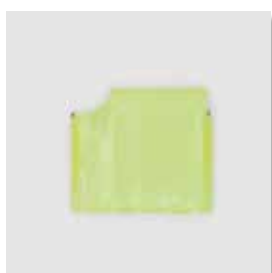

PH 08



PH 09

PH 08. BATTERY PACK

G0699	BATTERY PACK FOR "LIPOZERO EXCEL"
-------	-----------------------------------

F 09. BATTERY CHARGER

G2965	BATTERY CHARGER FOR "HT-922"
-------	------------------------------



PH 10

PH 10. GEL

G0484	260 ml GEL BOTTLE, 10 PCs
G0869	260 ml GEL BOTTLE, 1 PC
G4034	1000 ml GEL CAN
G3609	5000 ml GEL CAN

LASER THERAPY ACCESSORIES



PH 01

PH 01. HANDPIECES AND OPTICS FROM LEFT TO RIGHT

G5996	HANDPIECE 1064 nm		
G6090	Ø 11 mm OPTIC	G6092	Ø 22 mm OPTIC
G5993	HANDPIECE HP 808 nm		
G5995	HANDPIECE HP 980 nm		
G5994	MHANDPIECE MP 980 nm		
G5992	HANDPIECE MP 808 nm		



PH 02

PH 02. HANDPIECES FROM LEFT TO RIGHT

G5796	HANDPIECE GREEN 808 nm
G5748	HANDPIECE BLACK 808 nm
G5767	HANDPIECE WHITE 808 nm



PH 03

PH 03. LASER OPTIC

G5547	Ø 3 mm "SMALL" OPTIC ONLY FOR PODCARE 2.0
-------	---



PH 04



PH 05

PH 04. LASER GLASSES

G1462	OPERATOR LASER GLASSE
-------	-----------------------

PH 05. LASER GLASSES

G1461	PATIENT LASER GLASSES
-------	-----------------------



PH 06



PH 07

PH 06. POWER SUPPLY

G4146	LASER POWER SUPPLY FOR "PRO" VERSION
-------	--------------------------------------

PH 07. POWER SUPPLY

G4048	LASER POWER SUPPLY
-------	--------------------



PH 08



PH 09

PH 08. BATTERY PACK

G1127	BATTERY PACK FOR "PRO" VERSION
-------	--------------------------------

PH 09. BATTERY PACK

G0699	BATTERY PACK
-------	--------------

RADIOFREQUENCY ACCESSORIES



PH 01

PH 01. "RF CLINIC 100" HANDPIECE

G1420 | 2-SPHERE HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 02

PH 02. "RF CLINIC 300" HANDPIECE

G1512 | Ø 30 mm HEAD HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 03

PH 03. "RF CLINIC 500" HANDPIECE

G1511 | HEAD HANDPIECE Ø MM 50 HEAD HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 04

PH 04. CARRYING CASE

G1364 | ALUMINUM CASE, BLACK
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 05

PH 05. CARRYING CASE

G1365 | ALUMINUM CASE, RED FOR "RF BEAUTY MINI"



PH 06

PH 06. POWER SUPPLY

G4147 | POWER SUPPLY FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 07

PH 07. BATTERY CHARGER

G 1120 | BATTERY CHARGER FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 08

PH 08. BATTERY CHARGER

G3471 | BATTERY CHARGER FOR "HT-906"

PRESSOTHERAPY ACCESSORIES



PH 01

PH 01. PRESSCARE SLEEVE

G5642	PRESSCARE SLEEVE, 4 SECTIONS
-------	------------------------------



PH 02

PH 02. PRESSCARE BOOTS

G5635	"PRESSCARE", BOOT, 4 SECTIONS, SMALL
-------	--------------------------------------

G5631	"PRESSCARE", BOOT EXTENSION, SMALL
-------	------------------------------------

G5636	"PRESSCARE", BOOT, 4 SECTIONS, MEDIUM
-------	---------------------------------------

G5633	"PRESSCARE", BOOT EXTENSION, MEDIUM
-------	-------------------------------------

G5637	"PRESSCARE", BOOT, 4 SECTIONS, LARGE
-------	--------------------------------------

G5634	"PRESSCARE", BOOT EXTENSION, LARGE
-------	------------------------------------



PH 03

PH 03. PRESSCARE ABS/GLUTES BAND

G5641	"PRESSCARE", ABS/GLUTES BAND, 4 SECTIONS
-------	--



PH 04



PH 05

PH 04. Y CABLE

G5639	Y CABLE TO USE 3 ACCESSORIES AT THE SAME TIME
-------	---

F 05. CONNECTORS

G5640	SINGLE CONNECTOR FOR "PRESSCARE G300M" AND "PRESSCARE G200M"
-------	--

G5638	SINGLE CONNECTOR FOR "PRESSCARE G-SPORT 3"
-------	--



PH 06


PH 06. BAG

G5927	BAG FOR "PRESSCARE G-SPORT 3"
-------	-------------------------------


PRESSOTHERAPY ACCESSORIES

SIZE OF THE ACCESSORIES (4 SECTIONS)

BOOTS SIZE "S"	(A) LENGTH GROIN - HEEL 72 cm (B) THIGH CIRCUMFERENCE 60 cm with EXTENSION 70 cm (C) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "M"	(A) LENGTH GROIN - HEEL 82 cm (B) THIGH CIRCUMFERENCE 67 cm with EXTENSION 77 cm (C) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "L"	(A) LENGTH GROIN - HEEL 92 cm (B) THIGH CIRCUMFERENCE 74 cm with EXTENSION 84 cm (C) ANKLE CIRCUMFERENCE 40 cm with EXTENSION 45 cm
SLEEVE	(A) 84 cm - (B) 63 cm (C) SHOULDER CIRCUMFERENCE 43 cm with EXTENSION 53 cm (D) WRIST CIRCUMFERENCE 31 cm with EXTENSION 36 cm

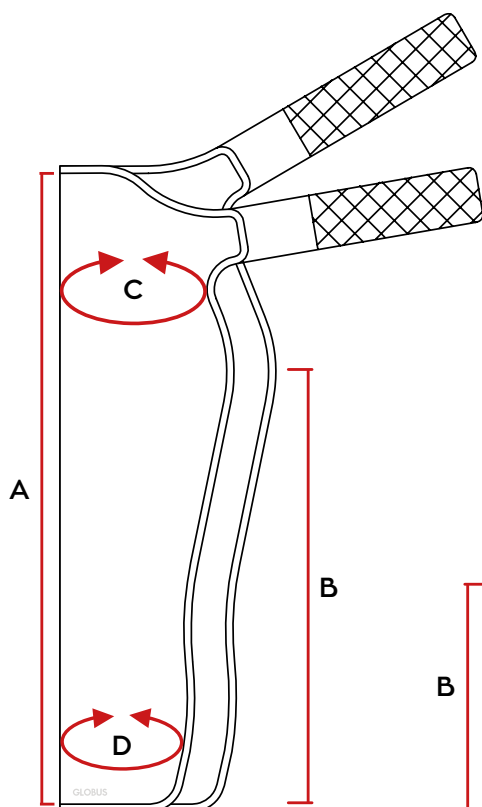

ABDOMINAL BAND	(A) 133 cm - (B) 38 cm ABDOMEN CIRCUMFERENCE 82 cm to 114 cm
----------------	---

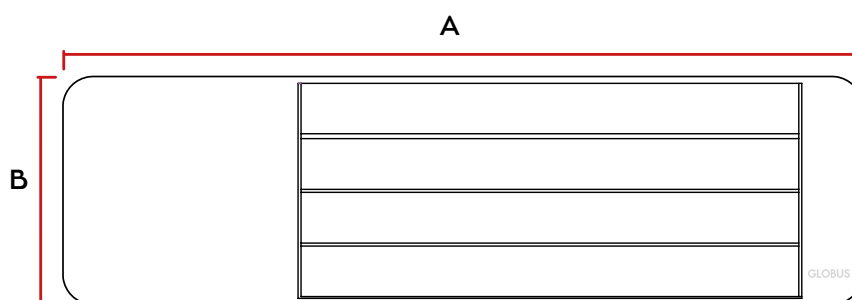

Y CABLE to use three accessories at the same time	SINGLE CABLE
---	--------------

NOTE All **EXTENSIONS** are to be considered an optional accessory available on request.

The **ACCESSORIES**, depending on the type of configuration, can be optional or included.

BALL-SHOOTING MACHINES ACCESSORIES



PH 01

PH 01. SPEED MEASURING GUN

G3048	SPEED MEASURING GUN
-------	---------------------



PH 02

PH 02. "EUROGOAL" LINE POWER SUPPLY

G1488	"EUROGOAL" LINE POWER SUPPLY
-------	------------------------------



PH 03

PH 03. "EUROGOAL" LINE BATTERY CHARGER

G1486	"EUROGOAL" LINE BATTERY CHARGER
-------	---------------------------------



PH 04

PH 04. BATTERY HOLDER AND BATTERY

G2946	BATTERY HOLDER
G3808	RECHARGEABLE BATTERY 12 V - 55 Ah



PH 05

PH 05. "WINSHOT" LINE BATTERY KIT

G6019	BATTERY KIT FOR "WINSHOT 1500" AND "WINSHOT 1000"
G2955	BATTERY KIT FOR "WINSHOT 800" AND "WINSHOT 500"

BALL-SHOOTING MACHINES ACCESSORIES



PH 06. "EUROGOAL" COVER

G3063 | COVER

PH 06



PH 07. "WINSHOT 1500" FRONT PROTECTION

G5073 | KIT FRONT PROTECTION

PH 07



VETERINARY

Line

CONTENTS

- TECARTHERAPY	168
- MAGNETOTHERAPY	170
- ELECTROTHERAPY	172
- ULTRASOUND THERAPY	176
- LASER THERAPY	178
- ACCESSORIES	182

TECARVET 4000 | TECARVET 4000 RE

TECARTHERAPY

TECARVET 4000 IS A DEVICE THAT USES DIATHERMY TO OBTAIN EXTRAORDINARY THERAPEUTIC EFFECTS ON SMALL ANIMALS SUCH AS DOGS AND CATS. THE DEVICE TRANSFERS ENERGY IN BOTH CAPACITIVE AND RESISTIVE MODE. IN THIS WAY THE EFFECT OF THE TREATMENT CAN BE DIRECTED WITH GOOD EFFECTIVENESS BOTH ON DEEP AND MORE SUPERFICIAL BIOLOGICAL TISSUES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Dimension 4.3" touch-screen

Power

160 W, $\pm 20\%$

Emission frequencies

470 KHz, $\pm 10\%$

Processing methods

Capacitive and resistive

Capacitive electrodes

3 \varnothing (30, 50, 70 mm)

Resistive electrodes

3 \varnothing (30, 50, 70 mm)

Mains power supply

100-240 V AC, 50-60 Hz, 1.8 A

Size


170x220x60 mm

Weight

1.04 Kg

STANDARD EQUIPMENT

- 1 Professional carrying case
- 1 TECARVET 4000 device
- 1 Counter electrode (plate)
- 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)
- 1 250 ml cream bottle
- 3 Elastic bands
- 1 Power supply
- 1 Operating manual

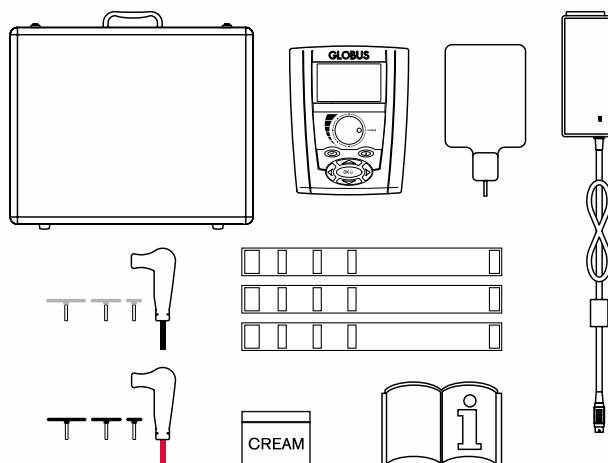

PURCHASE OPTIONS

- Code **G3865** TECARVET 4000
- Code **G3973** TECARVET 4000 RE *

* RECHARGEABLE TECARTHERAPY:

- 5 HOURS OF TREATMENT INCLUDED IN THE COST OF THE DEVICE
- CHEAP CHARGING
- GUARANTEED SUPPORT FOR 5 YEARS AND FORKLIFT
- CHARGING VIA INTERNET



STRONG POINTS



TECARVET 2000

TECARTHERAPY

THE TECARVET 2000 MODEL PROVIDES ALL THE QUALITY OF THE GLOBUS TECARTHERAPY IN A MACHINE OF SMALL SIZE AND WEIGHT, TO GUARANTEE MAXIMUM PORTABILITY AND KEEP ITS EFFECTIVENESS AND PROFESSIONALISM UNCHANGED. THE SPECIAL RESISTIVE BIPOLAR HANDPIECE ALLOWS TO SPEED UP THE PREPARATION FOR TREATMENT.

SPECIFIC PROGRAMS FOR DOGS AND CATS

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Power

50 W, $\pm 10\%$

Emission frequencies

470 KHz, $\pm 10\%$

Processing methods

Resistive bipolar

Power supply

- Mains
- Battery pack

Size


170x220x60 mm

Weight

1.04 Kg

STANDARD EQUIPMENT

- Code G5617
- 1 Professional carrying case
- 1 TECARVET 2000 device
- 1 Ø 50 mm bipolar handpiece
- 1 260 ml gel bottle
- 1 Power supply
- 1 Operating manual

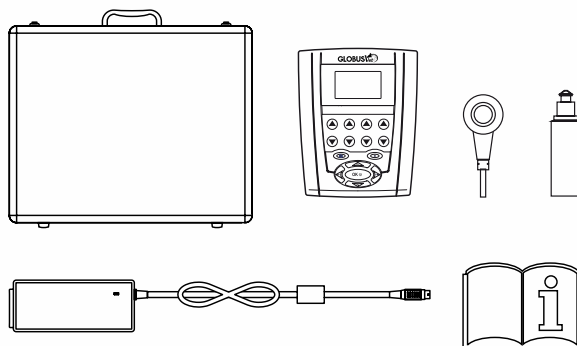

CE



STRONG POINTS



MAGNETOVET 4000

MAGNETOTHERAPY

MAGNETOVET 4000 IS A HIGH INTENSITY PULSED ELECTROMAGNETIC FIELD MAGNETOTHERAPY DEVICE. THE POSSIBILITY OF USING FOUR OUTPUT CHANNELS ALLOWS TO CONNECT MANY ACCESSORIES AND TO TREAT MORE BODY PARTS AT THE SAME TIME OR EVEN TWO ANIMALS SIMULTANEOUSLY.

MAGNETOTHERAPY HAS COUNTLESS AND PROVEN BENEFICIAL EFFECTS ON BONE AND ARTICULAR TISSUE. THE ELECTROMAGNETIC PULSES STIMULATE THE REGENERATION OF NEW CELLS AND THE REMODELING OF EXISTING STRUCTURES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Channels

4

Intensity

200 average Gauss per channel with flexible solenoids

800 total average Gauss (for 4 channels)

Frequency

Adjustable from 5 to 200 Hz

Duty Cycle

From 5% to 50% in steps of 5

Power supply

- Mains
- Rechargeable battery (not included)

Size


170x220x60 mm

Weight

820 gr


STANDARD EQUIPMENT

- Code G5628
- 1 Bag
- 1 MAGNETOVET 4000 device
- 4 Flexible solenoids 300x100 mm
- 1 Power supply
- 1 Operating manual

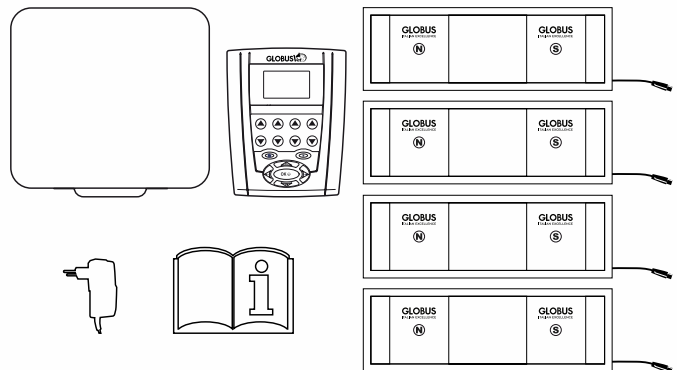

PURCHASE OPTION

- Code G6039
- Therapeutic bed in memory foam, 4 coils, 700x440 mm for small dogs and cats. For large dogs, buy 2 single beds to join with the velcro.

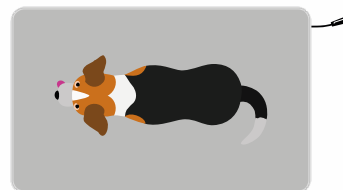

STRONG POINTS



MAGNETOVET 200 with THERAPEUTIC BED

MAGNETOTHERAPY

MAGNETOVET 200 IS A PORTABLE MAGNETOTHERAPY DEVICE THAT CAN PRODUCE A HIGH-INTENSITY ELECTROMAGNETIC FIELD. THE PULSES PRODUCED HAVE A THERAPEUTIC AND REGENERATIVE EFFECT ON THE MAIN OSTEOARTICULAR PATHOLOGIES THAT AFFECT ANIMALS.

BEING VERY EASY TO USE IN COMBINATION WITH THE THERAPEUTIC BED, MAGNETOVET 200 IS PARTICULARLY SUITABLE FOR THE HOME TREATMENT OF CHRONIC PAIN.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Channels

2

Intensity

200 average Gauss per channel with flexible solenoids

Frequency

Adjustable from 5 to 200 Hz

Duty Cycle

From 5% to 50% in steps of 5

Power supply

- Mains
- Rechargeable battery (not included)

Size


100x160x35 mm

Weight

280 gr

STANDARD EQUIPMENT

- Code G6115
- 1 Bag
- 1 MAGNETOVET 200 device
- 1 Therapeutic bed mm 700x440
- 2 Flexible solenoids 300x100 mm
- 1 Soft fleece blanket
- 1 Power supply
- 1 Operating manual

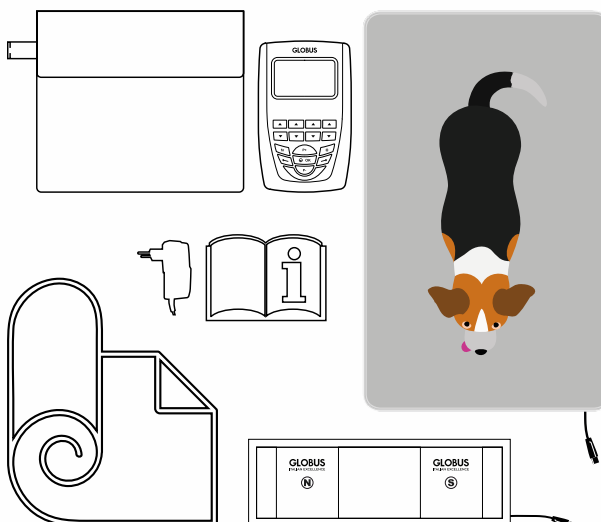

PURCHASE OPTIONS

- Code G5620 Therapeutic bed not included
- Code G6039 Memory foam therapeutic bed, 4 coils, 700x440 mm for small dogs and cats. For large dogs, buy 2 single beds to join with the velcro.



STRONG POINTS



STIMVET 4000

4-CHANNEL ELECTROTHERAPY

STIMVET 4000 IS A PROFESSIONAL ELECTROTHERAPY DEVICE SUITABLE FOR VETERINARY REHABILITATION. THE DEVICE INCLUDES CURRENTS FOR THE RECOVERY OF MUSCULAR ATROPHY AND TENS CURRENTS FOR ANTALGIC TREATMENT OF CHRONIC AND ACUTE PAIN. IN ADDITION, THERE ARE CURRENTS TO TREAT DENERVATED MUSCLES AND MICROCURRENTS INDICATED MAINLY IN CARTILAGE LESIONS AND WOUNDS. THE DEVICE IS FULLY EQUIPPED WITH PRE-SET PROGRAMS FOR THE TREATMENT OF DOGS, CATS AND HORSES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150 Hz

Pulse width

50-450 μ s

Intensity

100 mA per channel

Power supply

- Mains
- Rechargeable battery

Size

170x220x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

· Code G5500

· 1 Bag

· 1 STIMVET 4000 device

· 4 Elastic bands

· 2 Gray cables for microcurrents and ionophoresis

· 4 Colored cables for electrode

· 4 Silicon square electrodes


· 4 Silicon rectangular electrodes

· 4 Ionophoresis sponges

· 1 260 ml gel bottle

· 1 Power supply


· 1 Operating manual

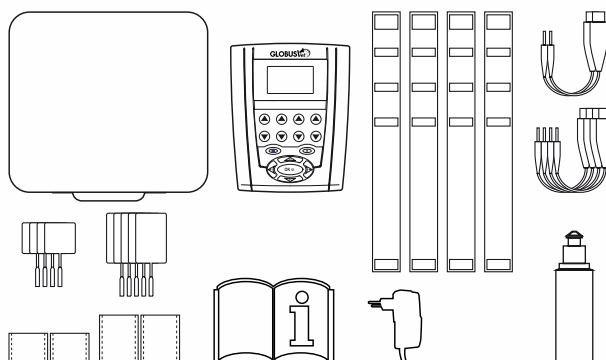

CE



STRONG POINTS

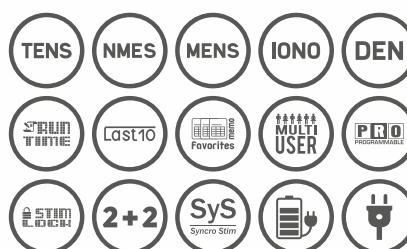

STIMVET 2000

4-CHANNEL ELECTROTHERAPY

STIMVET 2000 IS A PORTABLE ELECTROSTIMULATOR SUITABLE FOR THE MUSCLE REHABILITATION AND THE ANTALGIC TREATMENT OF DOGS, CATS AND HORSES. THE DEVICE IS ALSO SUITABLE FOR SPORTING OR WORKING DOGS WHO NEED TO INCREASE THEIR ATHLETIC SKILLS. STIMVET 2000 ALSO INCLUDES SOME MICROCURRENT PROGRAMS FOR WOUND CICATRIZATION.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150 Hz

Pulse width

50-450 μ s

Intensity

100 mA per channel

Power supply

- Mains
- Rechargeable battery

Size


100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT


- Code **G5614**
- 1 Bag
- 1 STIMVET 2000 device
- 4 Elastic bands
- 2 Gray cables for microcurrents
- 4 Colored cables for electrode
- 4 Silicon square electrodes
- 4 Silicon rectangular electrodes
- 1 260 ml gel bottle
- 1 Power supply
- 1 Operating manual

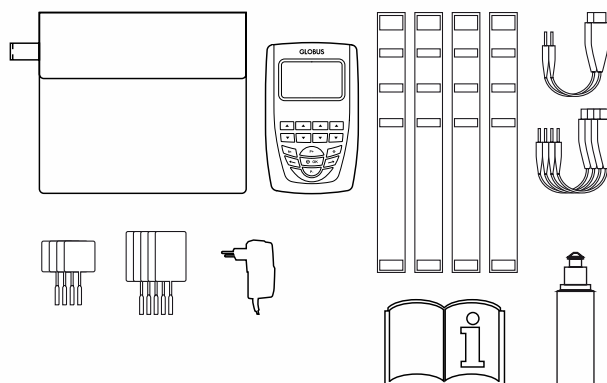

CE



STRONG POINTS

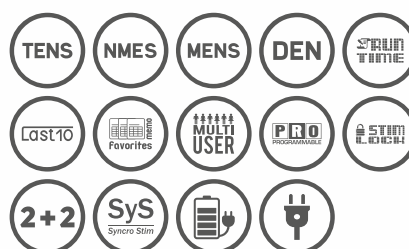

STIMVET 200

4-CHANNEL ELECTROTHERAPY

STIMVET 200 IS A PORTABLE ELECTROSTIMULATOR SUITABLE FOR THE MUSCLE REHABILITATION AND THE ANTALGIC TREATMENT OF DOGS AND CATS.

MUSCLE STRENGTHENING PROGRAMS ARE SUITABLE BOTH FOR ATROPHY RECOVERY AND FOR SPORTING OR WORKING DOGS WHO NEED TO INCREASE THEIR ATHLETIC SKILLS. TENS CURRENTS, INSTEAD, ARE SUITABLE FOR CHRONIC AND ACUTE PAIN AND MAKE THE PET FEEL BETTER WHILE REDUCING THE USE OF DRUGS.

SPECIFIC PROGRAMS FOR DOGS AND CATS

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150 Hz

Pulse width

50-450 μ s

Intensity

80 mA per channel

Power supply

Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code G5608

· 1 Bag

· 1 STIMVET 200 device

· 4 Elastic bands

· 2 Gray cables for microcurrents

· 4 Colored cables for electrode


· 4 Silicon square electrodes

· 4 Silicon rectangular electrodes


· 1 260 ml gel bottle

· 1 Power supply

· 1 Operating manual

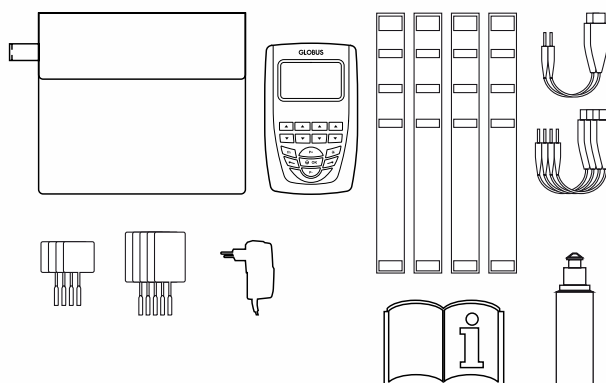

STRONG POINTS

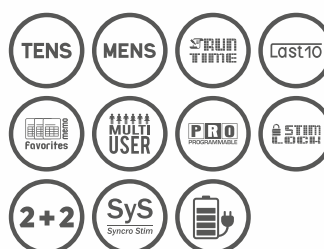



ULTRASOUNDVET 4000

ULTRASOUND THERAPY

THE THERAPY WITH ULTRASOUNDVET 4000 IS EFFECTIVE IN THE TREATMENT OF TENDON, MUSCLE AND LIGAMENT DISEASES OF DOGS, CATS AND HORSES.

IT IS ALSO USED ON SCAR ADHESION AND JOINT STIFFNESS. THIS INSTRUMENTAL PHYSICAL THERAPY PROVIDES THERMAL AND ATHERMAL MECHANICAL STIMULATION TO FACILITATE THE HEALING OF THE DAMAGED TISSUES OF MUSCLES, TENDONS, ARTICULAR CAPSULES, LIGAMENTS AND BONE.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Outputs

1 Ø 42 mm handpiece

Frequency

1 MHz and 3 MHz

Power

3 W/cm², ±20%

Duty Cycle

From 10% to 100% in steps of 10%

Emission

Continuous and pulsed

Mode

Direct contact and immersion

Power supply

- Mains
- Rechargeable battery (not included)

Size

170x220x60 mm

Weight

820 gr

STANDARD EQUIPMENT

· Code G3869

· 1 Bag


· 1 ULTRASOUNDVET 4000 device

· 1 Ø 42 mm handpiece - IP 68

· 1 260 ml gel bottle

· 1 Power supply

· 1 Operating manual

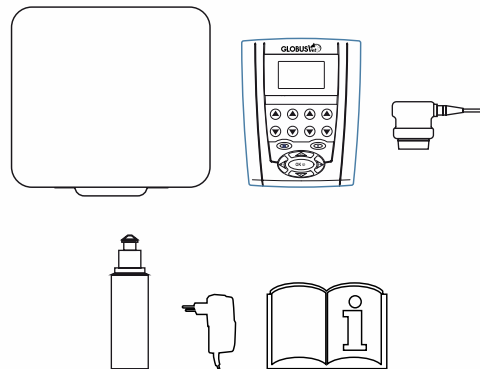

CE



STRONG POINTS



ULTRASOUNDVET 200

ULTRASOUND THERAPY

THE THERAPY WITH ULTRASOUNDVET 200 IS EFFECTIVE IN THE TREATMENT OF TENDON, MUSCLE AND LIGAMENT DISEASES OF PETS. IT IS ALSO USED ON SCAR ADHESION AND JOINT STIFFNESS.

SPECIFIC PROGRAMS FOR DOGS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Outputs

1 Ø 42 mm handpiece

Frequency

1 MHz

Power

1,4 W/cm², ±20%

Duty Cycle

From 10% to 100% in steps of 10%

Emission

Continuous and pulsed

Mode

Direct contact and immersion

Power supply

- Mains
- Rechargeable battery (not included)

Size


100x160x35 mm

Weight

280 gr

STANDARD EQUIPMENT

- Code G5626
- 1 Bag
- 1 ULTRASOUNDVET 200 device
- 1 Ø 42 mm handpiece - IP 68
- 1 260 ml gel bottle
- 1 Power supply
- 1 Operating manual

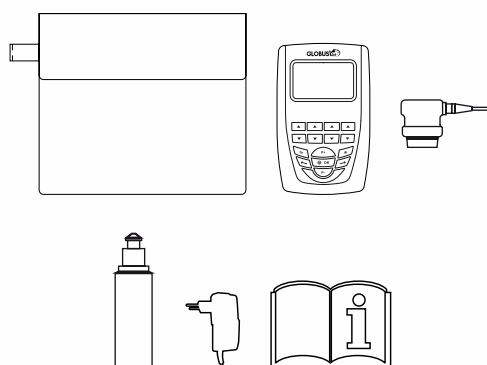

CE



STRONG POINTS



LASERVET 15.0

LASER THERAPY

LASERVET 15.0 IS A PORTABLE DEVICE FOR LASER THERAPY WITH HIGH PULSED POWER THAT ALLOWS TO PERFORM VETERINARY TREATMENTS QUICKLY. THE PATHOLOGIES TO TREAT ARE DIVIDED BY ANIMAL WITH GREAT ATTENTION TO EQUINES. THE ACTION OF THE LASER IS AIMED AT TISSUE REGENERATION, POST-TRAUMATIC AND ANTI-INFLAMMATORY TREATMENT AND THE REDUCTION OF BOTH CHRONIC AND ACUTE PAIN.

SPECIFIC PROGRAMS FOR DOG, CAT, HORSE, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

15 W $\pm 20\%$

Compatible with a wavelength of

808 nm

980 nm and 1064 nm (optional)

Emission

Continuous and pulsed

Available optics

\varnothing 11 mm and \varnothing 22 mm optics

\varnothing 60 mm optic (optional)

Frequency

Up to 10000 Hz

Power supply

- Mains

- Rechargeable battery

Size and weight

220x170x60 mm

Weight

1.12 Kg

STANDARD EQUIPMENT

· Code G6036

· 1 Bag


· 1 LASERVET 15.0 device

· 1 Handpiece HP 808 nm with \varnothing 11 mm and \varnothing 22 mm optics

· 2 Protection glasses

· 1 Power supply

· 1 Operating manual

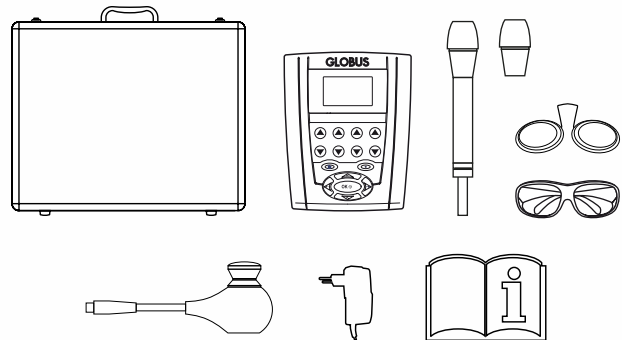

PURCHASE OPTIONS

· Code G6348 Handpiece HP 980 nm

· Code G5996 Handpiece 1064 nm

· Code G6089 \varnothing 60 mm optic



STRONG POINTS



LASERVET 8.0

LASER THERAPY

LASERVET 8.0 IS A PORTABLE DEVICE FOR LASER THERAPY WITH HIGH PULSED POWER THAT ALLOWS TO PERFORM VETERINARY TREATMENTS. THE ACTION OF THE LASER AIMS AT TISSUE REGENERATION, POST-TRAUMATIC AND ANTI-INFLAMMATORY TREATMENT AND THE REDUCTION OF CHRONIC AND ACUTE PAIN.

SPECIFIC PROGRAMS FOR DOG, CAT, HORSE, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

8 W $\pm 20\%$

Compatible with a wavelength of

808 nm

980 nm and 1064 nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm and Ø 60 mm optics (optional)

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code G6035

· 1 Bag


· 1 LASERVET 8.0 device

· 1 Handpiece MP 808 nm with Ø 11 mm optic

· 2 Protection glasses

· 1 Power supply

· 1 Operating manual

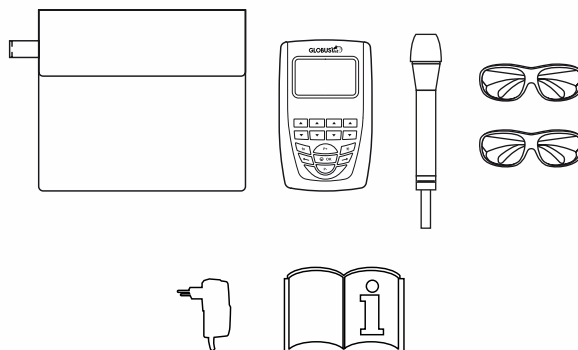

OPZIONI DI ACQUISTO

· Code G6347 Hanpiece MP 980 nm

· Code G5996 Hanpiece 1064 nm

· Code. G6089 Ø 60 mm optic

· Code G6092 Ø 22 mm optic



STRONG POINTS



LASERVET 3.0

LASER THERAPY

THE TREATMENT WITH LASERVET 3.0 ACCELERATES THE HEALING OF DAMAGED TISSUES AND IMPROVES CIRCULATION. IT IS IDEAL TO STIMULATE THE CICATRIZATION OF SKIN WOUNDS AND IS PARTICULARLY EFFECTIVE FOR THE TREATMENT OF DIFFERENT MUSCULOSKELETAL PATHOLOGIES CHARACTERIZED BY INFLAMMATION AND PAIN. LASER THERAPY ELIMINATES PAIN FASTER, BOTH USED ALONE AND WITH OTHER THERAPEUTIC METHODS..

SPECIFIC PROGRAMS FOR DOG, CAT, HORSE, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

3 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optic

\varnothing 11 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code G6192

· 1 Bag


· 1 LASERVET 3.0 device

· 1 Handpiece Black 808 nm with \varnothing 11 mm optic

· 2 Protection glasses

· 1 Power supply

· 1 Operating manual

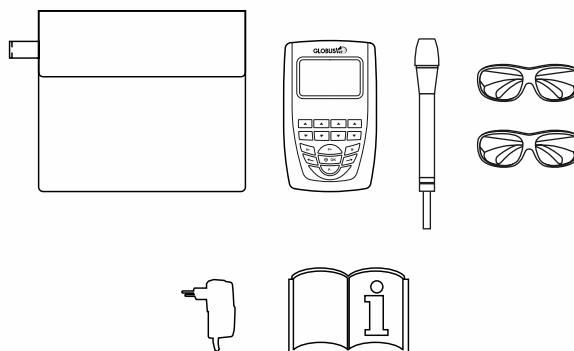

CE



STRONG POINTS



LASERVET 1.0

LASER THERAPY

LASERVET 1.0 IS THE ONLY CLASS 4 LASER IN THE GLOBUSVET CATALOG. THE DEVICE CAN REACH 1 WATT OF POWER AND IS FULLY EQUIPPED WITH PRE-SET PROGRAMS FOR DOGS, CATS AND HORSES. MOREOVER, IT GIVES THE POSSIBILITY TO PROGRAM AND STORE NEW PROGRAMS.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

1 W $\pm 20\%$

Compatible with a wavelength of

808 nm

Emission

Continuous and pulsed

Available optic

Ø 11 mm optic

Frequency

Up to 10000 Hz

Power supply

- Mains
- Rechargeable battery

Size

100x160x35 mm

Weight

450 gr

STANDARD EQUIPMENT

· Code G5629

· 1 Bag


· 1 LASERVET 1.0 device

· 1 Handpiece White 808 nm with Ø 11 mm optic

· 2 Protection glasses

· 1 Power supply

· 1 Operating manual

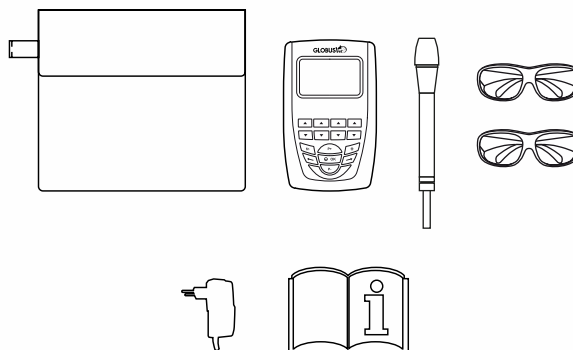

STRONG POINTS



ACCESSORIES

VET Line



THERAPEUTIC BED FOR MAGNETOTHERAPY

INDICATED FOR **ARTHROSIS** AND **ARTICULAR PAIN**

TECARTHERAPY

G3483	MONOPOLAR CAPACITIVE HANDPIECE
G5393	NEUTRAL FLEXIBLE PLATE 250X200 mm
G3485	MONOPOLAR RESISTIVE HANDPIECE
G3486	BIPOLAR HANDPIECE, Ø 30 mm
G3487	BIPOLAR HANDPIECE, Ø 50 mm
G3306	CAPACITIVE ELECTRODE, Ø 30 mm
G3308	ECAPACITIVE ELECTRODE, Ø 50 mm
G3310	CAPACITIVE ELECTRODE, Ø 70 mm
G5337	RESISTIVE ELECTRODE, Ø 30 mm
G5338	RESISTIVE ELECTRODE, Ø 50 mm
G5339	RESISTIVE ELECTRODE, Ø 70 mm
G3368	TECAR CREAM, 250 ml
G3369	TECAR CREAM, 1000 mm
G3503	NEUTRAL PLATE CABLE
G3517	POWER SUPPLY 8 V FOR "TECARVET 4000" AND "TECARVET 2000"
G1120	BATTERY PACK, 4000 mA, 12 V FOR "TECARVET 2000"
G1364	BLACK ALUMINUM CASE FOR "TECARVET 2000"
G3370	CREAM DISPENSER AND CAP
G3371	GRAY ALUMINUM CASE FOR "TECARVET 4000"
G1364	BLACK ALUMINUM CASE FOR "TECARVET 2000"
G6341	"TECARVET" LINE HOLDER FOR ACCESSORIES
G6342	"TECARVET" LINE HOLDER FOR DEVICES

MAGNETOTHERAPY

G6097	FLEXIBLE SOLENOIDE, 300X100 mm
G1408	SFLEXIBLE SOLENOID, 400X150 mm, 4 COILS
G4275	SINGLE SOFT SOLENOID 120X120 mm, 1 COIL
G4276	"SOFT 4" MAT 250X250 mm, 2 CHANNELS, 4 COILS
G6039	MEMORY FOAM THERAPEUTIC BED EXTENSION 440X700 mm, 4 COILS
G6134	GLOBUSVET SOFT FLEECE BLANKET
G5878	FLEXIBLE SOLENOID COVER KIT, 10 PCs
G5879	SOFT SOLENOID COVER KIT, 10 PCs
G5867	MAGNUM TEST, INDICATOR OF CORRECT EMISSION OF THE MAGNETIC FIELD
G4048	POWER SUPPLY FOR "MAGNETOVET 200"
G4147	POWER SUPPLY FOR "MAGNETOVET 4000"
G0699	BATTERY PACK FOR "MAGNETOVET 200"
G1127	BATTERY PACK FOR "MAGNETOVET 4000"
G4337	Y CABLE



PH 01



PH 02

FOR ALL DEVICES AND ACCESSORIES

PH 01. TROLLEY

G5521 | TROLLEY WITH MULTIPLE SPACES TO TRANSPORT THE DEVICES

PH 02. BACKPACK

G5927 | GLOBUS BACKPACK, 30x45x15 mm (WxHxD)

ELECTROTHERAPY

G0136	MYOTRODE PLATINUM ELECTRODES, 50X50 mm, 4 PCs
G0137	MYOTRODE PLATINUM ELECTRODES, 50X90 mm, 4 PCs
G0464	MYOTRODE PLUS ELECTRODES, 50X50 mm, 4 PCs
G0465	MYOTRODE PLUS ELECTRODES, 50X90 mm, 4 PCs
G0451	ELECTRODES, Ø 32 mm, WITH CABLE, 4 PCs
G0461	CONDUCTIVE SILICONE ELECTRODE PER METER
G0463	KIT 6 CONDUCTIVE SILICONE ELECTRODES
G0890	MEDIUM CARBON ELECTRODE FOR IONOPHORESIS, WITH POUCH, 50X50 mm
G0885	LARGE CARBON ELECTRODE FOR IONOPHORESIS, WITH POUCH, 60X85 mm
G0155	POUCH FOR LARGE IONOPHORESIS ELECTRODE, 60X85 mm
G1210	POUCH FOR MEDIUM IONOPHORESIS ELECTRODE, 50X50 mm
G1839	COLORLED CABLES KIT FOR 4-CHANNEL DEVICES, 4 PCs
G0437	SPARE CABLE, WHITE, FOR 4-CHANNEL DEVICES
G0439	KIT 2 Y CABLES, UNIVERSAL
G0441	KIT 2 Y CABLES FOR 4-CHANNEL DEVICES
G0484	GLOBUS GEL, 260 ml, 10 PCs
G0869	GLOBUS GEL, 260 ml - 1 PC
G4034	GLOBUS GEL, 1000 ml - 1 PC
G3609	GLOBUS GEL, 5000 ml - 1 PC
G1156	MOTOR POINT PEN
G1225	ELASTIC BANDS SIZE "S", 800X80 mm, 2 PCs
G1226	ELASTIC BANDS SIZE "M", 400X80 mm, 2 PCs
G4249	ADAPTER FOR ACUPUNCTURE, BLACK
G4250	ADAPTER FOR ACUPUNCTURE, RED
G1127	BATTERY PACK FOR "STIMVET 4000"
G0699	BATTERY PACK FOR "STIMVET 200"
G4154	BATTERY CHARGER FOR "STIMVET 200"
G4146	POWER SUPPLY FOR "STIMVET 4000"

ULTRASOUND THERAPY


G1402	HANDPIECE Ø 42 mm, FOR "ULTRASOUNDVET 4000" - SD4 BOARD
G5374	HANDPIECE Ø 42 mm, FOR "ULTRASOUNDVET 4000" - SD5 BOARD
G3152	HANDPIECE Ø 42 mm, FOR "ULTRASOUNDVET 200"
G5516	HANDPIECE Ø 16 mm, FOR "ULTRASOUNDVET 4000" - SD5 BOARD
G4060	HANDPIECE Ø 16 mm, FOR "ULTRASOUNDVET 4000" - SD4 BOARD
G0484	GLOBUS GEL, 260 ml, 10 PCs
G0869	GLOBUS GEL, 260 ml - 1 PC
G4034	GLOBUS GEL, 1000 ml - 1 PC
G3609	GLOBUS GEL, 5000 ml - 1 PC
G1120	BATTERY PACK FOR "ULTRASOUNDVET 4000"
G5216	POWER SUPPLY FOR "ULTRASOUNDVET 4000"
G4093	POWER SUPPLY FOR "ULTRASOUNDVET 200"

LASER THERAPY

G6348	HANDPIECE, HP 980 nm	G6089	Ø 60 mm OPTIC
G6347	HANDPIECE, MP 980 nm	G6092	Ø 22 mm OPTIC
G6194	HANDPIECE, HP 808 nm	G6090	Ø 11 mm OPTIC
G6282	HANDPIECE, MP 808 nm		
G5996	HANDPIECE, 1064 nm		
G6208	HANDPIECE BLACK 808 nm		
G6275	HANDPIECE WHITE 808 nm		
G1462	OPERATOR PROTECTION GLASSES - IPL		
G4048	POWER SUPPLY FOR LASERVET 1.0 / 3.0 / 8.0		
G4146	POWER SUPPLY FOR LASERVET 15.0		

The background of the entire page is a dark blue gradient. Overlaid on this are several abstract geometric patterns. These consist of thin, light blue lines that connect small, glowing blue dots. The dots vary in brightness, with some appearing as sharp points of light and others as softer glows. The lines form a complex web of triangles and other polygons, some of which are more prominent than others. The overall effect is one of a dynamic, interconnected network or a stylized representation of a celestial map.

Key **STRONG POINTS**

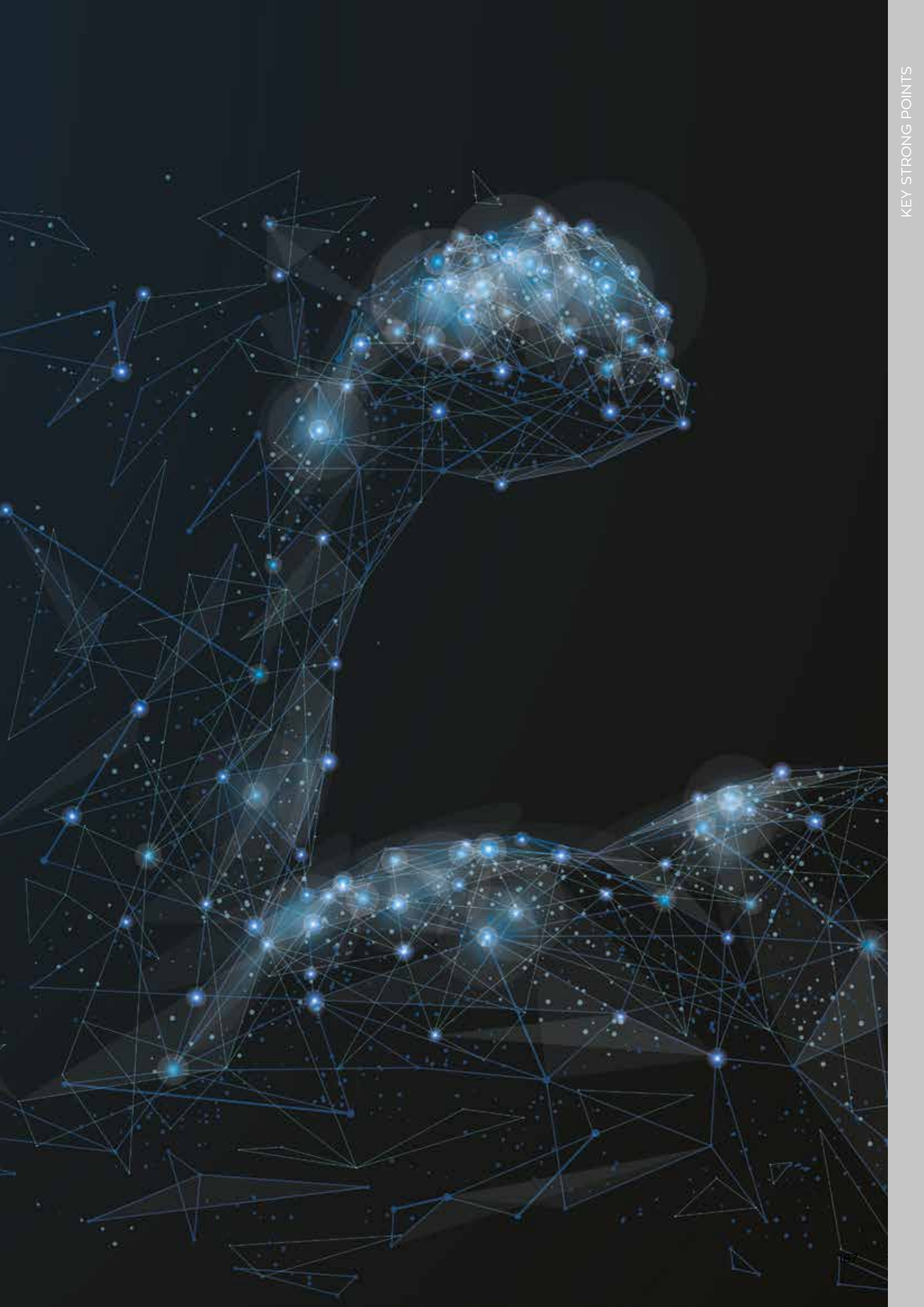

KINEO LINE



ISOTONIC NIT (NO INERTIA)

The NIT Isotonic method introduces an isotonic training which, by eliminating inertia, allows to set a constant load for the entire range of motion. In the training without inertia, the load is constant throughout the movement regardless of speed. Muscle activation is complete, from the beginning to the end of the movement, and the fixed load allows to perform every workout in the best possible way with maximum safety and never losing contact.



ISOTONIC WIT (WITH INERTIA)

In the innovative Kineo system, a robotic system allows to work with the Isotonic method, controlling the load and the work parameters completely and using different loads in the 2 phases of the movement, concentric and eccentric, for a more functional and effective treatment.



ELASTIC

The Elastic method adapts the resistance of the exercise according to your range of motion in a functional way: minimum load in the unfavorable points, higher load in the optimal ones. Setting an elastic with Kineo is simple: define the starting and final load. Your custom elastic will automatically be created.



SMART VISCOUS

The Kineo VISCOUS method is based on the principle of training in the water as a working method in common strengthening exercises. The value of the viscous load is determined by 2 parameters: the set viscosity level and the speed of movement during the exercise. Kineo allows to use 6 levels of viscosity that simulate workout in different fluids, such as water, oil, honey and more, with different levels of muscle commitment.



VARIABLE LOAD

In traditional Isotonic training (weights, barbells, machines), the workload consists of a weight that always remains the same during the exercise. Kineo, instead, allows to set a load that varies and follows a trend (curve) which reflects the real strength capacity of the subject. It also protects and safeguards the joints in the critical work angles, while producing the maximum strength demand and the maximum training effect in the most favorable work angles.



ISOMETRIC

This method consists in producing a muscle contraction, even maximal, without shortening the muscles and moving the affected joint. Kineo can program isometric exercises on any working angle for a more complete strength development. The Biofeedback system also allows to control the strength expressed (in real time) to reach the training Target. To overcome the limitations of the traditional isometric method, Kineo has also introduced the measurement of isometric strength in relation to the time taken to develop it. This method is mainly used in functional rehabilitation because it is safe, adjustable, it does not produce any joint movement and stress, and it improves the recruitment of muscle fibers.



ISOKINETIC

The isokinetic method is characterized by movement at a constant speed, regardless of the strength applied: this allows the rehabilitation process to be carried out in total safety. With Kineo, the functionalities of the isokinetic method were implemented with isokinetic exercises in open and closed kinetic chain and with differentiated load in the concentric and the eccentric phase.

DIACARE | TECAR BEAUTY LINE



TREATMENT ATLAS

Wide library of preset programs divided into anatomical areas with illustrative images.



HANDS FREE TREATMENT

The program facilitates the absorption of the active ingredients of a drug or cream.



MULTI-FREQUENCY

Different working frequencies can be selected to treat the patient accurately at different depths.



TRANSDERMAL DELIVERY

The program facilitates the absorption of the active ingredients of a drug or cream.



AUTOMATIC SCANNING OF FREQUENCIES

During the same treatment, the program can transfer energy evenly at different depths.



ATHERMAL TREATMENTS

Biostimulation of tissues without heat production to treat pathologies that cannot be treated at high temperatures.



PULSED TREATMENTS

Pulsed emission programs with different duty cycles, to keep the energy transferred to the tissues high while maintaining a good temperature control.



FASCIA TOOLS

Special ergonomic tools designed for the mobilization and the advanced treatment of all soft tissue disorders.



RESISTIVE METHOD

An unshielded steel handpiece is used to work on stronger and deeper tissues such as bones, very tight muscles and ligaments.



CAPACITIVE

A steel handpiece with shielding is used concentrating the ions under the handpiece to a maximum depth of 2/3 cm. The tissues that can be treated with the tecar capacitive handpiece are the most superficial and soft ones, such as muscle tissue, loose connective tissue or the lymphatic system.



BIPOLAR

A single handpiece with two built-in poles.



COLOR DISPLAY



RECHARGE

Function that allows to buy the device at an advantageous price, recharging the minutes according to the needs.



TOUCH SCREEN DISPLAY



MAINS POWER SUPPLY

REHAB | BEAUTY | SPORT | VET LINE

FOR **MAGNETOTHERAPY | ELECTROTHERAPY AND ELECTROSTIMULATION | ULTRASOUND THERAPY | LASER THERAPY**



TENS

Transcutaneous Electrical Nerve Stimulation (TENS) is performed by means of small electrical pulses to the painful area. The stimulation reaches the nerve endings, inhibiting the transmission of the painful sensation.



NMES

Neuromuscular Electrical Stimulation (NMES) is a current suitable for stimulating the muscles with the aim of improving muscle tone and trophism.



MENS

Microcurrent Electrical Nerve Stimulation (MENS) consists of tiny and imperceptible electrical pulses capable of relieving pain, stimulating tissues and promoting the healing of various types of injuries and osteoarticular pathologies.



IONO

Ionophoresis is a technique that allows the use of a direct current for the administration of a drug transcutaneously.



DENERVATED

Specific current for the stimulation of partially or totally innervated muscles.



3S

Serial Sequential Stimulation is a new stimulation method which, in the medical/physiotherapy field, allows the muscles involved to be stimulated in a sequential manner.



ACTION NOW

It allows to activate muscle contraction through an external control managed by an opera.



RUN TIME

This function allows to modify the time, the frequency and the pulse amplitude during the treatment.



LAST 10

The stimulator stores (as mobile phones) the last ten executed programs. In this way they are quickly available for the next execution.



FAVORITES

The therapists can store up to 30 programs for each registered user.



MULTI USER

This special function facilitates the use of the stimulator as it allows to store the personal settings of different users on the same device.



PROGRAMMABLE

The user can create some programs with personalized parameters.



STIM LOCK

To lock the device and carry out only the saved treatment.



1+1 FUNCTION

It allows you to use 2 different programs at the same time.



2+2 FUNCTION

It allows to use different programs and stimulate 2 parts of the body or 2 people at the same time.



MY TRAINER

This function, with some wave forms (EMS, TENS, MENS), shows the images of the electrode position directly on the display.



AUTO STIM

This function allows to execute an already executed program automatically, setting the intensity value used in the previous program.



SYNCHRO STIM

This function allows to increase or decrease the stimulation intensity of all the channels simultaneously using a single button.



WORK TIME

It indicates the total time the device has been used in terms of stimulation.



RENTAL OPTION

It allows you to lock the device keeping in memory only the number of treatments to be performed. Suitable for renting the device.



RECHARGEABLE BATTERY

The device is supplied with a rechargeable battery and charger to always have the device ready to use wherever you are.



POWER SUPPLY

The device can also operate on the mains.

GLOBUS

ITALIAN EXCELLENCE

DOMINO S.R.L.

Via Vittorio Veneto, 52
31013 Codognè (TV) - ITALY
Phone (+39) 0438 7933
globuscorporation.com

Information for sector operators only.

The images, features and data in this brochure are not binding and are subject to change by the manufacturer without prior notice.

For further information on our product range, treatments and our videos, visit the website:

www.globuscorporation.com

GLOBUS

ITALIAN EXCELLENCE

