



KINEO

MEDICAL

BEAUTY

SPORT

VET

CATALOGUE 2024



CONTENTS

ABOUT US	04
KINEO INTELLIGENT LOAD	06
MEDICAL	18
- TECAR THERAPY	20
- MAGNETOTHERAPY	28
- ELECTROTHERAPY	40
- ULTRASOUND THERAPY	50
- LASER THERAPY for physiotherapy	54
- LASER THERAPY for podiatry	64
- PRESSOTHERAPY	70
BEAUTY & WELLNESS	74
- TECAR BEAUTY	76
- ELECTROSTIMULATION	86
- RADIOFREQUENCY	98
- ULTRASOUND	102
SPORT	108
- ELECTROSTIMULATION	110
- PRESSOTHERAPY	118
- RELAX	120
- VIBRATION PLATFORM	124
- BALL SHOOTING MACHINE	126
ACCESSORIES	134
- TECAR THERAPY AND TECAR BEAUTY	136
- MAGNETOTHERAPY	144
- ELECTROTHERAPY AND ELECTROSTIMULATION	147
- ULTRASOUND THERAPY	151
- LASER THERAPY	152
- RADIOFREQUENCY	153
- PRESSOTHERAPY	154
- BALL SHOOTING MACHINE	156
VET	158
- TECAR THERAPY	160
- MAGNETOTHERAPY	162
- ELECTROTHERAPY	168
- ULTRASOUND THERAPY	172
- LASER THERAPY	174
- ACCESSORIES	178
KEY STRONG POINTS	182

ABOUT US

GLOBUS has been a world leading company in the production of portable electro-medical devices for over 30 years. Our wide range includes complete lines of products for electrotherapy, diathermy, ultrasound therapy, magnetotherapy, laser therapy, vibration platforms and a functional evaluation system.

OUR MISSION

The mission of **GLOBUS** is to contribute to promoting health and well-being in the world by providing high quality products to rehabilitation and physiotherapy professionals, as well as to private users who wish to improve their quality of life in total safety and benefit from the most modern technologies available on the market today.

OUR STORY

GLOBUS was founded in 1984 in Codognè, in the province of Treviso, by its current President and CEO, Pierpaolo Lucchetta, bronze in volleyball at the Los Angeles Olympics and award-winning professional until 1994.

The high-level experience of its President as an athlete has been fundamental for **GLOBUS** to design and develop different product lines for the treatment of several physical problems, muscle strengthening and the general improvement of sports performance. In order to develop safe, effective and highly technological products, **GLOBUS** collaborates with several universities and research institutes all over the world and scrupulously adheres to the most rigorous scientific guidelines.

Today, **GLOBUS** is a leading company in the rehabilitation and physiotherapy sector, exporting to more than 50 countries around the world. **GLOBUS** products are certified UNI CEI EN ISO 13485:2012 and CE of the Quality Assurance System in Europe, FDA in the United States, and IEx in Brazil.





KINEO INTELLIGENT LOAD

CONTENTS

- KINEO MULTISTATION	08
- KINEO PULLEY & SQUAT	09
- KINEO LEG PRO	10
- KINEO LEG PRESS	11
- KINEO CROSSOVER CABLE	12
- ACCESSORIES	14

KINEO MULTISTATION

TECHNICAL SPECIFICATIONS

· Code K0155M | Code K0155

Dimensions (WxDxH)

260x90x185 cm

Working area dimensions (WxDxH)

305x90x185 cm

Weight

342 kg

Electric motor

Rated power 1.6kW

Power supply

Voltage: 230V - Frequency: 50/60Hz

Maximum current: 8A

Maximum rope extraction

4 m

Min/Max Leg Extension Load

10-630Nm

Min/Max Pulley Load

1-110 kg

Min/Max Squat Load

1-220 kg (up to 440 kg with dedicated accessory)

Min/Max Isokinetic Speeds

0.1-4m/s

Isokinetic Speed Min/Max Leg Extension

12-778°/s

COMPLIANCE



FUNCTIONALITY

· Wide range of ISOKINETIC, ISOMETRIC AND DYNAMIC tests for continuous monitoring of the rehabilitation process

· Wi-Fi connection and remote control

· 2 15.6" touchscreens

· 2 USB connection sockets

· Integrated memory and automatic storage

· High-quality Italian leather

HARDWARE

· Electric motor. Provides variable and controlled effort torque

· Pulley workstation

· Leg Extension workstation

· Squat workstation

ACCESSORIES

· Leg Curl kit

· Leg Press kit

· Lat Machine

· Wooden Squat platform

TRAINING METHODS



WARM UP



ISOTONIC WITH AND WITHOUT INERTIA



ELASTIC



SMART VISCOUS



ISOMETRIC



ISOKINETIC



VARIABLE LOAD



KINEO PULLEY & SQUAT

TECHNICAL SPECIFICATIONS

· Code K0156M | · Code K0156

Dimensions (WxDxH)

185x90x185 cm

Working area dimensions (WxDxH)

185x90x185 cm

Weight

274 kg

Electric motor

Rated power 1.6kW

Power supply

Voltage: 230V - Frequency: 50/60Hz

Maximum current: 8A

Maximum rope extraction

4 m

Min/Max Pulley Load

1-110 kg

Min/Max Squat Load

1-220 kg (up to 440 kg with dedicated accessory)

Min/Max Isokinetic Speeds

0.1-4m/s

COMPLIANCE



CE
0476

CE

FUNCTIONALITY

- Wide range of ISOKINETIC, ISOMETRIC AND DYNAMIC tests for continuous monitoring of the rehabilitation process
- Wi-Fi connection and remote control
- 2 15.6" touchscreens
- 2 USB connection sockets
- Integrated memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. Provides variable and controlled effort torque
- Pulley workstation
- Squat workstation

ACCESSORIES

- Leg Curl kit
- Leg Press kit
- Lat Machine
- Wooden Squat platform

TRAINING METHODS



WARM UP



ELASTIC



ISOMETRIC



ISOTONIC WITH AND WITHOUT INERTIA



SMART VISCOUS



ISOKINETIC



VARIABLE LOAD



KINEO LEG PRO

TECHNICAL SPECIFICATIONS

· Code K0417M | · Code K0417

Dimensions (WxDxH)

260x90x185 cm

Working area dimensions (WxDxH)

310x90x185 cm

Weight

365 kg

Electric motor

Rated power 1.6kW

Power supply

Voltage: 230V - Frequency: 50/60Hz

Maximum current: 8A

Min/Max Leg Extension Load

10-630nm

Min/Max Leg Curl Load

1-110 kg

Min/Max Isokinetic Speeds

0.1-4m/s

Isokinetic Speed Min/Max Leg Extension

12-778°/s

COMPLIANCE



CE 0476

CE

FUNCTIONALITY

· Wide range of ISOKINETIC, ISOMETRIC AND DYNAMIC tests for continuous monitoring of the rehabilitation process

· Wi-Fi connection and remote control

· 2 15.6" touchscreens

· 2 USB connection sockets

· Integrated memory and automatic storage

· High-quality Italian leather

HARDWARE

· Electric motor. Provides variable and controlled effort torque

· Leg Extension workstation

· Leg Curl workstation

ACCESSORIES

· Leg Press kit

· Squat kit

· Wooden Squat platform

TRAINING METHODS



WARM UP



ISOTONIC WITH AND WITHOUT INERTIA



ELASTIC



SMART VISCOUS



ISOMETRIC



ISOKINETIC



VARIABLE LOAD



KINEO LEG PRESS

TECHNICAL SPECIFICATIONS

· Code K1350M | · Code K2004

Dimensions (WxDxH)

320x75x170 cm

Working area dimensions (WxDxH)

320x75x170 cm

Weight

500 kg

Electric motor

Rated power 1.6kW

Power supply

Voltage: 230V - Frequency: 50/60Hz

Maximum current: 8A

Min/max load

10-440 kg

Min/Max Isokinetic Speeds

0.1-4m/s

FUNCTIONALITY

- Wide range of ISOKINETIC, ISOMETRIC AND DYNAMIC tests for continuous monitoring of the rehabilitation process
- Wi-Fi connection and remote control
- 1 15.6" touchscreen
- Adjustable seat depth, 7-position adjustable backrest
- Integrated memory and automatic storage
- High-quality Italian leather

HARDWARE

- Electric motor. Provides variable and controlled effort torque
- Leg Press workstation

ACCESSORIES

- Proprioceptive cushion with easy magnetic hooking
- Chest Press kit

COMPLIANCE



CE
0476

CE

TRAINING METHODS



WARM UP



ELASTIC



ISOMETRIC



VARIABLE LOAD



ISOTONIC WITH AND WITHOUT INERTIA



SMART VISCOUS



ISOKINETIC



KINEO CROSSOVER CABLE



KINEO CROSSOVER CABLE IS THE UNION OF TWO KINEO DEVICES EQUIPPED WITH THE LAT MACHINE ACCESSORY, WHICH ARE CONNECTED BY A CENTRAL BRIDGE. THIS GIVES NEW POSSIBILITIES AND NEW GEOMETRIES TO PERFORM EXERCISES WITH A FREE CABLE, ALLOWING TO USE A CABLE FOR EACH LIMB. THE TWO KINEO DEVICES WORK INDEPENDENTLY OF EACH OTHER, WITH THE TRADITIONAL SOFTWARE FUNCTIONS INCLUDED.



ACCESSORIES

LEG PRESS KIT

TECHNICAL SPECIFICATIONS

· Code K1161

Dimensions (WxDxH)

220x75x170 cm

Working area dimensions (WxDxH)

220x75x170 cm

Weight

230 kg

Maximum seat weight

135 kg

Min/Max load*

10-440 kg

Min/Max Isokinetic Speeds*

0.1-4m/s

* The technical data refer to the KINEO device when connected.

KINEO LEG PRESS “Accessory” cannot work autonomously.

FUNCTIONALITY

· Adjustable seat depth, 7-position adjustable backrest

· High-quality Italian leather

HARDWARE

· Leg Press and Rowing workstation with Rowing bar

· Wheels with lock system

ACCESSORIES

· Proprioceptive cushion with easy magnetic hooking

· Chest Press kit

COMPLIANCE



ACCESSORIES



PH 01

PH 01. LAT MACHINE

K1640	LAT MACHINE KIT, COMPOSED OF LAT MACHINE ARM AND LAT MACHINE BENCH
-------	--

KO495	LAT MACHINE ARM
-------	-----------------



PH 02

PH 02. CONNECTION BRIDGE FOR 2 LAT MACHINES

KO329	BRIDGE FOR 2 LAT MACHINES
-------	---------------------------



PH 03

PH 03. LAT MACHINE BENCH

KO378	LAT MACHINE BENCH
-------	-------------------



PH 04

PH 04. STEP

KOO34	STEP
-------	------



PH 05

PH 05. LEG CURL KIT

KO133	LEG CURL KIT
-------	--------------

ACCESSORIES



PH 06

PH 06. SQUAT KIT

KO449 | SQUAT KIT



PH 07

PH 07. JUMP PLATFORM

KO393 | JUMP PLATFORM



PH 08

PH 08. CHEST PRESS KIT


K2201 | CHEST PRESS KIT



PH 09

PH 09. KINEO ACCESSORIES CART

KO762 | KINEO ACCESSORIES CART

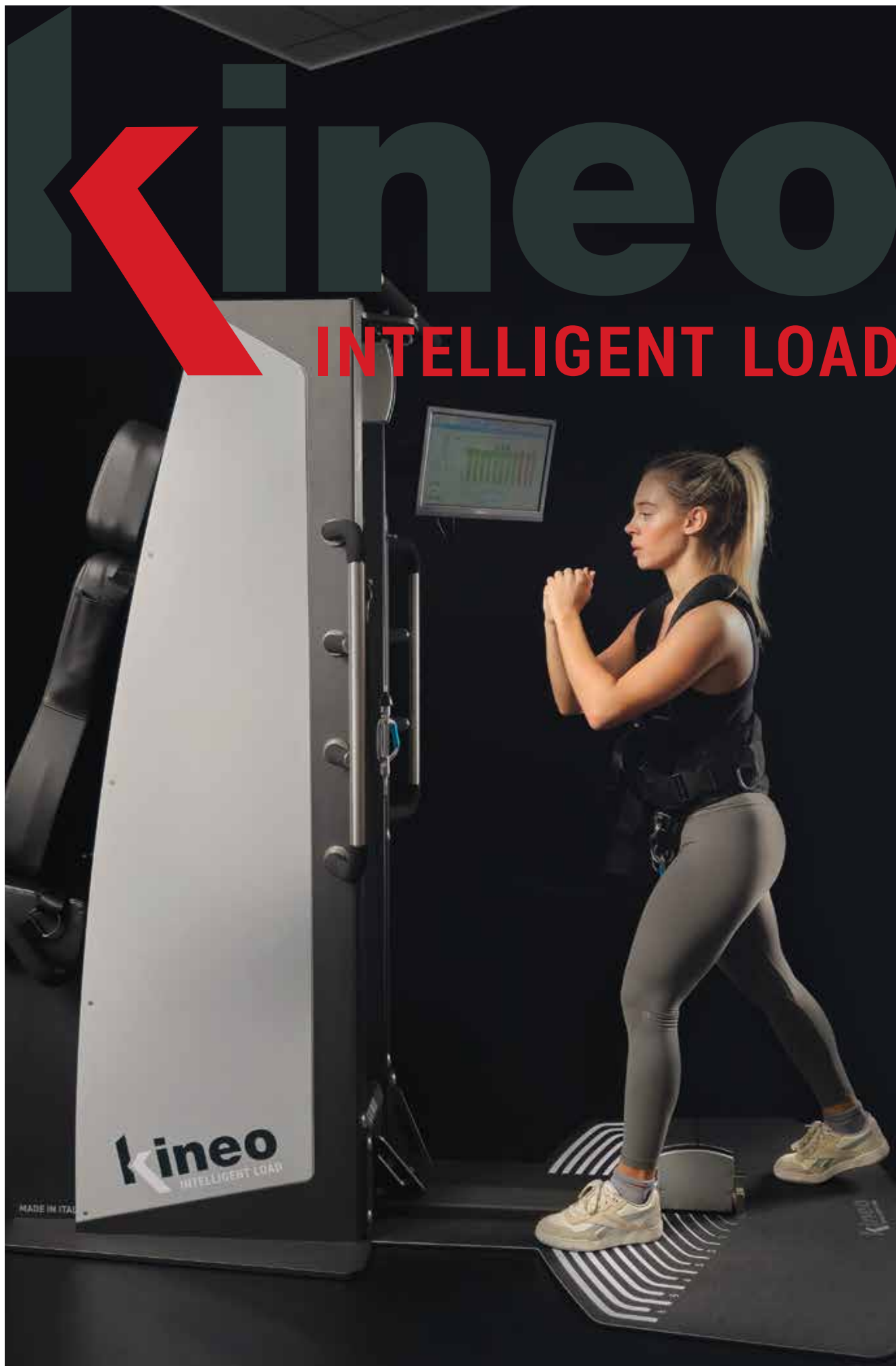


MEDICAL

CONTENTS

- TECAR THERAPY	20
• DIACARE 7000	22
• DIACARE 5000 DIACARE 5000 RE	23
• "FASCIA TOOLS" KIT	24
• "INTIMA" KIT	26
- MAGNETOTHERAPY	28
• MAGNUM 3500 PRO	30
• MAGNUM 3000 PRO	31
• MAGNUM 2 PRO DRIVE	32
• MAGNUM XL PRO	33
• MAGNUM 2500	34
• MAGNETOCEMP ELITE	35
• MAGNUM XL	36
• MAGNUM L	37
• PROGRAMS TABLE	38
- ELECTROTHERAPY	40
• GENESY 3000	42
• GENESY 1500	43
• GENESY 600	44
• GENESY 300 PRO	45
• GENESY S2	46
• DUO TENS	47
• DYSPHAGIA REHABILITATION	48
• AMI PRO AMI LIFE	49
- ULTRASOUND THERAPY	50
• MEDISOUND 3000	52
- LASER THERAPY for physiotherapy	54
• PHYSIOLASER 12.0 PRO	56
• PHYSIOLASER 12.0	57
• PHYSIOLASER 6.0 PRO	58
• PHYSIOLASER 6.0	59
• PHYSIOLASER 2.0	60
• PHYSIOLASER 1000	61
• PHYSIOLASER 500	62
- LASER THERAPY for podiatry	64
• PODCARE 12.0 PRO	66
• PODCARE 6.0 PRO	67
• PODCARE 6.0	68
• PODCARE 2.0	69
- PRESSOTHERAPY	70
• PRESSCARE G300M	72
• PRESSCARE G200M	73

TECAR THERAPY

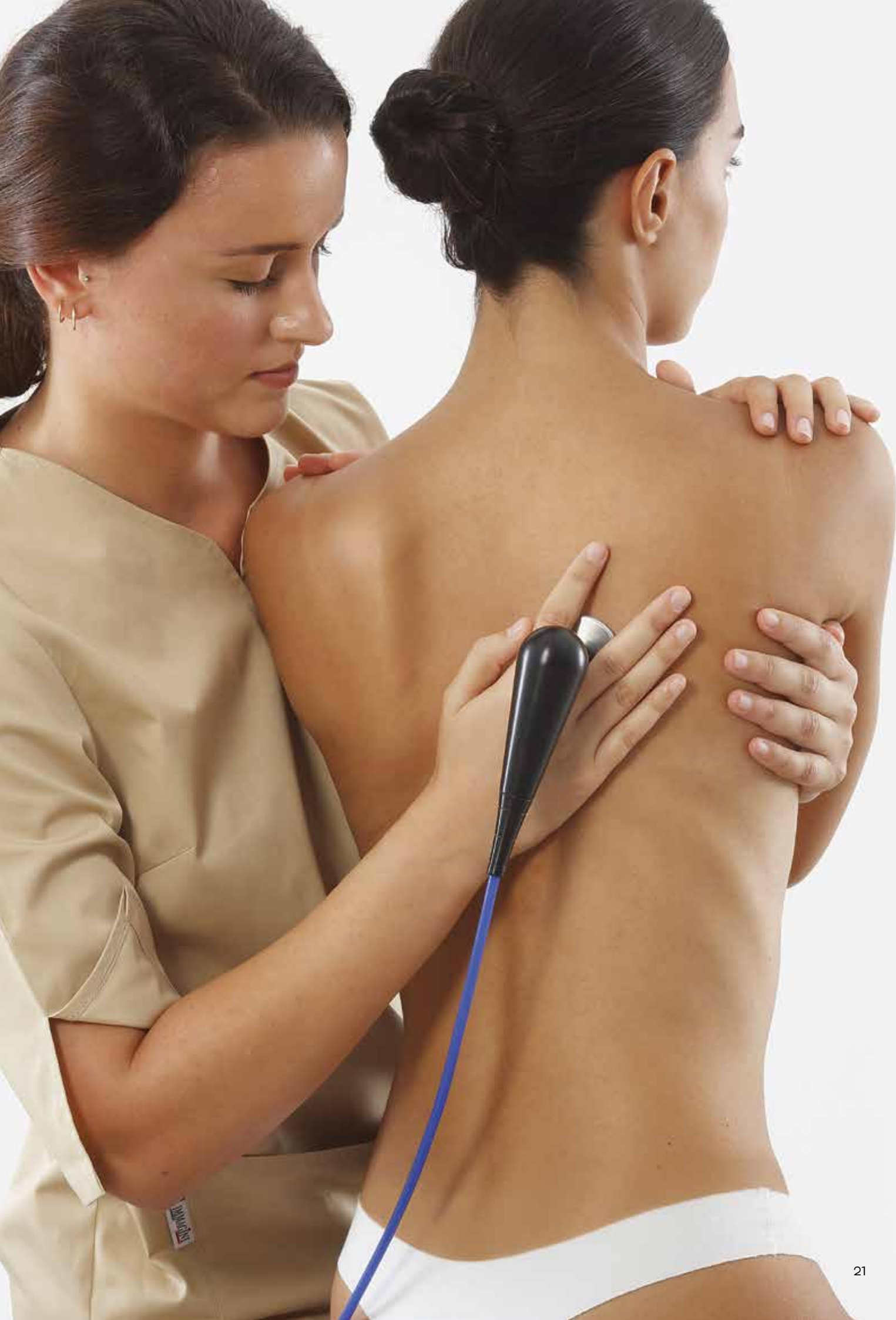

DIACARE 7000

TECAR THERAPY

DIACARE 7000 IS THE MOST POWERFUL, VERSATILE AND COMPLETE DEVICE IN THE GLOBUS TECAR THERAPY RANGE. IT USES BOTH RESISTIVE AND CAPACITIVE MODES FOR THE TREATMENT OF ACUTE AND CHRONIC CONDITIONS OF THE MUSCULOSKELETAL SYSTEM. THE DEVICE CAN BE USED BY ACTING ON THE PATIENT BOTH IN MANUAL MODE AND COMBINING PASSIVE MOVEMENTS OR, THANKS TO THE HANDS-FREE KIT, IT CAN ALSO BE USED DURING ACTIVE MOVEMENTS OF THE SUBJECT. DIACARE 7000 IS READY TO USE WITH APPROPRIATE ERGONOMIC TOOLS DESIGNED FOR THE MOBILISATION AND THE ADVANCED TREATMENT OF ALL SOFT TISSUE DISEASES (FASCIA TOOLS AND DYNAMIC HANDLE).

TECHNICAL SPECIFICATIONS

Display

Size 7" colour touchscreen display

Nominal Power (NP)

350W $\pm 10\%$

Emission frequencies

300kHz, 400kHz, 448kHz, 470kHz, 700kHz, 1000kHz, 1200kHz, $\pm 10\%$

Treatment methods

Capacitive and resistive

Capacitive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Resistive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

515x340x170 mm

Weight

11,5 kg



CERTIFIED
PROTOCOLS

STANDARD EQUIPMENT

· Code **G5210**

· 1 DIACARE 7000 Tecar therapy unit

· 1 Return electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes
 \emptyset (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes
 \emptyset (30, 50, 70 mm)


· 1 Bipolar handpiece \emptyset (50 mm)

· 1 Hands-free treatment kit

· 1 1000 ml jar of conductive cream

· 1 Power supply cable

· 1 Operating manual

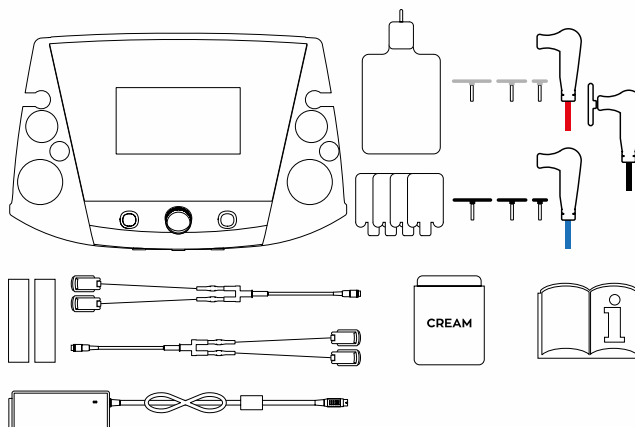

STRONG POINTS



DIACARE 5000 | DIACARE 5000 RE

TECAR THERAPY | 14 PROGRAMS

DIACARE 5000 HAS ALL THE QUALITY OF GLOBUS' TECAR THERAPY IN A COMPACT, LIGHTWEIGHT MACHINE TO ENSURE MAXIMUM PORTABILITY, WITHOUT COMPROMISING ON EFFECTIVENESS AND PROFESSIONALISM.

TECHNICAL SPECIFICATIONS

Display

Size 4.3" colour touchscreen display

Nominal Power (NP)

180W $\pm 10\%$

Emission frequencies

300kHz, 400kHz, 448kHz, 470kHz, 700kHz, 1000kHz, $\pm 10\%$

Treatment methods

Capacitive and resistive

Capacitive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Resistive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

170x220x60 mm

Weight

1.04 kg



CERTIFIED
PROTOCOLS

STANDARD EQUIPMENT

· Code **G3488**

· 1 Case

· 1 DIACARE 5000 Tecar therapy unit

· 1 Return electrode (plate)


· 1 Monopolar resistive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)

· 1 250 ml jar of conductive cream

· 1 Power supply unit

· 1 Operating manual

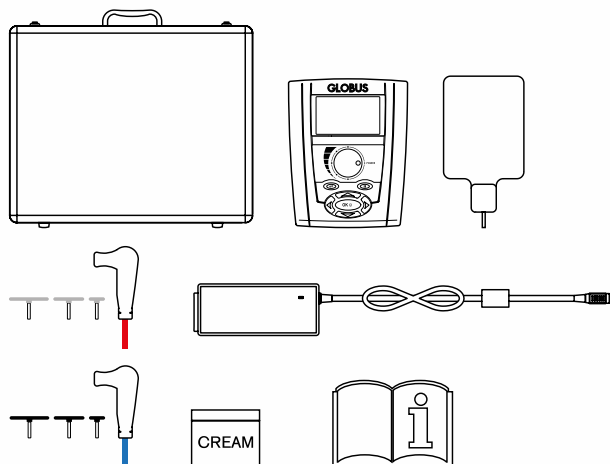

PURCHASE OPTIONS


· Code **G3996** DIACARE 5000 RE *

* RECHARGEABLE TECAR THERAPY :

- 5 HOURS OF TREATMENT INCLUDED IN THE COST OF THE DEVICE
- CHEAP CHARGING
- GUARANTEED SUPPORT FOR 5 YEARS AND FORKLIFT
- CHARGING VIA INTERNET



STRONG POINTS

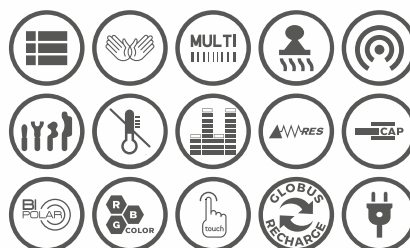

“FASCIA TOOLS” KIT

ADVANCED SOFT TISSUE MOBILISATION INSTRUMENTS FOR USE IN OUR LINE OF “DIACARE” DEVICES

Tecar therapy is a technology that is becoming increasingly popular in the world of physiotherapy as an indispensable instrument for any professional who wants to achieve excellent results

GLOBUS has successfully combined the TECAR technology with instrument assisted soft tissue mobilisation (IASTM), bringing together the advantages of these two extraordinary therapeutic possibilities.

The GLOBUS “FASCIA TOOLS” are special resistive handpieces that, due to their anatomical shape, are particularly effective for the mobilisation and the advanced treatment of all soft tissue diseases.

The combination of manual treatment with myofascial handpieces, together with the energy of GLOBUS DIACARE, offer a new and extraordinary treatment option and, with the advantages of the combined application of tecar therapy and the myofascial electrodes, it is possible to:

- **ELIMINATE FASCIAL RESTRICTIONS AND SCAR ADHESIONS**
- **INCREASE THE EXTENSIBILITY OF THE CONNECTIVE TISSUE**
- **IMPROVE TISSUE VASCULARISATION**
- **IMPROVE MUSCLE RECRUITMENT AND MOTOR CONTROL**



Large DiaBlade: tool to treat Lower Limb and Trunk disorders.




Small DiaBlade: tool to treat Upper Limb and Cervical Spine disorders.



Dia Trigger Point: Trigger Points treatment tool.



Dia Spine: Back treatment and spinal mobilisation tool.

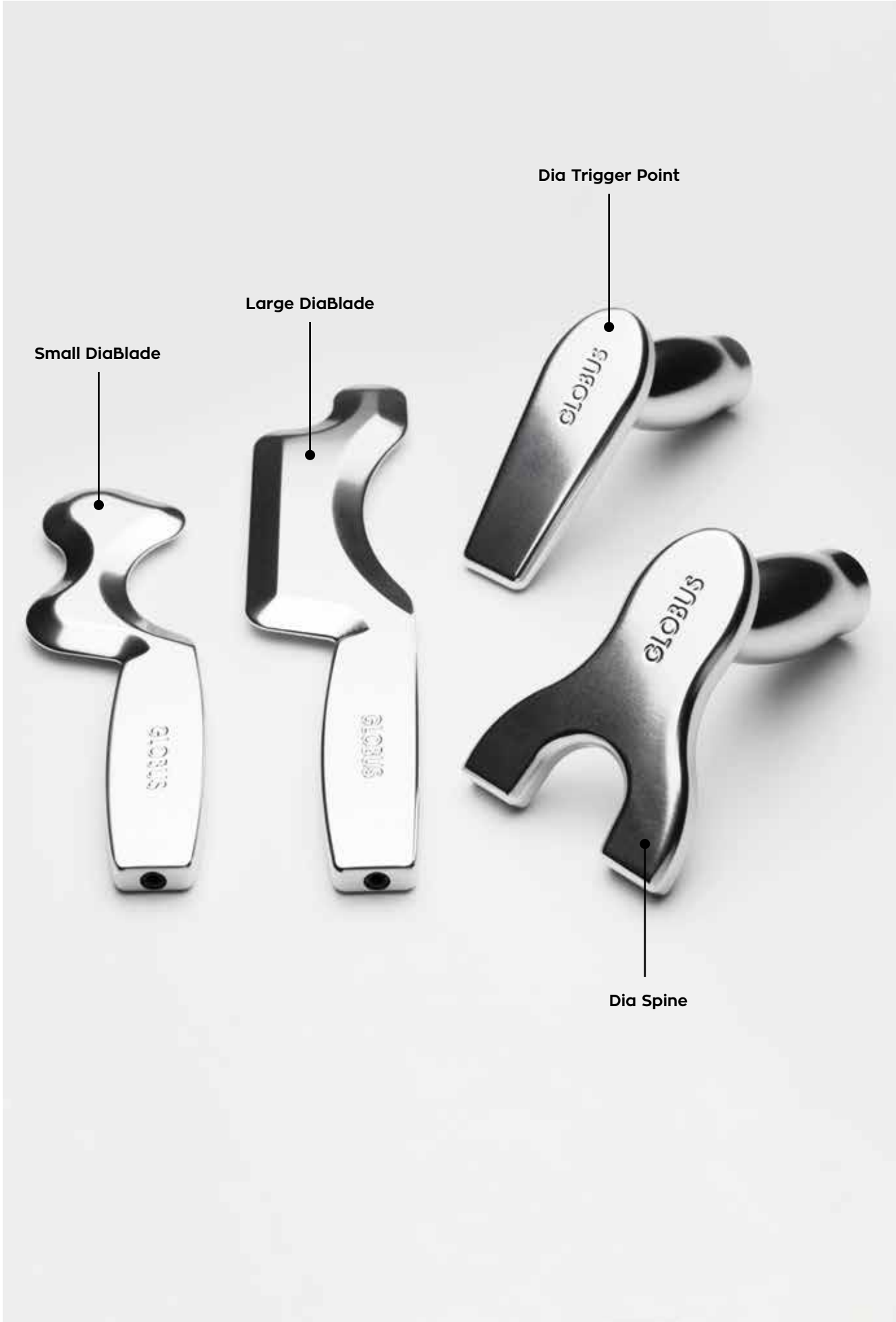

“INTIMA” KIT

ADVANCED TOOLS TO IMPROVE THE MUSCLE TONE OF THE PELVIC FLOOR FOR USE IN OUR LINE OF “DIACARE” DEVICES

The precision of our TECAR devices offers effective action in total safety and comfort, even for solving problems directly related to the genital area.

This is why GLOBUS has created the “INTIMA” KIT, with specific accessories for internal and external treatments to improve the health of the pelvic floor, both in men and women, and to restore vaginal tone and elasticity, and also to provide excellent results when treating childbirth scars.


The capacitive mode treatment of the “INTIMA” KIT enables deep action on the mucous membranes in the utmost comfort.

Furthermore, when combined with traditional handpieces and the “HANDS-FREE” KIT, it provides a complete system to increase and restore the tone of the pelvic floor.

The most frequent treatments with the “INTIMA” KIT, ideal for restoring intimate comfort, are:

- VAGINAL LAXITY
- CHILDBIRTH SCARS
- ABDOMINAL DIASTASIS
- ENDOPELVIC MYOFASCIAL SYNDROMES
- PERINEAL REHABILITATION FOR MEN AND WOMEN
- PELVIC FLOOR TONE IMPROVEMENT

Code **G6421** “INTIMA” KIT for DIACARE 7000 | Code **G6420** “INTIMA” KIT for DIACARE 5000

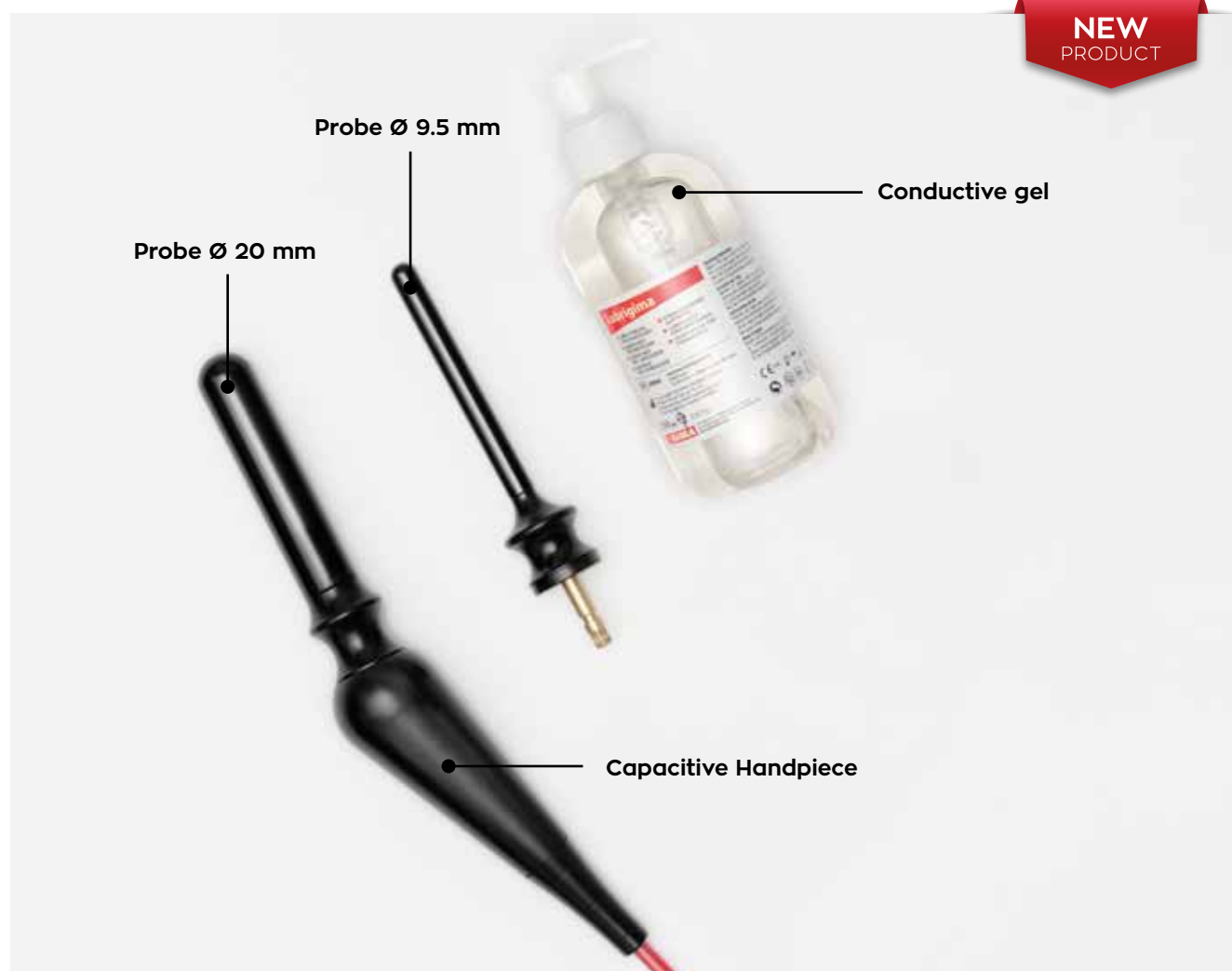



MAGNETOTHERAPY



MAGNUM 3500 PRO

PROFESSIONAL MAGNETOTHERAPY | 70 PROGRAMS

THE ONLY GLOBUS MAGNETOTHERAPY DEVICE WITH FOUR CHANNELS THAT ALLOWS TO TREAT DIFFERENT BODY PARTS OR TWO SUBJECTS AT THE SAME TIME. MAGNUM 3500 PRO IS EQUIPPED WITH A WIDE RANGE OF ACCESSORIES FOR ALL TREATMENT NEEDS.

26 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 6 INFLAMMATIONS | 6 ORTHOPAEDICS | 3 PAIN | 3 TRAUMAS/OEDEMAS/WOUNDS | 3 NIGHT PROGRAMS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4

Intensity control with flexible diffusers

- 200 average Gauss per channel
- 800 total average Gauss (for 4 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 10 custom programs

2+2 function

Possibility of running 2 different programs at the same time, using 2 channels each

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions | Weight

170x220x60 mm | 820 g

STANDARD EQUIPMENT

· Code G5434

· 1 Case

· 1 MAGNUM 3500 PRO magnetotherapy unit

· 2 Flexible solenoids 300x100 mm

· 4 "Pocket Pro" solenoids 155x95 mm


· 1 "Soft 4" mat 250x250 mm

· 4 Velcro elastic bands 800x80 mm

· 1 Magnum Test

· 1 Power supply unit

· 1 Operating manual

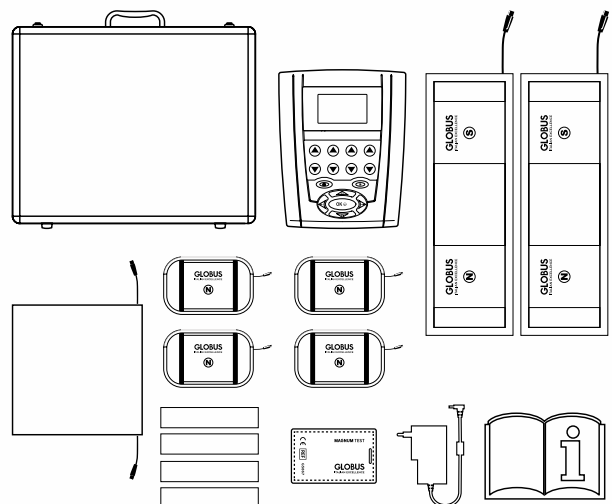

CE
0476

STRONG POINTS



MAGNUM 3000 PRO

PROFESSIONAL MAGNETOTHERAPY | 70 PROGRAMS

POWERFUL PROFESSIONAL MAGNETOTHERAPY WITH MANY SPECIAL PROGRAMS AND FEATURES THAT MAKE IT EFFECTIVE TO TREAT NUMEROUS OSTEOARTICULAR PATHOLOGIES AND FOR RECOVERY AFTER ORTHOPAEDIC OPERATIONS. OPTION OF LOCAL AND TOTAL-BODY TREATMENTS.

26 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 6 INFLAMMATIONS | 6 ORTHOPAEDICS | 3 PAIN | 3 TRAUMAS/OEDEMAS/WOUNDS | 3 NIGHT PROGRAMS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2

Intensity control with flexible diffusers

- 200 average Gauss per channel
- 400 total average Gauss (for 2 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 10 custom programs

1+1 function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply


- Mains
- Rechargeable battery (not included)

Dimensions | Weight

170x220x60 mm | 820 g

STANDARD EQUIPMENT

- Code G5335
- 1 Bag
- 1 MAGNUM 3000 PRO magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 2 Velcro elastic bands 800x80 mm
- 1 Magnum Test
- 1 Power supply unit
- 1 Operating manual

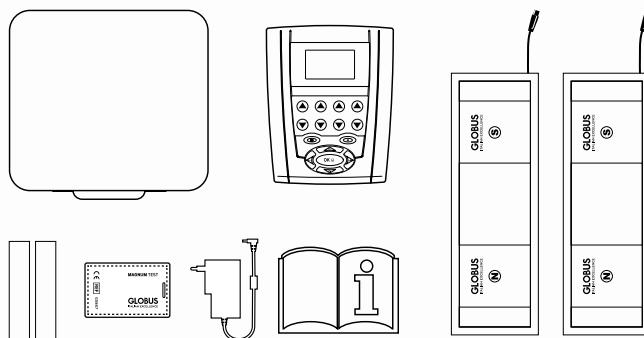

PURCHASE OPTION

- Code G6034 MAGNUM 3000 PRO, 2 "Pocket Pro" solenoids with 1 velcro elastic band.



CE
0476

STRONG POINTS



MAGNUM 2 PRO DRIVE

MAGNETOTHERAPY | 41 PROGRAMS

MAGNUM 2 PRO DRIVE IS DESIGNED TO OFFER A SOLUTION TO THOSE WHO SUFFER FROM PAIN IN THE BACK OR OTHER PARTS OF THE BODY DUE TO TOO MUCH TIME SPENT DRIVING. SEDENTARY LIFESTYLE AND MANY HOURS BEHIND THE WHEEL, SUBJECTED TO THE VIBRATIONS OF THE VEHICLE, ARE SURELY AMONG THE MOST FREQUENT CAUSES OF MUSCULOSKELETAL PAIN. MAGNUM 2 PRO DRIVE IS AN EXCELLENT INSTRUMENTAL THERAPY THAT CAN HELP TO REDUCE THE INFLAMMATION AND PAIN OF THIS CATEGORY OF PEOPLE.

15 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 6 ORTHOPAEDICS | 6 INFLAMMATIONS | 3 PAIN | 3 TRAUMAS/OEDEMAS/WOUNDS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2

Intensity control with flexible diffusers

- 200 average Gauss per channel
- 400 total average Gauss (for 2 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 5 custom programs

1+1 function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply

- Mains
- Rechargeable battery


Dimensions and Weight

100x160x35 mm - 450 g

STANDARD EQUIPMENT

Code **G5429**

- 1 Bag
- 1 MAGNUM 2 PRO DRIVE magnetotherapy unit
- 1 "Soft 4" mat 250x250 mm
- 2 Flexible solenoids 300x100 mm
- 2 Velcro elastic bands 800x80 mm
- 1 Magnum Test
- 1 USB cable with 5-12V converter
- 1 Power supply unit
- 1 Operating manual

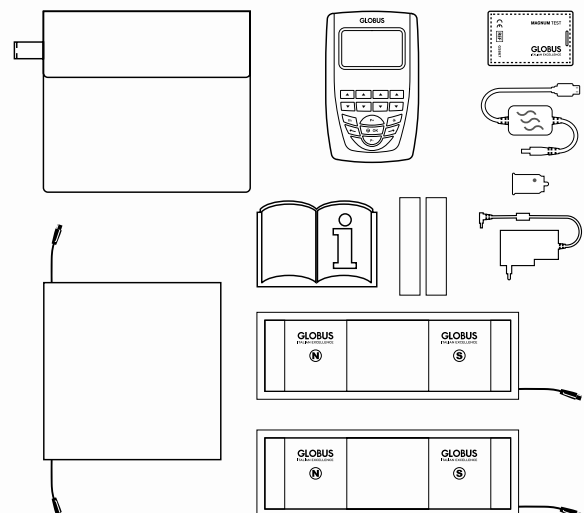

CE
0476

STRONG POINTS



MAGNUM XL PRO

MAGNETOTHERAPY | 46 PROGRAMS

A MAGNETOTHERAPY DEVICE WITH A FULL RANGE OF PRESET PROGRAMS TO TREAT THE MOST COMMON OSTEOARTICULAR DISORDERS. THE DEVICE IS AMONG THE MOST POWERFUL ON THE MARKET AND IS FITTED WITH TWO SOLENOIDS.

32 HEALTH AND REHAB | 14 AESTHETICS AND FITNESS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2

Intensity control with flexible diffusers

- 200 average Gauss per channel
- 400 total average Gauss (for 2 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 15 custom programs

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

100x160x35 mm


Weight

280 g

STANDARD EQUIPMENT

Code G3956

- 1 Bag
- 1 MAGNUM XL PRO magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 2 Velcro elastic bands 800x80 mm
- 1 Magnum Test
- 1 Power supply unit
- 1 Operating manual

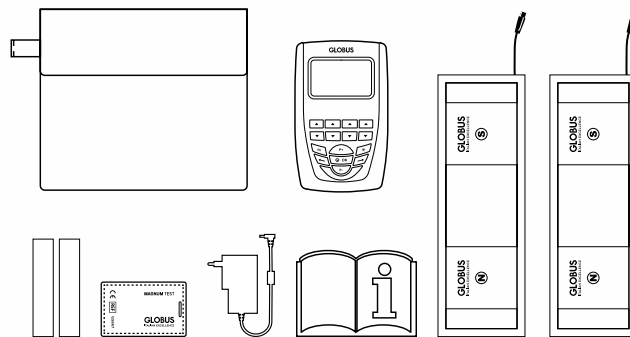

PURCHASE OPTION

- Code G6033 MAGNUM XL PRO, 2 "Pocket Pro" solenoids with 1 velcro elastic band.



CE
0476

STRONG POINTS



MAGNUM 2500

MAGNETOTHERAPY | 52 PROGRAMS

PORTABLE MAGNETOTHERAPY, IDEAL FOR HOME TREATMENTS. MAGNUM 2500 IS INDICATED FOR A BETTER RECOVERY FROM TRAUMAS AND FRACTURES, BUT ALSO TO PROVIDE RELIEF AND HELP ARTICULAR REGENERATION IN PATHOLOGIES LIKE ARTHROSIS AND ARTHRITIS.

16 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 14 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 5 INFLAMMATIONS 4 ORTHOPAEDICS | 3 PAIN | 2 TRAUMAS/OEDEMAS/WOUNDS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2

Intensity control with flexible diffusers

- 160 average Gauss per channel
- 320 total average Gauss (for 2 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 5 custom programs

1+1 function

Possibility of running 2 different programs at the same time, using 1 channel each

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

100x160x35 mm


Weight

280 g

STANDARD EQUIPMENT

Code G5438

- 1 Bag
- 1 MAGNUM 2500 magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 2 Velcro elastic bands 800x80 mm
- 1 Magnum Test
- 1 Power supply unit
- 1 Operating manual

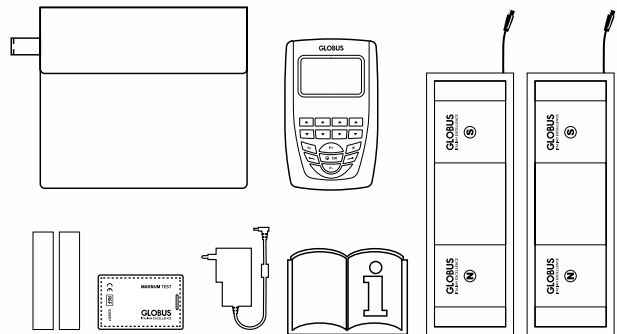

PURCHASE OPTION

- Code G6191 MAGNUM 2500, 2 "Pocket Pro" solenoids with 1 velcro elastic band.



CE
0476

STRONG POINTS



MAGNETOCEMP ELITE

MAGNETOTHERAPY | 52 PROGRAMS

PROFESSIONAL MAGNETOTHERAPY DEVICE WITH A FULL RANGE OF PRESET PROGRAMS FOR TREATING ARTHROSIS, JOINT INFLAMMATION, BONE OEDEMA AND OTHER CONDITIONS AFFECTING THE MUSCULOSKELETAL SYSTEM.

16 PRESET PROGRAMS WITH DIFFERENT FREQUENCY | 14 BONE PATHOLOGIES | 8 ARTHROSIS/ARTHRITIS | 5 INFLAMMATIONS 4 ORTHOPAEDICS | 3 PAIN | 2 TRAUMAS/OEDEMAS/WOUNDS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2

Intensity control with flexible diffusers

- 160 average Gauss per channel
- 320 total average Gauss (for 2 channels)

Frequency

5-200Hz

Free programming

Possibility of adding up to 5 custom programs

1+1 function

Possibility of running 2 different programs at the same time, using 1 channel each

Duty cycle

5-50%

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

100x160x35 mm


Weight

280 g

STANDARD EQUIPMENT

Code **G6408**

- 1 Case
- 1 MAGNETOCEMP ELITE magnetotherapy unit
- 2 Flexible solenoids 300x100 mm
- 2 Velcro elastic bands 800x80 mm
- 1 Magnum Test
- 1 Power supply unit
- 1 Operating manual

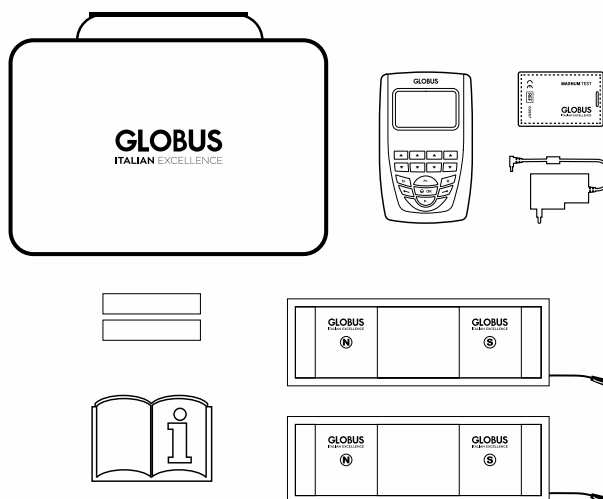

NOT AVAILABLE
ON LINE



CE
0476

STRONG POINTS



MAGNUM XL

MAGNETOTHERAPY | 26 PROGRAMS

DEVICE FOR HOME MAGNETOTHERAPY SUITABLE FOR TREATING THE MAIN OSTEOARTICULAR DISORDERS. EVEN EASIER TO USE WITH THE NEW EASY PROGRAM FUNCTION.

16 HEALTH AND REHAB | 10 AESTHETICS AND FITNESS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

2 independent

Intensity control with flexible diffusers

- 140 average Gauss per channel
- 280 total average Gauss (for 2 channels)

Frequency

5-200Hz

Duty cycle

5-50%

Free programming

Possibility of adding up to 1 custom program

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


100x160x35 mm

Weight

260 g

STANDARD EQUIPMENT

- Code G3216
- 1 Bag
- 1 MAGNUM XL magnetotherapy unit
- 1 Flexible solenoid 300x100 mm
- 1 Velcro elastic band 800x80 mm
- 1 Power supply unit
- 1 Operating manual

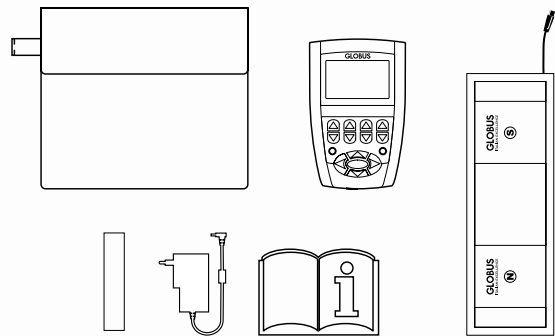

CE
0476

STRONG POINTS



MAGNUM L

MAGNETOTHERAPY | 8 PROGRAMS

VERY SIMPLE MAGNETOTHERAPY DEVICE WITH 8 PRESET PROGRAMS FOR TREATING MAINLY FRACTURES AND ARTHROSIS.

8 HEALTH AND REHAB

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

1

Intensity control with flexible diffuser

120 average Gauss

Frequency

5-200Hz

Duty cycle

5-50%

Power supply

Mains

Dimensions

100x160x35 mm

Weight

260 g

STANDARD EQUIPMENT

· Code **G3947**

· 1 Bag


· 1 MAGNUM L magnetotherapy unit

· 1 Flexible solenoid 300x100 mm

· 1 Velcro elastic band 800x80 mm

· 1 Power supply unit

· 1 Operating manual

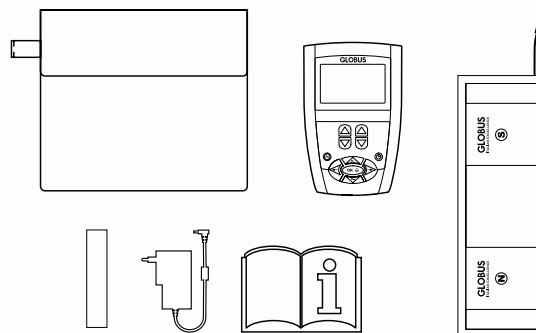

CE
0476

STRONG POINTS



PROGRAMS TABLE

PROGRAM	MODEL						
	3500 PRO 3000 PRO	2 PRO DRIVE	XL PRO	2500	MAGNETOCEMP ELITE	XL	L
Algodystrophy/Sudeck's Disease	✓	✓	/	✓	✓	/	/
Anti-inflammatory	✓	✓	/	✓	✓	/	/
Arthritis	✓	✓	✓	✓	✓	/	✓
Arthrosis - Osteoarthritis	✓	✓	✓	✓	✓	✓	✓
Ankle sprain	✓	✓	✓	✓	✓	/	/
Joint pain	✓	✓	✓	✓	✓	/	/
Bone oedemas	✓	✓	/	✓	✓	/	/
Oedemas-bruises	✓	✓	✓	✓	✓	/	/
Epicondylitis / Epitrochleitis	✓	✓	✓	✓	✓	/	✓
Fracture with plaster	✓	✓	*	✓	✓	*	*
Fractures	✓	✓	✓	✓	✓	✓	✓
Spinal fusion	✓	✓	/	/	/	/	/
Cartilage injury	✓	✓	✓	✓	✓	/	/
Backache	✓	✓	✓	✓	✓	✓	✓
Osteonecrosis	✓	✓	/	✓	✓	/	/
Osteoporosis	✓	✓	✓	✓	✓	✓	✓
Osteotomy	✓	✓	/	/	/	/	/
Post-ACL surgery	✓	✓	✓	✓	✓	/	/
Hip-knee prostheses	✓	✓	✓	✓	✓	/	/
Pseudoarthrosis	✓	✓	/	/	/	/	/
Delayed fracture consolidation	✓	✓	✓	✓	✓	✓	/
Shoulder tendonitis / rotator cuff	✓	✓	✓	✓	✓	✓	/
Tendinitis	✓	✓	✓	/	/	/	/
Carpal tunnel	✓	✓	✓	✓	✓	✓	/
Ulcers	✓	✓	/	/	/	/	/
Multifrequency program	✓	/	/	1	1	/	/
Preset frequency programs	25	/	/	15	15	/	/
Night programs	3	/	/	/	/	/	/
Cosmetic programs	/	/	14	/	/	10	/

* Use the "Fractures" program.

Some programs are arranged by body part. See the technical data sheets of the products for more information.

THE HEART OF GLOBUS PURE COPPER SOLENOIDS



More powerful more
effective higher
performance

An electromagnetic field is formed between the two pure copper components. The purer the copper, the more effectively and deeply the waves penetrate.

**This is Globus' secret:
the selection of the best
components and 100% Italian
technology.**



ELECTROTHERAPY



GENESY 3000

4-CHANNEL ELECTROTHERAPY | 423 PROGRAMS

TOP OF THE RANGE PROFESSIONAL FOUR-CHANNEL PRODUCT USED IN THE REHABILITATION SECTOR FOR THE TREATMENT OF THE MOST COMMON PATHOLOGIES OF THE MUSCULOSKELETAL SYSTEM.

236 ELECTROTHERAPY MENU (EMS/RUSSIAN) | 53 SPORT | 33 PAIN | 29 FITNESS | 23 MICROCURRENTS
16 BEAUTY | 11 REHAB | 9 DENERVATED | 9 INTERFERENTIALS | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

- Mains
- Rechargeable batteries

Dimensions

170x220x60 mm

Weight

1.12 kg

* Depending on current type

STANDARD EQUIPMENT

· Code **G1034**

· 1 Bag

· 1 GENESY 3000 stimulator

· 4 Electrode connection cables

· 2 Microcurrent and ionophoresis cables


· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 2 Ionophoresis electrodes + 2 bands

· 1 Power supply unit

· 1 Operating manual

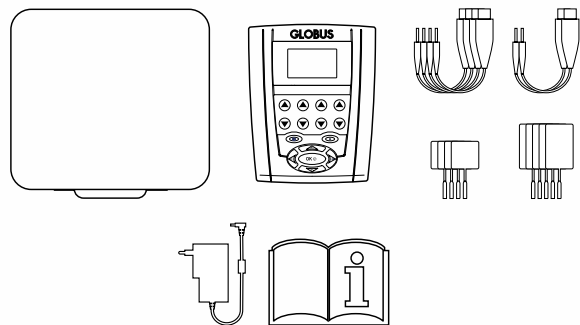


INCONTINENCE PROGRAM

DEVICE SUITABLE FOR THE REHABILITATION OF THE PELVIC FLOOR, BOTH MALE AND FEMALE, AND IDEAL FOR PROBLEMS RELATED TO URINARY OR FAECAL INCONTINENCE, PROLAPSES OF THE BLADDER, UTERUS, RECTUM OR FOR PELVIC PAIN. IT IS A TYPE OF INSTRUMENTAL PHYSICAL THERAPY PERFORMED WITH SPECIAL VAGINAL OR ANAL PROBES WHICH CONDUCT ELECTRICAL IMPULSES THAT STIMULATE THE MUSCLES TO CONTRACT, NORMALISING THE MUSCLE TONE OF THE PERINEUM.



CE
0476

STRONG POINTS

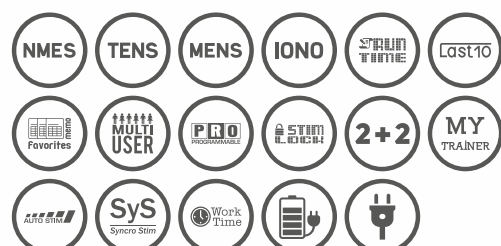

GENESY 1500

4-CHANNEL ELECTROTHERAPY | 414 PROGRAMS

PROFESSIONAL FOUR-CHANNEL PRODUCT USED IN THE MEDICAL AND REHABILITATION SECTOR. IT HAS ALL THE CURRENTS USED FOR REHABILITATION, PHYSIOTHERAPY AND MOTOR REHABILITATION THERAPY.

92 ELECTROTHERAPY MENU (EMS/RUSSIAN) | 84 ACTION NOW | 54 SERIAL SEQUENTIAL STIMULATION | 53 SPORT
33 PAIN | 29 FITNESS | 23 MICROCURRENTS | 16 BEAUTY | 11 REHAB | 9 DENERVATED | 6 INTERFERENTIALS
3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

- Mains
- Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code **G3554**

· 1 Bag

· 1 GENESY 1500 stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Power supply unit

· 1 Operating manual

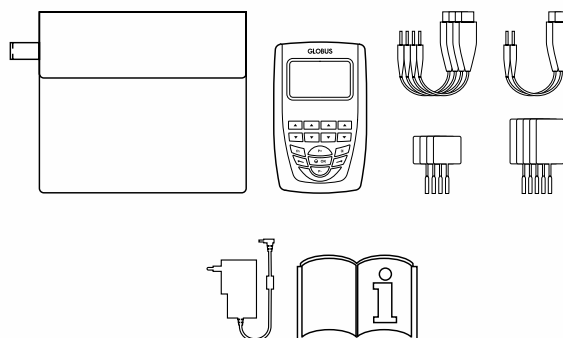


INCONTINENCE PROGRAM

DEVICE SUITABLE FOR THE REHABILITATION OF THE PELVIC FLOOR, BOTH MALE AND FEMALE, AND IDEAL FOR PROBLEMS RELATED TO URINARY OR FAECAL INCONTINENCE, PROLAPSES OF THE BLADDER, UTERUS, RECTUM OR FOR PELVIC PAIN. IT IS A TYPE OF INSTRUMENTAL PHYSICAL THERAPY PERFORMED WITH SPECIAL VAGINAL OR ANAL PROBES WHICH CONDUCT ELECTRICAL IMPULSES THAT STIMULATE THE MUSCLES TO CONTRACT, NORMALISING THE MUSCLE TONE OF THE PERINEUM.



CE
0476

STRONG POINTS

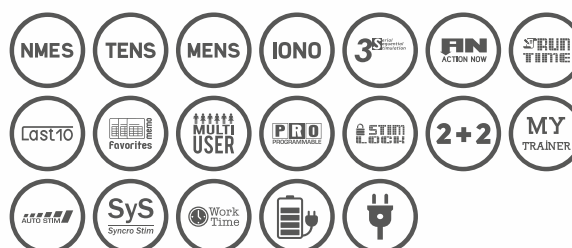

GENESY 600

4-CHANNEL ELECTROTHERAPY | 149 PROGRAMS

PROFESSIONAL FOUR-CHANNEL PRODUCT FOR DEMANDING USERS AND REHABILITATION PROFESSIONALS.

53 SPORT | 29 FITNESS | 20 PAIN | 12 MICROCURRENTS | 12 SERIAL SEQUENTIAL STIMULATION | 11 REHAB
7 ACTION NOW | 3 DENERVATED | 1 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G3553

· 1 Bag

· 1 GENESY 600 stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

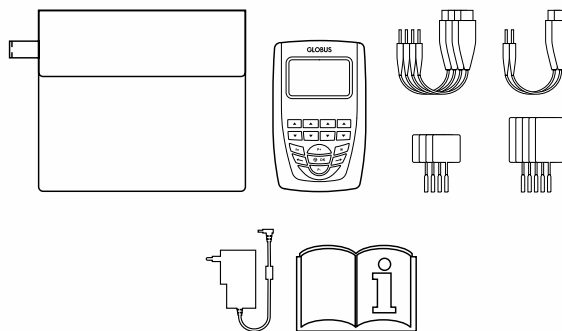


INCONTINENCE PROGRAM

DEVICE SUITABLE FOR THE REHABILITATION OF THE PELVIC FLOOR, BOTH MALE AND FEMALE, AND IDEAL FOR PROBLEMS RELATED TO URINARY OR FAECAL INCONTINENCE, PROLAPSES OF THE BLADDER, UTERUS, RECTUM OR FOR PELVIC PAIN. IT IS A TYPE OF INSTRUMENTAL PHYSICAL THERAPY PERFORMED WITH SPECIAL VAGINAL OR ANAL PROBES WHICH CONDUCT ELECTRICAL IMPULSES THAT STIMULATE THE MUSCLES TO CONTRACT, NORMALISING THE MUSCLE TONE OF THE PERINEUM.



CE
0476

STRONG POINTS

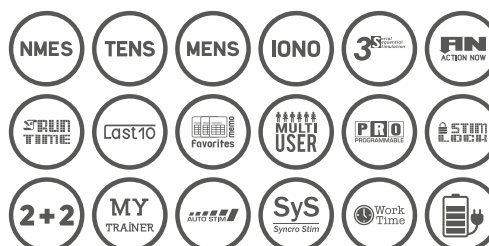

GENESY 300 PRO

4-CHANNEL ELECTROTHERAPY | 91 PROGRAMS

FOUR-CHANNEL PRODUCT WITH A WIDE RANGE OF PROGRAMS FOR PHYSIOTHERAPY.

29 FITNESS, BEAUTY, WELLNESS | 20 PAIN | 22 SPORT | 12 MICROCURRENTS | 6 REHAB
1 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

440 g

* Depending on current type

STANDARD EQUIPMENT

· Code G3222

· 1 Bag

· 1 GENESY 300 PRO stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

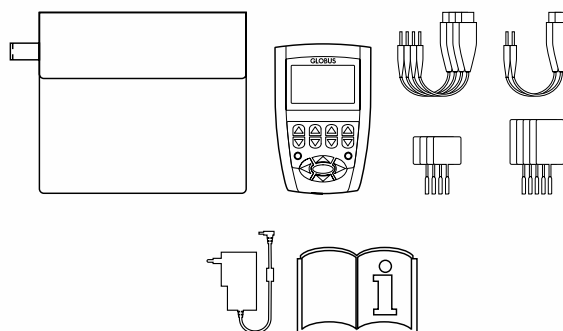


INCONTINENCE PROGRAM

DEVICE SUITABLE FOR THE REHABILITATION OF THE PELVIC FLOOR, BOTH MALE AND FEMALE, AND IDEAL FOR PROBLEMS RELATED TO URINARY OR FAECAL INCONTINENCE, PROLAPSES OF THE BLADDER, UTERUS, RECTUM OR FOR PELVIC PAIN. IT IS A TYPE OF INSTRUMENTAL PHYSICAL THERAPY PERFORMED WITH SPECIAL VAGINAL OR ANAL PROBES WHICH CONDUCT ELECTRICAL IMPULSES THAT STIMULATE THE MUSCLES TO CONTRACT, NORMALISING THE MUSCLE TONE OF THE PERINEUM.



CE
0476

STRONG POINTS

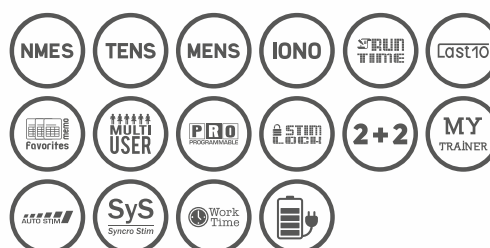

GENESY S2

2-CHANNEL ELECTROTHERAPY | 60 PROGRAMS

THE TWO-CHANNEL PRODUCT FOR PHYSIOTHERAPY.

16 TENS AND PAIN | 12 SPORT | 11 REHAB | 9 PREVENTION | 8 FITNESS, BEAUTY, WELLNESS
3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.5"

Channels

2 independent (4 electrodes)

Frequency

1-150Hz

Pulse width

50-400µs

Intensity

0-100mA per channel

Power supply

Rechargeable batteries

Dimensions

80x130x25 mm

Weight

220 g

STANDARD EQUIPMENT

· Code G3725

· 1 Bag

· 1 GENESY S2 stimulator


· 2 Electrode connection cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

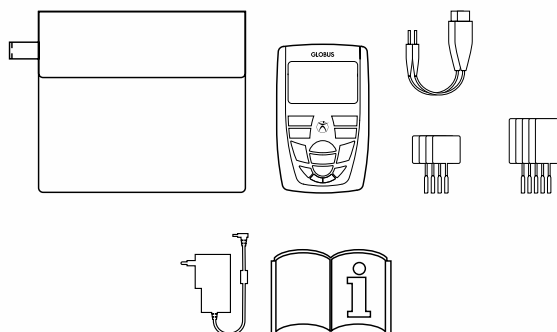

INCONTINENCE PROGRAM

DEVICE SUITABLE FOR THE REHABILITATION OF THE PELVIC FLOOR, BOTH MALE AND FEMALE, AND IDEAL FOR PROBLEMS RELATED TO URINARY OR FAECAL INCONTINENCE, PROLAPSES OF THE BLADDER, UTERUS, RECTUM OR FOR PELVIC PAIN. IT IS A TYPE OF INSTRUMENTAL PHYSICAL THERAPY PERFORMED WITH SPECIAL VAGINAL OR ANAL PROBES WHICH CONDUCT ELECTRICAL IMPULSES THAT STIMULATE THE MUSCLES TO CONTRACT, NORMALISING THE MUSCLE TONE OF THE PERINEUM.



CE
0476

STRONG POINTS



DUO TENS

2-CHANNEL ELECTROTHERAPY | 20 PROGRAMS

THE IDEAL PRODUCT TO TREAT THE MOST COMMON PAINFUL SYMPTOMS.

16 TENS AND PAIN | 3 SPORT | 1 FITNESS, BEAUTY, WELLNESS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.5"

Channels

2 independent (4 electrodes)

Frequency

1-150Hz

Pulse width

50-400µs

Intensity

0-100mA per channel

Power supply

Rechargeable batteries

Dimensions


80x130x25 mm

Weight

220 g

STANDARD EQUIPMENT

- Code **G3729**
- 1 Bag
- 1 DUO TENS stimulator
- 2 Electrode connection cables
- 4 Square self-adhesive electrodes
- 4 Rectangular self-adhesive electrodes
- 1 Charger
- 1 Operating manual

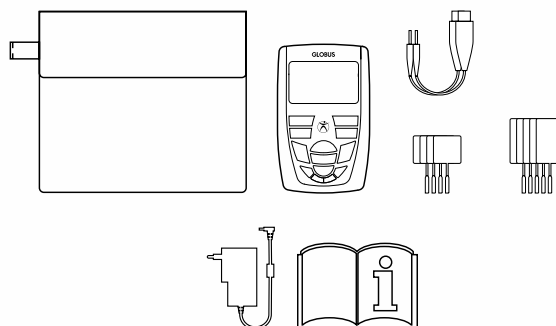

CE
0476

STRONG POINTS





DYSPHAGIA REHABILITATION

Dysphagia is a difficulty in carrying out the act of swallowing.

For those affected, it becomes very complicated, or even impossible, to swallow both saliva and liquid and/or solid food.

In most cases, dysphagia is due to neurological pathologies, but it frequently affects elderly people who suffer from a physiological weakening of the muscles involved in swallowing.

OUR DEVICES

AMI PRO and **AMI LIFE** are support devices for the rehabilitation of these patients and they act by promoting muscle trophism and improving proprioceptive neuromuscular sensitivity with the aim of re-educating the muscles involved in swallowing.

The devices intervene by inducing a non-invasive neuromuscular stimulation (NMES) that stimulates the muscles through electrodes positioned on the neck.

The protocols comply with the guidelines of the main scientific evidence currently available in the literature.

AMI PRO is the top-of-the-range device, fully equipped with different stimulation programs with motor, sensorimotor and sensory goals, with different parameters to enable the most appropriate choice of the work cycle. It also features various facial rehabilitation programs as well as programs that can be edited to suit the professional's therapeutic needs.

AMI PRO has the exclusive special **"ACTION NOW"** function, to perform the stimulation only at the request of the operator or the patient.

By pressing the **"ACTION NOW"** button on the device, it is possible to coordinate muscle activation with the patient's voluntary swallowing, creating an actual self-feedback system.

AMI LIFE has 3 stimulation programs with motor, sensorimotor and sensory goals and a facial rehabilitation program.

It is not possible to edit the preset programs and for this reason it is suitable for patients who need to continue the therapy at home (always under the supervision of the doctor).

AMI PRO | AMI LIFE

4-CHANNEL ELECTROTHERAPY FOR DYSPHAGIA

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent

Frequency

1-80Hz

Pulse width

80-260µs

Intensity

max 20mA per channel

Power supply

Rechargeable batteries

Dimensions | Weight


100x160x35 mm | 450 g



CE
0476

STANDARD EQUIPMENT

- 1 Bag
- 1 Stimulator
- 4 Electrode connection cables
- 4 Butterfly electrodes
- 1 Charger
- 1 Operating manual

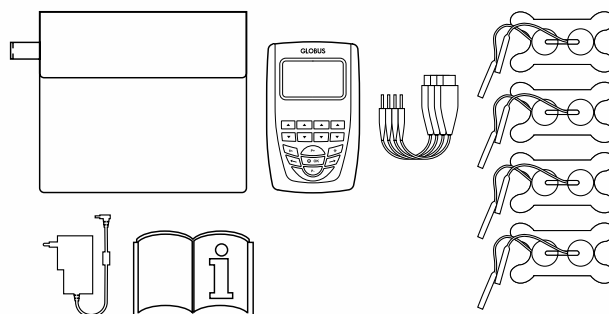

PURCHASE OPTIONS

- Code G5573 AMI PRO
- Code G5852 AMI LIFE

STRONG POINTS



STRONG POINTS



The use of these devices is only intended for doctors, authorised health professionals and their patients only under medical prescription and under the continuous supervision of the professional.

ULTRASOUND THERAPY



MEDISOUND 3000

MAGNETOTHERAPY | 49 PROGRAMS

TOP-OF-THE-RANGE PRODUCT SUITABLE FOR PHYSIOTHERAPY CENTRES AND HOSPITALS. BEING COMPACT, IT CAN BE EASILY TRANSPORTED IN ALL OPERATIONAL SITUATIONS.

31 AESTHETICS | 18 HEALTH

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

3W/cm² ± 20%

Contact management

Yes

Emitted frequency

3MHz and 1MHz

Emission

Continuous and pulsed

Duty cycle

10%-100%

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

170x220x60 mm

Weight

820 g

STANDARD EQUIPMENT

· Code G1033

· 1 Bag


· 1 MEDISOUND 3000 ultrasound therapy unit

· 1 Handpiece

· 1 Gel pack

· 1 Power supply unit

· 1 Operating manual

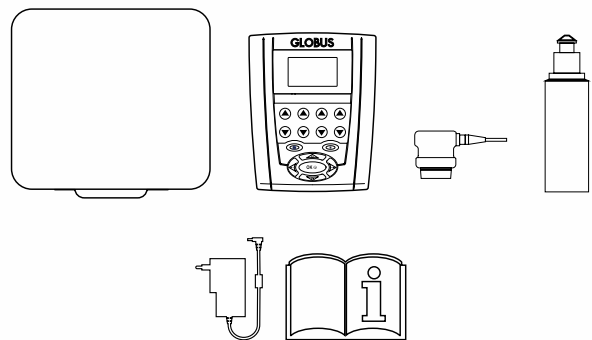

STRONG POINTS





LASER THERAPY for physiotherapy



PHYSIOLASER 12.0 PRO

LASER THERAPY

THE "PRO" VERSION OF THE PHYSIOLASER 12.0 IS DESIGNED FOR THE PROFESSIONAL WHO, FOR EASE OF USE, PREFERS A SLIGHTLY LARGER DEVICE THAT CAN ALWAYS BE KEPT READY FOR USE ON A CART. THE 12W OF MAXIMUM POWER WILL ALLOW QUICK AND EFFECTIVE TREATMENTS EVEN ON LARGER AREAS OF THE BODY BY USING THE 22MM OPTICS SUPPLIED WITH THE DEVICE. THE RICH LIBRARY OF PROGRAMS ALSO ENSURES THE PROFESSIONAL CAN TREAT ANY MUSCULOSKELETAL AILMENT AND SOME PROBLEMS RELATED TO WOUNDS AND SKIN ULCERS.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12W $\pm 20\%$

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

220x170x60 mm

Weight

1.12 kg

STANDARD EQUIPMENT

· Code **G5803**

· 1 Case

· 1 PHYSIOLASER 12.0 PRO laser therapy unit


· 1 Handpiece HP 808nm with Ø 11 mm and Ø 22 mm optics

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

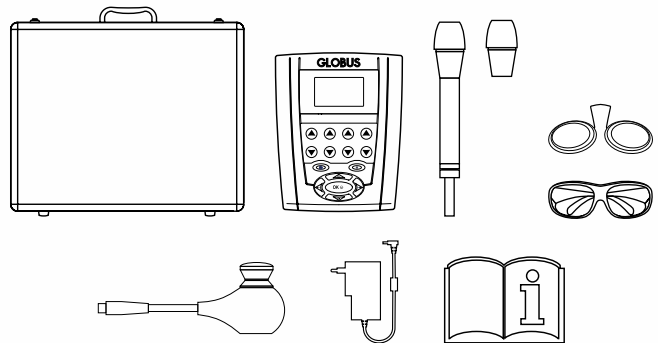

PURCHASE OPTIONS

· Code **G5995 HP Handpiece 980nm**

· Code **G5996 Handpiece 1064nm**



STRONG POINTS



PHYSIOLASER 12.0

LASER THERAPY

PHYSIOLASER 12.0 IS THE HANDIEST AND MOST POWERFUL LASER THERAPY DEVICE IN THE GLOBUS LINE. WITH ITS RICH LIBRARY OF PROGRAMS, PROFESSIONALS CAN TREAT ANY MUSCULOSKELETAL DISORDER, AS WELL AS SKIN WOUNDS AND ULCERS. DEFINITELY ONE OF THE DEVICES WITH THE BEST PRICE-QUALITY RATIO IN THE ENTIRE RANGE.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12W $\pm 20\%$

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains

- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code **G5741**

· 1 Case

· 1 PHYSIOLASER 12.0 laser therapy unit


· 1 Handpiece HP 808nm with Ø 11 mm and Ø 22 mm optics

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

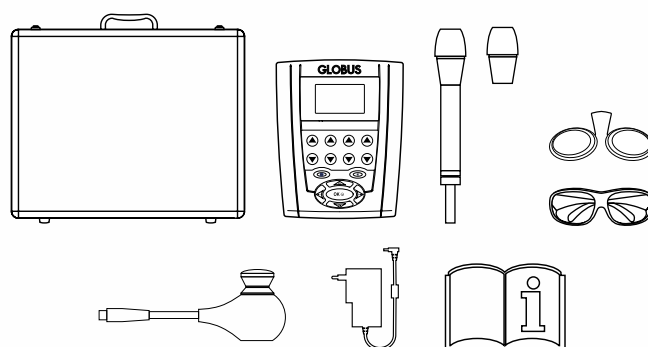

PURCHASE OPTIONS

· Code **G5995 HP Handpiece 980nm**

· Code **G5996 Handpiece 1064nm**



STRONG POINTS



PHYSIOLASER 6.0 PRO

LASER THERAPY

THE PHYSIOLASER 6.0 PRO IS A DEVICE OF LIMITED WEIGHT AND SIZE THAT CAN BE PLACED IN A PRACTICE TROLLEY OR TRANSPORTED TO OTHER ROOMS IF NECESSARY. WITH ITS MAXIMUM POWER OF 6W, IT IS IDEAL FOR TREATING OSTEOARTICULAR AND MUSCULAR DISORDERS, QUICKLY REDUCING PAIN AND OEDEMA AND ACCELERATING THE RECOVERY OF INJURED AND INFLAMED TISSUES.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6W $\pm 20\%$

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains

- Rechargeable battery

Dimensions

220x170x60 mm

Weight

1.12 kg

STANDARD EQUIPMENT

· Code G5764

· 1 Case

· 1 PHYSIOLASER 6.0 PRO laser therapy unit


· 1 Handpiece MP 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

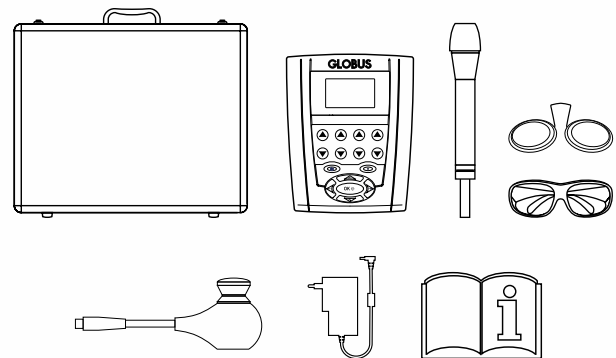

PURCHASE OPTIONS

· Code G5994 Handpiece MP 980nm

· Code G5996 Handpiece 1064nm

· Code G6092 Ø 22 mm optic



STRONG POINTS



PHYSIOLASER 6.0

LASER THERAPY

THE PHYSIOLASER 6.0 IS A LIGHTWEIGHT AND PORTABLE DEVICE, SUITABLE FOR CARRYING AROUND WITH YOU WHEN PERFORMING TREATMENTS IN MULTIPLE CLINICS. WITH ITS MAXIMUM POWER OF 6W, IT IS IDEAL FOR TREATING OSTEOARTICULAR AND MUSCULAR DISORDERS, QUICKLY REDUCING PAIN AND OEDEMA AND ACCELERATING THE RECOVERY OF INJURED AND INFLAMED TISSUES.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6W $\pm 20\%$

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains

- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G5740

· 1 Case

· 1 PHYSIOLASER 6.0 laser therapy unit


· 1 Handpiece MP 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

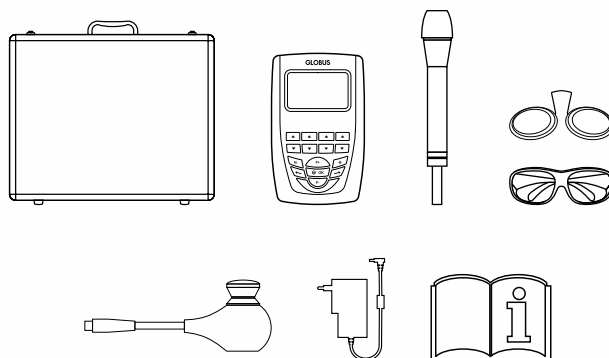

PURCHASE OPTIONS

· Code G5994 Handpiece MP 980nm

· Code G5996 Handpiece 1064nm

· Code G6092 Ø 22 mm optic



CE
0476

STRONG POINTS



PHYSIOLASER 2.0

LASER THERAPY

LIGHTWEIGHT AND PORTABLE DEVICE, CONTAINING ALL MAJOR PROGRAMS FOR THE TREATMENT OF CHRONIC AND ACUTE MUSCULOSKELETAL DISORDERS.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

2W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G5739

· 1 Case

· 1 PHYSIOLASER 2.0 laser therapy unit


· 1 Handpiece Black 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

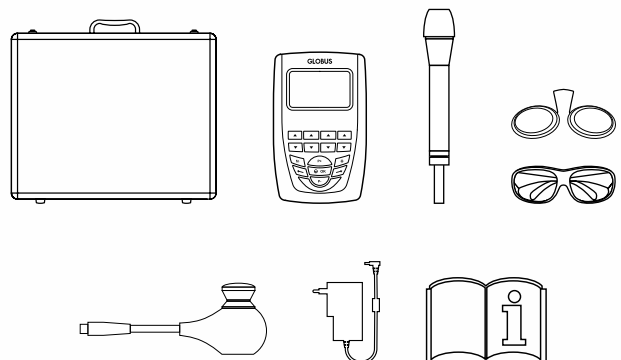

CE
0476



STRONG POINTS



PHYSIOLASER 1000

LASER THERAPY

PHYSIOLASER 1000 IS A PRODUCT OF THE PROFESSIONAL LINE. POWERFUL AND PORTABLE, THE PRODUCT IS DESIGNED FOR ACCELERATING THE RECOVERY OF INFLAMED TISSUES, REGENERATING CELLS AND IMPROVING MICROCIRCULATION.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

1W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code **G3579**

· 1 Case

· 1 PHYSIOLASER 1000 laser therapy unit


· 1 Handpiece Green 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

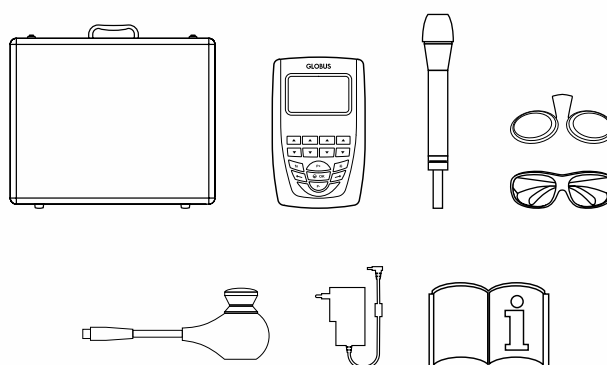

CE
0476



STRONG POINTS



PHYSIOLASER 500

LASER THERAPY

PRODUCT DESIGNED TO ACCELERATE THE RECOVERY OF INFLAMED TISSUE, REGENERATE CELLS AND IMPROVE MICROCIRCULATION. PHYSIOLASER 500 PROMOTES NATURAL HEALING PROCESSES.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

0.5W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code **G3786**

· 1 Case

· 1 PHYSIOLASER 500 laser therapy unit


· 1 Handpiece White 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button


· 1 Power supply unit

· 1 Operating manual

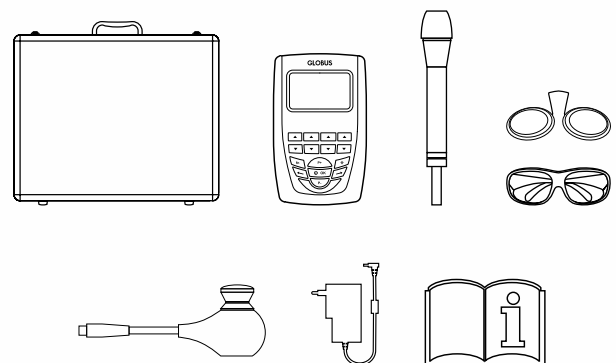

STRONG POINTS



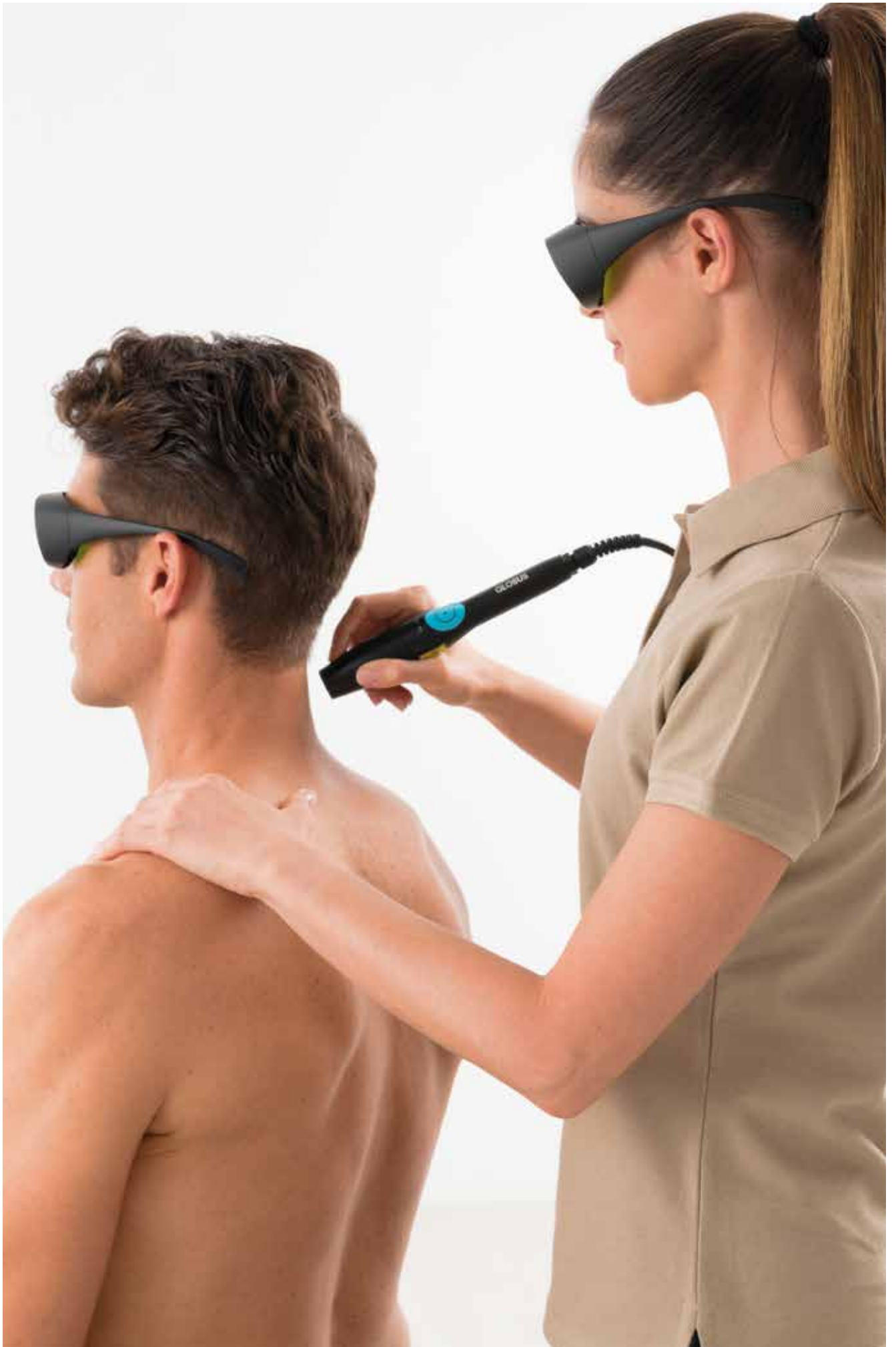

LASER THERAPY for podiatry



PODCARE 12.0 PRO

LASER THERAPY

TOP-OF-THE-RANGE DEVICE FOR PODIATRY LASERS. WITH ITS MAXIMUM POWER OF 12W, IT ALLOWS DIFFERENT PODIATRIC AND MUSCULOSKELETAL TREATMENTS TO BE PERFORMED, OPTIMISING TIME AND PROVIDING A GREATER DEPTH OF ACTION.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

12W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

220x170x60 mm

Weight

1.12 kg

STANDARD EQUIPMENT

· Code G5800

· 1 Case

· 1 PODCARE 12.0 PRO laser therapy unit


· 1 Handpiece HP 808nm with Ø 11 mm and Ø 22 mm optics

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

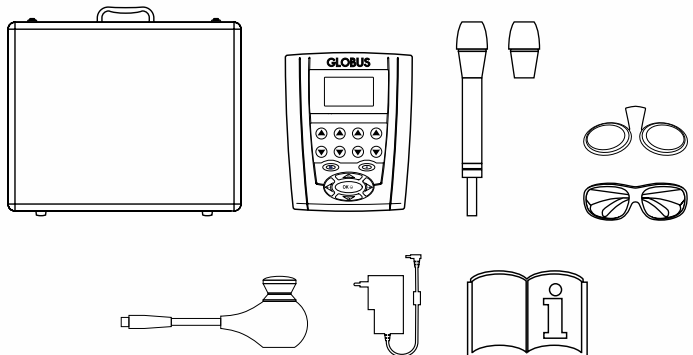

PURCHASE OPTIONS

· Code G5995 HP Handpiece 980nm

· Code G5996 Handpiece 1064nm



STRONG POINTS



PODCARE 6.0 PRO

LASER THERAPY

THE PODCARE 6.0 PRO IS A DEVICE OF LIMITED WEIGHT AND SIZE THAT CAN BE PLACED IN A PRACTICE TROLLEY OR TRANSPORTED TO OTHER ROOMS IF NECESSARY. THE DEVICE IS THE IDEAL DEVICE FOR FOOT AND ANKLE PODIATRY APPLICATIONS, BUT IS COMPLEMENTED BY A RANGE OF OTHER SPECIFIC PROGRAMS FOR ULCERS, WOUNDS AND MUSCULOSKELETAL DISORDERS.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions


220x170x60 mm

Weight

1.12 kg

STANDARD EQUIPMENT

- Code G5762
- 1 Case
- 1 PODCARE 6.0 PRO laser therapy unit
- 1 Handpiece MP 808nm with Ø 11 mm optic
- 2 Pairs of protective goggles (patient and operator)
- 1 Emergency button
- 1 Power supply unit
- 1 Operating manual

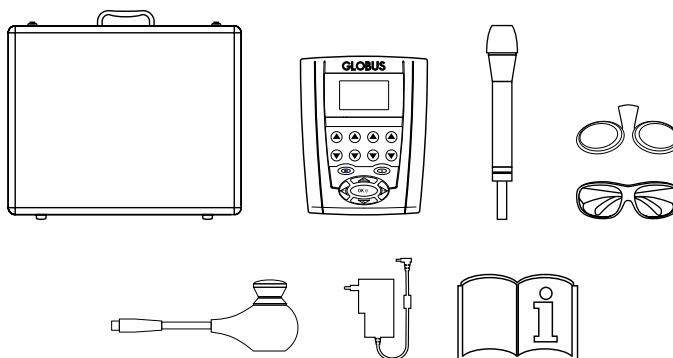

PURCHASE OPTIONS

- Code G5994 Handpiece MP 980nm
- Code G5996 Handpiece 1064nm
- Code G6092 Ø 22 mm optic



CE
0476

STRONG POINTS



PODCARE 6.0

LASER THERAPY

THE PODCARE 6.0 IS A LIGHTWEIGHT AND PORTABLE DEVICE, SUITABLE FOR CARRYING AROUND WITH YOU WHEN PERFORMING TREATMENTS IN MULTIPLE CLINICS. IT IS EQUIPPED WITH AN 808 nm DIODE WITH A MAXIMUM POWER OF 6W. THE DEVICE IS THE IDEAL DEVICE FOR FOOT AND ANKLE PODIATRY APPLICATIONS, BUT IS COMPLEMENTED BY A RANGE OF OTHER SPECIFIC PROGRAMS FOR ULCERS, WOUNDS AND MUSCULOSKELETAL DISORDERS.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

6W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G5737

· 1 Case

· 1 PODCARE 6.0 laser therapy unit


· 1 Handpiece MP 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles (patient and operator)

· 1 Emergency button

· 1 Power supply unit

· 1 Operating manual

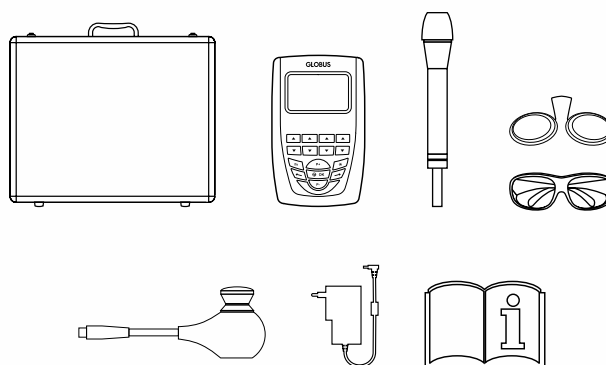

PURCHASE OPTIONS

· Code G5994 Handpiece MP 980nm

· Code G5996 Handpiece 1064nm

· Code G6092 Ø 22 mm optic



CE
0476

STRONG POINTS



PODCARE 2.0

LASER THERAPY

PODCARE 2.0 IS A PORTABLE DEVICE WITH PROTOCOLS TO TREAT DIFFERENT MUSCULOSKELETAL PATHOLOGIES AND WITH SPECIAL PODIATRY APPLICATIONS SUCH AS ONYCHOMYCOSIS, ONYCHOCRYPTOSIS AND OTHER TYPICAL INFLAMMATORY PATHOLOGIES OF THE FOOT AND THE LEG.

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Laser source

Class 4 GaAlAs diode

Maximum power

2W $\pm 20\%$

Compatible wavelength(s)

808 nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions


100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

- Code **G5736**
- 1 Case
- 1 PODCARE 2.0 laser therapy unit
- 1 Handpiece Black 808nm with Ø 11 mm optic
- 2 Pairs of protective goggles (patient and operator)
- 1 Emergency button
- 1 Power supply unit
- 1 Operating manual

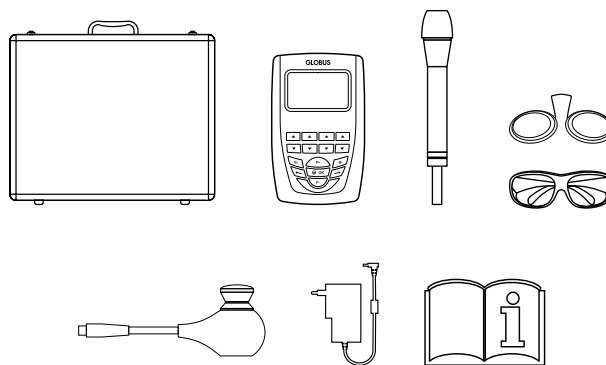

PURCHASE OPTION

- Code **G5547** Small Ø 3 mm optic



CE
0476



STRONG POINTS



PRESSOTHERAPY



PRESSCARE G300M

PRESSOTHERAPY | 2 OUTPUTS - 3 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, AESTHETICS AND RELAXATION NEEDS. SUITABLE FOR TREATING LYMPHOEDEMA, VENOUS CIRCULATION PROBLEMS, HEAVY LEGS, OEDEMA, SWOLLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, HAEMATOMA MASSAGE, VENOUS STASIS. IT CAN BE USED SIMULTANEOUSLY WITH THREE DIFFERENT ACCESSORIES. IT IS RECOMMENDED TO CONSULT YOUR DOCTOR BEFORE USE.

TECHNICAL SPECIFICATIONS

Outputs

2

Programs

3

Pressure

0-180mmHg $\pm 10\%$

Timer

Up to 99 minutes

Mains power supply

230V-AC, 50Hz

Dimensions

185x180x230 mm


Weight

3 kg




STANDARD EQUIPMENT

- Depending on the configuration
- 1 Power supply unit
- 1 Operating manual

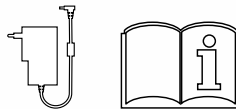

PURCHASING CONFIGURATIONS

- Code **G5411** PRESSCARE G300M-1, with 1 boot size "S"
- Code **G5412** PRESSCARE G300M-1, with 1 boot size "M"
- Code **G5413** PRESSCARE G300M-1, with 1 boot size "L"
- Code **G5414** PRESSCARE G300M-2, with 2 boots size "S"
- Code **G5415** PRESSCARE G300M-2, with 2 boots size "M"
- Code **G5416** PRESSCARE G300M-2, with 2 boots size "L"
- Code **G5417** PRESSCARE G300M-3, with 2 boots size "S", + abdominal/glutes band
- Code **G5418** PRESSCARE G300M-3, with 2 boots size "M", + 1 abdominal/glutes band
- Code **G5419** PRESSCARE G300M-3, with 2 boots size "L", + 1 abdominal/glutes band
- Code **G5420** PRESSCARE G300M-1B, with 1 sleeve

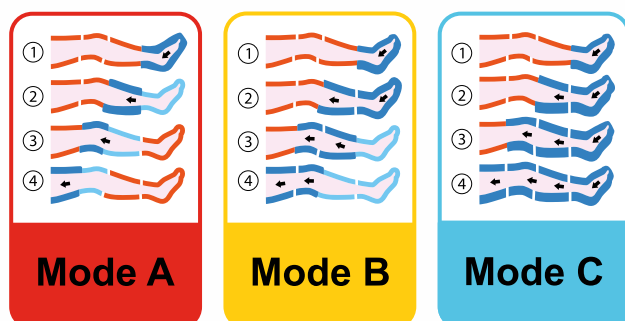

MODE A: the program sequence involves inflating and deflating the boot one sector at a time, starting from the most distal one (foot or hand) to the most proximal one.

MODE B: The program sequence involves inflating two chambers at a time, simulating the function of stacked-chamber boots.

MODE C: the program sequence applies pressure starting from the foot and keeping the previous chambers inflated until the entire part (leg, arm or abdomen) is inflated.

PRESSCARE G200M

PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, AESTHETICS AND RELAXATION NEEDS. SUITABLE FOR TREATING LYMPHOEDEMA, VENOUS CIRCULATION PROBLEMS, HEAVY LEGS, OEDEMA, SWOLLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, HAEMATOMA MASSAGE, VENOUS STASIS. IT CAN BE USED SIMULTANEOUSLY WITH THREE DIFFERENT ACCESSORIES. IT IS RECOMMENDED TO CONSULT YOUR DOCTOR BEFORE USE.

TECHNICAL SPECIFICATIONS


Outputs
2
Programs
2
Pressure
0-180mmHg $\pm 10\%$
Timer
Up to 99 minutes
Mains power supply
230V-AC, 50Hz
Dimensions
185x180x230 mm
Weight
3 kg



CE
0598


STANDARD EQUIPMENT

- Depending on the configuration
- 1 Power supply unit
- 1 Operating manual

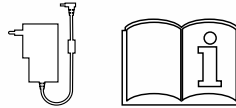

PURCHASING CONFIGURATIONS

- Code **G5401** PRESSCARE G200M-1, with 1 boot size "S"
- Code **G5402** PRESSCARE G200M-1, with 1 boot size "M"
- Code **G5403** PRESSCARE G200M-1, with 1 boot size "L"
- Code **G5404** PRESSCARE G200M-2, with 2 boots size "S"
- Code **G5405** PRESSCARE G200M-2, with 2 boots size "M"
- Code **G5406** PRESSCARE G200M-2, with 2 boots size "L"
- Code **G5407** PRESSCARE G200M-3, with 2 boots size "S", + 1 abdominal/glutes band
- Code **G5408** PRESSCARE G200M-3, with 2 boots size "M", + 1 abdominal/glutes band
- Code **G5409** PRESSCARE G200M-3, with 2 boots size "L", + 1 abdominal/glutes band
- Code **G5410** PRESSCARE G200M-1B, with 1 sleeve

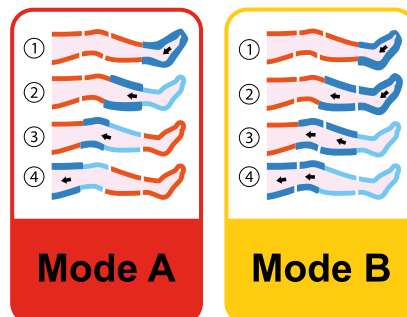

MODE A: the program sequence involves inflating and deflating the boot one sector at a time, starting from the most distal one (foot or hand) to the most proximal one.

MODE B: the program sequence involves inflating two chambers at a time, simulating the function of stacked-chamber boots.

BEAUTY & WELLNESS

CONTENTS

- TECAR BEAUTY	76
• TECAR BEAUTY 7000 MED	78
• TECAR BEAUTY 7000	79
• TECAR BEAUTY ELITE	80
• TECAR BEAUTY 6000 MED 6000 MED RE	81
• "FASCIA TOOLS" KIT	82
• "INTIMA" KIT	84
- ELECTROSTIMULATION	86
• SHAPE & BEAUTY	89
• ACTIVA 700	90
• PREMIUM 400	91
• ELITE 150	92
• ELITE	93
• ELITE S2	94
• FAST BAND	96
• FAST PAD	97
- RADIOFREQUENCY	98
• RF CLINIC PRO	100
• RF CLINIC BODY	101
- ULTRASOUND	102
• LIPOZERO G39	104
• LIPOZERO G150	105
• LIPOZERO G800	106
• LIPOZERO EXCEL	107

TECAR BEAUTY



TECAR BEAUTY 7000 MED

COSMETIC TECAR THERAPY

THE BEST MEDICALLY CERTIFIED DEVICE IN THE TECAR BEAUTY LINE FROM GLOBUS. THE DEVICE HAS AN EXTENSIVE LIST OF PROGRAMS THAT INCLUDES BOTH COSMETIC TREATMENTS AND PROTOCOLS FOR VARIOUS MUSCULOSKELETAL DISORDERS. THE TOUCHSCREEN NAVIGATION MENU ON A COLOUR SCREEN INCLUDES AN ATLAS WITH SEVERAL PICTURES AND A PRACTICAL DESCRIPTION OF HOW TO PERFORM THE TREATMENTS. MOREOVER, TECAR BEAUTY 7000 MED ALLOWS EXPERT USERS TO CREATE AND STORE NEW PROGRAMS IN JUST A FEW STEPS.

TECHNICAL SPECIFICATIONS

Display

Size 7" colour touchscreen display

Nominal Power (NP)

350W $\pm 10\%$

Emission frequencies

300kHz, 400kHz, 448kHz, 470kHz, 700kHz, 1000kHz, 1200kHz, $\pm 10\%$

Treatment methods

Capacitive and Resistive

Capacitive Electrodes

3 \emptyset sizes (30, 50, 70 mm)

Resistive Electrodes

3 \emptyset sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

515x340x170 mm

Weight

11.5 kg

STANDARD EQUIPMENT

· Code **G6244**

· 1 TECAR BEAUTY 7000 MED Tecar therapy unit

· 1 Return electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)


· 1 Bipolar handpiece \emptyset (30 mm)

· 1 Hands-free treatment kit

· 1 1000 ml jar of conductive cream

· 1 Power supply cable

· 1 Operating manual

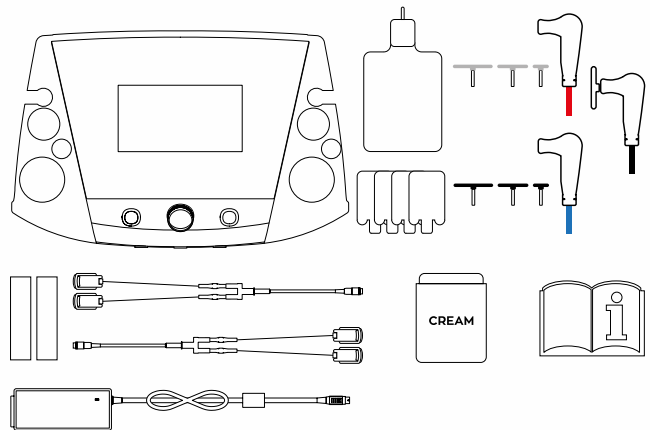

**CERTIFIED
PROTOCOLS**

STRONG POINTS



TECAR BEAUTY 7000

COSMETIC TECAR THERAPY

POWERFUL AND VERSATILE DEVICE FROM THE TECAR BEAUTY RANGE. NAVIGATION SYSTEM WITH COLOUR TOUCHSCREEN TO MAKE TREATMENT SEARCH AND PARAMETER MODIFICATION EASIER, AN ATLAS THAT DESCRIBES THE TREATMENTS DIVIDED INTO DIFFERENT PARTS OF THE BODY WITH COLOUR PHOTOS AND A SECTION WHERE THE OPERATOR CAN STRUCTURE AND STORE PERSONALISED TREATMENTS.

THE BEST EXPRESSION OF TECAR BEAUTY FOR BEAUTY TREATMENTS!

TECHNICAL SPECIFICATIONS

Display

Size 7" colour touchscreen display

Power

350W $\pm 10\%$

Emission frequencies

400kHz, 448kHz, 470kHz, 700kHz, 1000kHz, 1200kHz, $\pm 10\%$

Treatment methods

Capacitive and Resistive

Capacitive Electrodes

3 \varnothing sizes (30, 50, 70 mm)

Resistive Electrodes

3 \varnothing sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

515x340x170 mm

Weight

11.5 kg

STANDARD EQUIPMENT

· Code **G6248**

· 1 TECAR BEAUTY 7000 Tecar therapy unit

· 1 Return electrode (flexible plate)

· 1 Monopolar resistive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Monopolar capacitive handpiece with 3 electrodes \varnothing (30, 50, 70 mm)

· 1 Bipolar handpiece \varnothing (30 mm)


· 1 Hands-free treatment kit

· 1 500 ml pack of conductive hyaluronic acid face cream

· 1 1000 ml pack of conductive lymph drainage cream with aescin, fucus and caffeine

· 1 Power supply cable

· 1 Operating manual

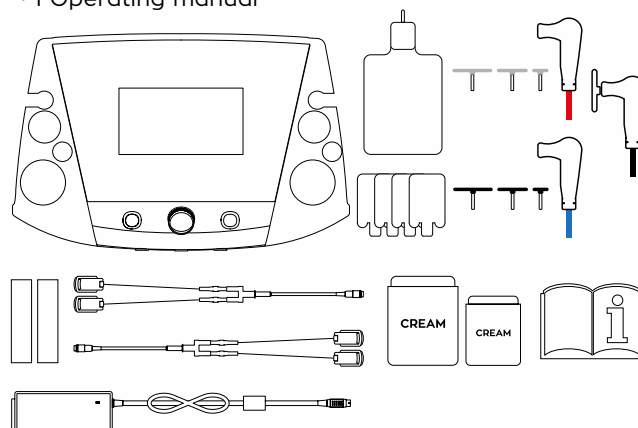

STRONG POINTS



TECAR BEAUTY ELITE

COSMETIC TECAR THERAPY

TECAR BEAUTY IS A NON-INVASIVE METHOD THAT FIGHTS AGEING AND SAGGING SKIN. THE ELECTROMAGNETIC WAVES PRODUCED BY THE DEVICE PENETRATE INTO THE DERMIS AND RELEASE HEAT IN THE TISSUES. THIS STIMULATES THE SYNTHESIS OF NEW COLLAGEN AND ELASTIN, VISIBLY IMPROVING THE COMPACTNESS OF THE SKIN. IN ADDITION, THE HEAT FACILITATES MICROCIRCULATION, PROMOTING A BETTER OXYGENATION OF THE TISSUES. PROFESSIONAL DEVICE SUITABLE FOR BEAUTY CENTRES.

TECHNICAL SPECIFICATIONS

Display

Size 7" colour touchscreen display

Power

50W $\pm 10\%$ Capacitive

25W $\pm 10\%$ Resistive

Emission frequencies

400kHz, 448kHz, 470kHz, 700kHz,
1000kHz, 1200kHz, $\pm 10\%$

Treatment methods

Capacitive and Resistive

Capacitive Electrodes

3 \emptyset sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

515x340x170 mm

Weight

11.5 kg

STANDARD EQUIPMENT

· Code **G6243**

· 1 TECAR BEAUTY ELITE Tecar therapy unit

· 1 Return electrode (flexible plate)

· 1 Monopolar capacitive handpiece with 3 electrodes
 \emptyset (30, 50, 70 mm)


· 1 Bipolar handpiece \emptyset (30 mm)

· 1 500 ml pack of conductive hyaluronic acid face
cream

· 1 1000 ml pack of conductive lymph drainage
cream with aescin, fucus and caffeine

· 1 Power supply cable

· 1 Operating manual

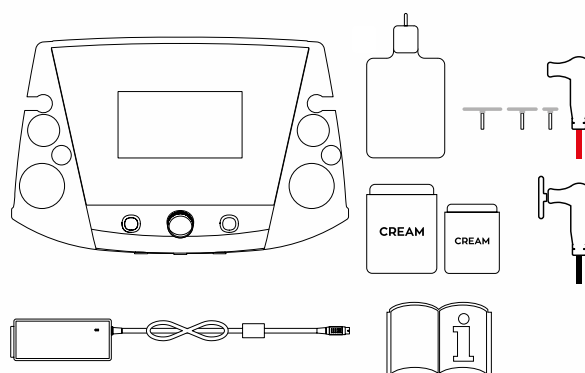

In compliance with the
Interministerial Decree no.
206/2015, the device is suitable
for use in beauty centres.

STRONG POINTS



“FASCIA TOOLS” KIT

ADVANCED SOFT TISSUE MOBILISATION INSTRUMENTS FOR USE IN OUR LINE OF “TECAR BEAUTY” DEVICES

The anatomical shape of the special resistive handpieces makes the treatment particularly effective on thighs, gluteus, hips and arms. The combination of the manual treatment with the myofascial handpieces, together with the energy of TECAR BEAUTY, offers an extraordinary and new treatment possibility.

TECAR BEAUTY with the special resistive handpieces makes it possible to:

- **MOBILISE CONNECTIVE DENSIFICATIONS**
- **VASCULARISE AND TREAT FIBROUS AND SCLEROTIC NODULES**
- **WORK WITH MAXIMUM ERGONOMICS AND TREATMENT PRECISION**



Large DiaBlade



Small DiaBlade



Large DiaBlade



Small DiaBlade





“INTIMA” KIT

ADVANCED TOOLS TO IMPROVE THE MUSCLE TONE OF THE PELVIC FLOOR FOR USE IN OUR LINE OF “TECAR BEAUTY” DEVICES

The precision of our TECAR devices offers effective action in total safety and comfort, even for solving problems directly related to the genital area.

This is why GLOBUS has created the “INTIMA” KIT, with specific accessories for internal and external treatments to improve the health of the pelvic floor, both in men and women, and to act on vaginal rejuvenation and elasticity, and also to provide excellent results when treating childbirth scars.

The most frequent treatments with the “INTIMA” KIT, ideal for restoring intimate comfort, are:

- VAGINAL LAXITY
- CHILDBIRTH SCARS
- ABDOMINAL DIASTASIS
- ENDOPELVIC MYOFASCIAL SYNDROMES
- PERINEAL REHABILITATION FOR MEN AND WOMEN
- PELVIC FLOOR REJUVENATION
- LIFTING

Code **G6421** “INTIMA” KIT for TECAR BEAUTY 7000 LINE | Code **G6420** “INTIMA” KIT for TECAR BEAUTY 6000 MED





ELECTROSTIMULATION





SHAPE & BEAUTY

4-CHANNEL ELECTROSTIMULATION

SHAPE & BEAUTY IS A PROFESSIONAL DEVICE FOR TONING AND FIRING THE BODY TO SHAPE THE CURVES AND IMPROVE THE SHAPE OF THE BODY IN JUST A FEW SESSIONS.

THE PRESET PROTOCOLS INSIDE THE DEVICE WILL GUIDE THE CHOICE OF THE MOST SUITABLE APPLICATION ACCORDING TO THE GOAL, FOLLOWING A WORK PROGRESSION THAT WILL OBTAIN VISIBLE RESULTS AFTER A FEW WEEKS.

THE DEVICE ALSO HAS RELAXING AND DECONTRACTING MASSAGE PROGRAMS AND VASCULARISATION AND CAPILLARISATION PROGRAMS TO IMPROVE THE DRAINAGE OF EXCESS LIQUIDS.

THE G-TRODE HANDPIECE INCLUDED HAS A VERY POWERFUL SKIN REGENERATING AND BIOSTIMULATING EFFECT, ACTING ON EXPRESSION WRINKLES, STRETCH MARKS AND SMALL ACNE SCARS. SHAPE & BEAUTY IS THE PERFECT ALLY FOR YOUR BEAUTY CLINIC!

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

170x220x60 mm

Weight

1.12 kg

* Depending on current type

STANDARD EQUIPMENT

· Code G6230

· 1 Case

· 1 SHAPE & BEAUTY stimulator

· 1 G-Trode handpiece - Skin regeneration

· 4 Electrode connection cables


· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 260 ml gel pack

· 1 Charger

· 1 Operating manual

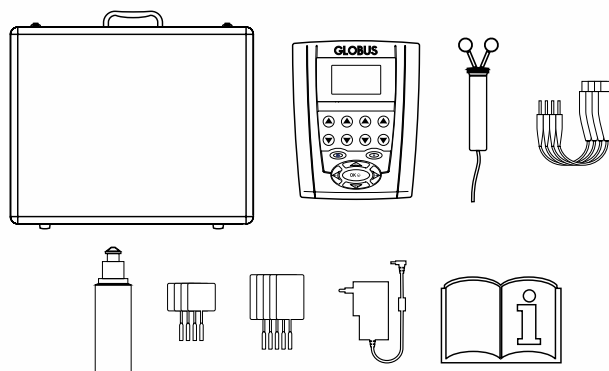

CE



STRONG POINTS



ACTIVA 700

4-CHANNEL ELECTROSTIMULATION | 271 PROGRAMS

PRODUCT FOR BEAUTY AND PHYSICAL SHAPE, DESIGNED WITH SPECIFIC PROGRAMS FOR FIRMING, TONING AND REJUVENATION. THE SPECIAL G-TRODE HANDPIECE (INCLUDED) MAKES ACTIVA 700 A TOP-OF-THE-RANGE IN THE ELECTROSTIMULATION LINE.

60 BEAUTY | 58 FITNESS | 53 SPORT | 36 SERIAL SEQUENTIAL STIMULATION | 21 SPECIAL SPORTS | 13 WRINKLES-SKIN
12 MICROCURRENTS | 12 TENS | 3 REHAB | 3 INCONTINENCE

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G3550

· 1 Case

· 1 ACTIVA 700 stimulator

· 1 G-Trode handpiece

· 4 Electrode connection cables

· 2 Microcurrent and ionophoresis cables


· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes


· 1 260 ml gel pack

· 1 Charger

· 1 Operating manual

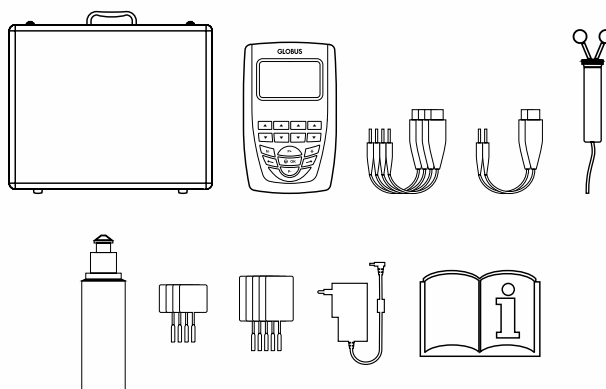

STRONG POINTS

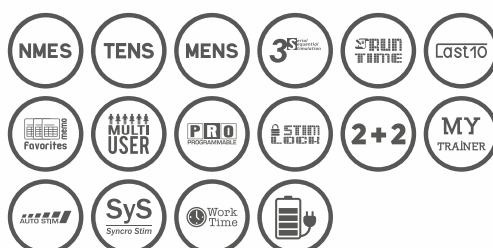

PREMIUM 400

4-CHANNEL ELECTROSTIMULATION | 258 PROGRAMS

MORE THAN 250 PROGRAMS ON THE PALM OF ONE HAND. PREMIUM 400 IS THE IDEAL ELECTROSTIMULATOR FOR SPORT AND FOR THE CARE OF YOUR BEAUTY, FOR HIM AND FOR HER.

58 FITNESS | 53 SPORT | 44 SPECIAL SPORTS | 36 BEAUTY | 23 MICROCURRENTS | 18 SERIAL SEQUENTIAL STIMULATION 12 TENS AND PAIN | 7 ACTION NOW | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G3551

· 1 Bag

· 1 PREMIUM 400 stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger


· 1 Operating manual

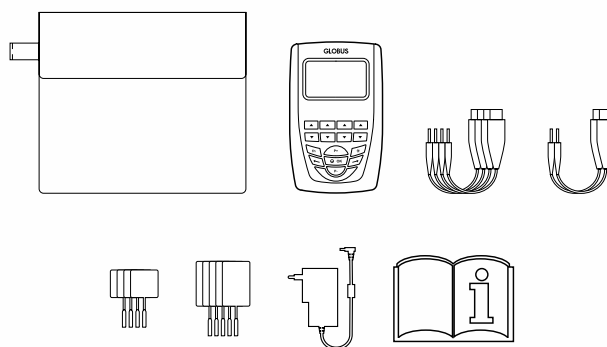

CE
0476



STRONG POINTS

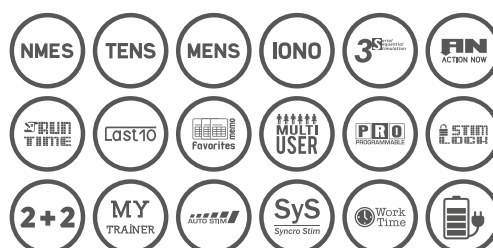

ELITE 150

4-CHANNEL ELECTROSTIMULATION | 150 PROGRAMS

PRODUCT DESIGNED FOR SPORT AND FITNESS, WITH SPECIAL ATTENTION TO BEAUTY TREATMENTS THANKS TO THE NEW SEQUENTIAL STIMULATION CURRENTS, EXCELLENT FOR THEIR DRAINING EFFECTS.

36 BEAUTY | 36 SERIAL SEQUENTIAL STIMULATION | 32 FITNESS | 31 SPORT | 12 TENS AND PAIN | 3 REHAB

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G1416

· 1 Bag

· 1 ELITE 150 stimulator


· 4 Electrode connection cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

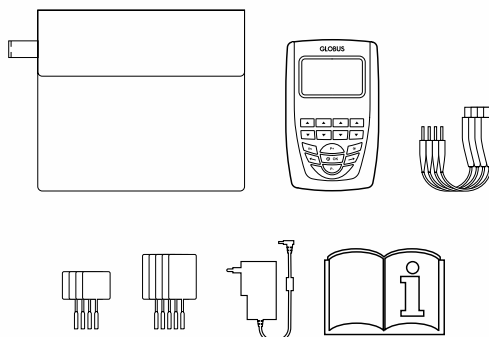

STRONG POINTS



ELITE

4-CHANNEL ELECTROSTIMULATION | 98 PROGRAMS

FOUR-CHANNEL PRODUCT FOR EVERY TREATMENT NEED.

32 SPORT | 29 BEAUTY | 26 FITNESS | 9 TENS | 2 REHAB

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

440 g

* Depending on current type

STANDARD EQUIPMENT

· Code G4300

· 1 Bag

· 1 ELITE stimulator


· 4 Electrode connection cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

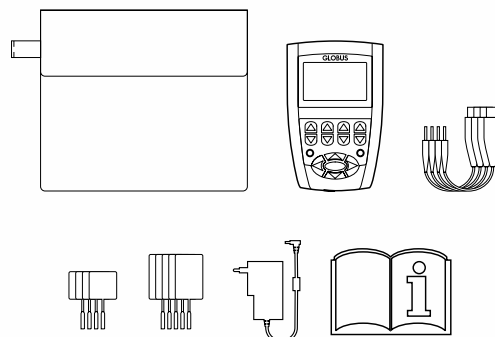

CE
0476



STRONG POINTS



ELITE S2

2-CHANNEL ELECTROSTIMULATION | 100 PROGRAMS

THE TWO-CHANNEL PRODUCT FOR EVERY TREATMENT NEED.

44 FITNESS, BEAUTY, WELLNESS | 30 SPORT | 16 TENS | 6 PREVENTION | 3 REHAB | 1 INCONTINENCE

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.5"

Channels

2 independent (4 electrodes)

Frequency

1-150Hz

Pulse width

50-400µs

Intensity

0-100mA per channel

Power supply

Rechargeable batteries

Dimensions

80x130x25 mm

Weight

220 g

STANDARD EQUIPMENT

· Code G3552

· 1 Bag

· 1 ELITE S2 stimulator


· 2 Electrode connection cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

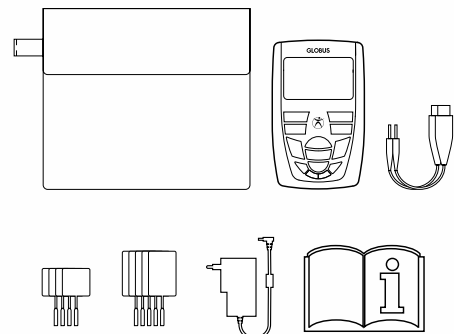

STRONG POINTS





FAST BAND

ABDOMINAL BAND

FAST BAND IS A BAND THAT IS EASILY CONNECTED TO ALL GLOBUS 2- AND 4-CHANNEL ELECTROSTIMULATORS. WITH THE FAST BAND, WITHOUT USING THE CLASSIC ELECTRODES, YOU CAN STIMULATE THE ABDOMINAL AND OBLIQUE MUSCLES, THE LUMBAR AREA AND THE GLUTES SIMPLY AND EFFECTIVELY.

TECHNICAL SPECIFICATIONS

Dimensions


980x150 mm

Weight

214 g


STANDARD EQUIPMENT

- Code **G0487**
- 1 FAST BAND
- 4 Clip adapters

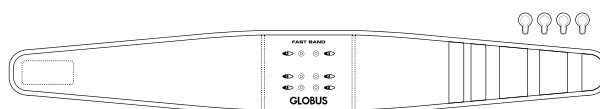

PURCHASE OPTIONS

- Code **G0488** FAST BODY KIT: "FAST BAND" + "FAST PAD"

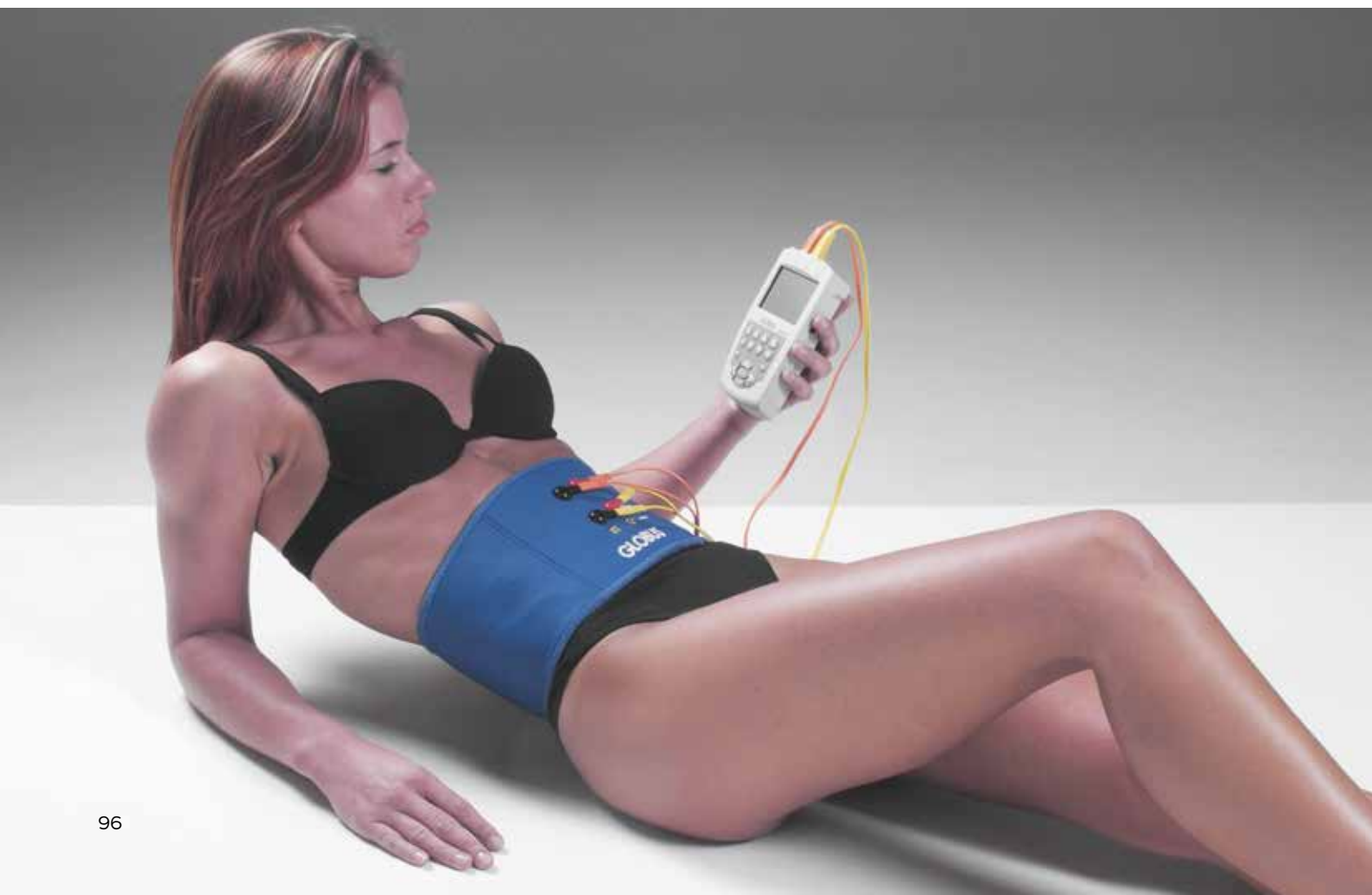

FAST PAD

LEG BANDS

FAST PADS ARE BANDS THAT ARE EASILY CONNECTED TO ALL GLOBUS 2- AND 4-CHANNEL ELECTROSTIMULATORS. THEY ARE IDEAL TO PERFORM SPECIFIC STIMULATION TREATMENTS FOR THE LEGS IN A COMFORTABLE AND EFFECTIVE WAY. THEY ARE ALSO CHARACTERISED BY AN OVAL SHAPE WITH LARGE CONTACT SURFACE AND ARE CONNECTED TO THE ELECTROSTIMULATOR THROUGH CLIP ADAPTERS.

TECHNICAL SPECIFICATIONS

Dimensions


250x150 mm

Weight

60 g (single band)

STANDARD EQUIPMENT

- Code **G0489**
- 2 FAST PADS
- 4 Clip adapters
- Gel

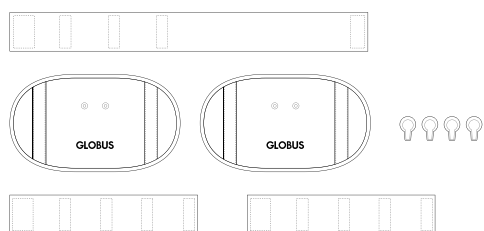

PURCHASE OPTIONS

- Code **G0488** FAST BODY KIT: "FAST BAND" + "FAST PAD"



RADIOFREQUENCY



RF CLINIC PRO

RADIOFREQUENCY | 30 PROGRAMS

RF CLINIC PRO IS AN INNOVATIVE PROFESSIONAL DEVICE THAT GENERATES RADIOFREQUENCY WAVES FOR AESTHETIC APPLICATIONS AGAINST SKIN AGEING ON THE FACE AND THE WHOLE BODY. THANKS TO THE 50 MM BIPOLAR HANDPIECE, IT IS POSSIBLE TO PERFORM FIRMING AND CELLULITE PROGRAMS.

13 SKIN | 10 CELLULITE | 7 FACE

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

- 1 RF Clinic 100 Ø 50 mm handpiece
- 1 RF Clinic 100 2-sphere handpiece

Power

200W (50 Ω)

Emitted frequency

470kHz to 2MHz

Treatment methods

Resistive

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

170x220x60 mm

Weight

920 g

STANDARD EQUIPMENT

· Code G1506

· 1 Case

· 1 RF CLINIC PRO radio frequency unit


· 1 RF Clinic 500 handpiece

· 1 RF Clinic 100 handpiece

· 1 250 ml jar of conductive cream

· 1 Power supply unit

· 1 Operating manual

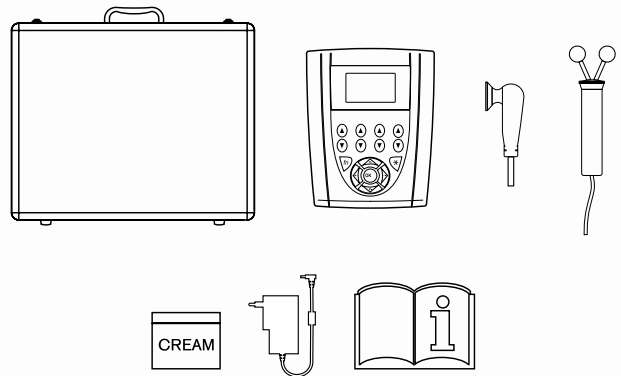

CE



STRONG POINTS



RF CLINIC BODY

RADIOFREQUENCY | 10 PROGRAMS

RF CLINIC BODY IS AN INNOVATIVE DEVICE THAT GENERATES RADIOFREQUENCY WAVES FOR AESTHETIC APPLICATIONS AGAINST SKIN AGEING ON THE FACE AND THE WHOLE BODY. IT ALSO BOOSTS VASCULARISATION OF THE SKIN TISSUE, ELIMINATING IMPERFECTIONS CAUSED BY EXPOSURE TO SUN RAYS AND SKIN AGEING.

6 SKIN | 4 FACE

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

1 RF Clinic 100 2-sphere handpiece

Power

200W (50 Ω)

Emitted frequency

470kHz to 2MHz

Treatment methods

Resistive

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


170x220x60 mm

Weight

920 g

STANDARD EQUIPMENT

- Code G1507
- 1 Case
- 1 RF CLINIC BODY radio frequency unit
- 1 RF Clinic 100 handpiece
- 1 250 ml jar of conductive cream
- 1 Power supply unit
- 1 Operating manual

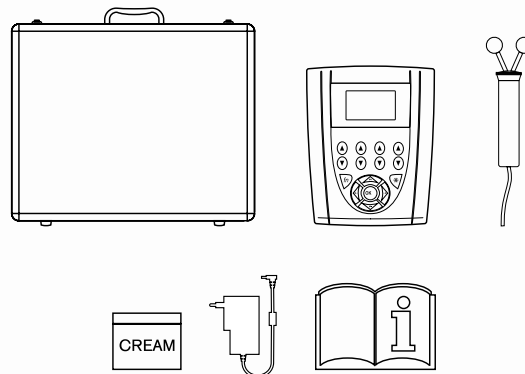

CE



STRONG POINTS



ULTRASOUND



LIPOZERO G39

AESTHETIC ULTRASOUND | 49 PROGRAMS

LIPOZERO G39 IS THE PROFESSIONAL LIPOLYTIC ULTRASOUND THAT FIGHTS CELLULITE AND ELIMINATES SUPERFLUOUS FAT THANKS TO PROGRAMS DESIGNED TO INTERVENE AT DIFFERENT FREQUENCIES AND ENSURE MAXIMUM EFFECTIVENESS IN EVERY BEAUTY TREATMENT.

10 MF-PLUS | 10 BEAUTY | 10 SKIN | 6 CAVIMAX 150 kHz CAVITATION
5 ADIPOSITY MOBILISATION NorAd-526 | 4 MALE ADIPOSITY | 4 ULTRACAV 39 kHz

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

- 1 Ø 50 mm - IP68 handpiece
- 1 Ø 42 mm - IP68 handpiece

Power

3W/cm² ± 20%

Emitted frequency

3MHz - 1MHz - 526kHz - 150kHz - 39kHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


170x220x60 mm

Weight

920 g

STANDARD EQUIPMENT

- Code G1404
- 1 Case
- 1 LIPOZERO G39 ultrasound
- 1 Multi-frequency handpiece Ø 42 mm - IP68
- 1 Handpiece for use with 39kHz, Ø 50 mm - IP68
- 1 250 ml gel pack
- 1 Power supply unit
- 1 Operating manual

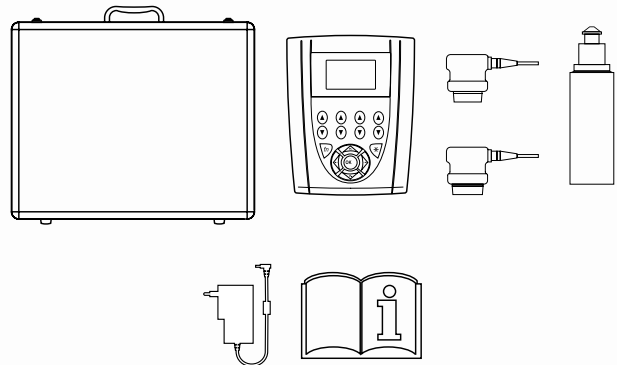

CE



STRONG POINTS



LIPOZERO G150

AESTHETIC ULTRASOUND | 42 PROGRAMS

ULTRASONIC CAVITATION IS THE MOST EFFECTIVE METHOD OF WORKING AGAINST CELLULITE AND ELIMINATING SUPERFLUOUS FAT. WITH LIPOZERO G150 YOU WILL BE ABLE TO BENEFIT FROM THE PROFESSIONAL TREATMENTS OF A BEAUTY CENTRE COMFORTABLY AT HOME.

10 MF-PLUS | 9 CAVIMAX 150 kHz CAVITATION | 8 BEAUTY | 8 SKIN
5 ADIPOSITY MOBILISATION NorAd-526 | 2 MALE ADIPOSITY

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

3W/cm² ± 20%

Emitted frequency

3MHz - 1MHz - 526kHz - 150kHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


170x220x60 mm

Weight

920 g

STANDARD EQUIPMENT

- Code G1398
- 1 Case
- 1 LIPOZERO G150 ultrasound
- 1 Multi-frequency handpiece Ø 42 mm - IP68
- 1 Gel pack
- 1 Power supply unit
- 1 Operating manual

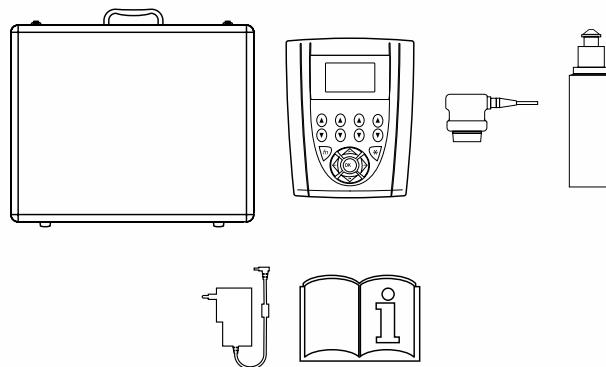

CE



STRONG POINTS



LIPOZERO G800

AESTHETIC ULTRASOUND | 30 PROGRAMS

A DEVICE FOR PROFESSIONAL AESTHETICS, SUITABLE FOR COMBATING CELLULITE AND LOCALISED FAT DEPOSITS. AND FOR IMPROVING THE APPEARANCE OF THE SKIN. THE POSSIBILITY OF DELIVERING WITH DIFFERENT FREQUENCIES GUARANTEES A SPECIFIC AND EFFECTIVE TREATMENT FOR EACH OBJECTIVE.

10 MF-PLUS | 10 BEAUTY | 10 SKIN

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

- 3W/cm² ± 20% at 1.3MHz and 3MHz
- 1.5W/cm² ± 20% at 1MHz and 800kHz

Emitted frequency

3MHz - 1.3MHz - 1MHz - 800kHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

170x220x60 mm

Weight

920 g

STANDARD EQUIPMENT

· Code G6232

· 1 Case


· 1 LIPOZERO G800 ultrasound

· 1 Multi-frequency handpiece Ø 42 mm - IP68

· 1 Gel pack

· 1 Power supply unit

· 1 Operating manual

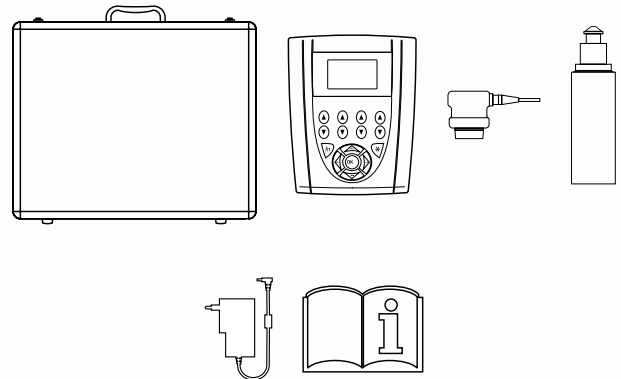

CE

In compliance with the Interministerial Decree no. 206/2015, the device is suitable for use in beauty centres.



STRONG POINTS



LIPOZERO EXCEL

AESTHETIC ULTRASOUND | 7 PROGRAMS

LIPOZERO EXCEL IS PART OF THE GLOBUS HOME CARE LINE AND HAS BEEN DESIGNED TO CARRY OUT THE PROFESSIONAL TREATMENTS OF AESTHETIC CENTRES AT HOME. THE PORTABLE DESIGN AND THE GRAPHIC MENU HAVE BEEN SPECIFICALLY DESIGNED FOR HOME USE AND MAKE THE DEVICE EASY TO USE EVEN BY INEXPERIENCED USERS. LIPOZERO EXCEL ALSO ALLOWS THE TREATMENT OF THE MOST COMMON IMPERFECTIONS RELATED TO CELLULITE OR FAT ACCUMULATION.

5 BEAUTY | 2 CELLULITE

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Outputs

1

Handpiece

1 Ø 42 mm - IP68 handpiece

Power

1.2W/cm² ± 20%

Emitted frequency

3MHz

Emission

Continuous and pulsed

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


100x160x35 mm

Weight

270 g

STANDARD EQUIPMENT

- Code G1297
- 1 Case
- 1 LIPOZERO EXCEL ultrasound
- 1 3MHz Ø 42 mm - IP68 handpiece
- 1 Gel pack
- 1 Power supply unit
- 1 Operating manual

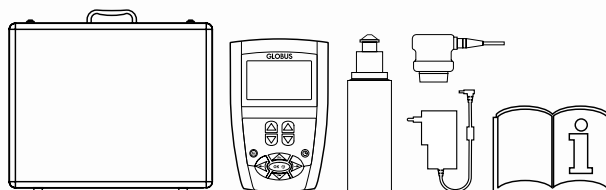

CE



STRONG POINTS



SPORT

CONTENTS

- ELECTROSTIMULATION	110
• THE CHAMPION	112
• TRIATHLON PRO	113
• CYCLING PRO	114
• SOCCER PRO	115
• RUNNER PRO	116
- PRESSOTHERAPY	118
• PRESSCARE G-SPORT ³	119
RELAX	120
• 4 VIBE MASSAGE GUN	122
- VIBRATION PLATFORM	124
• PHYSIOPLATE MY GOLD	125
- BALL SHOOTING MACHINE	126
• EUROGOAL 1500	128
• EUROGOAL 600	129
• WINSHOT 1500	130
• WINSHOT 800	131
• WINSHOT 500	132
• RUGBYKICK	133



ELECTROSTIMULATION



THE CHAMPION

4-CHANNEL ELECTROSTIMULATION | 441 PROGRAMS

IDEAL DEVICE FOR SPORTS LOVERS. THE DEVICE INCLUDES 12 SPORTS DISCIPLINES, EACH OF WHICH PROVIDES SPECIFIC PROGRAMS TO BUILD MUSCLES AND IMPROVE PERFORMANCE. THE PROGRAMS ARE SUITABLE BOTH FOR THOSE WHO PRACTICE SPORTS AT AMATEUR LEVEL AND FOR PROFESSIONALS WHO CARRY OUT STRENGTHENING ACTIVITIES THAT THEY WANT TO INTEGRATE WITH TARGETED ELECTROSTIMULATION.

90 SPECIAL SPORTS | 84 ACTION NOW | 70 BEAUTY | 58 FITNESS | 54 SERIAL SEQUENTIAL STIMULATION | 53 SPORT
23 MICROCURRENTS | 12 TENS AND PAIN | 3 REHAB | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G5015

· 1 Bag

· 1 THE CHAMPION stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger


· 1 Operating manual

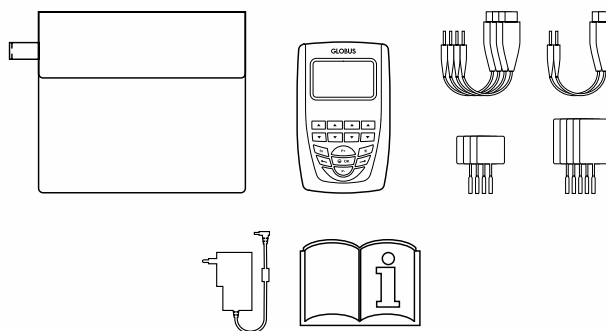

CE
0476



STRONG POINTS

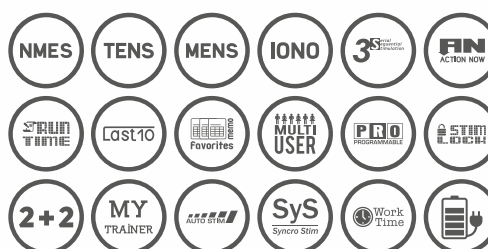

TRIATHLON PRO

4-CHANNEL ELECTROSTIMULATION | 424 PROGRAMS

DEVICE THAT, IN ADDITION TO A WIDE RANGE OF SPORTS, FITNESS, BEAUTY AND REHAB PROGRAMS, INCLUDES A SERIES OF SPECIFIC PROGRAMS AIMED AT IMPROVING THE PHYSICAL PERFORMANCE OF ATHLETES. THE SPECIFIC PROGRAMS FOR TRIATHLON PRO ARE DIVIDED INTO AREAS CORRESPONDING TO THE THREE SPORTS DISCIPLINES (SWIMMING, CYCLING, RUNNING). A PRACTICAL USER GUIDE WILL HELP YOU TO CHOOSE THE MOST SUITABLE PROGRAMS ACCORDING TO THE DISTANCES YOU PRACTICE.

84 ACTION NOW | 60 SPECIAL TRIATHLON | 60 BEAUTY | 58 FITNESS | 54 SERIAL SEQUENTIAL STIMULATION
53 SPORT | 23 MICROCURRENTS | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 3 REHAB | 3 INCONTINENCE | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm


Weight

450 g


* Depending on current type

STANDARD EQUIPMENT

- Code **G3666**
- 1 Bag
- 1 TRIATHLON PRO stimulator
- 4 Electrode connection cables
- 2 Microcurrent and ionophoresis cables
- 4 Square self-adhesive electrodes
- 4 Rectangular self-adhesive electrodes
- 1 Charger
- 1 Operating manual

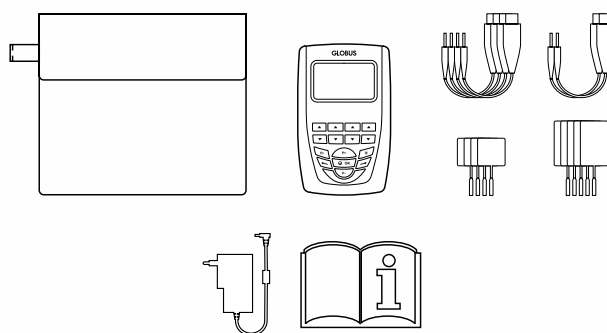

STRONG POINTS

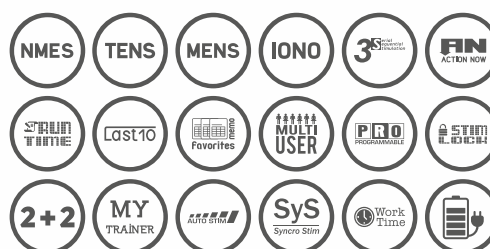

CYCLING PRO

4-CHANNEL ELECTROSTIMULATION | 270 PROGRAMS

DEVICE FOR CYCLISTS WHO WANT TO BE GOOD AND FEEL FIT DURING THEIR BICYCLE ACTIVITIES. THE SPECIFIC PROGRAMS FOR CYCLING ARE DIVIDED INTO MACRO-AREAS: PHYSICAL TRAINING, PREVENTION, PAIN-INJURY AND RECOVERY.

58 FITNESS | 53 SPORT | 46 SPECIAL CYCLING | 36 BEAUTY | 23 MICROCURRENTS
18 SERIAL SEQUENTIAL STIMULATION | 13 G-PULSE | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G4230

· 1 Bag

· 1 CYCLING PRO stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes


· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

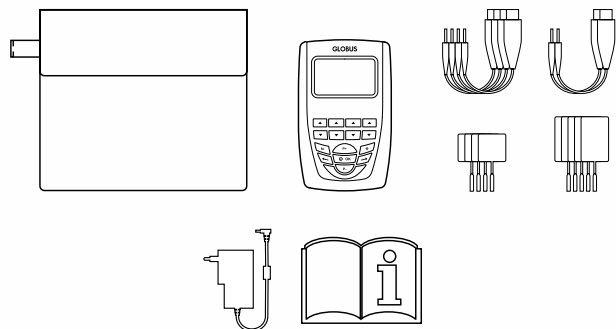

STRONG POINTS

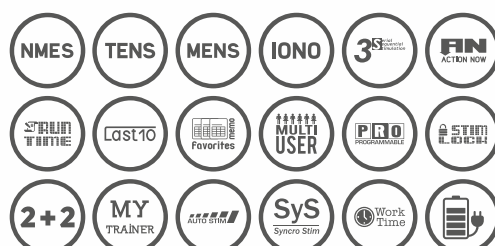

SOCCER PRO

4-CHANNEL ELECTROSTIMULATION | 258 PROGRAMS

SOCCER PRO IS THE ONLY ELECTROSTIMULATOR ON THE MARKET WITH SPECIFIC PROGRAMS TO IMPROVE THE PHYSICAL PERFORMANCE OF FOOTBALL PLAYERS. THE PROGRAMS ARE DIVIDED INTO FOUR LARGE AREAS CORRESPONDING TO THE MAIN FIELDS OF ATHLETIC TRAINING IN THE GAME OF FOOTBALL.

58 FITNESS | 53 SPORT | 36 BEAUTY | 34 SPECIAL SOCCER | 23 MICROCURRENTS

18 SERIAL SEQUENTIAL STIMULATION | 13 G-PULSE | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G4228

· 1 Bag

· 1 SOCCER PRO stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes


· 4 Rectangular self-adhesive electrodes

· 1 Charger

· 1 Operating manual

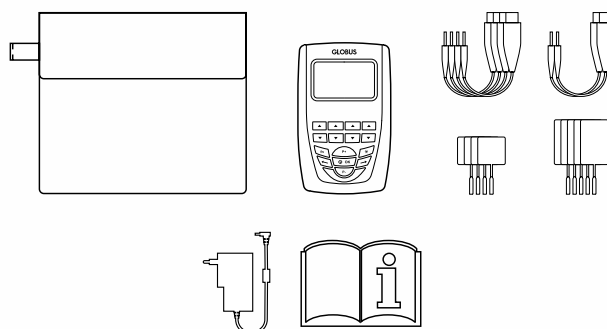

STRONG POINTS

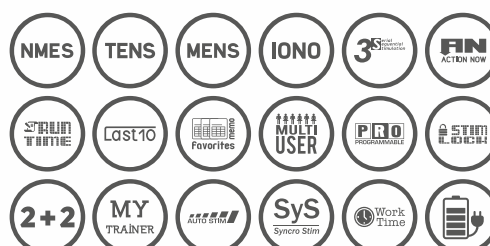

RUNNER PRO

4-CHANNEL ELECTROSTIMULATION | 254 PROGRAMS

RUNNER PRO IS A FUNDAMENTAL TOOL FOR THE MODERN RUNNER WHO WANTS TO OBTAIN IMPROVEMENTS IN PERFORMANCE. RUNNER PRO IS EQUIPPED WITH PROGRAMS SUITABLE TO STIMULATE THE MUSCLES OF THOSE WHO RUN AND EVERY RUNNER WILL BE ABLE TO PROGRAM AND CUSTOMISE THE TRAINING SESSIONS ACCORDING TO PERSONAL NEEDS. THE PROGRAMS, DIVIDED INTO FOUR MACRO-GOALS, ARE SUITABLE FOR BOTH AMATEUR AND PROFESSIONAL RUNNERS.

58 FITNESS | 53 SPORT | 36 BEAUTY | 30 SPECIAL RUNNING | 23 MICROCURRENTS
18 SERIAL SEQUENTIAL STIMULATION | 13 WRINKLES-SKIN | 12 TENS AND PAIN | 7 ACTION NOW | 3 REHAB | 1 IONOPHORESIS

TECHNICAL SPECIFICATIONS

Display

Visible area size 2.6"

Channels

4 independent (8 electrodes)

Frequency

0.3-150Hz *

Pulse width

40-450µs *

Intensity

0-120mA per channel

Power supply

Rechargeable batteries

Dimensions

100x160x35 mm

Weight

450 g

* Depending on current type

STANDARD EQUIPMENT

· Code G4229

· 1 Bag

· 1 RUNNER PRO stimulator

· 4 Electrode connection cables


· 2 Microcurrent and ionophoresis cables

· 4 Square self-adhesive electrodes

· 4 Rectangular self-adhesive electrodes

· 1 Charger


· 1 Operating manual

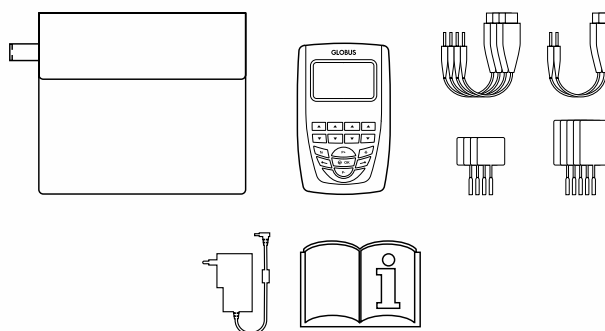

CE
0476



STRONG POINTS

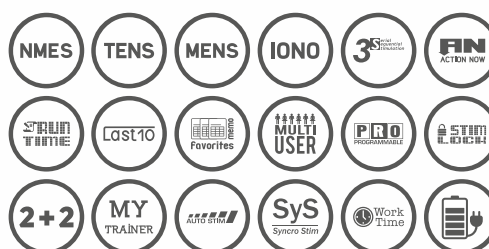



PRESSOTHERAPY

PRESSCARE G-SPORT³

PRESSOTHERAPY | 2 OUTPUTS - 3 INFLATION PROGRAMS

PNEUMATIC COMPRESSION DEVICE DESIGNED FOR ATHLETES, TO ACCELERATE MUSCLE RECOVERY FROM POST-TRAINING AND POST-COMPETITION FATIGUE. LIGHTWEIGHT AND TRANSPORTABLE WHEREVER YOU GO AND EQUIPPED WITH A POWER BANK THAT ENABLES TREATMENTS DIRECTLY ON THE FIELD AFTER A COMPETITION.

TECHNICAL SPECIFICATIONS


Outputs
2
Programs
3
Pressure
0-180mmHg $\pm 10\%$
Timer
Up to 99 minutes
Power supply
- Mains, 100-240V-AC, 50-60Hz - Rechargeable battery (Power Bank)
Dimensions
180x190x135 mm
Weight
3 kg

STANDARD EQUIPMENT

· Code G5380
· 1 Bag
· 1 G-SPORT ³ motor unit
· 2 Boots size "S"
· 2 Connection cables
· 1 Power supply unit
· 1 Power Bank
· 1 Operating manual

PURCHASE OPTIONS

· Code G5381 2 Boots size "M" included
· Code G5382 2 Boots size "L" included

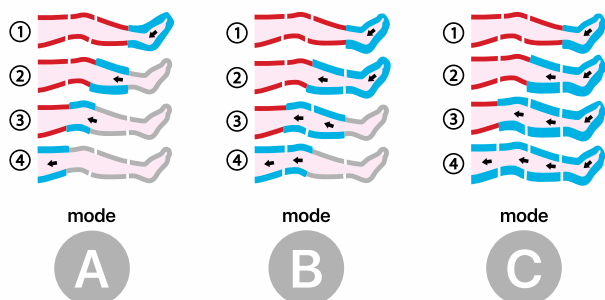


MODE A: the program sequence involves inflating and deflating the boot one sector at a time to carry out an activating massage, ideal before a competition or training.

MODE B: the program sequence consists in the inflation and deflation of two chambers at a time to avoid downward fluid reflux and make the treatment effective and very comfortable.

MODE C: the program sequence consists in the sequential inflation of all the chambers and their simultaneous deflation, ideal for a deeper massage.



RELAX

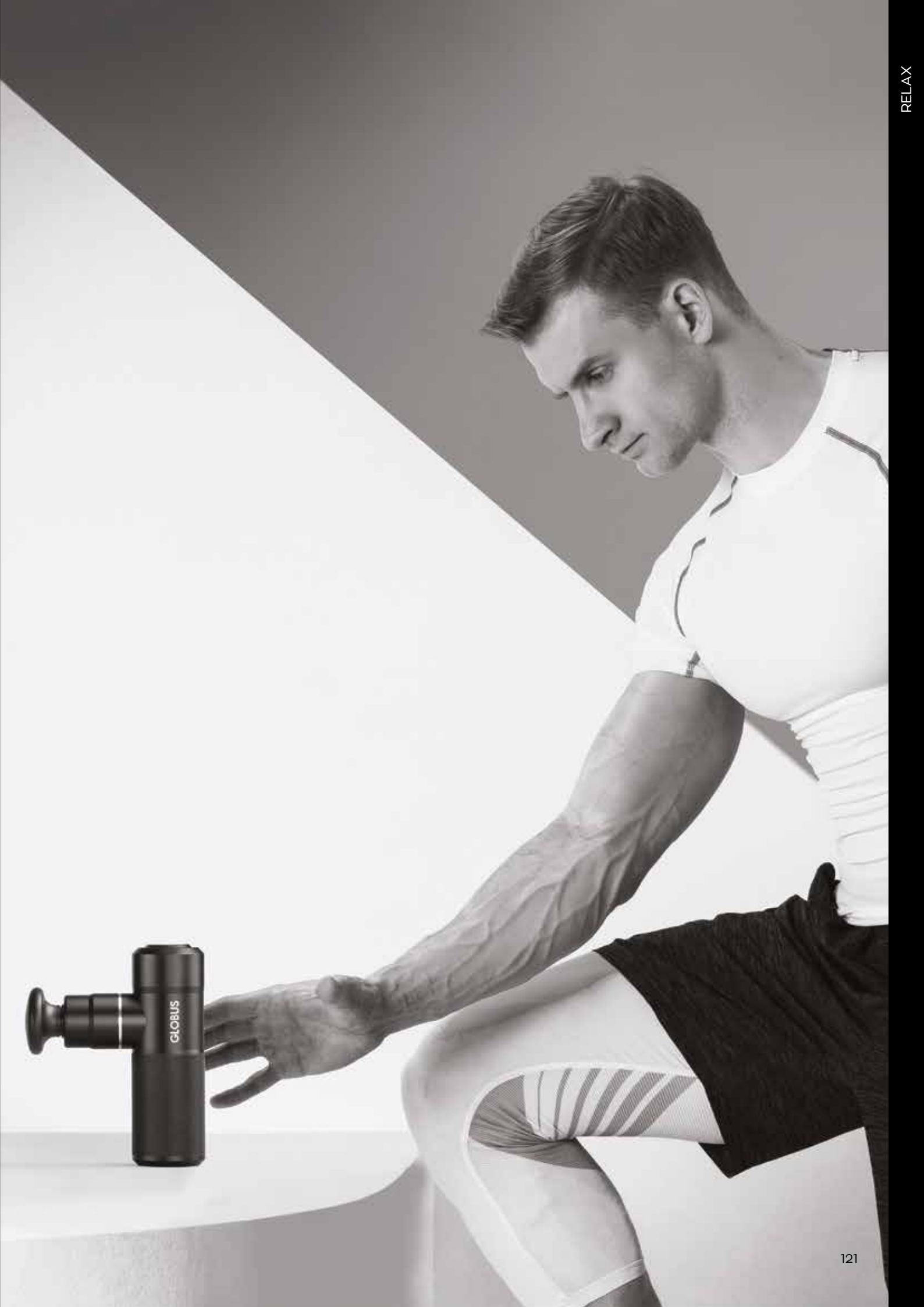

4 VIBE MASSAGE GUN

MASSAGE GUN | 4 SPEEDS

"4 VIBE MASSAGE GUN" IS THE MINI MASSAGE GUN MADE OF ALUMINIUM. IT HAS 4 DIFFERENT SPEED LEVELS AND IS EQUIPPED WITH VARIOUS HEADS TO ENSURE MAXIMUM EFFECTIVENESS ON THE DIFFERENT PARTS OF THE BODY TO BE TREATED AND ACCORDING TO NEED. IT IS IDEAL FOR MUSCLE RELAXATION, DEEP TISSUE MASSAGE AND POST-EXERCISE RECOVERY.

TECHNICAL SPECIFICATIONS

RPM

1800 - 2200 - 2600 - 3200

Timer

15 minutes

Movement range

6 mm

Maximum pressure force

14 kg

Control

1 on/off / speed change button

Power supply

- USB-C charging
- Lithium battery
- Charging time: 4 hours
- Operating time: 5 continuous hours at speed 1

Dimensions

92x142x47 mm

Weight

500 g

STANDARD EQUIPMENT

· Code **G6517**

· 1 Bag

· 1 4 VIBE MASSAGE GUN unit

· 1 Round head


· 1 Flat head

· 1 Trigger Point head

· 1 Spine head

· 1 USB-C power supply unit

· 1 Operating manual

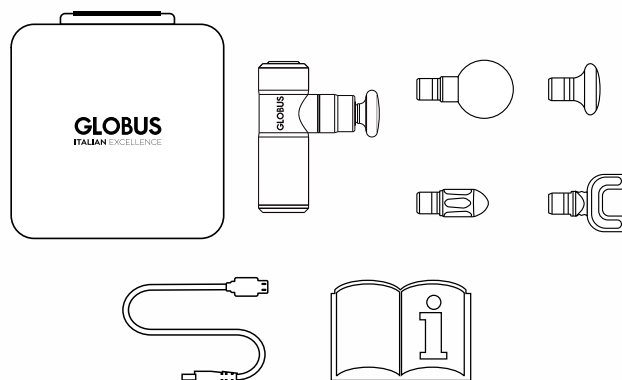



VIBRATION PLATFORM

PHYSIOPLATE MY GOLD

VIBRATION PLATFORM

PHYSIOPLATE MY GOLD IS THE IDEAL VIBRATION PLATFORM FOR ATHLETES AND FITNESS ENTHUSIASTS. IT HAS MORE THAN 100 PROTOCOLS DIVIDED INTO LEVELS AND WITH SPECIFIC EXERCISES FOR THE VARIOUS PARTS OF THE BODY.

TECHNICAL SPECIFICATIONS

Display

Size 7" colour touchscreen display

Power

18G

Type of vibration

Vertical

Frequency range

15-60Hz, with 1Hz variation

Vibration amplitude

Low/High

Voltage

230V

Mains power supply

200-230V-AC, 50Hz

Total overall dimensions

930x1550x760 mm

Plate dimensions

930x600 mm

Max. weight allowed

120 kg

Weight

170 kg

STANDARD EQUIPMENT

· Code **G1061**

· 1 PHYSIOPLATE MY GOLD vibration platform

· 1 Exercise mat

· 2 Exercise belts

· 1 Power supply unit

· 1 Operating manual

CE

ACCESSORIES



G1130

PHYSIOPLATE RED MAT

G1129

PHYSIOPLATE BLUE MAT




G1090

EXERCISE STRAP



BALL SHOOTING MACHINE

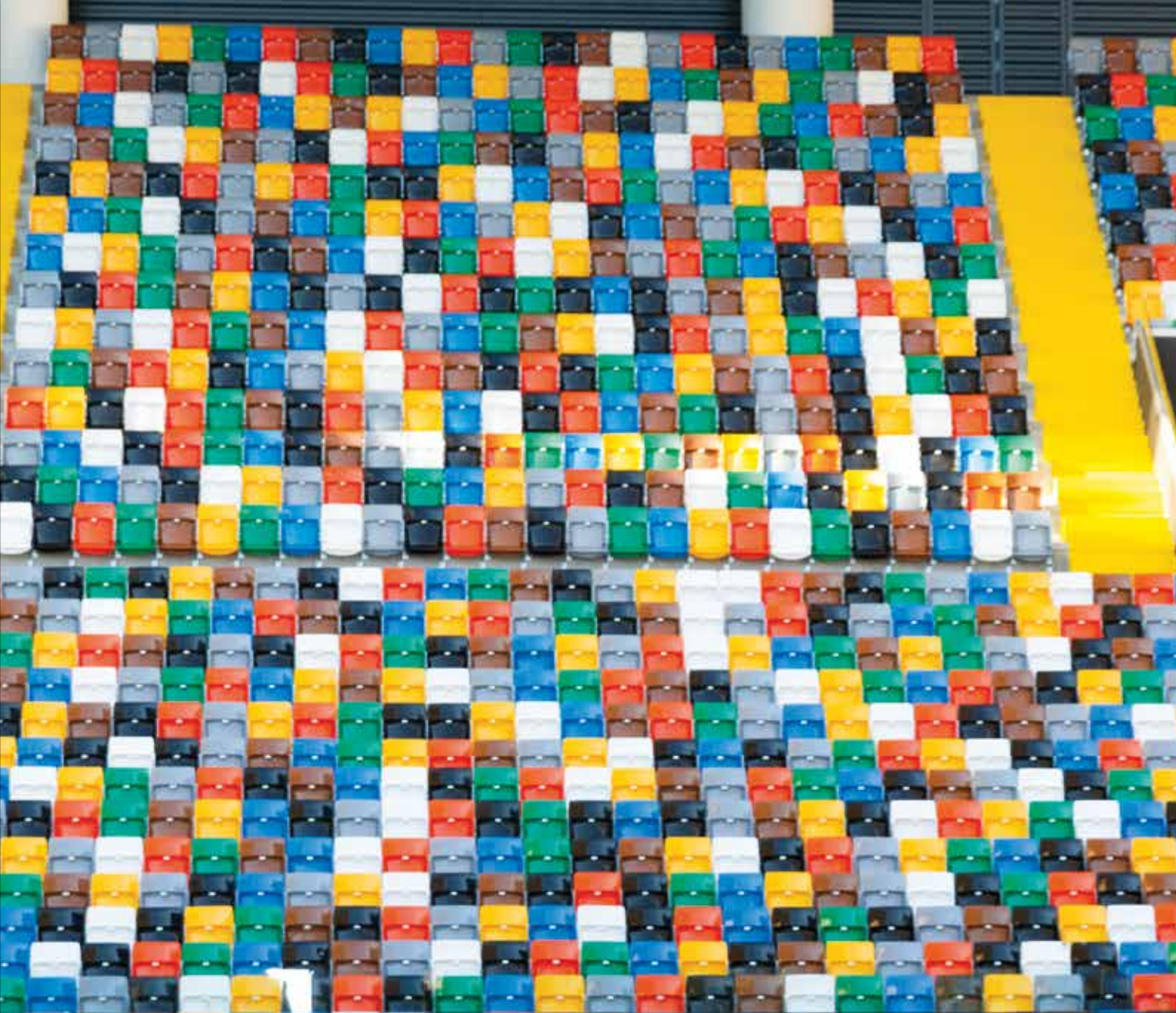

GLOBUS GLOBUS GLOBUS GLOBUS



EUROGOAL 1500

FOOTBALL BALL SHOOTING MACHINE

GLOBUS EUROGOAL 1500 IS A PROFESSIONAL BALL SHOOTING MACHINE FOR TRAINING THAT ENABLES CROSSES, THROWS AND SHOTS TYPICAL OF THE GAME OF FOOTBALL WITH EXTREME REALISM. ITS HIGH POWER ALLOWS THE BALL TO BE PROJECTED AT UP TO 140 KM/H AND IT ESPECIALLY OFFERS THE POSSIBILITY OF CREATING INFINITE TRAJECTORIES WITH ALL THE DESIRED EFFECTS. THANKS TO THE PATENTED GLOBUS TRIAXIAL SYSTEM, THE COACH CAN ADJUST THE SHOT ON THE VERTICAL PLANE (UPWARD OR DOWNWARD), ON THE HORIZONTAL PLANE (RIGHT OR LEFT) AND, ABOVE ALL, TILT THE DEVICE ON THE OBLIQUE PLANE TO ENHANCE THE CURVE OF THE TRAJECTORIES. GLOBUS EUROGOAL 1500 IS, WITHOUT A DOUBT, THE IDEAL PARTNER OF THE COACH FOR HUNDREDS OF TRAINING EXERCISES, WITH HIGH RHYTHM AND MAXIMUM SHOOTING PRECISION.

TECHNICAL SPECIFICATIONS

· Code **G3305**

Dimensions (WxDxH)

115x70x60 cm

Weight

100 kg

Launch wheels

2 in wear-resistant super grip rubber

Maximum ball shot height

35 m

Lateral movement of the head

+30°/-30°

Vertical inclination

+40°/-5°

Launching distance

100 m

Launch speed

From 10 to 140 km/h

Power supply

Double battery with high autonomy (5 effective hours or 1500 shots) with charger included

Batteries

2 12V - 55Ah

Shooting effects

Floating - Magnus - Top Spin - Back Spin

Electric actuation for vertical elevation

Oblique inclination of the head

Lateral shot range lock

Electronic control unit with 5 shooting memories

Simultaneous and digital electronic adjustment of the two wheels

STRONG POINTS

· Triaxial trajectory adjustment - GLOBUS patent

· Does not wear out the balls on the launch wheels - GLOBUS patent

· Powerful kicks and spin shots, both in-swinging and out-swinging, from the right and from the left

· Digital display

· With 5 buttons to store 5 types of shots with different speed and spin effects

· Ball speed control LED system indicating when both wheels have reached the set speed (Green LED)

· TOP SPIN: super trajectory curvature thanks to the oblique movement of the triaxial system

· Possibility of using balls of various sizes

· Easy to move on the wheels

· Ball insertion guide to ensure maximum launch repeatability

· Recommended for professional and semi-professional teams



EUROGOAL 600

FOOTBALL BALL SHOOTING MACHINE

THE GLOBUS EUROGOAL 600 IS A PROFESSIONAL TRAINING MACHINE. IT ENABLES CROSSES, THROWS AND SHOTS TYPICAL OF THE GAME OF FOOTBALL TO BE MADE WITH EXTREME PRECISION AND REALISM, WITH A POWER OF UP TO 110 KM/H. THANKS TO A SPECIAL DEVICE, THE COACH CAN ADJUST THE SHOT ON THE VERTICAL PLANE (UPWARDS OR DOWNWARDS), TO MAKE GROUND SHOTS, HALF-HEIGHT SHOTS OR TO SHOOT BETWEEN THE GOAL POSTS.

GLOBUS EUROGOAL 600 IS A TECHNICAL TRAINING TOOL THAT ENABLES THE PERFORMANCE OF MANY EXERCISES FOR GOALKEEPERS AND PLAYERS. IT CAN BE USED IN GLOBAL GAME SITUATIONS, SUCH AS CROSSES AND CORNER KICKS.

TECHNICAL SPECIFICATIONS

· Code **G1503**

Dimensions (WxDxH)

95x70x60 cm

Weight

90 kg

Launch wheels

2 in wear-resistant super grip rubber

Maximum ball shot height

25 m

Lateral movement of the head

+30°/-30°

Vertical inclination

+40°/-5°

Launching distance

60 m

Launch speed

10 to 120 km/h

Power supply

Double battery with high autonomy (5 effective hours or 1500 shots) with charger included

Batteries

2 12V - 55Ah

Shooting effects

Floating - Magnus

Electric actuation for vertical elevation

Electronic control unit

Manual adjustment by means of two independent potentiometers

STRONG POINTS

· Does not wear out the balls on the launch wheels - GLOBUS patent

· Powerful kicks and spin shots, both in-swinging and out-swinging, from the right and from the left

· Independent, micrometric wheel speed adjustment (999 levels)

· Possibility of using balls of various sizes

· Easy to move on the wheels

· Ball insertion guide to ensure maximum launch repeatability

· Also recommended for youth sections and football schools



WINSHOT 1500

VOLLEYBALL SHOOTING MACHINE

WINSHOT 1500 IS THE TOP OF THE RANGE GLOBUS BALL SHOOTING MACHINE LINE. IT IS PARTICULARLY SUITABLE FOR MEN'S AND WOMEN'S PROFESSIONAL TEAMS AS IT ENSURES HUNDREDS OF SERVES, SPIKES AND SETS WITH THE UTMOST PRECISION AND WITH A PAUSE BETWEEN SHOTS OF JUST 1.5 SECONDS. THIS MODEL CAN STORE 5 TYPES OF DIFFERENT SHOTS, BOTH IN TERMS OF POWER AND EFFECT, THANKS TO A SPECIAL ELECTRONIC BOARD. WINSHOT 1500 ALSO HAS THE "SPEED CONTROL SYSTEM", WHICH, BY MEANS OF TWO RED AND GREEN LEDS POSITIONED IN THE CENTRE OF THE CONSOLE, MAKES IT POSSIBLE TO CHECK THE ACTUAL SPEED AT WHICH THE BALL IS SHOT.

TECHNICAL SPECIFICATIONS

· Code G3504

Dimensions

Minimum height with head off 190 cm

Minimum height with head on 260 cm

Maximum total height 370 cm

Operator platform height 150 cm

Maximum width without baskets 80 cm

Maximum width with one basket 120 cm

Maximum width with two baskets 160 cm

Machine depth

150 cm

Weight

220 kg

Shooting length

230 cm to 340 cm

Launch speed

10 to 140 km/h

Power supply

2 rechargeable 12V - 55Ah batteries with charger included

STRONG POINTS

- Does not wear out the balls on the launch wheels - GLOBUS patent
- Digital display with simultaneous adjustment of the two pull wheels
- Speed Control System that indicates with a green LED when the ball speed reaches the set value
- Five different types of shots can be saved according to power and effect
- Articulated, lightweight and adjustable shooting head for variable trajectories
- Easy to move around the field in all areas for team and individual work
- Quiet and easy to use with targeted ball engagement for maximum shooting precision



WINSHOT 800

VOLLEYBALL SHOOTING MACHINE

WINSHOT 800 IS THE GLOBUS BALL SHOOTING MACHINE PARTICULARLY SUITABLE FOR PROFESSIONAL WOMEN'S TEAMS AND MID-LEVEL MEN'S TEAMS, AS IT ALLOWS FOR HUNDREDS OF SERVES, SPIKES AND SETS WITH MAXIMUM PRECISION. THIS MACHINE CAN LAUNCH BALLS, TO SERVE AND ATTACK, AT THE IMPRESSIVE SPEED OF 130 KM PER HOUR AT A HEIGHT OF 305 CM. THIS FEATURE MAKES WINSHOT 800 VERY USEFUL AND EFFECTIVE IN TRAINING, AS SERVES AND SPIKES IN COMPETITION HAVE REACHED A BALL SPEED THAT FORCES THE RECEPTION AND THE DEFENCE TO TRAIN WITH GREAT FREQUENCY AND ACCURACY, WITH SERVES HARDLY PERFORMED BY THE HUMAN ARM, SAFEGUARDING THE SHOULDERS OF ATHLETES AND COACHES.

TECHNICAL SPECIFICATIONS

· Code **G1334**

Dimensions

Minimum height with head off 175 cm

Minimum height with head on 255 cm

Maximum total height 350 cm

Operator platform height 110 cm

Maximum width without baskets 80 cm

Maximum width with one basket 120 cm

Maximum width with two baskets 160 cm

Machine depth

130 cm

Weight

- Without Battery kit 150 kg

- With Battery kit 190 kg

Shooting length

200 cm to 800 cm

Launch speed

10 to 130 km/h

Power supply

- Mains, with power supply included

2 rechargeable 12V - 55Ah batteries with charger included (not included)

STRONG POINTS

· Does not wear out the balls on the launch wheels - GLOBUS patent

· Separately regulated firing wheels with 999 speed levels

· Articulated, lightweight and adjustable shooting head for variable trajectories

· Mains power supply as standard or rechargeable batteries with a 10-hour operating time (optional)

· Easy to move around the field in all areas for team and individual work

· Quiet and easy to use with targeted ball engagement for maximum shooting precision



WINSHOT 500

VOLLEYBALL SHOOTING MACHINE

WINSHOT 500 IS THE GLOBUS BALL SHOOTING MACHINE PARTICULARLY SUITABLE FOR MEN'S AND WOMEN'S YOUTH TEAMS AND MID-LEVEL MEN'S TEAMS, AS IT ALLOWS FOR HUNDREDS OF SERVES, SPIKES AND SETS WITH MAXIMUM PRECISION. THIS MACHINE CAN LAUNCH BALLS, TO SERVE AND ATTACK, AT THE IMPRESSIVE SPEED OF 120 KM PER HOUR AT A MAXIMUM HEIGHT OF 270 CM. WINSHOT 500 IS EASY TO DISASSEMBLE AND MOVE TO DIFFERENT PLACES. IT IS VERY FLEXIBLE AND LIGHT IN THE GYM TO CARRY OUT SPECIFIC TECHNICAL SESSIONS FOR THE RECEPTION AND DEFENCE, BOTH INDIVIDUALLY AND IN A TEAM. THESE FEATURES MAKE WINSHOT 500 VERY USEFUL AND EFFECTIVE DURING TRAINING TO IMPROVE THE TECHNIQUE OF YOUNG VOLLEYBALLERS, THROUGH A SIGNIFICANT NUMBER OF HIGH-QUALITY REPETITIONS.

TECHNICAL SPECIFICATIONS

· Code **G1333**

Dimensions

Minimum height with head off 175 cm

Minimum height with head on 240 cm

Maximum total height 300 cm

Maximum width 80 cm

Machine depth

100 cm

Weight

- Without Battery kit 70 kg

- With Battery kit 110 kg

Shooting length

210 cm to 270 cm

Launch speed

10 to 120 km/h

Power supply

- Mains, with power supply included

2 rechargeable 12V - 55Ah batteries with charger included (not included)

STRONG POINTS

· Does not wear out the balls on the launch wheels - GLOBUS patent

· Separately regulated firing wheels with 999 speed levels

· Articulated, lightweight and adjustable shooting head for variable trajectories

· Mains power supply as standard or rechargeable batteries with a 10-hour operating time (optional)

· Easy to move around the field in all areas for team and individual work

· Quiet and easy to use with targeted ball engagement for maximum shooting precision

· Easy and quick to disassemble for transport by car to various training locations



RUGBYKICK

RUGBY BALL SHOOTING MACHINE

RUGBYKICK IS THE TOP OF THE RANGE OF GLOBUS TRAINING MACHINES. IT IS PARTICULARLY SUITABLE FOR PROFESSIONAL TEAMS, AS IT IS ABLE TO REPEAT HUNDREDS OF SERVICES, ATTACKING AND SETTING ACTIONS WITH MAXIMUM PRECISION, AT INTERVALS OF ONLY 1.5 SECONDS. THIS MODEL CAN STORE 5 TYPES OF DIFFERENT SHOTS, WITH THEIR POWER AND EFFECT PARAMETERS, THANKS TO A SPECIAL ELECTRONIC BOARD. RUGBYKICK'S SPEED CONTROL SYSTEM ALLOWS THE SPEED OF THE BALL TO BE CONTROLLED BY TWO LED LIGHTS (RED AND GREEN) LOCATED IN THE CENTRE OF THE CONSOLE.

TECHNICAL SPECIFICATIONS

· Code **G5371**

Dimensions

Maximum height 145 cm

Minimum height 130 cm

Width 70 cm

Length 100 cm

Weight

95 kg

Shooting length

230 cm to 80 cm

Launch speed

10 to 120 km/h

Power supply

2 rechargeable batteries 12V - 55Ah

STRONG POINTS

- Does not wear out the balls on the launch wheels - GLOBUS patent
- Digital display with simultaneous adjustment of the two pull wheels
- Speed Control System that indicates with a green LED when the ball speed reaches the set value
- Five different types of shots can be saved according to power and effect
- Articulated, lightweight and adjustable shooting head for variable trajectories
- Easy to move around the field in all areas for team and individual work
- Quiet and easy to use with targeted ball engagement for maximum shooting precision



ACCESSORIES

CONTENTS

- TECAR THERAPY AND TECAR BEAUTY	136
- MAGNETOTHERAPY	144
- ELECTROTHERAPY	147
- ELECTROSTIMULATION	147
- ULTRASOUND THERAPY	151
- LASER THERAPY	152
- RADIOFREQUENCY	153
- PRESSOTHERAPY	154
- BALL SHOOTING MACHINE	156

ACCESSORIES FOR DIACARE 7000 - TECAR BEAUTY 7000 MED TECAR BEAUTY 7000 - TECAR BEAUTY ELITE



PH 01

PH 01. BIPOLAR HANDPIECE FOR RESISTIVE TREATMENTS

G6650	BIPOLAR HANDPIECE Ø 50 mm
G6651	BIPOLAR HANDPIECE Ø 30 mm



PH 02

PH 02. MONOPOLAR RESISTIVE HANDPIECE AND ELECTRODES

G6652	MONOPOLAR RESISTIVE HANDPIECE
G6667	RESISTIVE ELECTRODE, Ø 70 mm HEAD
G6668	RESISTIVE ELECTRODE, Ø 50 mm HEAD
G6669	RESISTIVE ELECTRODE, Ø 30 mm HEAD



PH 03

PH 03. MONOPOLAR CAPACITIVE HANDPIECE AND ELECTRODES

G6653	MONOPOLAR CAPACITIVE HANDPIECE
G6670	CAPACITIVE ELECTRODE, Ø 70 mm HEAD
G6671	CAPACITIVE ELECTRODE, Ø 50 mm HEAD
G6672	CAPACITIVE ELECTRODE, Ø 30 mm HEAD



PH 04

PH 04. ACTIVE NEUTRAL HANDPIECE

G5133	RESISTIVE HANDPIECE TO USE INSTEAD OF THE NEUTRAL PLATE
-------	---



PH 05 (above) / PH 06 (below)

PH 05. SET OF NEUTRAL CABLES

G6656	NEUTRAL "A" CABLE, 4 mm, BLACK
G6697	ADAPTER FOR NEUTRAL CABLES, BLACK
G6692	NEUTRAL "C" BLACK CABLE WITH CLIP

PH 06. SET OF RESISTIVE CABLES

G6655	"A" - RESISTIVE CABLE, 4 mm, BLUE
G6698	ADAPTER FOR RESISTIVE CABLES, BLUE
G6693	"C" - RESISTIVE BLUE CABLE WITH CLIP

ACCESSORIES FOR DIACARE 5000 AND 5000 RE TECAR BEAUTY 6000 MED AND 6000 MED RE



PH 01

PH 01. BIPOLAR HANDPIECE FOR RESISTIVE TREATMENTS

G6658	BIPOLAR HANDPIECE Ø 50 mm
G6659	BIPOLAR HANDPIECE Ø 30 mm



PH 02

PH 02. MONOPOLAR RESISTIVE HANDPIECE AND ELECTRODES

G6660	MONOPOLAR RESISTIVE HANDPIECE
G6667	RESISTIVE ELECTRODE, Ø 70 mm HEAD
G6668	RESISTIVE ELECTRODE, Ø 50 mm HEAD
G6669	RESISTIVE ELECTRODE, Ø 30 mm HEAD



PH 03

PH 03. MONOPOLAR CAPACITIVE HANDPIECE AND ELECTRODES


G6661	MONOPOLAR CAPACITIVE HANDPIECE
G6670	CAPACITIVE ELECTRODE, Ø 70 mm HEAD
G6671	CAPACITIVE ELECTRODE, Ø 50 mm HEAD
G6672	CAPACITIVE ELECTRODE, Ø 30 mm HEAD



PH 04

PH 04. ACTIVE NEUTRAL HANDPIECE

G5035	RESISTIVE HANDPIECE TO USE INSTEAD OF THE NEUTRAL PLATE
-------	---

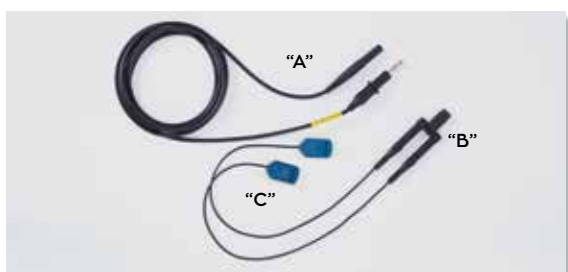

PH 05

PH 05. SET OF NEUTRAL CABLES

G6664	"A" - NEUTRAL CABLE, 4 mm, BLACK
G6697	ADAPTER FOR NEUTRAL CABLES, BLACK
G6692	"C" - NEUTRAL BLACK CABLE WITH CLIP



PH 06

PH 05. SET OF RESISTIVE CABLES

G6666	"A" - RESISTIVE CABLE, 4 mm, BLUE
G6698	ADAPTER FOR RESISTIVE CABLES, BLUE
G6693	"C" - RESISTIVE BLUE CABLE WITH CLIP

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



PH 01

PH 02



PH 03



PH 04



PH 05

PH 01. CONVEX RESISTIVE ELECTRODES

G6366	KIT WITH 3 CONVEX RESISTIVE ELECTRODES
G6675	CONVEX RESISTIVE ELECTRODE, Ø 70 mm
G6674	CONVEX RESISTIVE ELECTRODE, Ø 50 mm
G6673	CONVEX RESISTIVE ELECTRODE, Ø 30 mm

PH 02. CONVEX CAPACITIVE ELECTRODES

G6365	KIT WITH 3 CONVEX CAPACITIVE ELECTRODES
G6678	CONVEX CAPACITIVE ELECTRODE, Ø 70 mm
G6677	CONVEX CAPACITIVE ELECTRODE, Ø 50 mm
G6676	CONVEX CAPACITIVE ELECTRODE, Ø 30 mm

PH 03. CONCAVE CAPACITIVE ELECTRODES

G6540	KIT WITH 5 CONCAVE CAPACITIVE ELECTRODES
G6686	CONCAVE CAPACITIVE ELECTRODE, 30x100 mm
G6685	CONCAVE CAPACITIVE ELECTRODE, 30x80 mm
G6684	CONCAVE CAPACITIVE ELECTRODE, 30x60 mm
G6683	CONCAVE CAPACITIVE ELECTRODE, 30x40 mm
G6687	CONCAVE CAPACITIVE ELECTRODE, Ø 45 mm

PH 04. HANDS-FREE PLUS KIT FOR "DIACARE 7000" - "TECAR BEAUTY 7000" "TECAR BEAUTY 7000 MED" - "TECAR BEAUTY ELITE"

G5976	1 BAG 1 RESISTIVE CABLE 2 Y CABLES (1 BLUE, 1 BLACK) 4 ALLIGATOR CLIPS (2 BLUE, 2 BLACK) 4 MEDIUM ADHESIVE ELECTRODES 2 ELASTIC BANDS
-------	--

PH 05. FREE-HANDS PLUS KIT FOR "DIACARE 5000" AND "TECAR BEAUTY 6000 MED"

G5981	1 BAG 1 RESISTIVE CABLE 2 Y CABLES (1 BLUE, 1 BLACK) 4 ALLIGATOR CLIPS (2 BLUE, 2 BLACK) 4 MEDIUM ADHESIVE ELECTRODES 2 ELASTIC BANDS
-------	--

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



PH 06



PH 07

PH 06. FLEXIBLE NEUTRAL PLATES FOR ALL THE DEVICES IN THE "DIACARE" AND "TECAR BEAUTY" LINE

G6689	FLEXIBLE STEEL LARGE PLATE, 200X250 mm
-------	--

G6690	FLEXIBLE STEEL MEDIUM PLATE, 150X210 mm
-------	---

PH 07. FLEXIBLE NEUTRAL PLATES FOR "HANDS-FREE" TREATMENT FOR ALL THE DEVICES IN THE "DIACARE" AND "TECAR BEAUTY" LINE

G6691	FLEXIBLE STEEL SMALL PLATE, 75X130 mm, 1 PC
-------	---

G5983	KIT WITH FLEXIBLE STEEL SMALL PLATES, 75X130 mm, 2 PCS
-------	--



PH 08

PH 08. ADHESIVE "HANDS-FREE" TREATMENT ELECTRODES FOR ALL THE DEVICES IN THE "DIACARE" AND "TECAR BEAUTY" LINE

G6694	ADHESIVE TECAR ELECTRODE, LARGE, 101X202 mm, 1 PC
-------	---

G6695	ADHESIVE TECAR ELECTRODE, MEDIUM, 90X148 mm, 1 PC
-------	---

G6696	ADHESIVE TECAR ELECTRODE, SMALL, 74X88 mm, 1 PC
-------	---



PH 09

PH 09. RESISTIVE PLATE COVER KIT

G6607	RESISTIVE PLATE COVER KIT
-------	---------------------------

G6681	RESISTIVE PLATE COVER, Ø 70 mm
-------	--------------------------------

G6680	RESISTIVE PLATE COVER, Ø 50 mm
-------	--------------------------------

G6679	RESISTIVE PLATE COVER, Ø 30 mm
-------	--------------------------------



PH 10

PH 10. "DYNAMIC HANDLE", ACTIVE NEUTRAL HANDLE

G6588	DYNAMIC HANDLE
-------	----------------



PH 11



PH 12

PH 11. NEUTRAL PLATE CABLE FOR "DIACARE 7000" - "TECAR BEAUTY 7000" "TECAR BEAUTY 7000 MED" AND "TECAR BEAUTY ELITE"

G6656	NEUTRAL PLATE CABLE, BLACK
-------	----------------------------

PH 12. PLATE CABLE FOR "DIACARE 5000" AND "TECAR BEAUTY 6000 MED"

G6664	PLATE CABLE, BLACK
-------	--------------------



PH 13



PH 14

PH 13. POWER SUPPLY CABLE FOR "DIACARE 7000" - "TECAR BEAUTY 7000" "TECAR BEAUTY 7000 MED" AND "TECAR BEAUTY ELITE"

G4349	POWER SUPPLY CABLE, 230V- 50/60Hz, SCHUKO
-------	---

PH 14. POWER SUPPLY FOR THE FOLLOWING DEVICES: "DIACARE 5000" AND "TECAR BEAUTY 6000 MED"

G6663	POWER SUPPLY 48V
-------	------------------

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



PH 15

PH 15. "INTIMA" KIT

G6421	"INTIMA" KIT FOR "DIACARE 7000" - "TECARBEAUTY 7000" "TECARBEAUTY 7000 MED" - "TECARBEAUTY ELITE"
G6420	"INTIMA" KIT FOR "DIACARE 5000" - "TECARBEAUTY 6000 MED"
G6657	"INTIMA" KIT HANDPIECE, 7000 LINE
G6407	"INTIMA" KIT HANDPIECE, 5000 AND 6000 LINE
G6699	"INTIMA" KIT PROBE, Ø 20 mm
G6700	"INTIMA" KIT PROBE, Ø 9.5 mm
G6418	CONDUCTIVE GEL FOR "INTIMA" KIT
G6701	PROBE COVER FOR "INTIMA" KIT (144 PCS)



PH 16

PH 16. "FASCIA TOOLS" KIT, 4 PCS, FOR "DIACARE 7000" AND "DIACARE 5000"

G5651	1 LARGE DIAPHRAGM ELECTRODE 1 SMALL DIAPHRAGM ELECTRODE 1 DIA TRIGGER ELECTRODE 1 DIA SPINE ELECTRODE
G6703	LARGE DIAPHRAGM ELECTRODE
G6702	SMALL DIAPHRAGM ELECTRODE
G5655	DIA TRIGGER ELECTRODE
G5656	DIA SPINE ELECTRODE



PH 17

PH 17. "FASCIA TOOLS" BEAUTY KIT, 2 PCS, FOR "TECARBEAUTY 7000 MED" "TECARBEAUTY 7000", "TECARBEAUTY ELITE", "TECARBEAUTY 6000 MED"

G6041	1 LARGE DIAPHRAGM ELECTRODE 1 SMALL DIAPHRAGM ELECTRODE
G6703	LARGE DIAPHRAGM ELECTRODE
G6702	SMALL DIAPHRAGM ELECTRODE

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



G6488 PLEXIGLASS DEVICE AND ACCESSORY HOLDER ONLY FOR “DIACARE 5000” AND “TECAR BEAUTY 6000 MED”

G6489 PLEXIGLASS ACCESSORY HOLDER FOR “DIACARE 5000” AND “TECAR BEAUTY 6000 MED”



G6706 SINGLE PLEXIGLASS STAND FOR TECAR LINE ACCESSORIES

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



G3371 | GREY ALUMINIUM CASE FOR "DIACARE 5000" AND "TECAR BEAUTY 6000 MED"

G5591 | CART WITH WHEELS



G6510 | CART WITH "EASY GO" WHEELS

MORE ACCESSORIES FOR THE DIACARE AND TECAR BEAUTY LINE



G5521 | TROLLEY WITH MULTIPLE SPACES TO TRANSPORT THE DEVICES

G5520 | TROLLEY FOR DIACARE "7000" - "TECAR BEAUTY 7000" - "TECAR BEAUTY 7000 MED" AND "TECAR BEAUTY ELITE"



CONDUCTIVE TECAR CREAM FOR "DIACARE 7000" - "DIACARE 5000"

G3368 | CONDUCTIVE TECAR CREAM, 250 ml

G3369 | CONDUCTIVE TECAR CREAM, 1000 ml

G5396 | CONDUCTIVE TECAR CREAM WITH ARNICA, 1000 ml



CONDUCTIVE TECAR BEAUTY CREAM FOR "TECAR BEAUTY 7000" - "TECAR BEAUTY 7000 MED" - "TECAR BEAUTY ELITE" - "TECAR BEAUTY 6000 MED"

G6333 | CONDUCTIVE TECAR BEAUTY CREAM, 250 ml

G6332 | CONDUCTIVE TECAR BEAUTY CREAM, 1000 ml

G6226 | TECAR BEAUTY CREAM WITH HYALURONIC ACID, 500 ml

G6227 | TECAR BEAUTY CREAM FOR LYMPHATIC DRAINAGE, 1000 ml

MAGNETOTHERAPY ACCESSORIES



PH 01

PH 01. "SOFT 4" SOLENOID

G6631 | "SOFT 4" MAT, 250X250 mm, 4 COILS, 2 TOTAL CHANNELS



PH 02

PH 02. FLEXIBLE SOLENOID

G6629 | FLEXIBLE SOLENOID, 400X150 mm, 4 COILS



PH 03

PH 03. FLEXIBLE SOLENOID

G6627 | FLEXIBLE SOLENOID, 300X100 mm, 2 COILS



PH 04

PH 04. "POCKET PRO" SOLENOID

G6630 | "POCKET PRO" SOLENOID, 155x95 mm, 2 COILS.
THE CODE REFERS TO 1 SOLENOID BUT THE PURCHASE IS
MADE IN PAIRS

G6050 | "POCKET PRO" CASE



PH 05



PH 06

PH 05. MAGNUM TEST

G6713 | INDICATOR OF CORRECT EMISSION OF THE MAGNETIC FIELD,
85X55 mm

PH 06. VELCRO ELASTIC BAND

G6708 | VELCRO ELASTIC BAND, 80X800 mm, BLACK



PH 07

PH 07. SOLENOID COVER

G5878 | FLEXIBLE SOLENOID COVER KIT (10 PIECES)

MAGNETOTHERAPY ACCESSORIES



PH 08



PH 09

PH 08. POWER SUPPLY


G6711 | MAGNETO LINE POWER SUPPLY
FOR "MAGNUM XL" AND "MAGNUM L"

PH 09. POWER SUPPLY

G6712 | "MAGNUM 3000 PRO" POWER SUPPLY
"MAGNUM 2 PRO DRIVE" - "MAGNUM XL PRO"
"MAGNUM 2500" - "MAGNETOCEMP ELITE"



PH 10

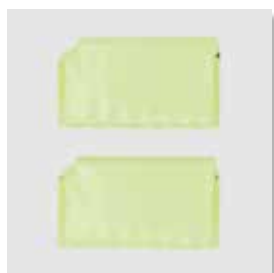

PH 11

PH 10. POWER SUPPLY

G6710 | "MAGNUM 3500 PRO" POWER SUPPLY

PH 11. BATTERY PACK

G6638 | BATTERY PACK - 12V AA NI-MH
FOR "MAGNUM 3500 PRO"



PH 12



PH 13

PH 12. BATTERY PACK

G6637 | BATTERY PACK - 7.2V NI-MH FOR
"MAGNUM 2 PRO DRIVE" - "MAGNUM XL PRO"
"MAGNUM 2500" - "MAGNETOCEMP ELITE" "MAGNUM XL" AND
"MAGNUM L"

PH 13. Y CABLE

G6709 | Y CABLE



PH 14

PH 14. EXTENSION

G6714 | EXTENSION FOR ACCESSORIES LENGTH 140 cm



PH 15

PH 15. USB ADAPTER

G6774 | USB ADAPTER WITH 5-12V CONVERTER


MAGNETOTHERAPY ACCESSORIES

MEMORY FOAM SINGLE MAT



DIMENSIONS: LENGTH 440 mm, WIDTH 700 mm

G5944	"MAT 100 MEMORY FOAM" MAT, 4 SOLENOIDS
G5945	"E-MAT 100" MAT, ACTIVE EXTENSION, MEMORY FOAM
G5946	"E-MAT 100" MAT, NEUTRAL EXTENSION, MEMORY FOAM
G6048	EXTERNAL COVER FOR "MAT 100 MEMORY FOAM"
G6049	INTERNAL COVER FOR "MAT 100 MEMORY FOAM"

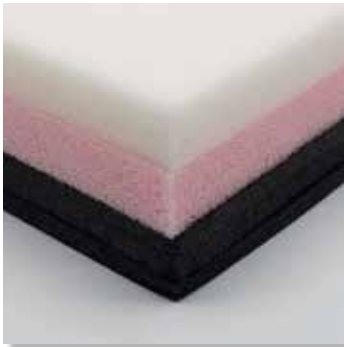

INTERNAL SECTION THAT HIGHLIGHTS THE DIFFERENT LAYERS OF THE MAT AND MATTRESS

"TOTAL BODY" MEMORY FOAM MATTRESS



DIMENSIONS: LENGTH 1760 mm, WIDTH 700 mm, TO CONNECT TO THE "GLOBUS" MAGNETHOTHERAPY DEVICES

G5947	"TOTAL BODY 200 MEMORY FOAM" MATTRESS, 2 ACTIVE MATS OUT OF 4, 8 TOTAL SOLENOIDS
G5948	"TOTAL BODY 400 MEMORY FOAM" MATTRESS, 4 ACTIVE MATS OUT OF 4, 16 TOTAL SOLENOIDS
G5949	"TOTAL BODY 400 XP MEMORY FOAM" MATTRESS, 4 ACTIVE MATS OUT OF 4, 16 TOTAL SOLENOIDS, ONLY COMPATIBLE WITH "MAGNUM 3500 PRO"

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 01

PH 01. PACK OF ELECTRODES - 4 PCS, 50X50 mm

GO136	MYOTRODE "PLATINUM"
GO464	MYOTRODE "PLUS"
GO466	MYOTRODE "PREMIUM"



PH 02

PH 02. PACK OF ELECTRODES - 4 PCS, 50X90 mm

GO137	MYOTRODE "PLATINUM"
GO465	MYOTRODE "PLUS"
GO467	MYOTRODE "PREMIUM"



PH 03

PH 03. PACK OF DOUBLE INLET ELECTRODES - 2 PCS, 50X100 mm

GO458	MYOTRODE "PLATINUM"
-------	---------------------



PH 04

PH 04. PACK OF 4 FACE ELECTRODES - Ø 32 mm WITH CLIP

GO451	4 FACE ELECTRODES
-------	-------------------



PH 05

PH 05. PACK OF 30 FACE ELECTRODES - Ø 26 mm WITH CLIP

GO871	30 FACE CLIP ELECTRODES
-------	-------------------------



PH 06



PH 07

PH 06. PACK OF CLIP SQUARE ELECTRODES - 4 PCS, 46X47 mm

GO650	4 CLIP ELECTRODES
-------	-------------------

PH 07. PACK OF CLIP RECTANGULAR ELECTRODES - 4 PCS, 45X98 mm

GO651	4 CLIP ELECTRODES
-------	-------------------

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 08


PH 08. BUCKSKIN ELECTRODES

G0452	MEDIUM BUCKSKIN ELECTRODES - 4 PCS, 80X60 mm
-------	--

G0453	LARGE BUCKSKIN ELECTRODES - 4 PCS, 80X120 mm
-------	--



PH 09

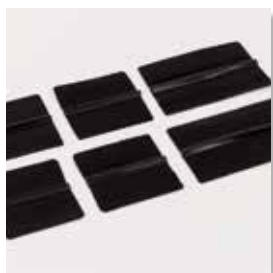

PH 10

PH 09. CONDUCTIVE SILICONE ELECTRODES

G6715	CONDUCTIVE SILICONE ELECTRODES PER METRE
-------	--

PH 10. CONDUCTIVE SILICONE ELECTRODES

G0463	KIT 6 CONDUCTIVE SILICONE ELECTRODES
-------	--------------------------------------



PH 11

PH 11. / PH 12. CONDUCTIVE ELASTIC BANDS

G0479	KIT 4 THIGH BANDS
-------	-------------------

G0480	KIT 8 "FITNESS TOP" LEG AND THIGH BANDS
-------	---

MORE ELASTIC BANDS

G1225	KIT 2 "S" ELASTIC BANDS, 80X400 mm, BLUE
-------	--

G1226	KIT 2 "M" ELASTIC BANDS, 80X800 mm, BLUE
-------	--



PH 12



PH 13



PH 14

PH 13. "EXTENSION BELT" BAND

G6723	"EXTENSION BELT" BAND, 450 mm (extension for "WINTEC SMART" and "FAST BAND")
-------	--

PH 14. IONOPHORESIS ELECTRODE

G0890	MEDIUM CARBON ELECTRODE + POUCH, 50x50 mm
-------	---

G0885	LARGE CARBON ELECTRODE + POUCH, 60x85 mm
-------	--

G6720	POUCH, 50x50 mm
-------	-----------------

G6719	POUCH, 60x85 mm
-------	-----------------



PH 15



PH 16

PH 15. VAGINAL PROBE

G6725	SINGLE-PATIENT VAGINAL PROBE
-------	------------------------------

PH 16. ANAL PROBE

G6724	SINGLE-PATIENT ANAL PROBE
-------	---------------------------

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 17

PH 17. G-TRODE HANDPIECE

G6726 | "G-TRODE" BIPOLAR HANDPIECE FOR "SHAPE & BEAUTY" AND "ACTIVA 700"



PH 18

PH 18. COLOURED CABLES FOR 4-CHANNEL DEVICES

G1839 | KIT 4 COLOURED CABLES



PH 19

PH 19. Y CABLES FOR 4-CHANNEL DEVICES

GO441 | KIT 2 Y CABLES FOR 4 CHANNELS



PH 20

PH 20. CABLES FOR 4-CHANNEL DEVICES

G6734 | GLOBUS SPARE WHITE CABLE FOR 4 CHANNELS



PH 21

PH 21. MICROCURRENT CABLES

GO440 | KIT 2 MICROCURRENT CABLES



PH 22

PH 22. CABLES FOR 2-CHANNEL DEVICES, ROUND PLUG

G6735 | "WINTEC" SPARE GREY CABLE



PH 23

PH 23. CABLE FOR 2-CHANNEL DEVICES, RECTANGULAR PLUG

G6736 | "WINTEC" BLUE SPARE CABLE



PH 24

PH 24. ADAPTERS

GO727 | BAG OF 4 ADAPTERS FOR CLIP ELECTRODES, Ø 2 mm



PH 25

PH 25. CABLES FOR 2-CHANNEL DEVICES

GO439 | KIT 2 Y CABLES (universal)

ELECTROTHERAPY AND ELECTROSTIMULATION ACCESSORIES



PH 26



PH 27

PH 26. BATTERY PACK

G6727 | BATTERY PACK FOR 2-CHANNEL DEVICES

PH 27. BATTERY PACK

G6640 | 1600mA BATTERY PACK FOR "ELITE"



PH 28



PH 29

PH 28. BATTERY PACK

G6637 | 1800mA BATTERY PACK FOR 4-CHANNEL DEVICES

PH 29. BATTERY PACK

G6639 | 2000mA BATTERY PACK FOR "GENESY 3000"



PH 30



PH 31

PH 30. POWER SUPPLY

G6728 | "GENESY 1500" POWER SUPPLY

PH 31. POWER SUPPLY

G6729 | "GENESY 3000" POWER SUPPLY



PH 32



PH 33

PH 32. BATTERY CHARGER

G6732 | BATTERY CHARGER FOR 2-CHANNEL DEVICES

PH 33. BATTERY CHARGER

G6731 | BATTERY CHARGER FOR 4-CHANNEL DEVICES



PH 34



PH 35

PH 34. BATTERY CHARGER

G6730 | BATTERY CHARGER FOR "GENESY 300 PRO" AND "ELITE"

PH 35. POWER SUPPLY

G3915 | CHARGER FOR AA/AAA BATTERIES ONLY FOR "DUO PRO"



PH 36



PH 37

PH 36. MOTOR POINT PEN

G6682 | MOTOR POINT PEN

PH 37. GEL

G0484 | 260 ml GEL BOTTLE, 10 PCS

G0869 | 260 ml GEL BOTTLE, 1 PC

G4034 | 1000 ml GEL CAN

G3609 | 5000 ml GEL CAN

ULTRASOUND THERAPY ACCESSORIES



PH 01

PH 01. Ø 16 mm HANDPIECES

G6741	Ø 16 mm HEAD HANDPIECE FOR MEDISOUND 3000" - SD5 BOARD
G6742	Ø 16 mm HEAD HANDPIECE FOR MEDISOUND 3000" - SD4 BOARD



PH 02



PH 03

PH 02. Ø 42 mm HANDPIECES

G6737	Ø 42 mm HEAD HANDPIECE FOR MEDISOUND 3000" AND "LIPOZERO G150"
G6738	Ø 42 mm HANDPIECE, WHITE HEAD FOR "MEDISOUND PRO 2" AND "LIPOZERO EXCEL"
G6739	Ø 42 mm HANDPIECE FOR "LIPOZERO G39" - SD4 BOARD

PH 03. Ø 50 mm HANDPIECE

G6740	Ø 50 mm HEAD HANDPIECE, 39kHz
-------	-------------------------------



PH 04



PH 05

PH 04. POWER SUPPLY

G6745	"LIPOZERO G39" POWER SUPPLY
-------	-----------------------------

PH 05. POWER SUPPLY

G6743	"MEDISOUND 3000", "LIPOZERO G800" AND "LIPOZERO G150" POWER SUPPLY
-------	---




PH 06

PH 06. POWER SUPPLY

G6711	"LIPOZERO EXCEL" POWER SUPPLY
-------	-------------------------------



PH 07

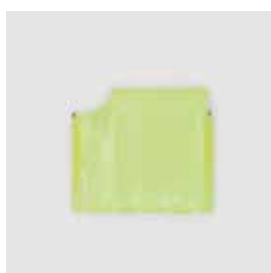

PH 08

PH 07. BATTERY PACK

G6638	BATTERY PACK FOR "MEDISOUND 3000", "LIPOZERO G39" "LIPOZERO G150" AND "LIPOZERO G800"
-------	---

PH 08. BATTERY PACK

G6637	BATTERY PACK FOR "LIPOZERO EXCEL"
-------	-----------------------------------



PH 09

PH 09. GEL

G0484	260 ml GEL BOTTLE, 10 PCS
G0869	260 ml GEL BOTTLE, 1 PC
G4034	1000 ml GEL CAN
G3609	5000 ml GEL CAN

LASERTERAPIA ACCESSORIES



PH 02

PH 01. HANDPIECES AND OPTICS

G5996	HANDPIECE 1064nm (3W)
G5993	HANDPIECE HP 808nm (12W)
G5995	HANDPIECE HP 980nm (12W)
G5994	HANDPIECE MP 980nm (6W)
G5992	HANDPIECE MP 808nm (6W)
G5796	HANDPIECE GREEN, 808nm (1W) Ø 11 mm OPTIC INCLUDED
G5748	HANDPIECE BLACK, 808nm (2W) Ø 11 mm OPTIC INCLUDED
G5767	HANDPIECE WHITE, 808nm (1W) Ø 11 mm OPTIC INCLUDED
G6090	Ø 11 mm OPTIC
G6092	Ø 22 mm OPTIC



PH 03

PH 03. SMALL OPTIC

G5547	Ø 3 mm "SMALL" OPTIC ONLY FOR PODCARE 2.0
-------	---



PH 04



PH 05

PH 04. LASER GLASSES

G6746	OPERATOR LASER GLASSES
-------	------------------------

PH 05. LASER GLASSES

G6747	PATIENT LASER GLASSES
-------	-----------------------



PH 06



PH 07

PH 06. POWER SUPPLY

G6729	LASER POWER SUPPLY FOR "PRO" VERSION
-------	--------------------------------------

PH 07. POWER SUPPLY

G6728	LASER POWER SUPPLY
-------	--------------------



PH 08



PH 09

PH 08. BATTERY PACK

G6639	2000mA BATTERY PACK FOR "PRO" VERSION
-------	---------------------------------------

PH 09. BATTERY PACK

G6637	1800mA BATTERY PACK
-------	---------------------

RADIOFREQUENCY ACCESSORIES



PH 01

PH 01. "RF CLINIC 100" HANDPIECE

G6765 | 2-SPHERE HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 02

PH 02. "RF CLINIC 300" HANDPIECE

G6764 | Ø 30 mm HEAD HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 03

PH 03. "RF CLINIC 500" HANDPIECE

G6763 | Ø 50 mm HEAD HANDPIECE
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 04

PH 04. CARRYING CASE

G1364 | ALUMINIUM CASE, BLACK
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 05

PH 05. POWER SUPPLY

G6710 | POWER SUPPLY
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"



PH 06

PH 06. BATTERY PACK

G6638 | BATTERY PACK
FOR "RF CLINIC PRO" AND "RF CLINIC BODY"

PRESSOTHERAPY ACCESSORIES



PH 01

PH 01. PRESSCARE SLEEVE

G6752	PRESSCARE SLEEVE, 4 SECTIONS
-------	------------------------------



PH 02

PH 02. PRESSCARE BOOTS

G6750	"PRESSCARE" BOOT, 4 SECTIONS, SMALL
-------	-------------------------------------

G6753	"PRESSCARE" BOOT EXTENSION, SMALL
-------	-----------------------------------

G6749	"PRESSCARE" BOOT, 4 SECTIONS, MEDIUM
-------	--------------------------------------

G6754	"PRESSCARE" BOOT EXTENSION, MEDIUM
-------	------------------------------------

G6748	"PRESSCARE" BOOT, 4 SECTIONS, LARGE
-------	-------------------------------------

G6755	"PRESSCARE" BOOT EXTENSION, LARGE
-------	-----------------------------------



PH 03

PH 03. PRESSCARE ABS/GLUTES BAND

G6751	"PRESSCARE" ABS/GLUTES BAND, 4 SECTIONS
-------	---



PH 04



PH 05

PH 04. Y CABLE

G6757	Y CABLE TO USE 3 ACCESSORIES AT THE SAME TIME
-------	--

PH 05. CONNECTORS

G6758	SINGLE CONNECTOR FOR "PRESSCARE G300M" AND "PRESSCARE G200M"
-------	---

G6759	SINGLE CONNECTOR FOR "PRESSCARE G-SPORT 3"
-------	--



PH 06



PH 06. BAG

G6762	BAG FOR "PRESSCARE G-SPORT 3"
-------	-------------------------------


PRESSOTHERAPY ACCESSORIES

SIZE OF ACCESSORIES (4 SECTIONS)

BOOTS SIZE "S"	(A) LENGTH GROIN - HEEL 72 cm (B) THIGH CIRCUMFERENCE 60 cm with EXTENSION 70 cm (B) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "M"	(A) LENGTH GROIN - HEEL 82 cm (B) THIGH CIRCUMFERENCE 67 cm with EXTENSION 77 cm (B) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "L"	(A) LENGTH GROIN - HEEL 92 cm (B) THIGH CIRCUMFERENCE 74 cm with EXTENSION 84 cm (B) ANKLE CIRCUMFERENCE 40 cm with EXTENSION 45 cm
SLEEVE	(A) 84 cm - (B) 63 cm (C) SHOULDER CIRCUMFERENCE 43 cm with EXTENSION 53 cm (D) WRIST CIRCUMFERENCE 31 cm with EXTENSION 36 cm

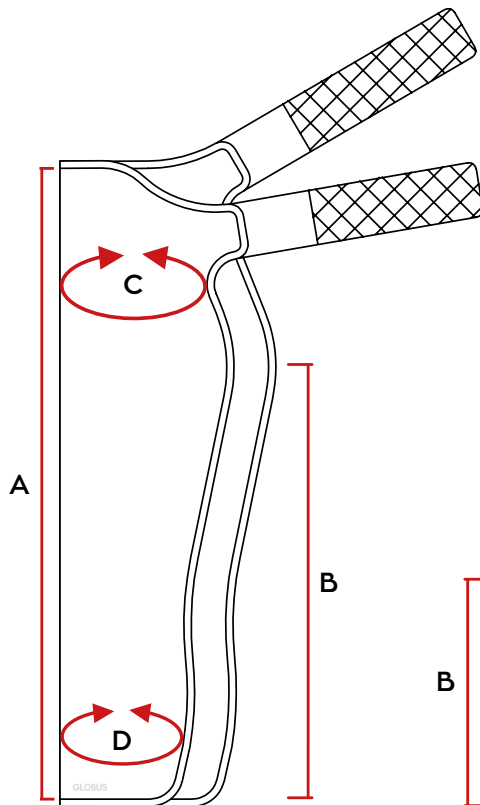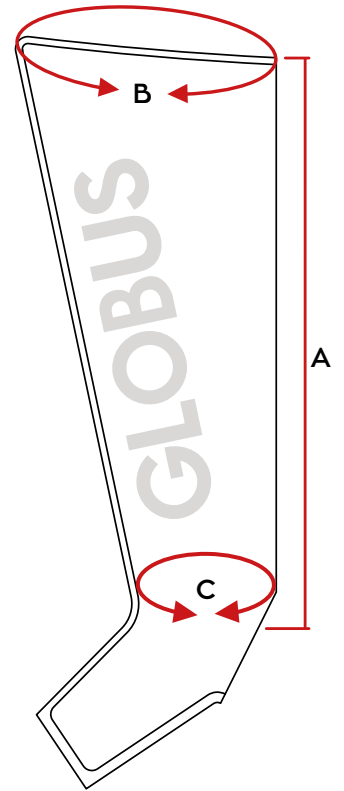

ABDOMINAL BAND	(A) 133 cm - (B) 38 cm ABDOMEN CIRCUMFERENCE 82 cm to 114 cm
----------------	---

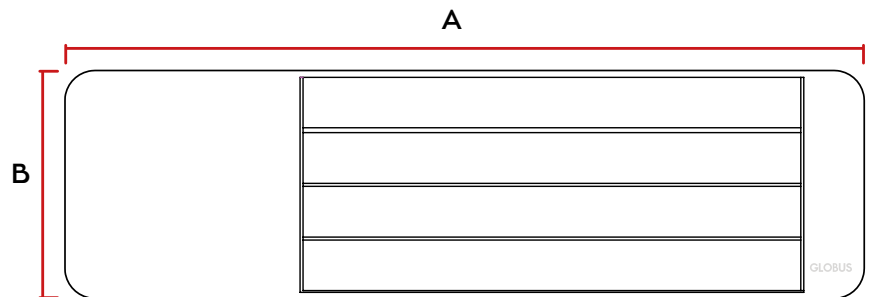

Y CABLE to use 3 accessories at the same time	SINGLE CABLE
---	--------------

NOTE: All **EXTENSIONS** are to be considered an optional accessory available on request.

The **ACCESSORIES**, depending on the type of configuration, can be optional or included.

BALL SHOOTING MACHINE ACCESSORIES



PH 01

PH 01. SPEED MEASURING GUN

G3048	SPEED MEASURING GUN
-------	---------------------



PH 02

PH 02. "EUROGOAL" LINE POWER SUPPLY

G1488	"EUROGOAL" LINE POWER SUPPLY
-------	------------------------------



PH 03

PH 03. "EUROGOAL" LINE BATTERY CHARGER

G1486	"EUROGOAL" LINE BATTERY CHARGER
-------	---------------------------------



PH 04

PH 04. BATTERY HOLDER AND BATTERY

G2946	BATTERY HOLDER
G3808	RECHARGEABLE BATTERY 12V - 55Ah



PH 05

PH 05. "WINSHOT" LINE BATTERY KIT

G6019	BATTERY KIT FOR "WINSHOT 1500" AND "WINSHOT 1000"
G2955	BATTERY KIT FOR "WINSHOT 800" AND "WINSHOT 500"

BALL SHOOTING MACHINE ACCESSORIES



PH 06. "EUROGOAL" COVER

G3063 | COVER

PH 06



PH 07. "WINSHOT 1500" FRONT PROTECTION

G5073 | KIT FRONT PROTECTION

PH 07



VET

CONTENTS

- TECAR THERAPY	160
• TECARVET 4000 TECARVET 4000 RE	160
• TECARVET 2000	161
- MAGNETOTHERAPY	162
• KIT MAGNETOVET 4000 HORSE PRO	162
• MAGNETOVET 4000	166
• MAGNETOVET 200 + THERAPEUTIC BED	167
- ELECTROTHERAPY	168
• STIMVET 4000	168
• STIMVET 2000	169
• STIMVET 200	170
- ULTRASOUND THERAPY	172
• ULTRASOUNDEVET 4000	172
- LASER THERAPY	174
• LASERVET 15.0	174
• LASERVET 8.0	175
• LASERVET 3.0	176
• LASERVET 1.0	177
- ACCESSORIES	178

TECARVET 4000 | TECARVET 4000 RE

TECAR THERAPY

TECARVET 4000 IS A DEVICE THAT USES DIATHERMY TO OBTAIN EXTRAORDINARY THERAPEUTIC EFFECTS ON SMALL ANIMALS SUCH AS DOGS AND CATS. THE DEVICE TRANSFERS ENERGY IN BOTH CAPACITIVE AND RESISTIVE MODES. IN THIS WAY THE EFFECT OF THE TREATMENT CAN BE DIRECTED WITH GOOD EFFECTIVENESS BOTH ON DEEP AND MORE SUPERFICIAL BIOLOGICAL TISSUES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 4.3" touchscreen

Power

160W, $\pm 20\%$

Emission frequencies

400kHz, 470kHz, 700kHz, 1000kHz $\pm 10\%$

Treatment methods

Capacitive and resistive

Capacitive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Resistive electrodes

3 \emptyset sizes (30, 50, 70 mm)

Mains power supply

100-240V-AC, 50-60Hz, 1,8A

Dimensions

170x220x60 mm

Weight

1.04 kg

STANDARD EQUIPMENT

· Code **G3865**

· 1 Professional carrying case

· 1 TECARVET 4000 device

· 1 Return electrode (plate)

· 1 Monopolar resistive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)


· 1 Monopolar capacitive handpiece with 3 electrodes \emptyset (30, 50, 70 mm)

· 1 Jar of cream, 250 ml

· 3 Elastic bands

· 1 Power supply unit

· 1 Operating manual

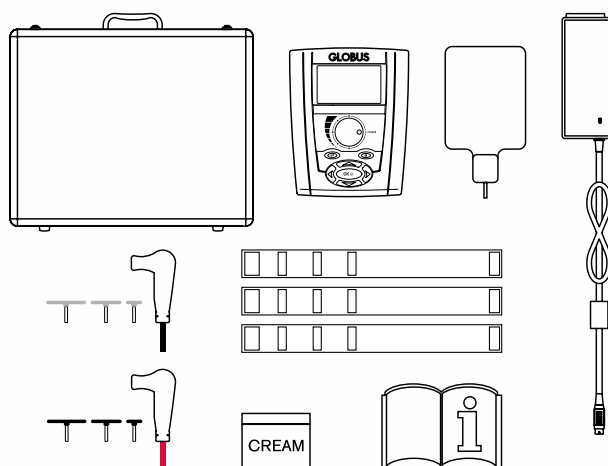

PURCHASE OPTION

· Code **G3973** TECARVET 4000 RE *

* RECHARGEABLE TECAR THERAPY :

- 5 HOURS OF TREATMENT INCLUDED IN THE COST OF THE DEVICE
- CHEAP CHARGING
- GUARANTEED SUPPORT FOR 5 YEARS AND FORKLIFT
- CHARGING VIA INTERNET



STRONG POINTS



TECARVET 2000

TECAR THERAPY

THE TECARVET 2000 MODEL PROVIDES ALL THE QUALITY OF GLOBUS' TECAR THERAPY IN A COMPACT, LIGHTWEIGHT MACHINE TO ENSURE MAXIMUM PORTABILITY, WITHOUT COMPROMISING ON EFFECTIVENESS AND PROFESSIONALISM. THE SPECIAL RESISTIVE BIPOLAR HANDPIECE ALLOWS FOR FASTER TREATMENT PREPARATION.

SPECIFIC PROGRAMS FOR DOGS AND CATS

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Power

50W, $\pm 10\%$

Emission frequencies

470kHz, $\pm 10\%$

Treatment methods

Bipolar resistive

Power supply

- Mains
- Battery pack

Dimensions

170x220x60 mm

Weight

1.04 kg

STANDARD EQUIPMENT

· Code G5617

· 1 Professional carrying case


· 1 TECARVET 2000 device

· 1 Bipolar handpiece Ø 50 mm

· 1 260 ml gel pack

· 1 Power supply unit

· 1 Operating manual

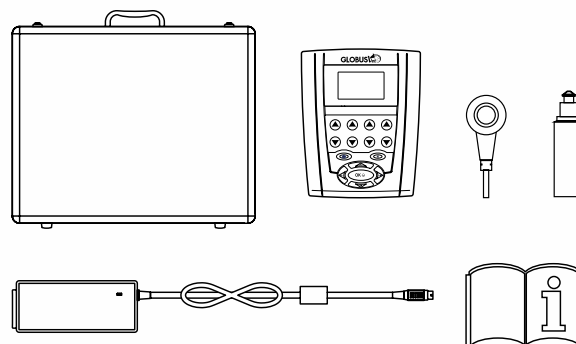

CE



STRONG POINTS




KIT MAGNETOVET 4000 HORSE PRO






MAGNETOTHERAPY FOR HORSES

THE NEW MAGNETOVET 4000 HORSE PRO KIT, IN ADDITION TO THE DEVICE FOR MAGNETOTHERAPY AND ITS NEW HORSE-DEDICATED PROGRAMS, INCLUDES ACCESSORIES DESIGNED FOR ITS TREATMENT: A THERAPEUTIC BLANKET AND NECK COVER, FLEXIBLE SOLENOIDS WITH RELATED Y CABLES, A PRACTICAL DEVICE CASE AND AN ADJUSTABLE STRAP FOR ATTACHING THE DEVICE TO THE HORSE. A CONVENIENT CARRYING BAG AND A DOUBLE BATTERY PACK COMPLETE THE KIT. MAGNETOVET 4000 HORSE IS THE NEW MAGNETOTHERAPY DEVICE WITH HIGH-INTENSITY ELECTROMAGNETIC FIELDS, SPECIFICALLY DESIGNED TO TREAT HORSES. THE FOUR OUTPUT CHANNELS ALLOW THE SOLENOIDS IN THE THERAPEUTIC BLANKET AND NECK COVER TO BE CONNECTED, SO THAT THE HORSE'S ENTIRE SPINAL COLUMN CAN BE TREATED.

THERAPEUTIC BLANKET

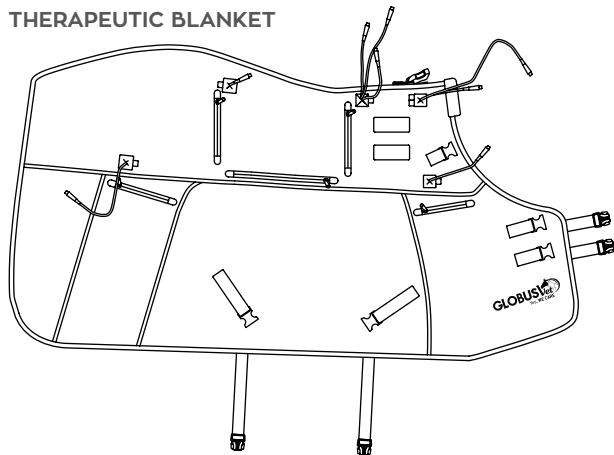

THERAPEUTIC NECK COVER

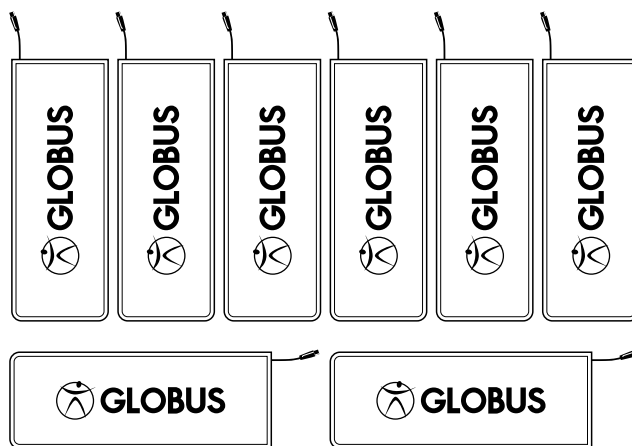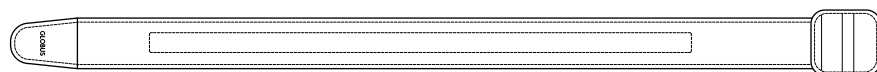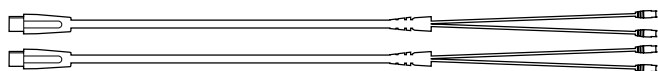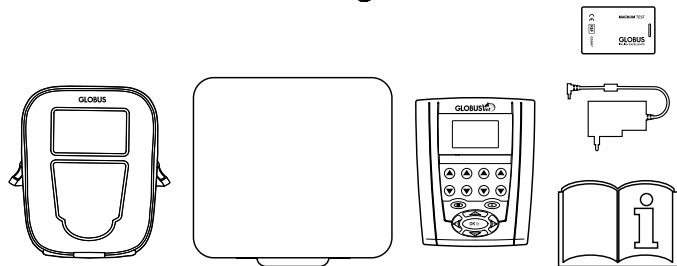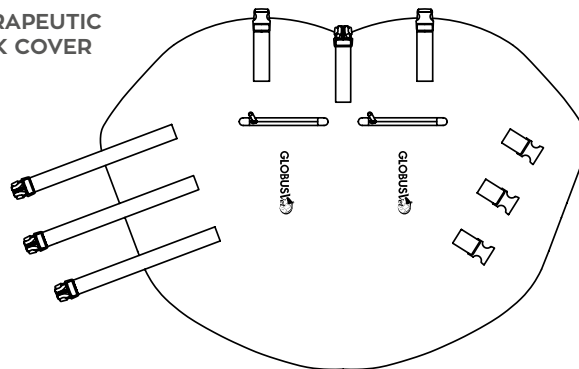

STANDARD EQUIPMENT

Code **G6593**

- 1 Therapeutic blanket
- 1 Therapeutic neck cover
- 1 Case with shoulder strap for MAGNETOVET 4000 HORSE
- 1 Anti-shock carrying bag
- 1 MAGNETOVET 4000 HORSE device
- 1 Battery pack (inside the device)
- 1 Magnum Test
- 1 Power supply unit
- 1 Operating manual
- 8 Flexible solenoids 400x150 mm for blanket and/or neck cover
- 2 Y cables
- 1 Velcro strap to fasten case
- 1 Kit carry bag



A WIDE RANGE OF ACCESSORIES COMPATIBLE WITH MAGNETOVET 4000 HORSE, IN ADDITION TO THOSE IN THE KIT, ARE AVAILABLE ON REQUEST

KIT MAGNETOVET 4000 HORSE PRO

MAGNETOTHERAPY FOR HORSES

MAGNETOVET 4000 HORSE DEVICE

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

4

Deliverable intensity

200 Gauss average per channel with flex solenoids

800 total average Gauss (for 4 channels)

The powers indicated refer to accessories with two coils. For accessories with more than two coils, the power is evenly distributed.

Frequency

Adjustable programming from 5 to 200Hz

Duty cycle

5% to 50% in steps of 5

Power supply

- Mains
- 4000mAh battery pack included

Dimensions

170x220x60 mm

Weight

1.37 kg

LIST OF HORSE PROGRAMS

THE PROGRAMS ARE DIVIDED BY THE ANATOMICAL AREAS OF THE NECK, BACK AND LIMBS

Arthrosis

Arthritis

Fracture

Delayed Fracture Consolidation

Discopathy

Wound/Ulcer

Anti-inflammatory

Injury/Contusion

Muscular Relaxation

PRE competition warm-up

POST competition recovery

Muscle injury

LIST OF DOG/CAT PROGRAMS

THE PROGRAMS INVOLVE THE USE OF THE THERAPEUTIC BED OR FLEXIBLE SOLENOIDS WITH 2 COILS, 300x100 mm

Fracture

Osteoarthritis

Discopathy

Cicatrisation

Inflammatory pain



EASY PROGRAM



2+2 FUNCTION



LAST 10



MULTI USER



PROGRAMMABLE



BATTERY PACK



MAINS POWER SUPPLY

KIT MAGNETOVET 4000 HORSE PRO

MAGNETOTHERAPY FOR HORSES

THERAPEUTIC BLANKET

TECHNICAL SPECIFICATIONS

HARD-WEARING AND WASHABLE

Dimensions

142 cm length from withers (A) to tail (B)

Weight

1.80 kg (without solenoids)

Size

Universal

Belts

5 for adjustable fastening

3 for connecting to the therapeutic neck cover

Pockets

8 for insertion of flexible solenoids with 4 coils

2 central for cable storage

Handles

3 for blanket positioning

Cables and extensions

Included

Material

96% polyester, 4% polypropylene



THERAPEUTIC NECK COVER

TECHNICAL SPECIFICATIONS

HARD-WEARING AND WASHABLE

Dimensions

121 cm length from ears (A) to withers (B)

Weight

560 g (without solenoids)

Size

Universal

Belts

3 for adjustable fastening

3 for connecting to the therapeutic blanket

Pockets

2 for insertion of flexible solenoids with 4 coils

Cables and extensions

Included

Material

94% polyester, 3% polypropylene,
3% expanded polyethylene





MAGNETOVET 4000

MAGNETOTHERAPY

MAGNETOVET 4000 IS A HIGH-INTENSITY PULSED ELECTROMAGNETIC FIELD MAGNETOTHERAPY DEVICE. THE POSSIBILITY OF USING FOUR OUTPUT CHANNELS ALLOWS SEVERAL ACCESSORIES TO BE CONNECTED, SO THAT SEVERAL BODY PARTS OR EVEN TWO ANIMALS CAN BE TREATED AT THE SAME TIME. MAGNETOTHERAPY HAS COUNTLESS AND PROVEN BENEFICIAL EFFECTS ON BONE AND ARTICULAR TISSUE. THE ELECTROMAGNETIC PULSES PRODUCED STIMULATE THE REGENERATION OF NEW CELLS AND REMODELLING OF EXISTING STRUCTURES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

4

Deliverable intensity

200 Gauss average per channel with flex solenoids

800 total average Gauss (for 4 channels)

Frequency

Adjustable programming from 5 to 200Hz

Duty cycle

5% to 50% in steps of 5

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


170x220x60 mm

Weight

820 g

STANDARD EQUIPMENT


- Code **G5628**
- 1 Bag
- 1 MAGNETOVET 4000 device
- 4 Flexible solenoids 300x100 mm
- 1 Power supply unit
- 1 Operating manual

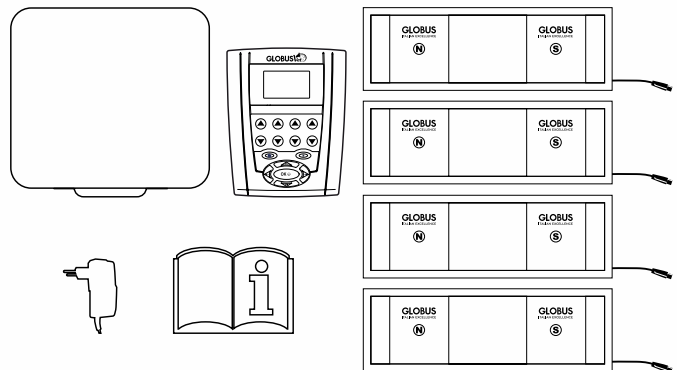

PURCHASE OPTION

- Code **G6039**

THERAPEUTIC BED IN MEMORY FOAM, 4 COILS, 70X44 cm FOR SMALL DOGS AND CATS.
FOR LARGE DOGS, BUY 2 SINGLE BEDS TO JOIN WITH THE VELCRO.

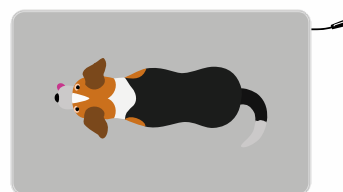

CE

STRONG POINTS



MAGNETOVET 200 WITH THERAPEUTIC BED

MAGNETOTHERAPY

MAGNETOVET 200 IS A PORTABLE MAGNETOTHERAPY DEVICE THAT CAN PRODUCE A HIGH-INTENSITY ELECTROMAGNETIC FIELD. THE PULSES PRODUCED HAVE A THERAPEUTIC AND REGENERATIVE EFFECT ON THE MAIN OSTEOARTICULAR PATHOLOGIES THAT AFFECT ANIMALS. BECAUSE OF ITS EASE OF USE IN COMBINATION WITH THE THERAPEUTIC BED, MAGNETOVET 200 IS ALSO PARTICULARLY SUITABLE FOR HOME TREATMENT OF CHRONIC PAIN.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

2

Deliverable intensity

200 Gauss average per channel with flex solenoids

Frequency

Adjustable programming from 5 to 200Hz

Duty cycle

5% to 50% in steps of 5

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions


100x160x35 mm

Weight

280 g

STANDARD EQUIPMENT

- Code G6115
- 1 Bag
- 1 MAGNETOVET 200 device
- 1 Therapeutic bed 700x440 mm
- 2 Flexible solenoids 300x100 mm
- 1 Soft fleece blanket
- 1 Power supply unit
- 1 Operating manual

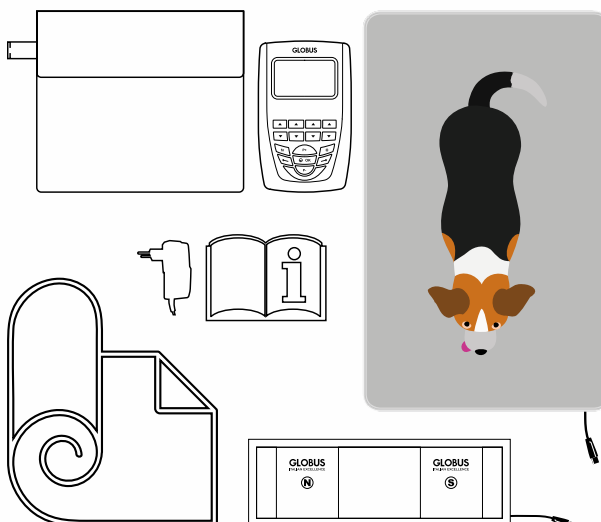

PURCHASE OPTION

- Code G5620 MAGNETOVET 200 SUPPLIED WITHOUT THERAPEUTIC BED AND WITHOUT FLEECE BLANKET
- Code G6039 THERAPEUTIC BED IN MEMORY FOAM, 4 COILS, 70X44 cm FOR SMALL DOGS AND CATS.
FOR LARGE DOGS, BUY 2 SINGLE BEDS TO JOIN WITH THE VELCRO.



STRONG POINTS



STIMVET 4000

4-CHANNEL ELECTROTHERAPY

THE STIMVET 4000 IS A PROFESSIONAL ELECTROTHERAPY DEVICE SUITABLE FOR VETERINARY REHABILITATION. THE DEVICE INCLUDES CURRENTS FOR THE RECOVERY OF MUSCULAR ATROPHY AND TENS CURRENTS FOR ANALGESIC TREATMENT OF CHRONIC AND ACUTE PAIN. IN ADDITION, THERE ARE CURRENTS TO TREAT DENERVATED MUSCLES AND MICROCURRENTS INDICATED MAINLY IN CARTILAGE LESIONS AND WOUNDS. THE DEVICE IS FULLY EQUIPPED WITH PRESET PROGRAMS FOR THE TREATMENT OF DOGS, CATS AND HORSES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150Hz

Pulse width

50-450µs

Intensity

100mA per channel

Power supply

- Mains
- Rechargeable battery

Dimensions

170x220x60 mm

Weight

1.37 kg

STANDARD EQUIPMENT

· Code G5500

· 1 Bag

· 1 STIMVET 4000 device

· 4 Elastic bands

· 2 Grey microcurrent and ionophoresis cables

· 4 Coloured cables to connect the electrodes

· 4 Square silicone electrodes


· 4 Rectangular silicone electrodes

· 4 Ionophoresis sponges

· 1 260 ml gel pack

· 1 Power supply unit


· 1 Operating manual

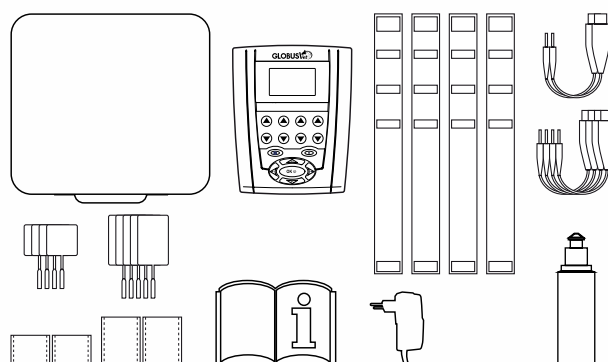

CE



STRONG POINTS

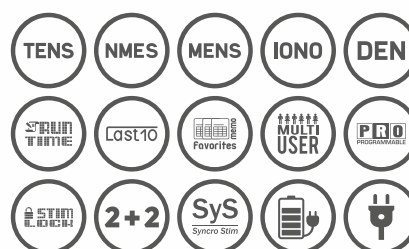

STIMVET 2000

4-CHANNEL ELECTROTHERAPY

STIMVET 2000 IS THE PORTABLE ELECTROSTIMULATOR CREATED FOR MUSCLE REHABILITATION AND ANALGESIC TREATMENT OF DOGS, CATS AND HORSES. THE DEVICE IS ALSO SUITABLE FOR SPORTING OR WORKING DOGS THAT NEED TO INCREASE THEIR ATHLETIC ABILITIES. STIMVET 2000 ALSO INCLUDES SOME MICROCURRENT PROGRAMS FOR WOUND CICATRISATION.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150Hz

Pulse width

50-450µs

Intensity

100mA per channel

Power supply

- Mains
- Rechargeable battery

Dimensions


100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT


- Code **G5614**
- 1 Bag
- 1 STIMVET 2000 device
- 4 Elastic bands
- 2 Grey microcurrent cables
- 4 Coloured cables to connect the electrodes
- 4 Square silicone electrodes
- 4 Rectangular silicone electrodes
- 1 260 ml gel pack
- 1 Power supply unit
- 1 Operating manual

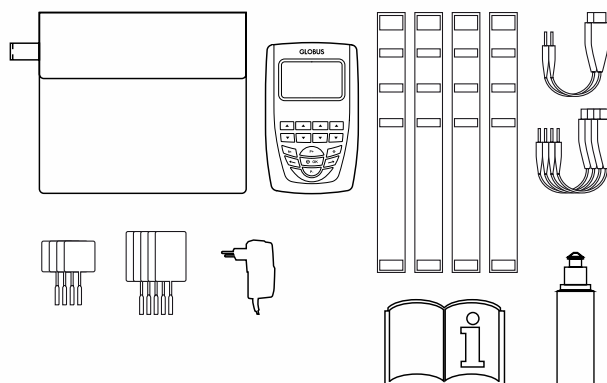

CE



STRONG POINTS

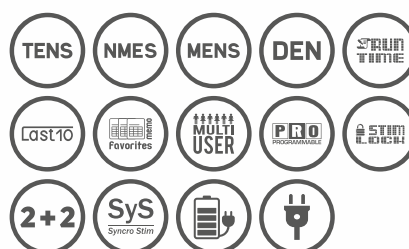

STIMVET 200

4-CHANNEL ELECTROTHERAPY

STIMVET 200 IS THE PORTABLE ELECTROSTIMULATOR CREATED FOR MUSCLE REHABILITATION AND ANALGESIC TREATMENT OF DOGS AND CATS. THE MUSCLE-STRENGTHENING PROGRAMS ARE SUITABLE BOTH FOR ANIMALS THAT NEED TO RECOVER FROM AN ATROPHIC SITUATION AND FOR SPORTING OR WORKING DOGS THAT NEED TO INCREASE THEIR ATHLETIC ABILITIES. ALTERNATIVELY, TENS CURRENTS ARE PARTICULARLY SUITABLE FOR CHRONIC AND ACUTE PAIN STATES, PROVIDING MORE COMFORTABLE RELIEF FOR THE ANIMAL AND HELPING TO REDUCE THE NEED FOR MEDICATION.

SPECIFIC PROGRAMS FOR DOGS AND CATS

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Channels

4 independent (8 electrodes)

Frequency

1-150Hz

Pulse width

50-450µs

Intensity

80mA per channel

Power supply

Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G5608

· 1 Bag

· 1 STIMVET 200 device

· 4 Elastic bands

· 2 Grey microcurrent cables

· 4 Coloured cables to connect the electrodes


· 4 Square silicone electrodes

· 4 Rectangular silicone electrodes


· 1 260 ml gel pack

· 1 Power supply unit

· 1 Operating manual

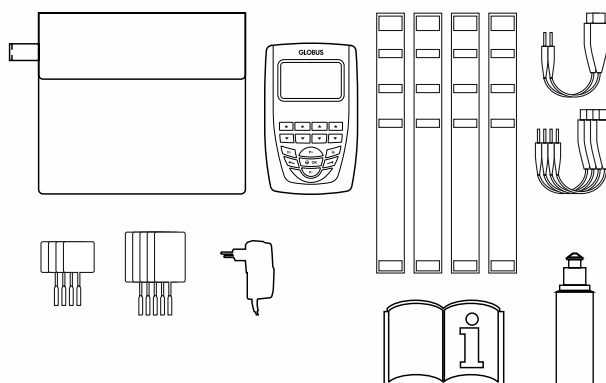

STRONG POINTS

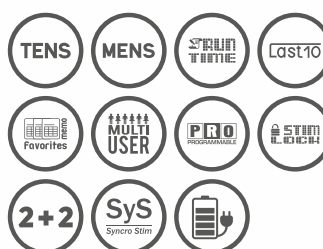



ULTRASOUNDVET 4000

ULTRASOUND THERAPY

THERAPY WITH ULTRASOUNDVET 4000 IS EFFECTIVE IN THE TREATMENT OF TENDON, MUSCLE AND LIGAMENT DISORDERS IN DOGS, CATS AND HORSES. IT IS ALSO USED TO TREAT SCAR ADHESIONS AND JOINT STIFFNESS. THIS INSTRUMENTAL PHYSICAL THERAPY PROVIDES MECHANICAL THERMAL AND NON-THERMAL STIMULATION TO PROMOTE HEALING OF DAMAGED TISSUE IN MUSCLES, TENDONS, JOINT CAPSULES, LIGAMENTS AND BONES.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Outputs

1 with Ø 42 mm handpiece

Frequency

1MHz and 3MHz

Power

3W/cm², ±20%

Duty cycle

From 10% to 100% in 10% steps

Emission

Continuous and pulsed

Mode

Direct contact and immersion

Power supply

- Mains
- Rechargeable battery (not included)

Dimensions

170x220x60 mm

Weight

820 g

STANDARD EQUIPMENT

· Code G3869

· 1 Bag


· 1 ULTRASOUNDVET 4000 device

· 1 Ø 42 mm - IP 68 handpiece

· 1 260 ml gel pack

· 1 Power supply unit

· 1 Operating manual

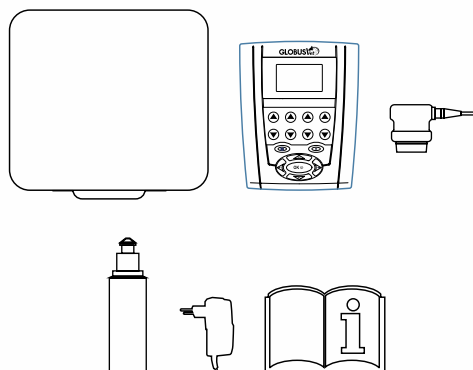

CE



STRONG POINTS





LASERVET 15.0

LASER THERAPY

LASERVET 15.0 IS A PORTABLE DEVICE FOR LASER THERAPY WITH HIGH PULSED POWER THAT ENABLES FAST VETERINARY TREATMENTS. THE PATHOLOGIES TREATED ARE DIVIDED BY ANIMAL, WITH A FOCUS ON EQUINES. THE LASER ACTION IS AIMED AT TISSUE REGENERATION, POST-TRAUMATIC AND ANTI-INFLAMMATORY TREATMENT, AND THE REDUCTION OF BOTH CHRONIC AND ACUTE PAIN.

SPECIFIC PROGRAMS FOR DOGS, CATS, HORSES, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

15W $\pm 20\%$

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

\varnothing 11 mm and \varnothing 22 mm optics

\varnothing 60 mm optic (optional)

Frequency

Up to 10000Hz

Power supply

- Mains

- Rechargeable battery

Dimensions and weight

220x170x60 mm

Weight

1.12 kg

STANDARD EQUIPMENT

· Code G6036

· 1 Bag


· 1 LASERVET 15.0 device

· 1 Handpiece HP 808nm with \varnothing 11 mm and \varnothing 22 mm optics

· 2 Pairs of protective goggles

· 1 Power supply unit

· 1 Operating manual

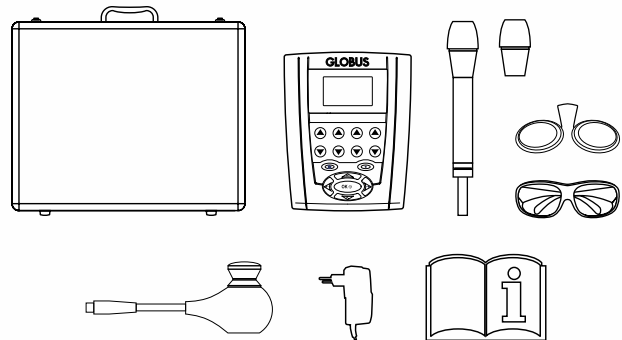

PURCHASE OPTIONS

· Code G6348 Handpiece HP 980nm

· Code G5996 Handpiece 1064nm

· Code G6089 \varnothing 60 mm optic



STRONG POINTS



LASERVET 8.0

LASER THERAPY

LASERVET 8.0 IS A HANDHELD, PORTABLE HIGH-POWER PULSED LASER THERAPY DEVICE FOR VETERINARY TREATMENTS. THE LASER ACTION IS AIMED AT TISSUE REGENERATION, POST-TRAUMATIC AND ANTI-INFLAMMATORY TREATMENT, AND THE REDUCTION OF BOTH CHRONIC AND ACUTE PAIN.

SPECIFIC PROGRAMS FOR DOGS, CATS, HORSES, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

8W \pm 20%

Compatible wavelength(s)

808nm

980nm and 1064nm (optional)

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Ø 22 mm and Ø 60 mm optics (optional)

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G6035

· 1 Bag


· 1 LASERVET 8.0 device

· 1 Handpiece MP 808nm with Ø 11 mm optic

· 2 Pairs of protective goggles

· 1 Power supply unit

· 1 Operating manual

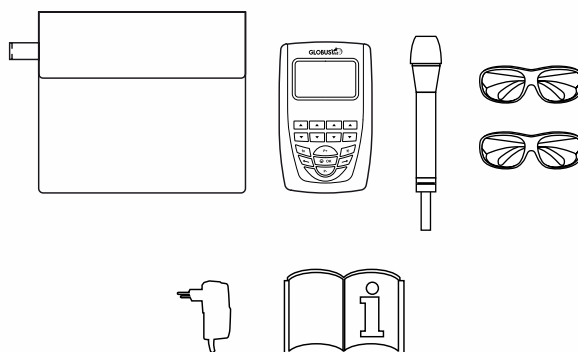

PURCHASE OPTIONS

· Code G6347 Handpiece MP 980nm

· Code G5996 Handpiece 1064nm

· Code G6089 Ø 60 mm optic

· Code G6092 Ø 22 mm optic



STRONG POINTS



LASERVET 3.0

LASER THERAPY

THE TREATMENT WITH LASERVET 3.0 ACCELERATES THE HEALING TIME OF INJURED TISSUE AND IMPROVES CIRCULATION. IT IS IDEAL FOR STIMULATING THE HEALING OF SKIN WOUNDS AND IS PARTICULARLY EFFECTIVE IN THE TREATMENT OF VARIOUS MUSCULOSKELETAL CONDITIONS CHARACTERISED BY INFLAMMATION AND PAIN. LASER THERAPY ELIMINATES PAIN FASTER, WHETHER USED ALONE OR WITH OTHER THERAPEUTIC METHODS.

SPECIFIC PROGRAMS FOR DOGS, CATS, HORSES, SMALL MAMMALS, REPTILES AND BIRDS

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

3W $\pm 20\%$

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

\varnothing 11 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions

100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT

· Code G6192

· 1 Bag


· 1 LASERVET 3.0 device

· 1 Handpiece Black 808nm with \varnothing 11 mm optic

· 2 Pairs of protective goggles

· 1 Power supply unit

· 1 Operating manual

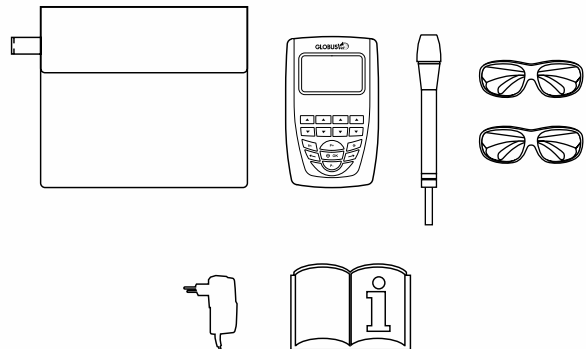

CE



STRONG POINTS



LASERVET 1.0

LASER THERAPY

LASERVET 1.0 IS THE ONLY CLASS 4 LASER IN THE GLOBUSVET CATALOGUE. THE DEVICE CAN REACH 1 WATT OF POWER AND IS FULLY EQUIPPED WITH PRESET PROGRAMS FOR DOGS, CATS AND HORSES. MOREOVER, IT CAN PROGRAM AND STORE NEW PROGRAMS.

SPECIFIC PROGRAMS FOR DOGS, CATS AND HORSES

TECHNICAL SPECIFICATIONS

Display

Size 2.6", backlit

Laser source

Class 4 GaAlAs diode

Peak power

1W \pm 20%

Compatible wavelength(s)

808nm

Emission

Continuous and pulsed

Available optics

Ø 11 mm optic

Frequency

Up to 10000Hz

Power supply

- Mains
- Rechargeable battery

Dimensions


100x160x35 mm

Weight

450 g

STANDARD EQUIPMENT


- Code G5629
- 1 Bag
- 1 LASERVET 1.0 device
- 1 Handpiece White 808nm with Ø 11 mm optic
- 2 Pairs of protective goggles
- 1 Power supply unit
- 1 Operating manual

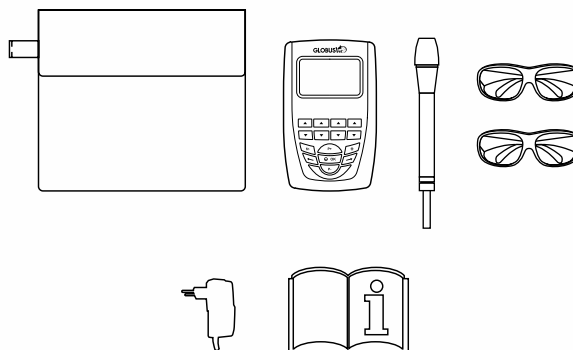

STRONG POINTS



VET ACCESSORIES

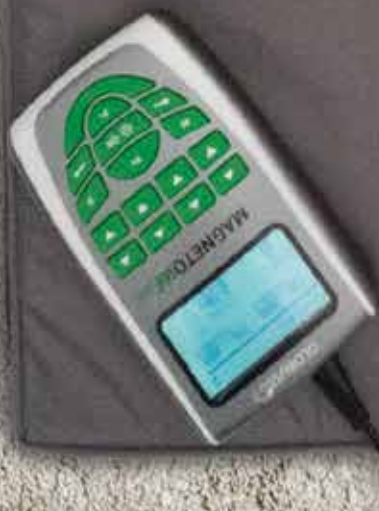

THERAPEUTIC BED FOR MAGNETOTHERAPY

INDICATED FOR **ARTHROSIS** AND **ARTICULAR PAIN**

TECAR THERAPY

G6661	MONOPOLAR CAPACITIVE HANDPIECE
G6689	NEUTRAL FLEXIBLE PLATE 250X200 mm
G6660	MONOPOLAR RESISTIVE HANDPIECE
G6659	BIPOLAR HANDPIECE, Ø 30 mm
G6658	BIPOLAR HANDPIECE, Ø 50 mm
G6672	CAPACITIVE ELECTRODE, Ø 30 mm
G6671	CAPACITIVE ELECTRODE, Ø 50 mm
G6670	CAPACITIVE ELECTRODE, Ø 70 mm
G6669	RESISTIVE ELECTRODE, Ø 30 mm
G6668	RESISTIVE ELECTRODE, Ø 50 mm
G6667	RESISTIVE ELECTRODE, Ø 70 mm
G3368	TECAR CREAM, 250 ml
G3369	TECAR CREAM, 1000 mm
G6664	NEUTRAL PLATE CABLE
G6663	POWER SUPPLY FOR "TECARVET 4000"
G6638	BATTERY PACK FOR "TECARVET 2000"
G1364	BLACK ALUMINIUM CASE FOR "TECARVET 2000"
G3370	CREAM DISPENSER AND CAP
G3371	GREY ALUMINIUM CASE FOR "TECARVET 4000"
G6489	"TECARVET" LINE HOLDER FOR ACCESSORIES
G6488	"TECARVET" LINE HOLDER FOR DEVICES AND ACCESSORIES

MAGNETOTHERAPY

G6589	MAGNETOVET 4000 HORSE
G6619	THERAPEUTIC BLANKET WITH SOLENOIDS
G6530	THERAPEUTIC BLANKET WITHOUT SOLENOIDS
G6620	THERAPEUTIC NECK COVER WITH SOLENOIDS
G6531	THERAPEUTIC NECK COVER WITHOUT SOLENOIDS
G6591	THERAPEUTIC FRONT LEG WRAP WITH SOLENOIDS
G6256	THERAPEUTIC FRONT LEG WRAP WITHOUT SOLENOIDS
G6592	THERAPEUTIC BACK LEG WRAP WITH SOLENOIDS
G6255	THERAPEUTIC BACK LEG WRAP WITHOUT SOLENOIDS
G6632	FLEXIBLE SOLENOID, 300X100 mm, 2 COILS
G6629	FLEXIBLE SOLENOID, 400X150 mm, 4 COILS
G6631	"SOFT 4" MAT, 250X250 mm, 2 CHANNELS, 4 COILS
G6039	THERAPEUTIC BED EXTENSION, 440X700 mm, 4 COILS
G6134	GLOBUSVET SOFT FLEECE BLANKET
G5878	FLEXIBLE SOLENOID COVER KIT, 10 PCS
G6533	CASE FOR MAGNETOVET 4000 HORSE
G6532	STRAP FOR MAGNETOVET 4000 HORSE CASE
G6713	MAGNUM TEST, INDICATOR OF CORRECT EMISSION OF THE MAGNETIC FIELD
G6728	POWER SUPPLY FOR "MAGNETOVET 200"
G6710	POWER SUPPLY FOR "MAGNETOVET 4000"
G6637	BATTERY PACK FOR "MAGNETOVET 200"
G6638	BATTERY PACK FOR "MAGNETOVET 4000"
G6709	Y CABLE

MORE

G5927	GLOBUS BACKPACK, 30x45x15 mm (WxHxD)
G5521	TROLLEY WITH MULTIPLE SPACES TO TRANSPORT THE DEVICES

ELECTROTHERAPY

G0136	MYOTRODE PLATINUM ELECTRODES, 50X50 mm, PACK OF 4 PCS
G0137	MYOTRODE PLATINUM ELECTRODES, 50X90 mm, PACK OF 4 PCS
G0464	MYOTRODE PLUS ELECTRODES, 50X50 mm, PACK OF 4 PCS
G0465	MYOTRODE PLUS ELECTRODES, 50X90 mm, PACK OF 4 PCS
G0451	ELECTRODES, Ø 32 mm, WITH CABLE, PACK OF 4 PCS
G6715	CONDUCTIVE SILICONE ELECTRODES PER METRE
G0463	KIT 6 CONDUCTIVE SILICONE ELECTRODES
G0890	MEDIUM CARBON ELECTRODE FOR IONOPHORESIS, WITH POUCH, 50X50 mm
G0885	LARGE CARBON ELECTRODE FOR IONOPHORESIS, WITH POUCH, 60X85 mm
G6719	POUCH FOR LARGE IONOPHORESIS ELECTRODE, 60X85 mm
G6720	POUCH FOR MEDIUM IONOPHORESIS ELECTRODE, 50X50 mm
G1839	COLOURED CABLES KIT FOR 4-CHANNEL DEVICES, PACK OF 4 PCS
G6734	SPARE CABLE, WHITE, FOR 4-CHANNEL DEVICES
G0439	KIT 2 Y CABLES, UNIVERSAL
G0441	KIT 2 Y CABLES FOR 4-CHANNEL DEVICES
G0484	GLOBUS GEL, 260 ml, PACK OF 10 PCS
G0869	GLOBUS GEL, 260 ml - 1 PC
G4034	GLOBUS GEL, 1000 ml - 1 PC
G3609	GLOBUS GEL, 5000 ml - 1 PC
G6682	MOTOR POINT PEN
G1225	ELASTIC BANDS SIZE "S", 400X80 mm, PACK OF 2 PCS
G1226	ELASTIC BANDS SIZE "M", 800X80 mm, PACK OF 2 PCS
G4249	ADAPTER FOR ACUPUNCTURE, BLACK
G4250	ADAPTER FOR ACUPUNCTURE, RED
G6639	BATTERY PACK FOR "STIMVET 4000"
G6637	BATTERY PACK FOR "STIMVET 200"
G6731	BATTERY CHARGER FOR "STIMVET 200"
G6729	POWER SUPPLY FOR "STIMVET 4000"

ULTRASOUND THERAPY

G6739	HANDPIECE Ø 42 mm, FOR "ULTRASOUNDVET 4000" - SD4 BOARD
G6737	HANDPIECE Ø 42 mm, FOR "ULTRASOUNDVET 4000" - SD5 BOARD
G6741	HANDPIECE Ø 16 mm, FOR "ULTRASOUNDVET 4000" - SD5 BOARD
G6742	HANDPIECE Ø 16 mm, FOR "ULTRASOUNDVET 4000" - SD4 BOARD
G0484	GLOBUS GEL, 260 ml, PACK OF 10 PCS
G0869	GLOBUS GEL, 260 ml - 1 PC
G4034	GLOBUS GEL, 1000 ml - 1 PC
G3609	GLOBUS GEL, 5000 ml - 1 PC
G6638	BATTERY PACK FOR "ULTRASOUNDVET 4000"
G6743	POWER SUPPLY FOR "ULTRASOUNDVET 4000"


LASER THERAPY

G6348	HANDPIECE, HP 980nm
G6347	HANDPIECE, MP 980nm
G6194	HANDPIECE, HP 808nm
G6282	HANDPIECE, MP 808nm
G5996	HANDPIECE, 1064nm
G6208	HANDPIECE BLACK 808nm
G6275	HANDPIECE WHITE 808nm
G6746	PROTECTIVE GOGGLES FOR OPERATOR
G6728	POWER SUPPLY FOR LASERVET 1.0 / 3.0 / 8.0
G6729	POWER SUPPLY FOR LASERVET 15.0
G6637	BATTERY PACK FOR LASERVET 1.0 / 3.0 / 8.0
G6639	BATTERY PACK FOR LASERVET 15.0

G6089	Ø 60 mm OPTIC
G6092	Ø 22 mm OPTIC
G6090	Ø 11 mm OPTIC

The background of the slide is a dark blue gradient. It features an abstract geometric pattern of thin, light blue lines and small dots. These lines and dots form various geometric shapes, including triangles, quadrilaterals, and larger, more complex polygons. Some dots are larger and more prominent than others, and some lines are thicker than others, creating a sense of depth and complexity. The overall effect is a modern, tech-oriented aesthetic.

Key **STRONG POINTS**

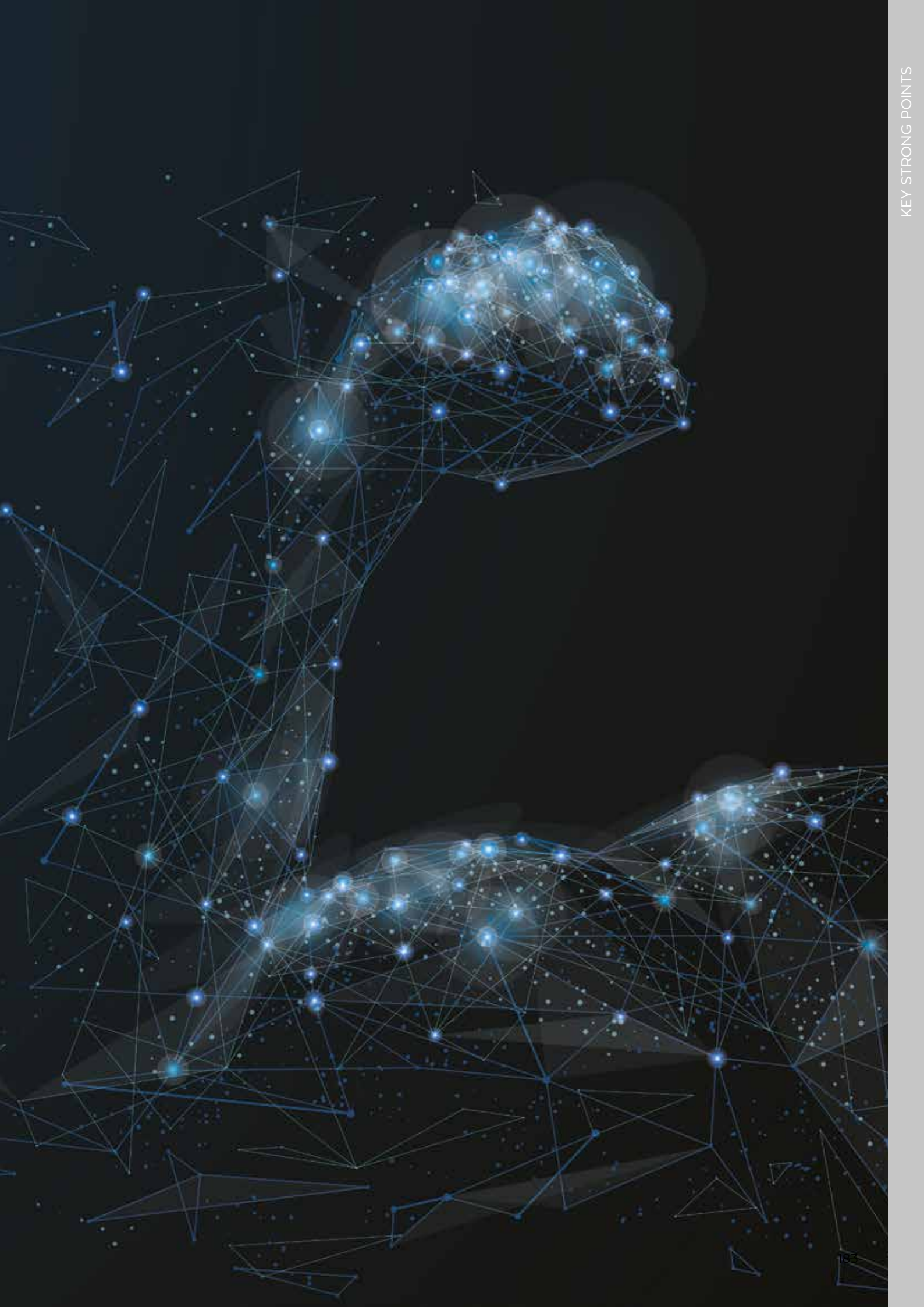

KINEO LINE



ISOTONICO NIT (NO INERTIA)

The NIT Isotonic method introduces an isotonic training which, by eliminating inertia, allows to set a constant load for the entire range of motion. In the training without inertia, the load is constant throughout the movement regardless of speed. Muscle activation is complete, from the beginning to the end of the movement, and the fixed load allows to perform every workout in the best possible way with maximum safety and never losing contact.



ISOTONIC WIT (WITH INERTIA)

In the innovative Kineo system, a robotic system allows to work with the Isotonic method, controlling the load and the work parameters completely and using different loads in the 2 phases of the movement, concentric and eccentric, for a more functional and effective treatment.



ELASTIC

The Elastic method adapts the resistance of the exercise according to your range of motion in a functional way: minimum load in the unfavourable points, higher load in the optimal ones. Setting an elastic with Kineo is simple: define the starting and final load. Your custom elastic will automatically be created.



SMART VISCOUS

The Kineo VISCOUS method is based on the principle of training in the water as a working method in common strengthening exercises. The value of the viscous load is determined by 2 parameters: the set viscosity level and the speed of movement during the exercise. Kineo allows to use 6 levels of viscosity that simulate workout in different fluids, such as water, oil, honey and more, with different levels of muscle commitment.



VARIABLE LOAD

In traditional Isotonic training (weights, barbells, machines), the workload consists of a weight that always remains the same during the exercise. Kineo, instead, allows to set a load that varies and follows a trend (curve) which reflects the real strength capacity of the subject. It also protects and safeguards the joints in the critical work angles, while producing the maximum strength demand and the maximum training effect in the most favourable work angles.



ISOMETRIC

This method consists in producing a muscle contraction, even maximal, without shortening the muscles and moving the affected joint. Kineo can program isometric exercises on any working angle for a more complete strength development. The Biofeedback system also allows to control the strength expressed (in real time) to reach the training Target. To overcome the limitations of the traditional isometric method, Kineo has also introduced the measurement of isometric strength in relation to the time taken to develop it. This method is mainly used in functional rehabilitation because it is safe, adjustable, it does not produce any joint movement and stress, and it improves the recruitment of muscle fibres.



ISOKINETIC

The Isokinetic method is characterised by movement at a constant speed, regardless of the strength applied: this allows the rehabilitation process to be carried out in total safety. With Kineo, the functionalities of the isokinetic method were implemented with isokinetic exercises in open and closed kinetic chain and with differentiated load in the concentric and the eccentric phase.

DIACARE | TECAR BEAUTY LINE



TREATMENT ATLAS

Extensive library of preset programs divided into anatomical areas with illustrative images.



HANDS-FREE TREATMENT

The program facilitates the absorption of the active ingredients of a drug or cream.



MULTI-FREQUENCY

Different working frequencies can be selected to treat the patient accurately at different depths.



TRANSDERMAL DELIVERY

The program facilitates the absorption of the active ingredients of a drug or cream.



AUTOMATIC SCANNING OF FREQUENCIES

During the same treatment, the program can transfer energy evenly at different depths.



ATHERMAL TREATMENTS

Biostimulation of tissues without heat production to treat pathologies that cannot be treated at high temperatures.



PULSED TREATMENTS

Pulsed emission programs with different duty cycles, to keep the energy transferred to the tissues high while maintaining a good temperature control.



FASCIA TOOLS

Special ergonomic tools designed for the mobilisation and the advanced treatment of all soft tissue disorders.



RESISTIVE METHOD

An unshielded steel handpiece is used to work on stronger and deeper tissues such as bones, very tight muscles and ligaments.



CAPACITIVE

A steel handpiece with shielding is used concentrating the ions under the handpiece to a maximum depth of 2/3 cm. The tissues that can be treated with the tecar capacitive handpiece are the most superficial and soft ones, such as muscle tissue, loose connective tissue or the lymphatic system.



BIPOLAR

A single handpiece with two built-in poles.



COLOUR DISPLAY



RECHARGE

Function that allows to buy the device at an advantageous price, recharging the minutes according to the needs.



TOUCH SCREEN DISPLAY



MAINS POWER SUPPLY

REHAB | BEAUTY | SPORT | VET LINE

FOR **MAGNETOTHERAPY | ELECTROTHERAPY AND ELECTROSTIMULATION | ULTRASOUND THERAPY | LASER THERAPY**



TENS

Transcutaneous Electrical Nerve Stimulation (TENS) is performed by means of small electrical pulses to the painful area. The stimulation reaches the nerve endings, inhibiting the transmission of the painful sensation.



PROGRAMMABLE

The user can create some programs with personalised parameters.



NMES

Neuromuscular Electrical Stimulation (NMES) is a current suitable for stimulating the muscles with the aim of improving muscle tone and trophism.



STIM LOCK

To lock the device and carry out only the saved treatment.



MENS

Microcurrent Electrical Nerve Stimulation (MENS) consists of tiny and imperceptible electrical pulses capable of relieving pain, stimulating tissues and promoting the healing of various types of injuries and osteoarticular pathologies.



1+1 FUNCTION

It allows you to use 2 different programs at the same time.



IONO

Ionophoresis is a technique that allows the use of a direct current for the administration of a drug transcutaneously.



2+2 FUNCTION

It allows to use different programs and stimulate 2 parts of the body or 2 people at the same time.



3S

Serial Sequential Stimulation is a new stimulation method which, in the medical/physiotherapy field, allows the muscles involved to be stimulated in a sequential manner.



MY TRAINER

This function, with some wave forms (EMS, TENS, MENS), shows the images of the electrode position directly on the display.



ACTION NOW

It allows to activate muscle contraction through an external control managed by an operator.



AUTO STIM

This function allows to execute an already executed program automatically, setting the intensity value used in the previous program.



RUN TIME

This function allows to modify the time, the frequency and the pulse amplitude during the treatment.



SYNCHRO STIM

This function allows to increase or decrease the stimulation intensity of all the channels simultaneously using a single button.



LAST 10

The stimulator stores (as mobile phones) the last ten executed programs. In this way they are quickly available for the next execution.



WORK TIME

It indicates the total time the device has been used in terms of stimulation.



FAVOURITES

The therapists can store up to 30 programs for each registered user.



RENTAL OPTION

It allows you to lock the device keeping in memory only the number of treatments to be performed. Suitable for renting the device.



MULTI USER

This special function facilitates the use of the stimulator as it allows to store the personal settings of different users on the same device.



RECHARGEABLE BATTERY

The device is supplied with a rechargeable battery and charger to always have the device ready to use wherever you are.



EASY PROGRAM

This function allows to start a magnetotherapy treatment, changing its parameters, with a few simple clicks.



POWER SUPPLY

The device can also operate on the mains.



DOMINO S.R.L.

Via Vittorio Veneto, 52
31013 Codognè - (TV) - ITALY
Phone (+39) 0438 7933
globuscorporation.com

Information for sector
operators only.

The images, features and
data in this brochure are not
binding and are subject to
change by the manufacturer
without prior notice.

For further information
on our product range,
treatments and our videos,
visit the website:

www.globuscorporation.com

GLOBUS

ITALIAN EXCELLENCE

KINEO

MEDICAL

BEAUTY

SPORT

VET

