

PRESSCARE G250

PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

A PRACTICAL AND LIGHTWEIGHT PRESSURE THERAPY DEVICE, DESIGNED TO RELIEVE LEGS, PROMOTE CIRCULATION, AND AID RECOVERY AFTER INTENSE DAYS. THANKS TO ITS COMPACT DESIGN, PRESSCARE G250 IS THE PERFECT COMPANION FOR DAILY WELL-BEING AND RELAXATION.

TECHNICAL SPECIFICATIONS

Outputs

2

Programs

- 2 inflation modes

Pressure

40-160 mmHg (in steps of 10)

Timer

Adjustable from 10 to 30 minutes (in steps of 5)

Power supply

Mains

Dimensions (LxWxH)

172x133x71 mm

Weight

645 gr

STANDARD EQUIPMENT

- 1 PRESSCARE G250 device
- 2 Leg sleeves, sizes S/M/L
- 1 Power supply
- 1 User manual



PRODUCT VARIANTS

- G7951 PRESSCARE G250 - S
- G7952 PRESSCARE G250 - M
- G7953 PRESSCARE G250 - L



KEY STRENGTHS



Compact



5 Chambers



2 Preset programs



13 Intensity levels



Pressure range
40-160 mmHg



Personalized timing



Mains power



GLOBUS
ITALIAN EXCELLENCE

PRESSCARE G250

PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

Rev. En_05_2026

PROGRAMS DESCRIPTION

MODE A - LIGHT ACTIVATION (Warm Up)
Progressive inflation of one chamber at a time for a gentle massage that promotes leg relaxation.
Ideal for: daily comfort, regaining a feeling of lightness between close training sessions, and a revitalizing massage before a competition.

MODE B - DEEP MASSAGE
Sequential inflation of all chambers followed by simultaneous deflation. The compression is more intense and enveloping, designed for a deeper action and a reshaping effect.
Ideal for: a more intense and deep massage.

CHAMBERS INFLATION SEQUENCE




ACCESSORIES DIMENSIONS

LEG SLEEVES: 5-chamber accessories, ambidextrous, with integrated connector, available in three sizes S/M/L.

- A - GROIN-TO-HEEL length
- B - THIGH circumference
- C - ANKLE circumference

SIZE	A (cm)	B (cm)	C (cm)
“S” LEG SLEEVES	75	76	47
“M” LEG SLEEVES	86	80	48
“L” LEG SLEEVES	97	84	49

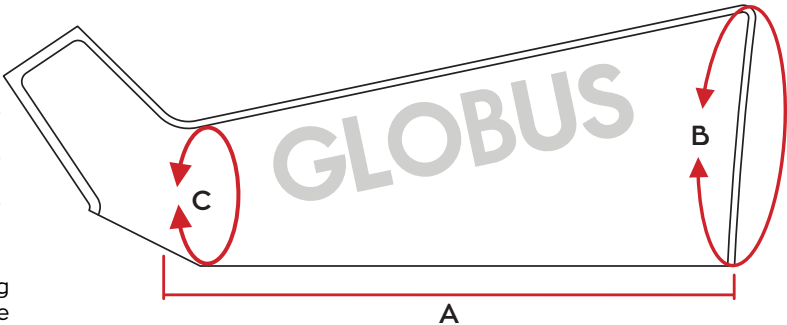

LEG SLEEVE EXTENSION: Extensions for the leg sleeves are available in three sizes (S/M/L) to increase the fit by 7 cm at the ankle and 9 cm at the thigh.

ABDOMINAL BELT: 5-chamber accessory for treatments on the abdomen, hips, and buttocks.

ONE SIZE	A (cm)	B (cm)
	125	35

ABDOMINAL BAND EXTENSION: Thanks to the Velcro, it can be attached to the abdominal band to improve its fit. Dimensions 40x35 cm.

